

ÖSSUR PROSTHETIC SOLUTIONS





Össur Prosthetic Solutions Catalogue 2014

Packed with the latest information on all our products, the 2014 Össur catalogue is a comprehensive aid to product selection. Working closely with prosthetists and technicians, we aim to maintain the highest industry standards, delivering effective clinical outcomes for amputees of every age and ability.

In these pages you can find technical information on product dimensions, weights and configurations, as well as recommendations relating to user weight and activity levels. Current copies of Instructions For Use (IFUs) for all Össur products can be found on our website. Visit the site regularly for news and views on all our products and services. www.ossur.co.uk

ACTIVITY / IMPACT LEVELS

- | | |
|------------|---|
| LOW - | Daily activities involving limited and steady actions at low levels of impact on the prosthesis.
Examples: Community ambulator. Restricted walking with ability or potential to demonstrate varied cadence. Low impact recreational activities - bowling, gardening, shopping. |
| MODERATE - | Daily activities involving typically normal and repetitive actions at moderate levels of impact on the prosthesis.
Examples: Prolonged walking with rapid and varied cadence. Occasional recreational sporting activities - jogging, tennis, hiking, golf. |
| HIGH - | Daily activities involving rigorous and repetitive actions at high levels of impact on the prosthesis.
Examples: Lifting, jumping, climbing. Frequent activities - basketball, volleyball, running, athletics, skiing, football, karate, landscaping, construction. |
| SPORTS - | Activities with repeated extreme level of impact on the foot.
Examples: track and field, sprinting and long-distance running. |

Table of Contents

Introduction	3
--------------------	---

Bionic Technology by Össur

PROPRIO FOOT™	7
RHEO® KNEE3	12
POWER KNEE™	16
SYMBIONIC® LEG	20

Interfaces

Iceross Seal-In® (TT)	35
Iceross Seal-In® V	36
Iceross Seal-In® X5	38
Iceross Dermo Seal-In® HSM	40
Iceross Synergy® Cushion	42
Iceross Dermo® Cushion	43
Iceross Comfort® Cushion	44
Iceross® Activa with TibiaGuard™	45
Iceross Synergy® Locking	47
Iceross Dermo® Locking	48
Iceross Comfort® Locking	50
Iceross® Sport	51
Iceross® Original Locking	52
Iceross Seal-In® X5 (TF)	53
Iceross® Transfemoral Seal-In®	54
Iceross® Transfemoral Locking	58
Iceross® Upper-X	59
Iceross® Junior	60
Iceross® Sleeve	61
Iceflex® Endurance	62
Iceflex® Balance	63
Distal Cup™	64
Iceross® Pads	64
Clean & Simple	65

Unity®

Unity® for PROPRIO FOOT®	69
Unity® for Re-Flex®	72
Unity® for LP Rotate™	76
Unity® for Vari-Flex®	80
Unity® for LP Vari-Flex®	83

Sockets

Icecast® Anatomy	89
Icecast® Compact	90
Modular Socket System	91
Pressure Pads	93

Locks

Icelock® 100 Series	97
Icelock® 200 Series	99
Icelock® 300 Series	103
Icelock® 500 Series	105
Icelock® 600 Series	106
Icelock® 700 Series	112

Adapters

Standard Adapters, Rated to 100kg	117
Standard Adapters, Rated to 166kg	120
Junior Adapters, Rated to 45kg	123
Total Shock™	124

Knees

Mauch® Knee	127
Mauch® Knee Plus	129
Balance® Knee	131
Locking Knee	133
Total Knee® 1900	134
Total Knee® 2000	137
Total Knee® 2100	140
Total Knee® Junior	143

Cosmetic Finishing

Cosmetic Finishing	146
--------------------------	-----

Feet


Flex-Foot® Balance	150
Balance™ Foot J	152
Flex-Foot® Assure	154
Talux® ST	158
Vari-Flex®	164
LP Vari-Flex®	170
Vari-Flex® Junior	174
Vari-Flex® XC	177
Vari-Flex® XC Rotate™	182
Re-Flex Rotate™	185
LP Rotate™	191
Vari-Flex® Modular	196
Re-Flex Shock™	202
Elation®	207
Flex-Symes™	211
Flex-Foot® Junior	217
Chopart™	220
Flex-Run™ with Nike Sole	222
Flex-Run™ Junior	225
Flex-Sprint™	228
Cheetah®	230
Cheetah® Xtreme with Nike Spike Pad	233
Cheetah® Xtend with Nike Spike Pad	236
Cheetah® Xplore Junior	239
Cheetah® Junior	242

Rehabilitation

Össur Rigid Dressing	247
Iceross Post-Op	248
Femurett	249

Additional Information

Torque Values	251
Terms of Sales/Policies	253
Warranty	254
Order form Special made	256
Returns/Complaints Form	257

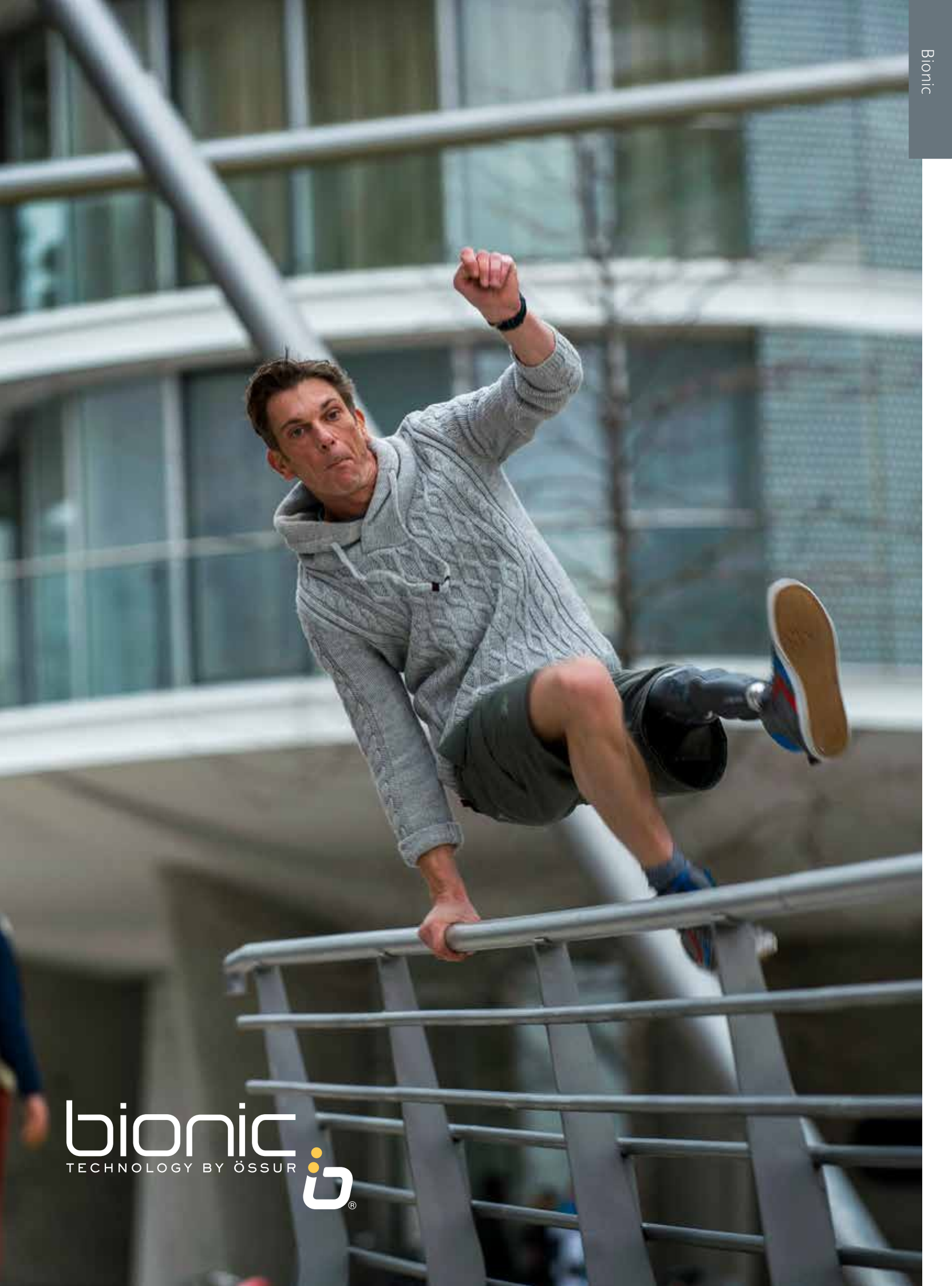

Bionic Technology by Össur

Introduction


Bionics = bi(o) + (electr)onics: the application of biological principles to the study and design of engineering systems, especially electronic systems.

The American heritage dictionary of the English language, 2000

Bionic Technology by Össur is beginning to restore the true anatomical function displaced by amputation. Our intelligent structures contain sensors, artificial intelligence and actuator technology which can respond in a human-like way to counteract loads, reduce vibration, change shape and enhance fluid and natural motion.

 SENSE	 THINK	 ACT
Sophisticated sensor technology mimics the body's mechanoreceptors, providing artificial proprioception (a sense of where the limb is in space) by tracking factors that include load and angular displacement up to a rate of 1600 times per second.	Patented forms of artificial intelligence (AI) process information from the sensors and instigate the most appropriate response for the next step.	The AI transmits a constant stream of signals which instruct high-precision actuator technology to act and deliver optimal function. The addition of powered motion represents a significant step in human energy conservation, and the ability to deliver the required response. A motor and battery generates appropriate prosthetic function, providing the precise amount of power necessary.

The human versus the bionic model

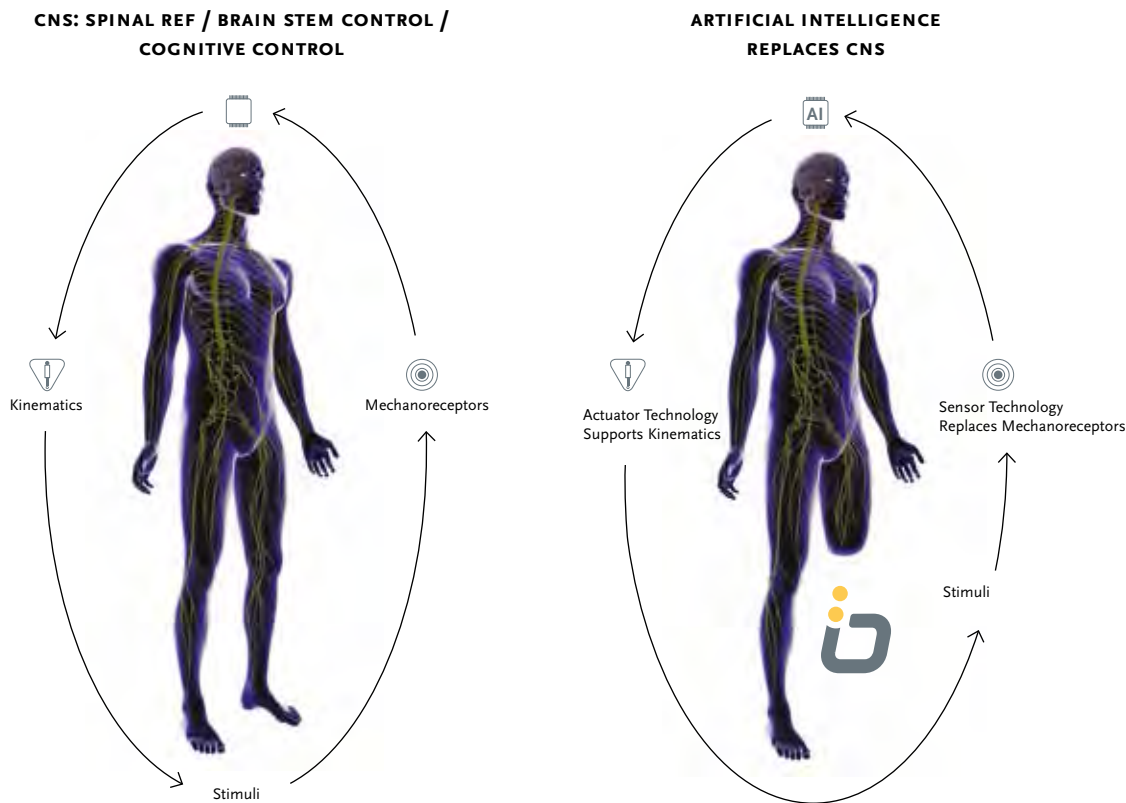





Introduction

The PROPRIO FOOT provides unparalleled physiological benefits for transtibial amputees. Automatic ankle flexion means function is as close as you can get today to the human foot. This device thinks for itself, responding beautifully to changing terrain and transforming the approach to stairs and slopes. Overall, the effect is a more balanced and symmetric gait, preserving the rest of the body.

 SENSE	 THINK	 ACT
Accelerometer technology measures real-time motion at a rate of 1600 times per second. Tracking the path of the ankle through space, it times gait characteristics and events, including heel strike and toe off. For each stride, the device constructs it's path by continuously analysing horizontal and vertical movement, 'tracing' the foot as it moves through space. This trace varies according to the terrain so gait pattern recognition algorithms in the foot can detect and identify when a user is walking on flat or sloped surfaces and up or down stairs, as well as when stationary or in a 'relaxed' position (seated).	As fast as it is measured, the data is fed into the 'brain' of the device, a specialised form of artificial intelligence called Terrain Logic™. Together with the motion analysis sensors Terrain Logic helps to rebuild the sense of proprioception lost after amputation and provide the most appropriate response for the next step. As toe-off is detected, for example, it directs a toe lift in order to clear the ground. Having sensed an inclined surface, it directs ankle flexion proportional to the slope. Learning from the user's previous stride, ankle movement parameters are adapted for the next step.	A control board receives a constant stream of signals from the artificial intelligence. The controls then instruct a high precision linear actuator to deliver optimal function – dorsiflexion (toe lift) in swing phase; ankle angle adjustment on varying terrain, in chair exit mode and relaxed mode; and heel height adjustment when changing shoes.

PROPRIO FOOT bionic model



PROPRIO FOOT®

Information



USER INFORMATION

Maximum patient weight	125kg
Minimum patient weight	45kg
Impact level	low and moderate


FOOT INFORMATION

Categories	1-8
Heel height adjustability up to	50mm
Build height	168 - 176mm
Sizes	22 - 30cm
Weight of PROPRIO FOOT (Size 27, Cat 5)	1468g

Category Selection Guide

	Weight (kg)							
	45-52	53-59	60-68	69-77	78-88	89-100	101-116	117-125
Low impact level	1	1	2	3	4	5	6	7
Moderate impact level	1	2	3	4	5	6	7	8


Dimensions

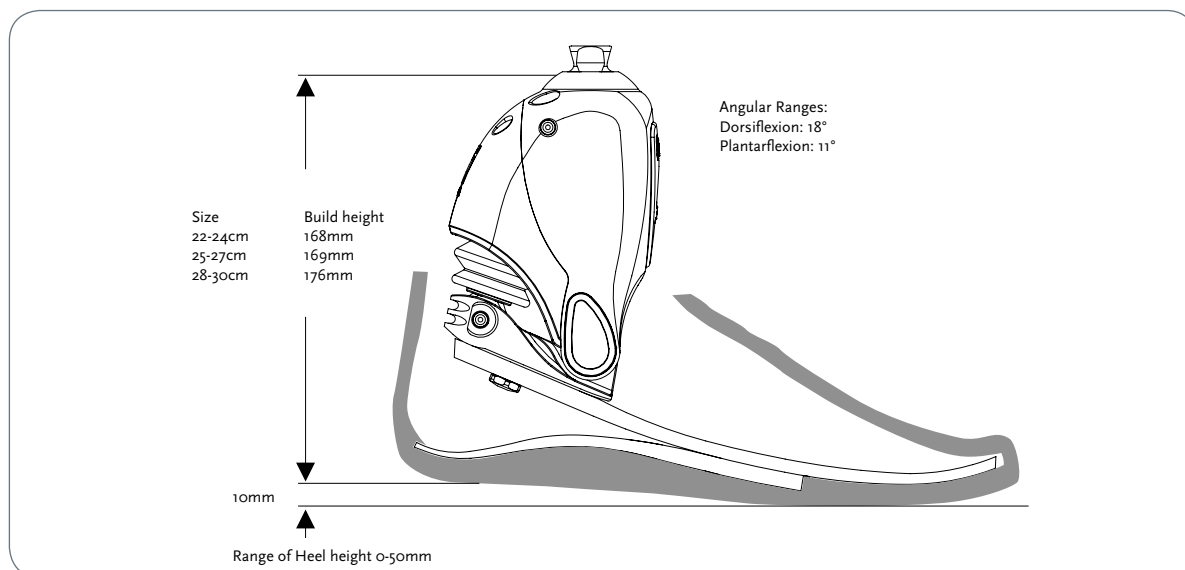

PROPRIO FOOT®

Kits

Part#

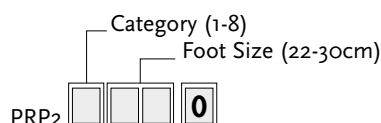

Includes:

- PROPRIO FOOT Ankle
- Foot module
- PROPRIO FOOT Battery dummy
- PROPRIO FOOT Power Supply
- Heel Wedge Kit
- Terms of Warranty
- Quick Reference Cards
- Technical Manual
- Fitting Guide
- Instructions for Use
- Troubleshooting Guide



Components

Picture	Part#	Description
1	PRXoo431	Battery Dummy
2	PRXoo421	Power supply Kit




Part#	Description	Category	Size
CM15oo11	Heel Wedge Kit	1-8	22-30



Available with Unity (see page 63)

Cosmetic Finishing

Part#

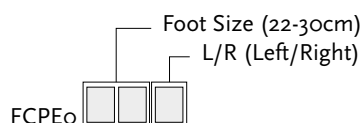

Please specify when ordering brown covers. Add "BR" to the part#.

Includes


- Foot Cover
- Flex-Foot EVO Sock
- Foot Cover Attachment

Part#	Description	Size
FCE63oo6	Flex-Foot EVO Sock	22-25cm
FCE63oo7	Flex-Foot EVO Sock	26-30cm



PROPRIO FOOT®

Diagram

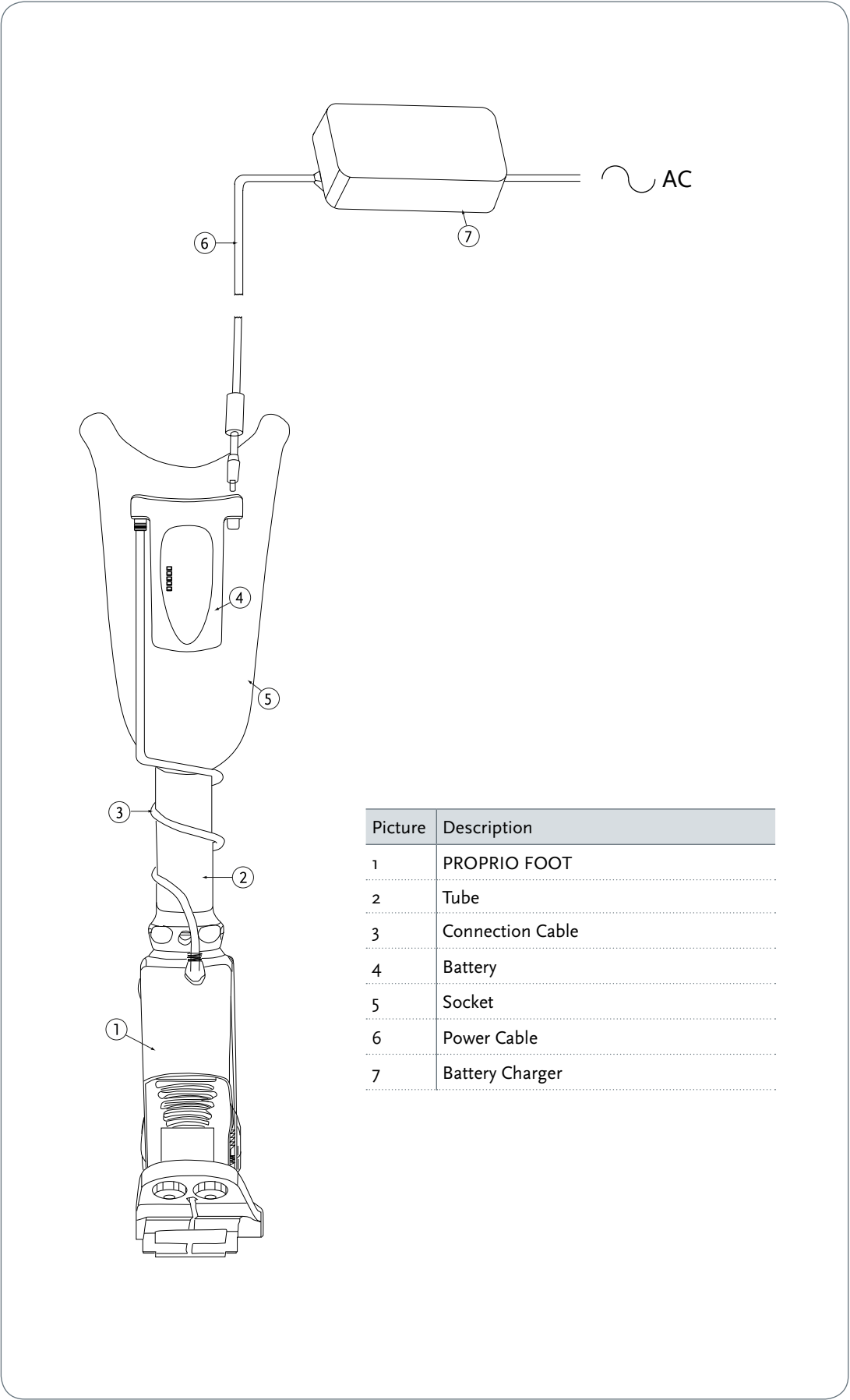



Introduction

The RHEO KNEE 3 provides the most natural knee function among all microprocessor knees because it continuously adapts to the user and the environment while providing excellent stability and safety.

Key Features

Improved Actuator design

- Unmatched response to gait and speed changes
- Complete and constant responsiveness to different and changing circumstances
- Increased maximum torque on stairs providing improved stability
- Easy, natural swing motion

Improved extension assist

- Faster and more natural swing extension
- Improved toe clearance translates to less risk of falling

Upgraded Software Technology

- Advanced control features continuously monitor and adapt resistance levels

Kinematic sensor

- Added stability in non-gait movements
- Accurate stand-to-walk transitions

Expanded User Profile

- User weight limit increased to 136kg (300lbs)
- Improved battery life with up to 3 days use per battery charge



RHEO® KNEE3

Information



USER INFORMATION

Maximum patient weight	136kg
Impact level	low and moderate


KNEE INFORMATION

Net weight	1610g
Flexion	120°
Build Height	236mm

Dimensions

RHEO KNEE 3

RHEO KNEE 3 WITH THREADED TOP ADAPTER (A - 145310)

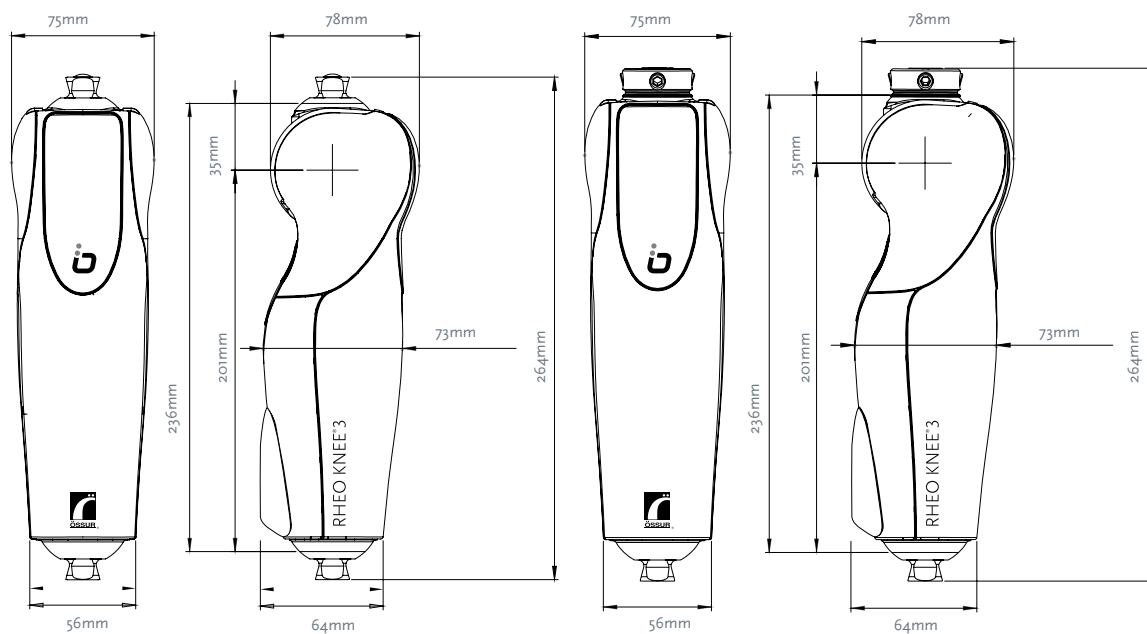

RHEO® KNEE3

Kits

RHEO KNEE 3 KIT

Part#	Description
RKN130002	2 Years Warranty (No Service Check)
RKN130003	3 Years Warranty (No Service Check)
RKN130005	5 Years Warranty (Incl. Service Check At 40 Months)

Includes:

- RHEO KNEE 3
- RHEO KNEE 3 Charger
- Warranty Certificate
- User Manual
- Technical Manual
- Quick Reference Card



RHEO KNEE 3 KIT WITH THREADED TOP ADAPTER:

RHEO KNEE 3 with threaded top adapter is the solution for knee disarticulation amputees and those with long residual limbs. To order RHEO KNEE 3 with threaded top adapter, order the RHEO KNEE 3 Kit and add Össur's Female Pyramid Insert to the order.

Part#	Description	Height	Weight
A-145310	Female Pyramid Insert / Threaded Top Adapter	9mm	33g



Accessories

Part#	Description
RKNLPC401	RHEOLOGIC Software Version 1.2
RKNLPCWB1	RHEOLOGIC Workbench Software Version 1.2

Part#	Description
BIC00110	Bionic Protector for RHEO KNEE



Components

Part#	Description
PRX00421	Power Supply w/Plug



Extended Warranty

Part#	Description
RKN3GE1YR	Additional 1 Year Warranty
RKN3GE2YR	Additional 2 Year Warranty
RKN3GE3YR	Additional 3 Year Warranty
RKN3SC20M	Voluntary 20 Month Service Check
RKN3SC40M	Mandatory 40 Month Service Check*

* Maximum total warranty period is 5 years, must be purchased within 20 months after original purchase date.

** If total warranty time acquired adds up to 4 or 5 years, then 40 month service check will be mandatory.



POWER KNEE®

Introduction

The world's first active prosthesis for above-knee amputees, offering unprecedented function and performance. The POWER KNEE:

- provides the power to walk further and for longer
- assists with the upward motion required for stairs and inclines
- learns and responds to unique patterns of gait for natural, efficient motion

The POWER KNEE works as an integrated extension of its user, restoring the power of lost muscles and the symmetry of movement so vital for those with limb loss. This highly intelligent, second-generation product is now lighter, smaller and more efficient than before.

 SENSE	 THINK	 ACT
Enhanced sensor technology now located solely within the knee combines accelerometers, gyroscopes, a torque sensor and a load cell, all of which work together to monitor the position, velocity and movement of the knee. They also measure with great accuracy the forces and load applied to the knee at any given moment. The sum of this data is fed through in a constant stream to the onboard computer within the knee, which, by utilising the power of artificial intelligence, can then work out exactly what type of activity the user is engaged in and successfully predict the next intended movement.	Artificial intelligence (AI) located in the central microprocessor acts as the 'brain' of the POWER KNEE, replacing the function that the body's central nervous system would normally play in processing information about changes occurring inside and outside the limb if it was still intact. AI allows the knee to think for itself, to automatically process the data gleaned from the sensors and to anticipate and initiate the required function. The AI in the POWER KNEE also assesses gait symmetry in order to stimulate highly efficient gait patterns, which help to minimise strain on the back, hips and sound side.	The electromechanical actuator in the POWER KNEE provides the power for active motion. As fast as information about the user's walking pattern is being processed, the microprocessor transmits a stream of signals which instruct the actuator on precisely how to act and the exact amount of power required. The actuator delivers a range of responses, including active lifting power, to lift the prosthetic leg up as the user ascends stairs or slopes, resistance for downhill slopes, and gentle forward propulsion for improved gait dynamics on level ground. Additionally, an integrated torque meter helps recreate the natural feeling of bounce when the leg comes into contact with the ground, absorbing shock and helping to protect back and hips.

POWER KNEE Bionic Model



POWER KNEE®

Information



USER INFORMATION


Amputation Level	Transfemoral
Maximum patient weight	125kg
Impact level	low to moderate

KNEE INFORMATION

Net weight	2700g
Battery weight	490g
Flexion	120°
Build Height	275mm
Minimum Clearance	380mm

- Aluminium Frame Construction

Dimensions

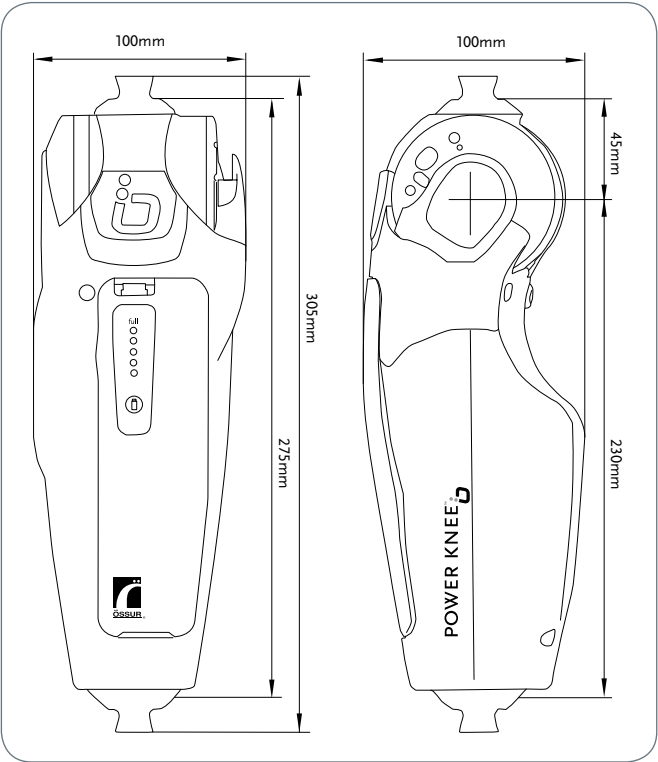

POWER KNEE®

Kits

POWER KNEE KIT

Part#	Description
PKN120003	POWER KNEE KIT EU


Includes:

- POWER KNEE
- 2 x Battery Packs
- Battery Charger
- Power Supply
- Warranty Certificate
- User Manual & Information DVD
- Technical Manual
- Bluetooth Dongle



Components

Picture	Part#	Description
1	PK220000	Battery Charger
2	PK230000	Battery Pack

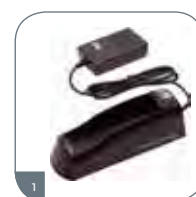




SYMBIONIC® LEG

Introduction

The SYMBIONIC LEG is the first and only available complete bionic leg. Combining a powered ankle and an adaptive microprocessor knee joint, the SYMBIONIC LEG provides unmatched benefits for transfemoral amputees. An active toe lift in swing clears the ground and minimizes the risk of tripping. Should a stumble occur, an advanced stumble recovery feature protects the user from a fall by providing instant stance support. An automatic terrain adaptation adjusts the ankle angle according to the gradient of the slope, generating better ground contact and optimal knee response when ascending or descending ramps. Using the heel height adjustment, the user can change shoes or go barefoot without compromising the alignment of the prosthesis

 SENSE	 THINK	 ACT
The SYMBIONIC LEG uses advanced sensor technology to gather information on terrain, the trajectory of the prosthesis through space and the load exerted by the user. This allows the device to analyze vertical and horizontal movements and gait events in order to trace the prosthesis during use. This information varies on different terrain allowing the device to detect which activity the user is engaging in and adjusting the ankle position and knee resistance accordingly in order to provide a natural motion that adapts to different walking speeds and terrain.	Information is collected from the sensors and analyzed in the SYMBIONIC LEG “brain” at the rate of 1000 times per second. The artificial intelligence uses the sensor data to determine which gait event the user is in and what type of motion to execute next. The leg adapts to the individual’s walking style, providing the most appropriate response for the next step. The leg immediately adapts to change in walking speed and increased activity. It also adapts to changing terrain. Having sensed a decline, the knee resistance and ankle flexion are adjusted proportionally to the slope.	The SYMBIONIC LEG provides a smooth and fast response whether it is in fast walking, moving in confined spaces or descending a slope. The two actuators in the knee module and the ankle module receive instructions from the “brain” on what motion to implement. As toe-off is detected, for example, the motor is instructed to lift the toe in order to clear the ground and to reduce the knee resistance to accommodate for a more natural swing phase. Similarly the motor adjusts the ankle angle and knee resistance on varying terrain and allows heel height adjustments when changing shoes. In case of a stumble, a stumble recovery function in the actuator is immediately activated in order to prevent a fall.

POWER KNEE Bionic Model

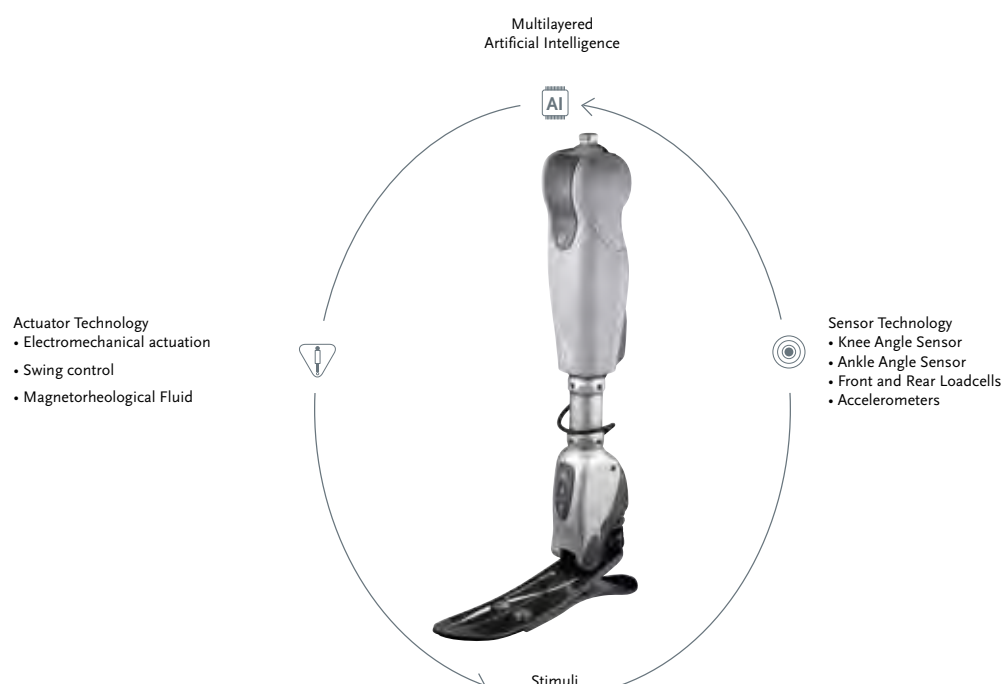


SYMBIONIC® LEG

Information



USER INFORMATION	
Amputation Level	Transfemoral
Maximum patient weight	125kg
Impact level	low to moderate
KNEE INFORMATION	
Categories:	Size 22-24; Cat. 1-5 Size 25-30; Cat. 1-8
Size:	22-30
Weight (Size 27, Cat 5):	3000g with Foot Cover
Heel Height Adjustability:	Up to 50mm
Build Height:	405mm
Knee Flexion:	120°

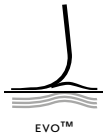


Category Selection Guide

	Weight (kg)							
	45-52	53-59	60-68	69-77	78-88	89-100	101-116	117-125
Low impact level	1	1	2	3	4	5	6	7
Moderate impact level	1	2	3	4	5	6	7	8

SYMBIONIC® LEG


Dimensions

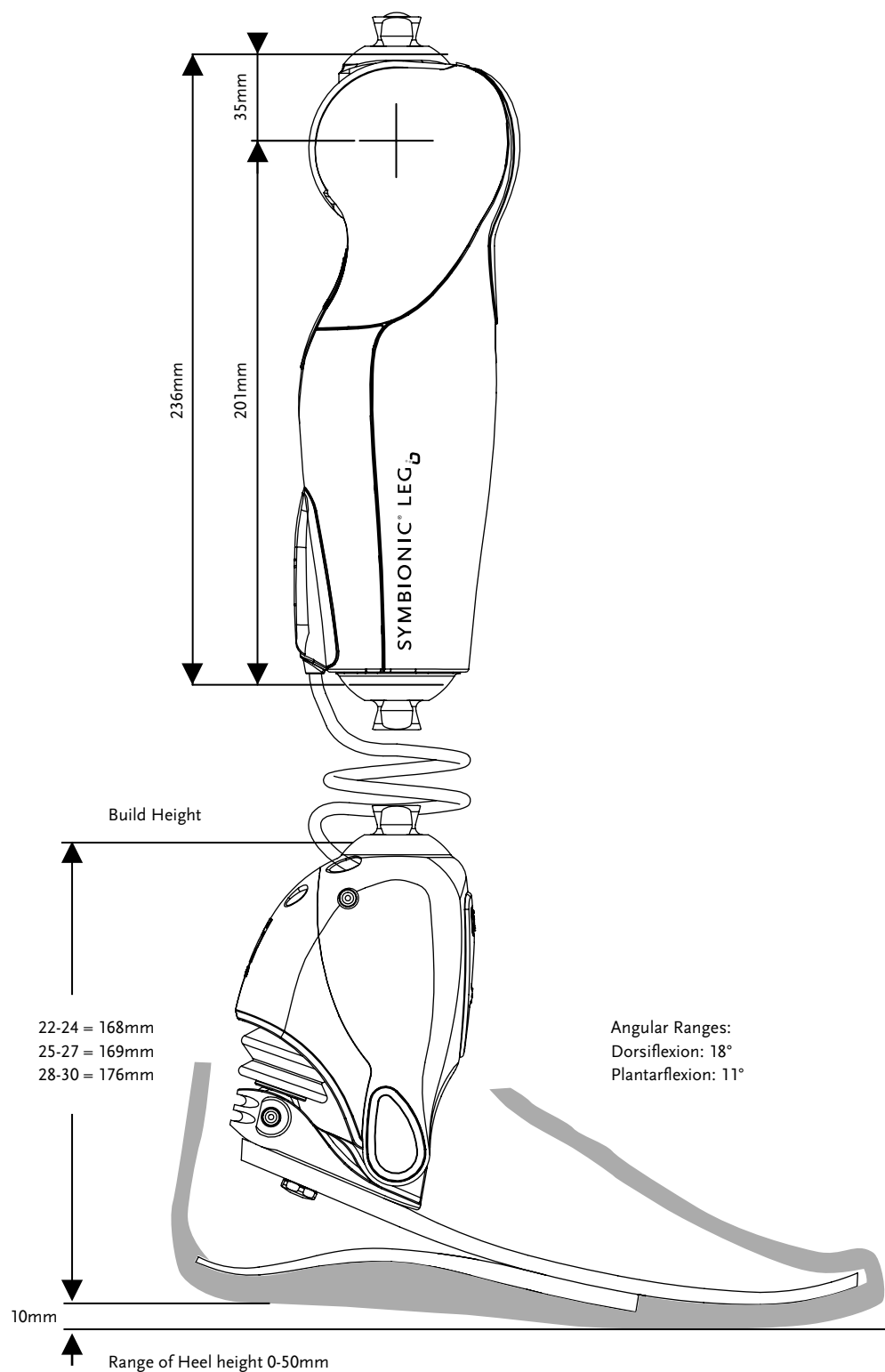

SYMBIONIC® LEG

Kits

Part#

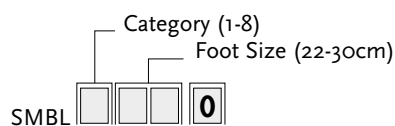

Includes:

- SYMBIONIC LEG Ankle
- Foot module
- Power Supply
- Heel Wedge Kit
- Quick Reference Cards
- Technical Manual
- Instructions for Use
- Terms of Warranty



Components

Part#	Description
SBLGE1Y	SYMBIONIC LEG 1 year guarantee extension **
SBLSC4oM	SYMBIONIC LEG 40 month periodic service *, including new LP Vari-Flex foot

* If more than 1 year guarantee extension is ordered, the 40 month periodic service becomes mandatory.

** Maximum guarantee extension is 3 years to a total of 5.

SYMBIONIC LEG KIT WITH THREADED TOP ADAPTER:

SYMBIONIC LEG with threaded top adapter is the solution for knee disarticulation amputees and those with long residual limbs. To order SYMBIONIC LEG with threaded top adapter, order the SYMBIONIC LEG Kit and add Össur's Female Pyramid Insert to the order.




Part#	Description	Height	Build Height	Weight
A-145310	Female Pyramid Insert / Threaded Top Adapter	9mm	9 mm	33g
A-554700	Male-Female Double Adapter	40mm	24 mm	153g




Components

Part#	Description	
RKNLPC401	RHEOLOGIC 4.0.1 Software	
RKNLPCWB1	RHEOLOGIC Workbench Software	
Picture	Part#	Description
1	PRX00421	Power Supply Kit

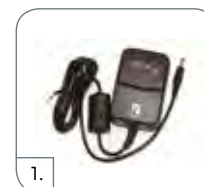


Part#	Description	Category	Size
CM150011	Heel Wedge Kit	1-8	22-30

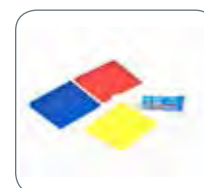

SYMBIONIC® LEG

Cosmetic Finishing

Part#

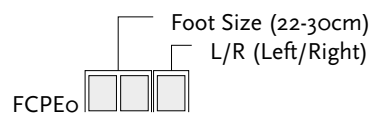

Please specify when ordering brown covers. Add "BR" to the part#.

Includes

- Foot Cover
- Flex-Foot EVO Sock

Part#	Description	Size
FCE63006	Flex-Foot EVO Sock	22-25cm
FCE63007	Flex-Foot EVO Sock	26-30cm





Interfaces

Introduction

Iceross stands for Icelandic Roll-On Silicone Socket. It was developed by Össur Kristinsson, a prosthetist and prosthetic user, as an alternative suspension system to more conventional methods. It functions as the interface between the skin and the inner socket wall to protect the residual limb and to provide greater comfort to transtibial and transfemoral amputees while wearing their prostheses. Iceross is available in a range of designs to suit different residual limb conditions as well as activity and impact levels.

Features and Benefits

COMFORT - all Iceross Interface systems are made exclusively from medical grade silicone. Silicone has a proven history of use in the medical field due to its unique material properties. With Interfaces, the properties provide benefits such as cushioning for bony prominences and sensitive areas, high elasticity to enable total contact and ease of conformity to most residual limb shapes, and high tear strength for increased durability.

PROTECTION - Iceross has been designed to function in a way similar to a compression stocking, i.e. by improving blood circulation while controlling swelling. It also protects skin and newly formed or sensitive scar tissue by reducing potentially damaging shear forces.

SUSPENSION - the Roll-On method ensures that Iceross grips the skin and displaces it distally as it relaxes, thereby stabilising the soft tissues of the limb. The patented matrix and size-specific distal attachment work together in Iceross locking liners to:

- Evenly distribute forces over the entire surface of the residual limb during swing phase.
- Minimise distal end pressures.
- Resist axial loading during the swing phase of gait.



Contents

Iceross Seal-In® (TT)	35
Iceross Seal-In® V	36
Iceross Seal-In® X5	38
Iceross Dermo Seal-In® HSM	40
Iceross Synergy® Cushion	42
Iceross Dermo® Cushion	43
Iceross Comfort® Cushion	44
Iceross® Activa with TibiaGuard™	45
Iceross Synergy® Locking	47
Iceross Dermo® Locking	48
Iceross Comfort® Locking	50
Iceross® Sport	51
Iceross® Original Locking	52
Iceross Seal-In® X5 (TF)	53
Iceross® Transfemoral Seal-In®	54
Iceross® Transfemoral Locking	58
Iceross® Upper-X	59
Iceross® Junior	60
Iceross® Sleeve	61
Iceflex® Endurance	62
Iceflex® Balance	63
Distal Cup™	64
Iceross® Pads	64
Clean & Simple	65

Interfaces

Product Overview

Transtibial Solutions

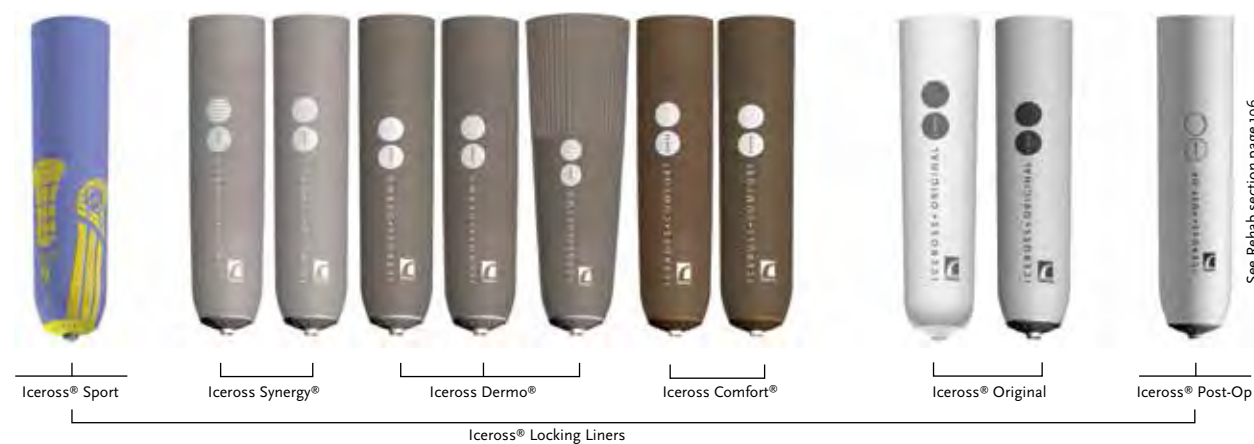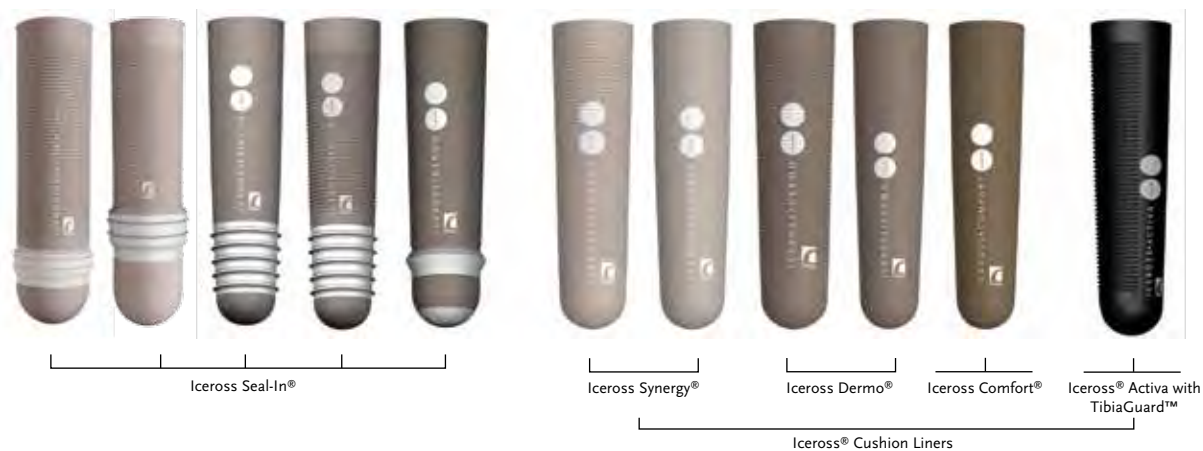


Transfemoral Solutions



Upper Extremity Solution Junior Solution



Sleeves

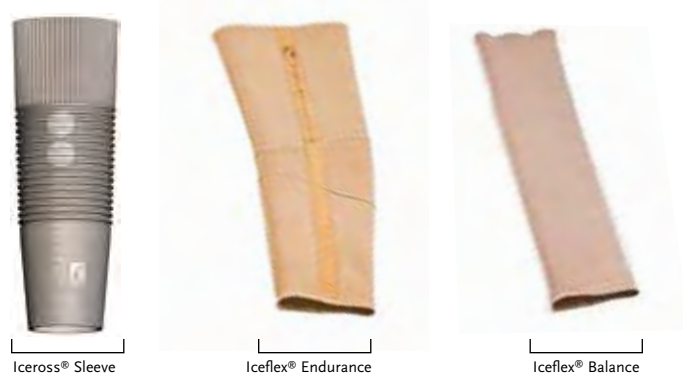

Product Characteristics

LOCKING

Available with either a locking or lanyard attachment, a size specific umbrella and stabilising matrix. Optimal function of both the umbrella and matrix is dependent on correct size selection of the Iceross liner. The system is designed to distribute forces evenly over the entire surface of the residual limb during swing phase, to minimise distal end pressure and to resist axial loading during the swing phase of gait.

**CUSHION**

A soft, extra thick distal pad improves user comfort and skin protection. With up to 16mm of soft silicone, the distal pad protects sensitive distal ends and conforms to irregular distal shapes.

**WAVE**

The Wave feature makes the liner easier to flex at the knee, allowing it to adapt more to the user's movements. And less resistance means that less energy is required to flex the knee. It is designed so that the material stretches effectively at the front and folds in a controlled way behind the knee.

**SEAL-IN®**

The Iceross Seal-In liner optimises suspension while offering ultimate user convenience for daily use.

- Transtibial hypobaric suspension without a knee sleeve
- Improved rotation control and ease of donning

**TIBIAGUARD™**

An integrated softer silicone pad provides cushioning for the tibial crest while a thinner stabilising posterior layer adds control and proprioception. The well known wave feature covers the entire anterior part of this liner adding flexibility at the knee. Optimal cushioning in combination with utmost flexibility is given by these two types of silicone in a non-uniform profile.

**STABILISING MATRIX**

Reinforced fabric integrated into the distal end of the Iceross liners minimises longitudinal stretching. Transfemoral and Upper-X liners have a full-length matrix cover.

- Minimises longitudinal stretching and pistoning
- Protects sensitive distal tissue
- Allows radial stretch

**ACTIVE SKIN CARE**

Active Skin Care ingredients incorporated into the silicone migrate to the outer layers of the skin helping to maintain a healthy, strong skin.

- Aids in controlling dermatological problems
- Vaseline® helps with skin hydration
- Aloe Vera helps to reduce skin irritation

**MENTHOL**

Menthol offers a pleasant, fresh smell and has a cooling and soothing effect on the skin

**DOUBLE DUROMETER**

A stiffer outer layer of silicone provides outstanding stability, while a softer inner layer nurtures the skin and provides shock absorption and comfort.

**SILKEN INNER SURFACE**

Microengineered surface for superior skin conformance.

**FABRIC COVER**

For durability and easier donning and doffing.

**CONICAL**

This latest feature is designed to provide a much more comfortable and intimate liner fit for transtibial amputees with conical-shaped residual limbs. The conical shape incorporates a series of vertical waves intended to reduce pressure and help prevent the proximal edge from rolling down.




Transtibial Selection

Careful selection of the Iceross liner is required to produce maximum benefit to the user. When selecting the appropriate transtibial Iceross liner the following must be considered:


1. Size of the residual limb
2. Matrix length
3. Activity/Impact level of the user
4. Use of fitting aids

1. SIZE

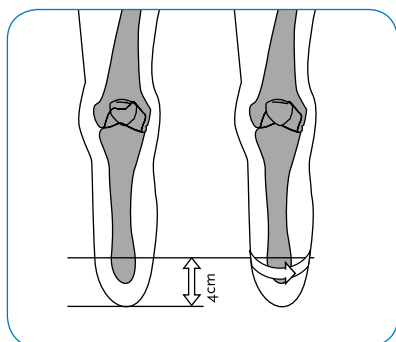

With the tissue of the residual limb hanging down, measure 4cm proximal to the distal end of the limb and mark this point on the skin.

Measure the circumference of the limb at this point.

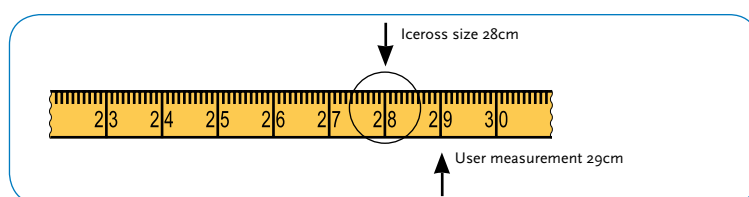

The measurement taken corresponds to the size of Iceross to be used. Select the liner size matching your measurement or, if between sizes select the next size below. Do not select a size larger than your measurement or more than one size below the measurement. e.g. if the circumference of the limb measures 29cm, select Iceross size 28.

Note locking liners:

The incorporated size-specific umbrella works in conjunction with the matrix. Optimal function of both the umbrella and matrix is dependent on correct size selection of the Iceross liner. Incorrect sizing may cause discomfort to the user, perspiration and reduced control of the prosthesis.

Note transtibial and transfemoral HSM Seal-In® liners:

For Iceross Seal-In liners, it is essential that the Hypobaric Sealing Membrane (HSM™) sits correctly against the liner around its full circumference. Once you have selected the liner according to our instructions, assess the membrane according to the pictures below, visually and manually. Picture 2 is the optimal situation. In some cases, the next size down may need to be selected to achieve this fitting (this is unique to Seal-In liners and should not be applied to other Iceross models). Please note that this does not apply for the Seal-In X5 liners.



1. Liner is too small

2. Liner fits correctly


3. Liner is too big

The sizing recommendations are intended as a guide. There may be occasions when alternative sizes need to be selected.

Transtibial Selection

2. MATRIX LENGTH

Most Iceross liners are available with 5cm and 10cm and custom matrix options (2cm-14cm). The length of the matrix within the Iceross liner should be considered on the residual limb.

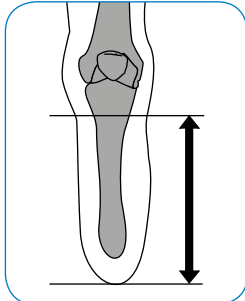

On transtibial applications the tibia end should be contained within the stabilising matrix. However, to prevent restriction of the knee during flexion, the matrix should not extend above the level of the fibula head

3. ACTIVITY AND IMPACT LEVELS

The activity and impact levels of the user are important considerations when selecting a liner. Activity levels relate to function whereas impact levels relate to the amount of loading placed through the Iceross liner, which may be associated with the activity level or the weight of the user. Assess the activity and impact level of the amputee according to the statements on page 3.

4. FITTING AIDS

The profile of the Iceross is cylindrical with a rounded distal end. This design accommodates most residual limb shapes. There are occasions, however, when a patient will exhibit an unusual or irregular residual limb shape that does not allow a total contact fit between the skin and the Iceross. The fitting aids are designed to work in conjunction with the Iceross to change the profile of unusual or irregular residual limb shapes to match that of the Iceross. This process is known as profile matching.



Assess the residual limb shape:

- If the shape is cylindrical, then no fitting aids are required.
- If the residual limb presents a square or flattened distal end, then an Iceross Pad should be installed distally to ensure profile matching between the end of the residual limb and the rounded end of the Iceross liner.
- If the residual limb is conical in shape, then a Distal Cup could be used in conjunction with the Iceross to help ensure correct profile matching. The Distal Cup will assist in providing a more cylindrical shape prior to fitting the Iceross liner.
- Iceross liner size is selected from a measurement taken over the Distal Cup on the residual limb.
- If the residual limb presents with a bulbous distal end (ankle disarticulation, Symes, Lisfranc etc), Iceross use is contraindicated and you should consult your local Össur representative to find an appropriate solution.


Now you can select the appropriate Iceross product for the user.

Transfemoral Selection

The transfemoral residual limb offers slightly different challenges to the prosthetist than the transtibial residual limb in terms of achieving optimal results from a prosthetic fitting. In recognising these challenges, Össur was the first to design and introduce a specific Iceross Transfemoral liner to offer comfort, protection and optimal suspension of the transfemoral prosthesis.


1. Size
2. Limb length
3. Profile
4. Use of fitting aids

1. SIZE

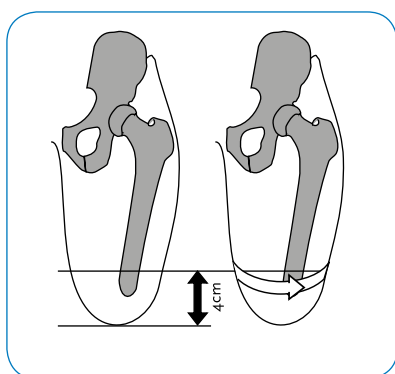

With the tissue of the residual limb hanging down, measure 4cm proximal to the distal end of the limb and mark this point on the skin.

Measure the circumference of the limb at this point with medium tension on the tape measure.

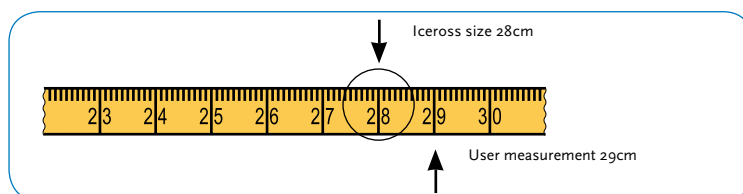

The measurement taken corresponds to the size of Iceross to be used. Select the liner size matching your measurement or, if between sizes select the next size below. Do not select a size larger than your measurement or more than one size below it. e.g. if the circumference of the limb measures 29cm, select Iceross size 28.

Note locking liners:

To achieve optimal results from the Iceross liner, it is important that the correct size of liner is selected. Incorrect sizing may cause discomfort to the user, perspiration and reduced control of the prosthesis.

Note HSM Seal-In liners:

Please refer to the note about size confirmation of the seal on page 25.

Record the measurement and refer to the vertical axis of the selection chart.

2. RESIDUAL LIMB LENGTH


Compare the residual limb length to the sound limb by dividing the sound limb femur into thirds.

- The proximal third is considered **SHORT** length - refer to selection chart part A.
- The middle third is considered **MEDIUM** length - refer to selection chart part B.
- The distal third is considered **LONG** length - refer to selection chart part C.

Transfemoral Selection

3. PROFILE

Measure and record the circumference of the limb at the level of the perineum, refer to the horizontal axis of the selection chart. If the intersection of the horizontal and vertical axes lies in the yellow area of the selection chart then a standard profile is recommended; if it lies in the blue area then a conical profile is recommended.

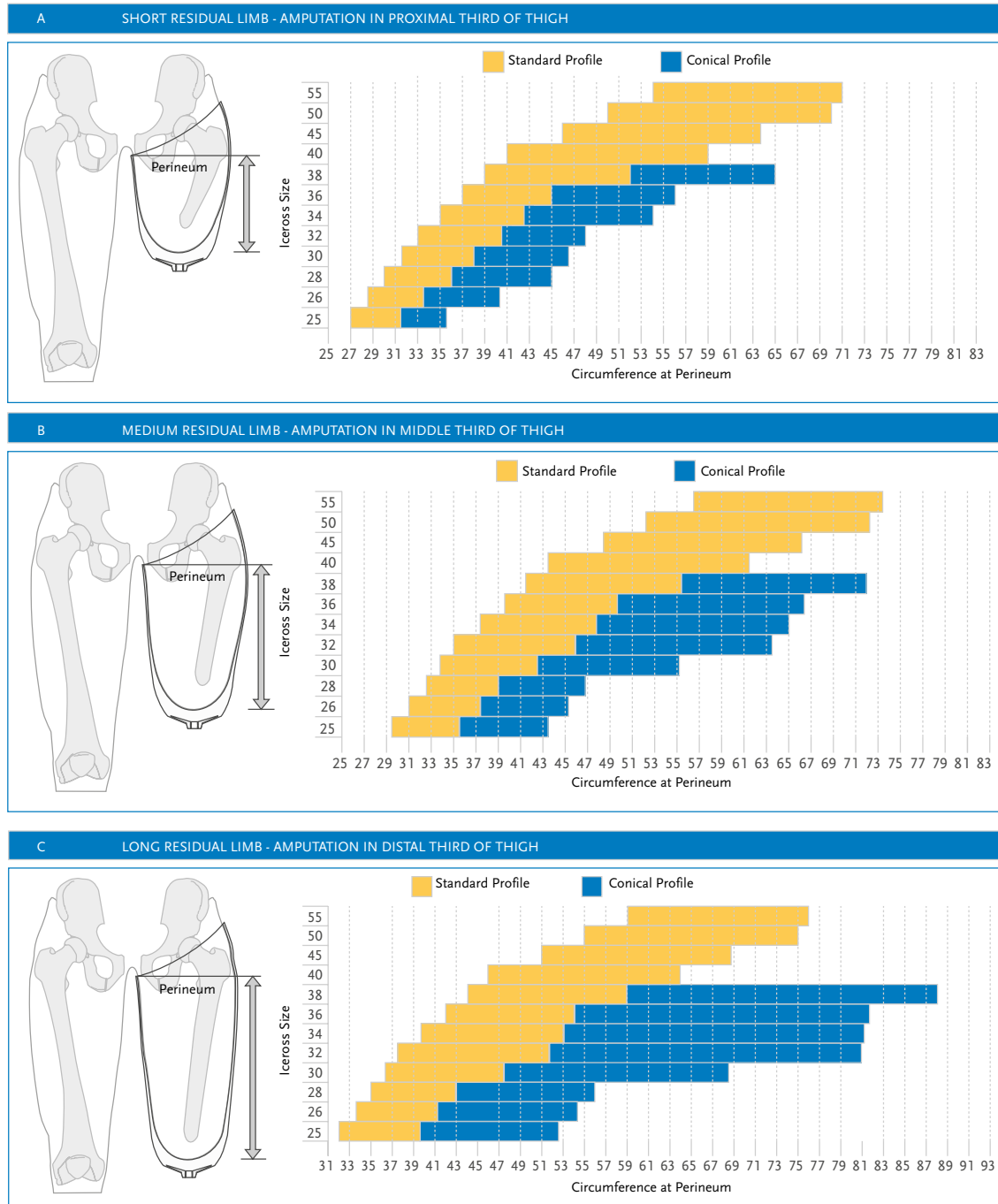

4. FITTING AIDS

Use the principles of profile matching as described in transtibial section. If there are any deep scars or voids fill with Iceross Pad or silicone putty.

Iceross Seal-In® (TT)

Introduction

The Iceross Seal-In range provides a number of options for suspension liners for low, medium and high activity users. Comfortable suspension is achieved without an external sleeve. The liberation from the restriction of a sleeve offers substantially improved freedom of movement, increased comfort and simplified maintenance.



User indications:

- Applicable to most residual limb shapes > 12cm with controlled volume fluctuation
- Recommended for all activity levels
- Ideal for vacuum socket users looking for sleeve-free, partial vacuum suspension

Medical benefits:

- Reduced pistoning in comparison to other suspension methods
- Excellent rotational stability and proprioceptive control
- Comfortable user experience supports active lifestyles

USER INFORMATION

Amputation Level	Transtibial
Impact Level	Low Moderate High

LINER INFORMATION

Available Suspension Methods



Seal-In

Features



Fabric cover



Matrix



Aloe Vera




Wave



Menthol

Iceross Seal-In® V

It is recommended that Iceross Seal-In V is used in conjunction with Icelock® Expulsion Valve L-551002. See page 105 for more details.

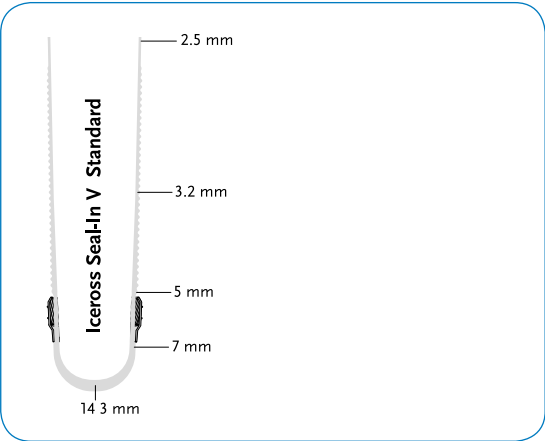

LINER INFORMATION

Size	18, 20, 22, 23, 25, 26, 28, 30, 32, 34, 36 (Size 23 = 23.5 and size 26 = 26.5)
------	---

For sizing and selection advice please refer to pages 24 - 28.

Profile	3mm
---------	-----

Suspension	
------------	--



Seal-In V

Features	
----------	--



Fabric cover



Wave



Active Skin
Care



Silken Inner
Surface

Ordering Information

ICEROSS SEAL-IN V

Part#	Profile	Suspension method
I-4713XX	3mm	Seal-In V

Iceross Seal-In® V High Profile (HP)

It is recommended that Iceross Seal-In V HP is used in conjunction with Icelock® Expulsion Valve L-551002. See page 105 for more details.



LINER INFORMATION

Size 18, 20, 22, 23, 25, 26, 28, 30, 32, 34, 36
(Size 23 = 23.5 and size 26 = 26.5)

For sizing and selection advice please refer to pages 24 - 28.

Profile 3mm

Suspension



Seal-In V

Features



Fabric cover




Wave



Active Skin
Care



Silken Inner
Surface

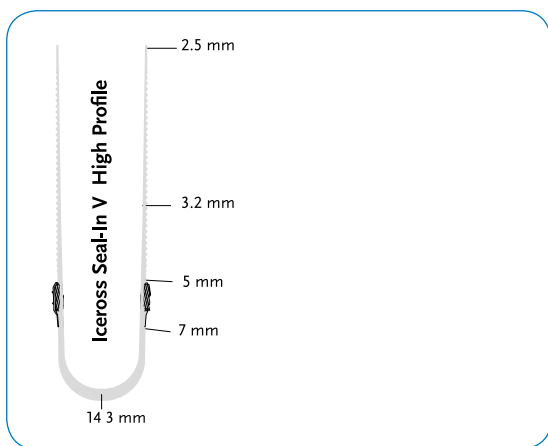


Ordering Information

ICEROSS SEAL-IN V HP

Part#	Profile	Suspension method
I-4723XX	3mm	Seal-In V

Iceross Seal-In® X5 Wave

It is recommended that Iceross Seal-In X5 is used in conjunction with Icelock® Expulsion Valve L-551002. See page 105 for more details.

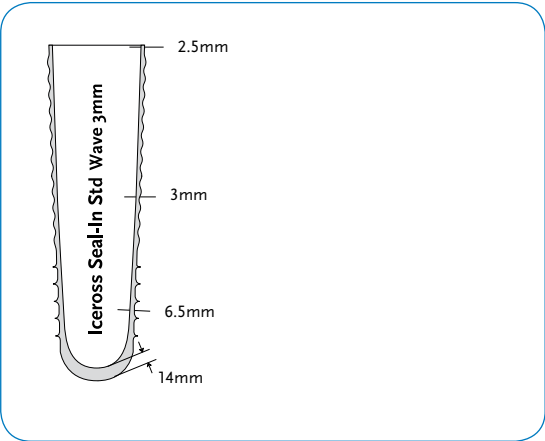

LINER INFORMATION

Size	18, 20, 22, 23, 25, 26, 28, 30, 32, 34, 36 (Size 23 = 23.5 and size 26 = 26.5)
------	---

For sizing and selection advice please refer to pages 24 - 28.

Profile	3mm
---------	-----

Matrix	10-13cm, size dependent
--------	-------------------------

Suspension	
------------	--



Seal-In

Features



Fabric cover



Matrix



Aloe Vera



Menthol



Wave


Ordering Information

ICEROSS SEAL-IN X5 WAVE

Part#	Profile	Suspension method	Matrix
I-3664XX	3mm	Seal-In	10-13cm, size dependent

Iceross Seal-In® X5

It is recommended that Iceross Seal-In X5 is used in conjunction with Icelock® Expulsion Valve L-551002. See page 105 for more details.

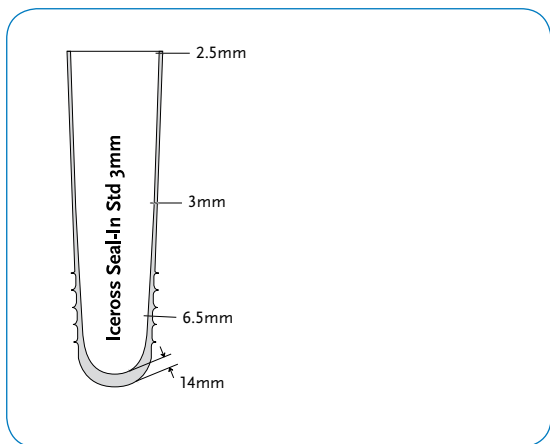

LINER INFORMATION

Size 18, 20, 22, 23, 25, 26, 28, 30, 32, 34, 36
(Size 23 = 23.5 and size 26 = 26.5)

For sizing and selection advice please refer to pages 24 - 28.

Profile 3mm

Matrix 10-13cm, size dependent

Suspension



Seal-In

Features



Fabric cover



Matrix



Aloe Vera



Menthol

Ordering Information

ICEROSS SEAL-IN X5

Part#	Profile	Suspension method	Matrix
I-3663XX	3mm	Seal-In	10-13cm, size dependent

Iceross Dermo Seal-In® HSM

It is recommended that Iceross Dermo Seal-In is used in conjunction with Icelock® Expulsion Valve L-551002. See page 105 for more details.



LINER INFORMATION

Size 18, 20, 22, 23, 25, 26, 28, 30, 32, 34, 36
(Size 23 = 23.5 and size 26 = 26.5)

For sizing and selection advice please refer to pages 24 - 28.

Profile 3mm

Matrix 10cm

*Please note that deviation in the stabilising matrix length can be +/- 1.5cm.

Suspension



Seal-In

Features




Fabric cover



Matrix



Aloe Vera

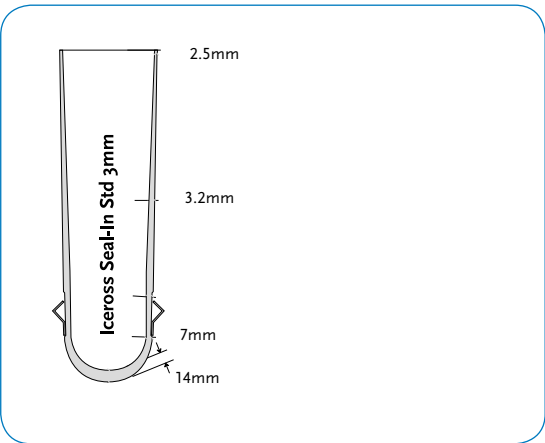

Part Number Information

ICEROSS DERMO SEAL-IN

Part#	Profile	Suspension method	Matrix
I-4613XX	3mm	Seal-In	10cm

Iceross Cushion liners

Introduction

User indications:

- Applicable to all shapes and sizes of residual limbs requiring gentle cushioning
- Recommended for all activity levels
- Ideal for transtibial amputees seeking the reassurance of a suspension sleeve

Medical benefits:

- Pressure relief, which reduces the risk of sensitive or diabetes-impaired skin breakdown
- Good control due to optimal surface contact with the residual limb
- Comfortable user experience supports active lifestyles

The wave feature makes the liner easier to flex at the knee, allowing it to adapt more to the user's movements. Less resistance means that less energy is required to flex the knee. It stretches easily over the patella and effectively controls folding of the liner behind the knee.



Iceross Synergy®

A firm silicone outer layer provides outstanding stability for active use, while a softer shock-absorbing inner layer nurtures the skin, offering superb suspension and comfort. The synergy of the different durometers is what makes this Iceross liner such a comfortable socket interface for daily life.

Iceross Dermo®

Superb softness, gentle skin contact and Össur's unique Active Skin Care ingredients. The combination of DermoGel® silicone with an ultra-strong and elastic Supplex® outer cover offers excellent durability and an intimate fit. Ideal for vascular patients and those with sensitive skin.

Iceross Comfort®


Gentle but strong, the comfortable option. This durable and affordable liner gently conforms to the user's profile, producing a comfortable interface.

Iceross® Activa with Tibiaguard™

Being well acquainted with the world of active people, Össur knows that freedom of movement and shock absorption are on top of the wish list. The Iceross Activa is the first silicone liner of its kind to accommodate these demanding user requirements in one liner. A unique new feature makes it possible - TibiaGuard™.

Iceross Synergy® Cushion

It is recommended that Iceross Synergy Cushion liners are used in conjunction with the Iceross Sleeve and Icelock® Expulsion Valve L-551002. See pages 61/105 for more details.

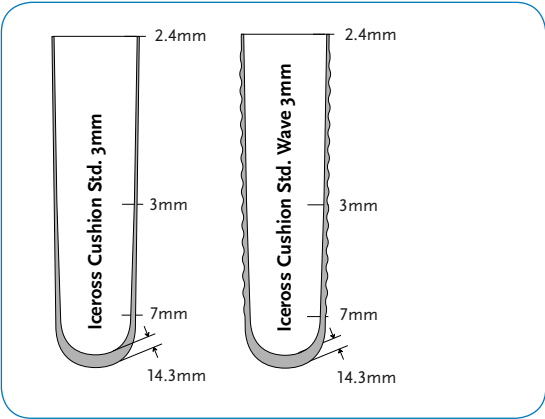

USER INFORMATION

Amputation Level:	Transtibial
Impact Level:	Moderate High Sport

LINER INFORMATION

Size*:	16, 18, 20, 22, 23, 25, 26, 28, 30, 32, 34, 36, 40 (Size 23 = 23.5 and size 26 = 26.5)
--------	---

*For sizing and selection advice please refer to pages 24 - 28.

Profile:	3mm
----------	-----

Suspension:



Cushion

Features:



Fabric cover



Silken inner surface



Double durometer



Active Skin Care



Menthol



Wave


Ordering Information

ICEROSS SYNERGY CUSHION

Part#	Profile	Suspension method
I-CW63XX	3mm Wave	Cushion
I-CL63XX	3mm	Cushion

Iceross Dermo® Cushion

It is recommended that Iceross Dermo Cushion liners are used in conjunction with the Iceross Sleeve and Icelock® Expulsion Valve L-551002. See pages 61/105 for more details.

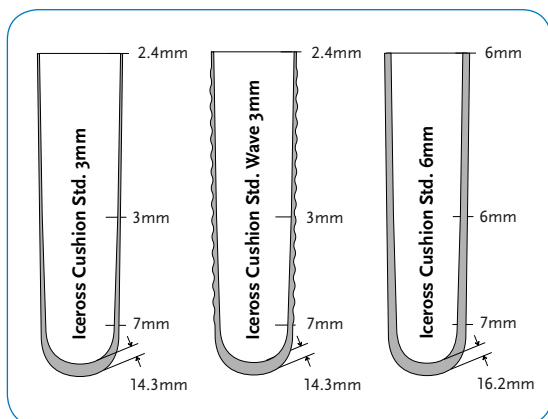

USER INFORMATION

Amputation Level:	Transtibial
Impact Level:	Low Moderate High (only for very sensitive skin or bony residual limb)
Conditions:	Sensitive skin Diabetic Vascular

LINER INFORMATION

Size*:	16, 18, 20, 22, 23, 25, 26, 28, 30, 32, 34, 36, 40 (Size 23 = 23.5 and size 26 = 26.5)
--------	---

*For sizing and selection advice please refer to pages 24 - 28.

Profile:	3mm and 6mm
----------	-------------

Suspension:



Cushion

Features:



Fabric cover

Silken inner surface

Active Skin Care

Wave

Ordering Information


ICEROSS DERMO CUSHION

Part#	Profile	Suspension method
I-CW43XX	3mm Wave	Cushion
I-CL43XX	3mm	Cushion
I-CL46XX	6mm*	Cushion

*Not available in size 40.

Iceross Comfort® Cushion

It is recommended that Iceross Comfort Cushion liners are used in conjunction with the Iceross Sleeve and Icelock® Expulsion Valve L-551002. See pages 61/105 for more details.

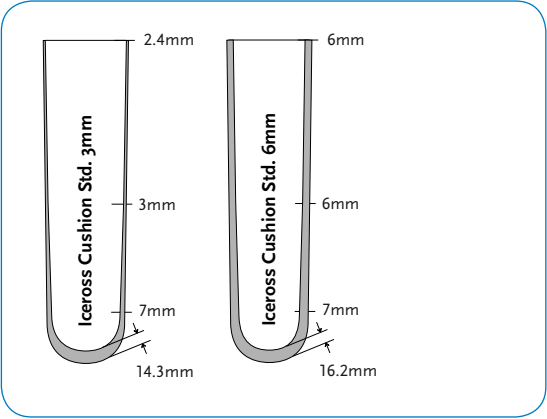

USER INFORMATION

Amputation Level:	Transtibial
Impact Level:	Low Moderate
Conditions:	Sensitive areas Bony prominences

LINER INFORMATION

Size*:	16, 18, 20, 22, 23, 25, 26, 28, 30, 32, 34, 36, 40 (Size 23 = 23.5 and size 26 = 26.5)
--------	---

*For sizing and selection advice please refer to pages 24 - 28.

Profile:	3mm and 6mm
----------	-------------

Suspension:



Cushion

Features:



Fabric cover

Ordering Information


ICEROSS COMFORT CUSHION

Part#	Profile	Suspension method
I-CL53XX	3mm	Cushion
I-CL56XX	6mm*	Cushion

*Not available in size 40.

Iceross® Activa with TibiaGuard™

It is recommended that the Iceross Activa with Tibiaguard liner is used in conjunction with the Iceross Sleeve and Icelock® Expulsion Valve L-551002. See pages 61/105 for more details.

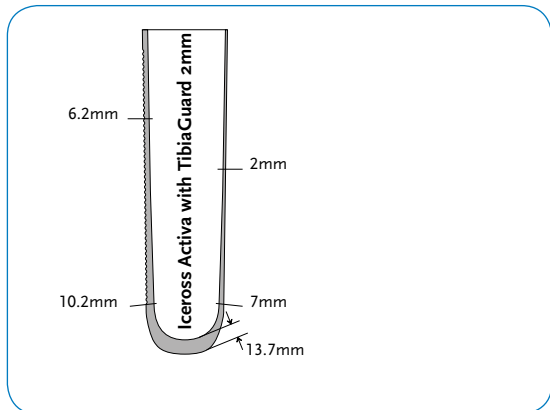

LINER INFORMATION

Amputation Level:	Transtibial
Impact Level:	Moderate, High and Sport

LINER INFORMATION

Size*:	18, 20, 22, 23, 25, 26, 28, 30, 32, 34 (Size 23 = 23.5 and size 26 = 26.5)
--------	---

*For sizing and selection advice please refer to Product Catalogue, pages 24 - 28.

Profile:	2mm
----------	-----

Suspension:



Cushion

Features:



TibiaGuard



Fabric Cover



Active Skin Care



Wave



Double
Durometer



Silken Inner
Surface

Ordering Information

ICEROSS ACTIVA

Part#	Profile	Suspension method
I-CT62XX	2mm Wave	Cushion

Iceross Locking Liners

Introduction

User indications:

- Especially good for those requiring soft tissue stabilisation
- Applicable to medium and shorter limb lengths (<12cm)

Medical benefits:

- Increased control due to soft tissue stabilisation
- Reduced pistoning in mid-swing
- Good compliance – easy to use for patients with compromised hand-agility

The wave feature makes the liner easier to flex at the knee, allowing it to adapt more to the user's movements. Less resistance means that less energy is required to flex the knee. It stretches easily over the patella and effectively controls folding of the liner behind the knee.



Iceross Synergy®

A firm silicone outer layer provides outstanding stability for active use, while a softer shock-absorbing inner layer nurtures the skin, offering superb suspension and comfort. The synergy of the different durometers is what makes this Iceross liner such a comfortable socket interface for daily life.

Iceross Dermo®


Superb softness, gentle skin contact and Össur's unique Active Skin Care ingredients. The combination of DermoGel® silicone with an ultra-strong and elastic Supplex® outer cover offers excellent durability and an intimate fit. Ideal for vascular patients and those with sensitive skin.

Iceross Comfort®

Gentle but strong, the comfortable option. This durable and affordable liner gently conforms to the user's profile, producing a comfortable interface.

Iceross Synergy® Locking

It is recommended that Iceross locking liners are used in conjunction with the Icelock® 200 or 600 Series. See pages 99/106 for more details.

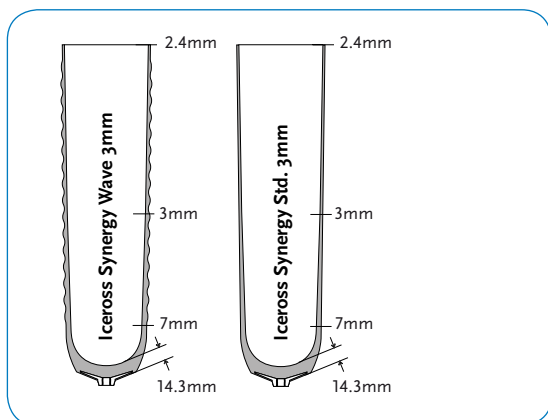

USER INFORMATION

Amputation Level:	Transtibial
Impact Level:	Moderate High Sport

LINER INFORMATION

Size*:	16, 18, 20, 22, 23, 25, 26, 28, 30, 32, 34, 36, 40, 45 (Size 23 = 23.5 and size 26 = 26.5)
--------	--

*For sizing and selection advice please refer to pages 24 - 28.

Profile:	3mm
----------	-----

Matrix:	10cm or a custom matrix** (2cm-14cm)
---------	--------------------------------------

**For custom matrix, please specify the required matrix length on your order. Custom and special order products are not returnable. Please note that deviation in the stabilising matrix length can be +/- 1.5cm.

Suspension:



Locking

Features:



Fabric cover



Silken inner surface



Stabilising matrix



Double durometer



Active Skin Care



Menthol



Wave


Ordering Information

ICEROSS SYNERGY

Part#	Profile	Suspension method	Matrix
I-6303XX	3mm Wave	Locking	10cm
I-6206XX	3mm Wave	Locking	Custom 2-14cm
I-6003XX	3mm	Locking	10cm
I-6203XX	3mm	Locking	Custom 2-14cm

Iceross Dermo® Locking

It is recommended that Iceross locking liners are used in conjunction with the Icelock® 200 or 600 Series. See pages 99/106 for more details.

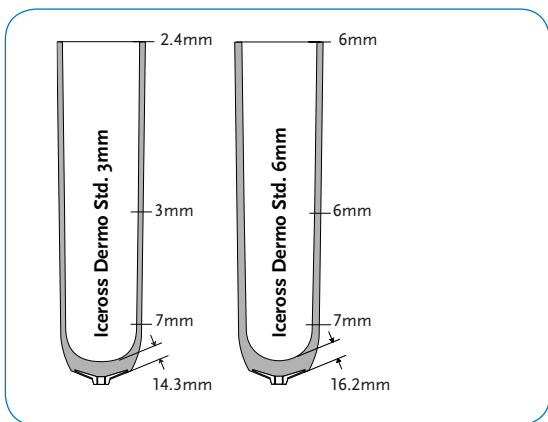

USER INFORMATION

Amputation Level:	Transtibial
Impact Level:	Low Moderate High (only for very sensitive skin or bony residual limb)
Conditions:	Sensitive skin Diabetic Vascular

LINER INFORMATION

Size*:	16, 18, 20, 22, 23, 25, 26, 28, 30, 32, 34, 36, 40, 45 (Size 23 = 23.5 and size 26 = 26.5)
--------	---

*For sizing and selection advice please refer to pages 24 - 28.


Profile:	3mm
Matrix:	10cm, 5cm or a custom matrix** (2cm-14cm)

**For custom matrix, please specify the required matrix length on your order. Custom and special order products are not returnable. Please note that deviation in the stabilising matrix length can be +/- 1.5cm.

Suspension:



Features:

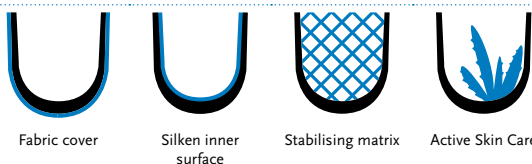


Ordering Information

ICEROSS DERMO

Part#	Profile	Suspension method	Matrix
I-4013XX	3mm	Locking	10cm
I-4513XX	3mm	Locking	5cm
I-4213XX	3mm	Locking	Custom 2-14cm
I-4016XX	6mm	Locking	10cm
I-4516XX	6mm	Locking	5cm
I-4216XX	6mm	Locking	Custom 2-14cm

Iceross Dermo® Locking with Wave

It is recommended that Iceross locking liners are used in conjunction with the Icelock® 200 or 600 Series. See pages 99/106 for more details.

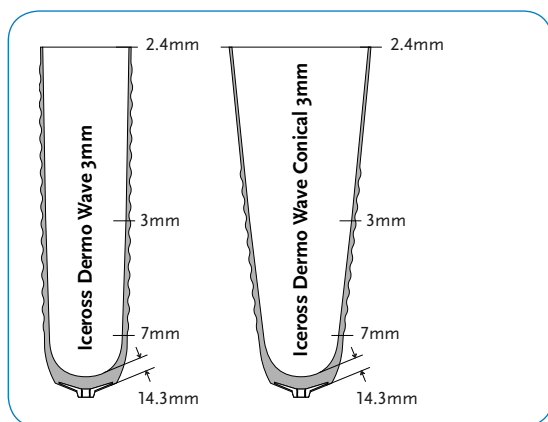

USER INFORMATION

Amputation Level:	Transtibial
Impact Level:	Low Moderate High (only for very sensitive skin or bony residual limb)
Conditions:	Sensitive skin Diabetic Vascular

LINER INFORMATION

Size*:	16, 18, 20, 22, 23, 25, 26, 28, 30, 32, 34, 36, 40, 45 (Conical sizes available 18-36) (Size 23 = 23.5 and size 26 = 26.5)
--------	---

*For sizing and selection advice please refer to pages 24 - 28.


Profile:	3mm
Matrix:	10cm or a custom matrix** (2cm-14cm)

**For custom matrix, please specify the required matrix length on your order. Custom and special order products are not returnable. Please note that deviation in the stabilising matrix length can be +/- 1.5cm.

Suspension:



Features:

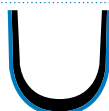

Fabric cover



Silken inner surface



Stabilising matrix



Active Skin Care



Wave



Conical


Ordering Information

ICEROSS DERMO

Part#	Profile	Suspension method	Matrix
I-4313XX	3mm Wave	Locking	10cm
I-4223XX	3mm Wave	Locking	Custom 2-14cm
I-4913XX	Conical 3mm Wave	Locking	10cm

Iceross Comfort® Locking

It is recommended that Iceross locking liners are used in conjunction with the Icelock® 200 or 600 Series. See pages 99/106 for more details.

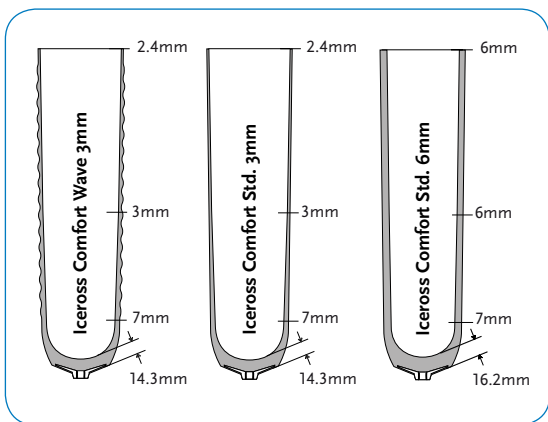

USER INFORMATION

Amputation Level:	Transtibial
Impact Level:	Low Moderate
Conditions:	Sensitive areas Bony prominences

LINER INFORMATION

Size*:	16, 18, 20, 22, 23, 25, 26, 28, 30, 32, 34, 36, 40, 45 (Size 23 = 23.5 and size 26 = 26.5)
--------	--

*For sizing and selection advice please refer to pages 24 - 28.

Profile:	3mm and 6mm
Matrix:	10cm, 5cm or a custom matrix** (2cm-14cm)

**For custom matrix, please specify the required matrix length on your order. Custom and special order products are not returnable. Please note that deviation in the stabilising matrix length can be +/- 1.5cm.

Suspension:



Locking

Features:



Fabric cover



Stabilising matrix



Wave

Ordering Information


ICEROSS COMFORT

Part#	Profile	Suspension method	Matrix
I-5406XX	3mm	Locking	10cm
I-5806XX	3mm	Locking	5cm
I-5506XX	3mm	Locking	Custom 2-14cm
I-5006XX	6mm	Locking	10cm
I-5706XX	6mm	Locking	5cm
I-5106XX	6mm	Locking	Custom 2-14cm
I-5306XX	3mm Wave	Locking	10cm

Iceross® Sport

Heavy-duty rotation control and an aluminium umbrella ensure this liner's performance enhances that of its user, whatever the sport. Interchangeable distal attachment.

It is recommended that Iceross locking liners are used in conjunction with the Icelock® 200 or 600 Series. See pages 99/106 for more details.

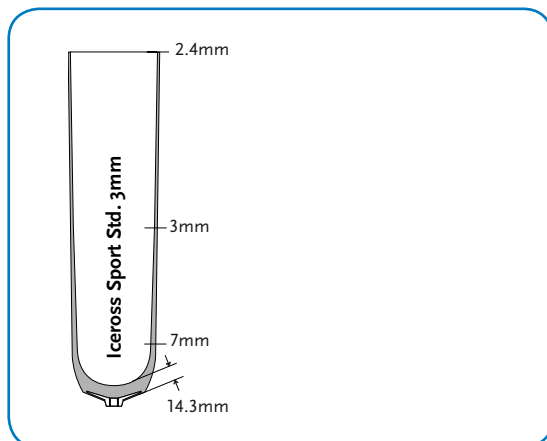

USER INFORMATION

Amputation Level	Transtibial
Impact Level	High & Sport
Conditions	Soft tissue control

LINER INFORMATION

Size	16, 18, 20, 22, 23, 25, 26, 28, 30, 32, 34 (Size 23 = 23.5 and size 26 = 26.5)
------	---

For sizing and selection advice please refer to pages 24 - 28.

Profile	3mm
Matrix	5cm, 10cm or custom matrix (2cm-14cm)

* For custom matrix, please specify the required matrix length on your order. Custom and special order products are not returnable. Please note that deviation in the stabilising matrix length can be +/- 1.5cm.

Suspension



Locking

Features



Fabric cover



Matrix



Double
durometer



Aloe Vera



Menthol

Ordering Information


ICEROSS SPORT

Part#	Profile	Suspension method	Matrix
I-1533XX	3mm	Locking	5cm
I-1033XX	3mm	Locking	10cm
I-1233XX	3mm	Locking	Custom 2-14cm*

Iceross® Original Locking

Iceross Original with cover

It is recommended that Iceross locking liners are used in conjunction with the Icelock® 200 or 600 Series. See pages 99/106 for more details.

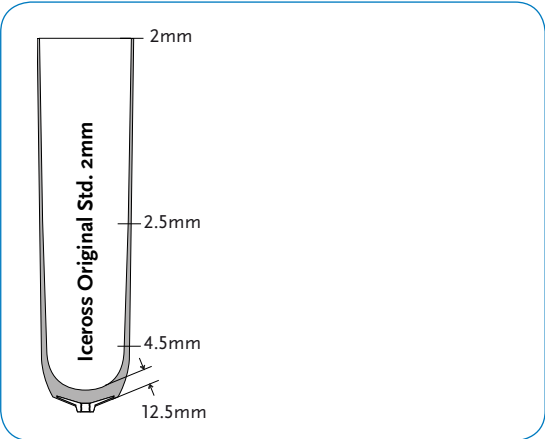

USER INFORMATION

Amputation Level	Transtibial
Impact Level	All

LINER INFORMATION

Size	16, 18, 20, 22, 23, 25, 26, 28, 30, 32, 34, 36, 38, 40, 42, 45 (Size 23 = 23.5 and size 26 = 26.5)
------	---

For sizing and selection advice please refer to pages 24 - 28.

Profile	2mm
Matrix	5cm, 10cm or a custom matrix (2cm-14cm)

*For custom matrix, please specify the required matrix length on your order. Custom and special order products are not returnable. Please note that deviation in the stabilising matrix length can be +/- 1.5cm.

Suspension



Locking

Features



Fabric cover



Matrix

Ordering Information


ICEROSS ORIGINAL WITH COVER

Part#	Profile	Suspension method	Matrix
I-2124XX	2mm	Locking	5cm
I-0124XX	2mm	Locking	10cm
I-2024XX	2mm	Locking	Custom 2-14cm*

Iceross® Original Locking

Iceross Original without cover

It is recommended that Iceross locking liners are used in conjunction with the Icelock® 200 or 600 Series. See pages 99/106 for more details.

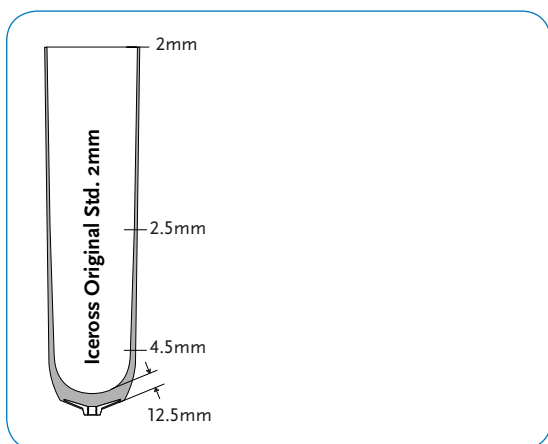

LINER INFORMATION

Size	12, 14, 16, 18, 20, 22, 23, 25, 26, 28, 30, 32, 34, 36, 38, 40, 42, 45 (Size 23 = 23.5 and size 26 = 26.5)
------	---

For sizing and selection advice please refer to pages 24 - 28.

Profile	2mm
Matrix	5cm, 10cm or a custom matrix (2cm-14cm)

*For custom matrix, please specify the required matrix length on your order. Custom and special order products are not returnable. Please note that deviation in the stabilising matrix length can be +/- 1.5cm.

Suspension



Locking

Features



Matrix

Ordering Information

ICEROSS ORIGINAL WITHOUT COVER

Part#	Profile	Suspension method	Matrix
I-2104XX	2mm	Locking	5cm
I-0104XX	2mm	Locking	10cm
I-2004XX	2mm	Locking	Custom 2-14cm*

Iceross Transfemoral Seal-In® Liners

Introduction

Iceross Transfemoral Seal-In Liners are specifically designed to provide secure suspension for low, moderate and high activity transfemoral users. The incorporated seals provide a firm, comfortable suspension controlling rotation within the socket. It has a full length stabilising matrix for added stability control.



USER INFORMATION

Amputation Level:	Transfemoral
Impact Level:	All

LINER INFORMATION

Available Suspension Methods:



Seal-In X5



Seal-In

Features:



Fabric cover



Silken inner surface



Stabilising matrix



Active Skin Care



Menthol



Conical

Iceross Seal-In® X5 (TF)

It is recommended that Iceross Seal-In X5 is used in conjunction with Icelock® Expulsion Valve L-551002 or L-552000. See page 105 for more details.




LINER INFORMATION

Size Standard*:	25, 26, 28, 30, 32, 34, 36, 38, 40, 45, 50, 55 (size 26 = 26.5)
Size Conical*:	25, 26, 28, 30, 32, 34, 36, 38 (size 26 = 26.5)

*For sizing and selection advice please refer to pages 24 - 28.

Profile: Standard and Conical 3mm

Suspension:

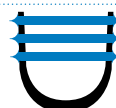

Seal-In X5

Features:



Fabric cover



Silken inner surface



Stabilising matrix




Active Skin Care



Menthol



Conical

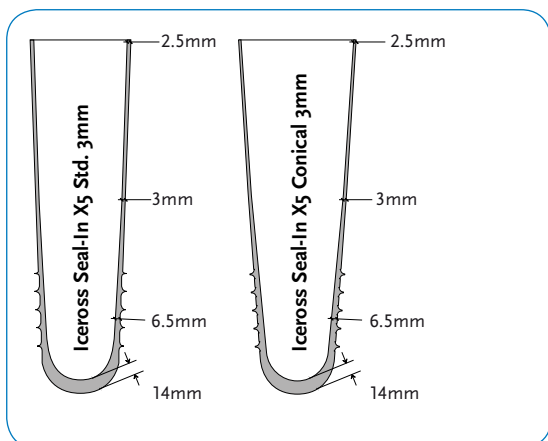

Ordering Information

ICEROSS SEAL-IN X5

Part#	Profile	Suspension method
I-TF673XX	Standard 3mm	Seal-In X5
I-TF678XX	Conical 3mm	Seal-In X5

Iceross Transfemoral Seal-In®

It is recommended that Iceross TF HSM Seal-In is used in conjunction with Icelock® Expulsion Valve L-551002 or L-552000. See page 105 for more details.



LINER INFORMATION

Size Standard*:	25, 26, 28, 30, 32, 34, 36, 38, 40, 45, 50, 55 (Size 26 = 26.5)
Size Conical*:	25, 26, 28, 30, 32, 34, 36, 38 (Size 26 = 26.5)




*For sizing and selection advice please refer to pages 24 - 28.

Profile:	Standard and Conical 3mm
----------	--------------------------

Suspension:



Features:

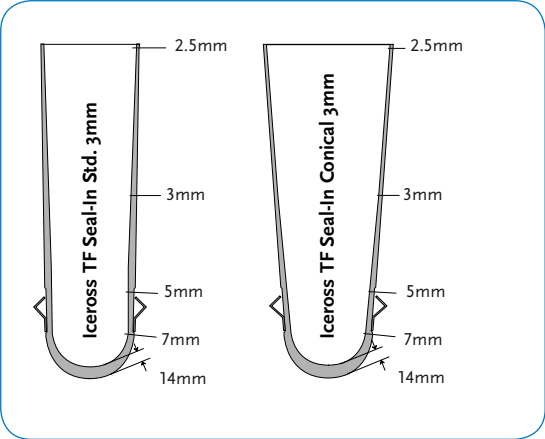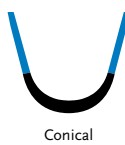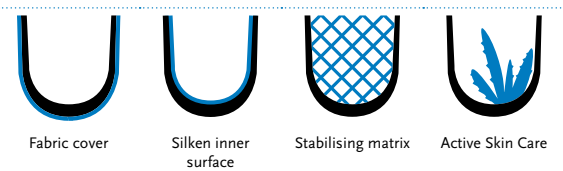

Ordering Information

ICEROSS TRANSFEMORAL SEAL-IN

Part#	Profile	Suspension method
I-7532XX	Standard 3mm	Seal-In
I-7632XX	Conical 3mm	Seal-In

Iceross® Transfemoral Locking Liners

Introduction

The first ever liner designed specifically for transfemoral applications, combines gentle skin care with a full length matrix for excellent suspension and stability.



USER INFORMATION

Amputation Level:	Transfemoral
Impact Level:	All

LINER INFORMATION

Available Suspension Methods:



Locking

Features:



Fabric cover



Silken inner surface



Stabilising matrix



Active Skin Care



Conical

Iceross® Transfemoral Locking

It is recommended that Iceross locking liners are used in conjunction with the Icelock® 200 or 600 Series. See pages 99/106 for more details.



LINER INFORMATION

Size Standard*:	25, 26, 28, 30, 32, 34, 36, 38, 40, 45, 50, 55 (Size 26 = 26.5)
Size Conical*:	25, 26, 28, 30, 32, 34, 36, 38 (Size 26 = 26.5)

*For sizing and selection advice please refer to pages 24 - 28.

Profile:	Standard 2mm and conical 3mm
----------	------------------------------

Suspension:



Locking

Features:



Fabric cover



Silken inner surface




Stabilising matrix



Active Skin Care



Conical

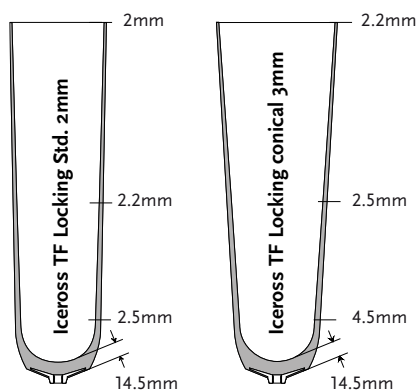

Ordering Information

ICEROSS TRANSFEMORAL


Part#	Profile	Suspension method
I-7032XX	Standard 2mm	Locking
I-7132XX	Conical 3mm	Locking

Iceross® Upper-X

Iceross Upper-X - Locking

Stability and anatomical fit for upper extremity amputees.

Upper-X has an M6 threaded distal attachment. It is recommended that Iceross Upper-X is used with the Icelock® 700 Series. See page 112 for more details.

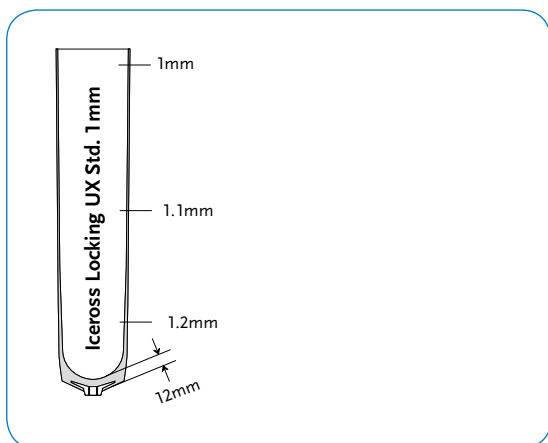

USER INFORMATION


Amputation Level	Transradial and Transhumeral
Impact Level	All
Conditions	Soft tissue control

LINER INFORMATION

Size	12, 14, 16, 18, 20, 22.
------	-------------------------

Measure the circumference 4cm from the distal end and select the corresponding size, or the next size down.

Note: Downsizing the liner is not recommended as it may compromise suspension and cause discomfort.

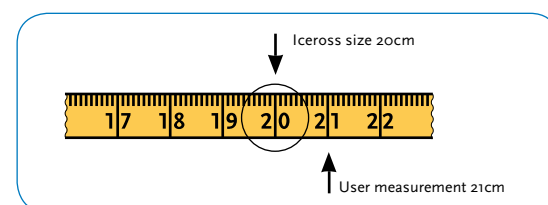

Profile	1mm
---------	-----

Suspension



Locking

Features



Fabric cover



Matrix



Aloe Vera

Ordering Information

ICEROSS UPPER-X - LOCKING

Part#	Profile	Suspension method
I-8131XX	1mm	Locking
L-792031	Pin Adapter, M10 - M6 *	

*To convert alternative liners with M10 thread to M6 for use with 700 series lock

Iceross® Junior NEW

Introduction

Available in a choice of two silicone blends to suit activity levels and limb conditions. A blank cover and a set of magic markers mean that children can now create their own designer liner, helping them to come to terms with their prosthesis and have some fun! It is recommended that Iceross Junior liners are used in conjunction with the Icelock® 700 series. See page 112 for more details.



USER INFORMATION

Amputation Level: Transtibial

LINER INFORMATION

Size*: 12, 14, 16, 18

*For sizing and selection advice please refer to pages 24 - 28.

Profile: 3mm

Matrix: 5mm

Suspension:



Locking

Features:



Fabric cover




Silken inner surface



Active Skin Care



Stabilising matrix

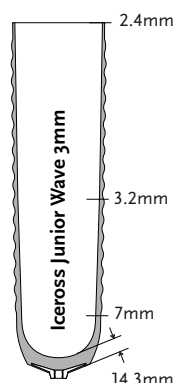

Ordering Information

ICEROSS JUNIOR

Part#	Profile	Material	Suspension method
I-8301XX	3mm	Stabilo - DermoSil	Locking
I-8302XX	3mm	Dermo - DermoGel	Locking
L-792031	Pin Adapter, M10 - M6 *		

*To convert alternative liners with M10 thread to M6 for use with 700 series lock

Iceross® Sleeve

Introduction

The Iceross Sleeve is designed for longer lasting suction, more security and less resistance at the knee.

Össur took Wave technology to the next level when developing the Iceross Sleeve. The Wave feature makes the liner easier to flex at the knee, allowing it to adapt more to the user's movements. Less resistance means that less energy is required to flex the knee. It stretches effectively at the front and folds in a controlled way behind the knee, reducing vaso-constriction and positively influencing blood-circulation in this area.

The profile of the sleeve is thicker around the socket brim, focusing on strength at selected areas to reduce wear and tear.

The Iceross Sleeve is for all activity levels.

It is recommended that Iceross Sleeves are used in conjunction with the Iceross Cushion Liner range and the Icelock® 500 Series. See pages 41 and 105 for more information.




SLEEVE INFORMATION

Size	35, 41, 49, 55, 62
Measure circumference 15cm above mid-patella and choose the correct size.	
Length	35cm

Ordering Information

ICEROSS SLEEVE

Part#	Size	Measurement
F-201135	35cm	35 - 41cm
F-201141	41cm	41 - 49cm
F-201149	49cm	49 - 55cm
F-201155	55cm	55 - 62cm
F-201162	62cm	62 - 76cm

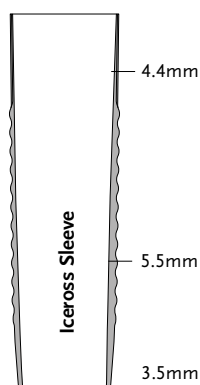

Iceflex® Endurance

Introduction

- Pre-flexed to facilitate knee flexion
- Radial stretch in distal and proximal areas, longitudinal stretch in knee area
- Vertical reinforcements for reduced pistoning
- Reinforced seal at proximal and distal ends
- For moderate and high impact levels



Ordering Information

Measure circumference 15cm above mid-patella and choose the correct size

Part#	Size	Length	Measurement
F-010735	35	36cm	35-41cm
F-010741	41	36cm	41-49cm
F-010749	49	36cm	49-55cm
F-010755	55	36cm	55-62cm
F-010762	62	36cm	62-69cm

Note: Correct alignment is important. Place the pre-flexed knee section over the brim of the socket and the vertical reinforcement ribbons along the sides in the M/L direction. Allow sleeve to extend approximately 5cm above the proximal edge of the interface to maintain an airtight seal.

Iceflex® Balance

Introduction

- Easy to don
- Soft, elastic cover
- Soft silicone inner layer for gentle skin contact
- For low and moderate impact levels



Ordering Information

Measure circumference 15cm above mid-patella and choose the correct size

Part#	Size	Length	Measurement
F-101035	35	36cm	35-49cm
F-101049	49	36cm	49-62cm
F-101062	62	36cm	62-76cm

Note: Allow sleeve to extend approximately 5cm above the proximal edge of the interface to maintain an airtight seal.

Distal Cup™

The Distal Cup is intended for use with Iceross liners. It is made from silicone and is designed for conical or sensitive residual limbs that require volume build-up.



PRODUCT INFORMATION

Size	06, 07, 09, 11, 12, 14, 16, 17, 19, 20, 22, 24, 26.
------	---

To select the correct Distal Cup size, the circumference is measured approximately 4cm above the distal end of the residual limb. Select the Distal Cup two or three sizes smaller than the measurement and roll on. To select the correct Iceross liner size to fit over the Distal Cup, the circumference is measured approximately 4cm above the end of the Distal Cup. Select Iceross liner size as measured or the next size down.

Ordering Information

Part#	Description
I-3000XX	Distal Cup

Iceross® Pads

Iceross Pads are ideal for pressure distribution in highly loaded areas, padding over bony prominences and as a filler for deep voids and scarred residual limbs.



Ordering Information

Part#	Size	Fits Iceross Liner Size
I-610338	38mm	12-20
I-610348	48mm	21-23,5
I-610358	58mm	25-28
I-610368	68mm	30-32
I-610378	78mm	34-36
I-610388	88mm	38-42
I-610398	98mm	45

Iceross® Accessories

Introduction

Iceross Clean & Simple skin care products have been carefully developed to help ensure that the use and maintenance of Iceross liners is convenient, comfortable and hygienic.

Ordering Information

W Also available on www.ossurwebshop.co.uk

CLEAN & SIMPLE SET 1

Part#	Description
S-010001	Soap, Lotion & Lubricant Spray W

This set includes one bottle each of soap, moisturising lotion and lubricant spray (250ml). For Iceross liners without fabric covers and Iceross Seal-In liners.

CLEAN & SIMPLE SET 2

Part#	Description
S-010002	2 Soaps & 1 Lotion W

This set includes two bottles of soap and one bottle of moisturising lotion (250ml). It is ideal for those using Iceross liners with covers when a lubricant spray is not required.

CLEAN & SIMPLE SOAP

Part#	Description
S-210250	Soap W

Iceross Clean & Simple soap is a mild liquid soap for daily hygienic treatment of the skin and the silicone liner. Free from fragrances and dyes, it is designed for sensitive skin and contains natural ingredients that nourish and soften the skin. It has a pH of 5.5 and is hypoallergenic. (250ml)

CLEAN & SIMPLE MOISTURISING LOTION

Part#	Description
S-410250	Moisturising Lotion W

Iceross Clean & Simple moisturising lotion is a non-greasy restorative lotion that nourishes and softens the skin. It is free from fragrance and dyes, has a pH of 5.5 and is hypoallergenic. (250ml)

CLEAN & SIMPLE LUBRICANT SPRAY

Part#	Description
S-610250	Lubricant Spray W

Iceross Clean & Simple lubricant spray aids the donning and doffing of Iceross Original without fabric cover and Iceross Seal-In liners. (250ml)

SILICONE ADHESIVE

Part#	Description
I-400310	For repair of minor damage to Iceross liners (8g)

Helpful for stabilisation of fabric cover on Iceross Transfemoral liners after trimming.

ICEROSS CURVEMASTER

Part#	Description
K-220000	Iceross Curvemaster
K-200150*	Replacement Cutting Wheels Set

*Includes

Upper cutting wheel (K-200126)

Lower cutting wheel (K-200127)

Friction ring (K-200118)



Iceross® Accessories

Ordering Information

Part#	Description	Length
K-190200	Interface Sheath	28cm

The Interface Sheath makes donning and doffing of the socket easier. Designed with a reinforced distal opening for durability.



Part#	Description	Length
K-520002	Iceross Sock 2-Ply – Two Socks in a Pack	32cm
K-520004	Iceross Sock 4-Ply – Two Socks in a Pack	33cm

The socks are thin distally to prevent the residual limb from being pushed out of the socket as more socks are added. This will ensure good function from the socket lock. Designed with a fully knitted hole for the distal connector.




ICEROSS SEAL-IN SOCK

Part#	Size	Iceross Size	Length	Ply
K-531110	10	16-25	35cm	1
K-531113	13	26.5-34	35cm	1
K-531116	16	36-45	35cm	1
K-531120	20	50-55	35cm	1
K-531310	10	16-25	35cm	3
K-531313	13	26.5-34	35cm	3
K-531316	16	36-45	35cm	3
K-531320	20	50-55	35cm	3



The Seal-In sock is designed to work in combination with the Iceross Seal-In liner to obtain an improved fit if volume reduction is experienced in the proximal area of the residual limb. The Seal-In sock is intended for use above the Hypobaric Sealing Membrane (HSM™) only and should never cross it. For X5 liners the sock should only cover the top seal.

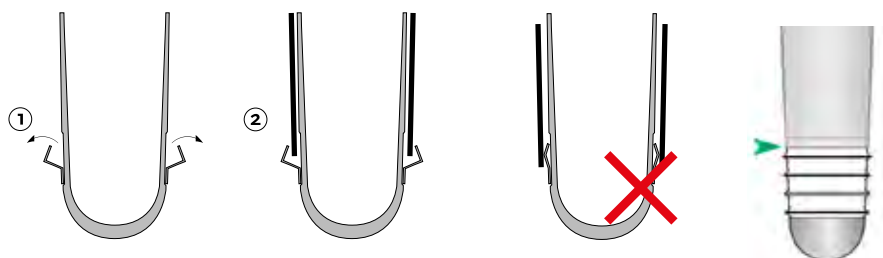

Note: Loss of suspension will be encountered if a sock is placed over the HSM or all 5 seals of the X5.

ICEROSS SEAL-IN SPRAY PEN

Part#	Description
I-990100	Iceross Seal-In Spray Pen ^W

A convenient 9ml spray for handbag or pocket.



ICEROSS PROSTHESIS BAGS

Part#	Description
P-960110	Iceross Below Knee Prosthesis Bag ^W
P-960120	Iceross Above Knee Prosthesis Bag ^W



^W Also available on www.ossurwebshop.co.uk



Unity® Vacuum Suspension System

Introduction

Building on renowned Iceross Seal-In® technology, Össur designers and engineers have developed a new form of advanced vacuum suspension that is also compatible with Flex-Foot® technology.

The design challenge was to apply elevated vacuum with a Seal-In liner, using a simple pump mechanism which had minimal impact on the function of the foot. By applying vacuum where it is needed – around the distal part of the limb – Unity helps to stabilise volume while creating highly effective and sleeveless vacuum suspension.

Without a sleeve, range-of-motion is not restricted and vacuum can be released quickly and easily at the touch of a button. The design also minimises the risk of proximal vacuum leakage and the puncture issues associated with more traditional 'above-knee' vacuum methods.

Key Advantages

Unity offers the security and freedom of vacuum suspension, but without the sensation and restrictions of a sleeve. This advanced system from Össur effectively addresses volume fluctuation and is totally compatible with Flex-Foot technology.

- No sleeve required: Greater mobility and user acceptance
- Distal vacuum: Effective volume stabilisation
- Simple design: Quick and easy to achieve/release vacuum
- Independent pump: Shock function uncompromised
- Light weight: Complete system weighs only 130g

Product Overview

UNITY® FOR BIONIC



Unity® for PROPRIO FOOT®

UNITY® FOR RE-FLEX™



Unity® for Re-Flex Rotate™



Unity® for Re-Flex Shock™



Unity® for LP Rotate™

UNITY® FOR VARI-FLEX®



Unity® for Vari-Flex®



Unity® for LP Vari-Flex®

Unity® for PROPRIO FOOT®

Introduction



Unity for PROPRIO FOOT combines motor powered ankle motion and sleeveless vacuum suspension in one seamless solution. For the low to moderately active user, Unity and Seal-In V® provide secure suspension and volume management while PROPRIO FOOT increases toe clearance and adapts to the ground for safe and effortless ambulation. Unity adds minimal weight to the prosthesis and does not interfere with the PROPRIO FOOT function.

USER INFORMATION

Amputation Level	Transtibial
Impact Level	All


SYSTEM INFORMATION

Foot Sizes	25-30
Maximum Vacuum	-23inHg (600mmHg)
Normal Operating Vacuum:	-16-20 inHg
Unity System Weight	130g w/ Pump, Socket Valve and Tube

Ordering Information

PROPRIO FOOT can be ordered with Unity in sizes 25 and above. Please select the appropriate Unity Pump Module kit according to foot module size when placing an order and Unity will come preassembled on the foot. For extra convenience and system compatibility, Unity TT Valve Kit is highly recommended as a means of connecting to the socket. For liner selection charts and item numbers please refer to Iceross Seal-In V pages.

For complete foot with male pyramid

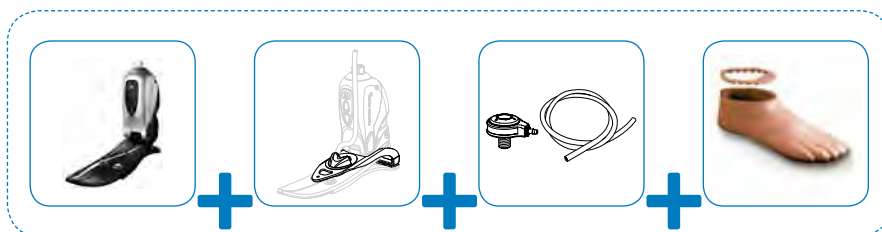

Unity® for PROPRIO FOOT®

Category Selection Guide



Weight Kg

Impact Level	45-52	53-59	60-68	69-77	78-88	89-100	101-116	117-125
Low	1	1	2	3	4	5	6	7
Moderate	1	2	3	4	5	6	7	8

Kits

UNITY KITS

Picture	Part#	Description	Applicable Foot Kits
1	UPR0027	Unity Pump Module PRX 25-27	PROPRIO FOOT size 25-27
1	UPR0030	Unity Pump Module PRX 28-30	PROPRIO FOOT size 28-30
2	UTT0001	Unity TT Valve Kit	All
3	UTT0002	Unity Distal Vacuum Plate Kit	All

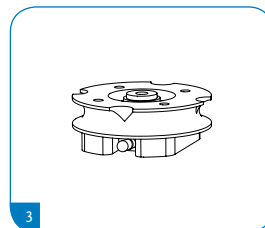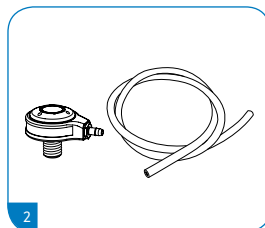


Spare Parts

SPARE COMPONENTS

Part#	Description
UVL00022	Unity VL Components 25-30

Includes

- Membrane
- Check valves
- Barb fittings
- Tube
- Tube fastener
- Heel pad

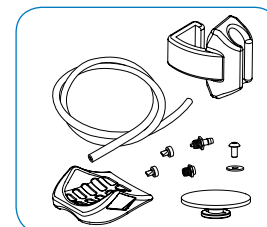

Unity® for PROPRIO FOOT®

PROPRIO Kit With Unity

PROPRIO FOOT KIT




Includes:

- PROPRIO FOOT Ankle
- Foot module
- PROPRIO FOOT Battery dummy
- PROPRIO FOOT Power Supply
- Heel Wedge Kit
- Terms of Warranty
- Quick Reference Cards
- Technical Manual
- Fitting Guide
- Instructions for Use
- Troubleshooting Guide



Cosmetic Finishing

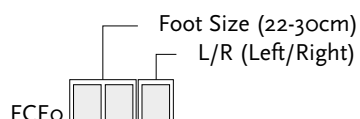

Please specify when ordering brown covers. Add "BR" to the part#.

Includes


- Foot Cover
- Flex-Foot EVO Sock
- Foot Cover Attachment



Part#	Description	Size
FCE63006	Flex-Foot EVO Sock	22-25cm
FCE63007	Flex-Foot EVO Sock	26-30cm

Unity® for Re-Flex™

Introduction

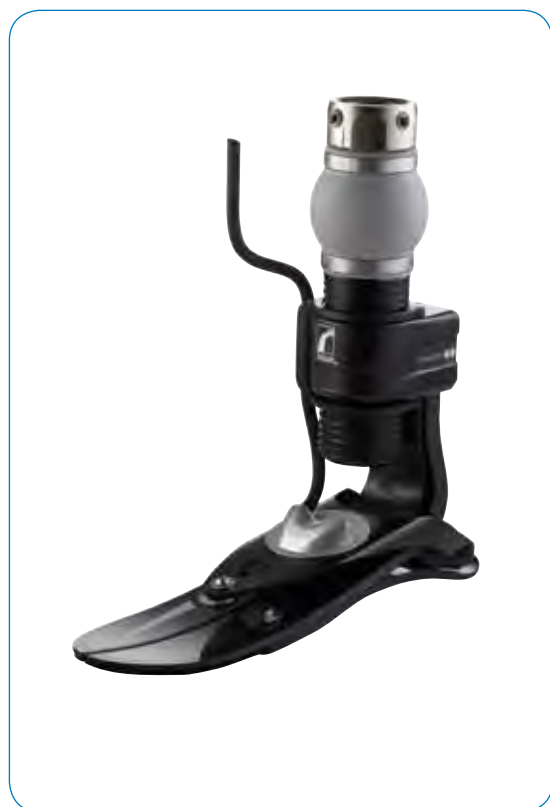

The Unity vacuum suspension system is available as an addition to the Re-Flex Rotate™ and Re-Flex Shock™ foot modules. Unity in combination with Iceross Seal-In V® provides elevated vacuum in the socket during ambulation, without the restricting effect of an external sealing sleeve. Unity adds minimal weight to the existing prosthesis and does not interfere with the superior function of Re-Flex foot and shock modules.

USER INFORMATION

Amputation Level	Transtibial
Impact Level	All


SYSTEM INFORMATION

Foot Sizes	22-30
Maximum Vacuum	-23inHg (600mmHg)
Normal Operating Vacuum:	-16-20 inHg
Unity System Weight	130g w/ Pump, Socket Valve and Tube


Ordering Information

Unity can only be ordered with Re-Flex Rotate with EVO and Re-Flex Shock with EVO in the standard foot module height. Please select the appropriate Unity Pump Module kit according to foot module size when placing an order and Unity will come preassembled on the foot. For extra convenience and system compatibility, Unity TT Valve Kit is highly recommended as a means of connecting to the socket.

Complete foot: Unity for Re-flex Rotate with male/female pyramid

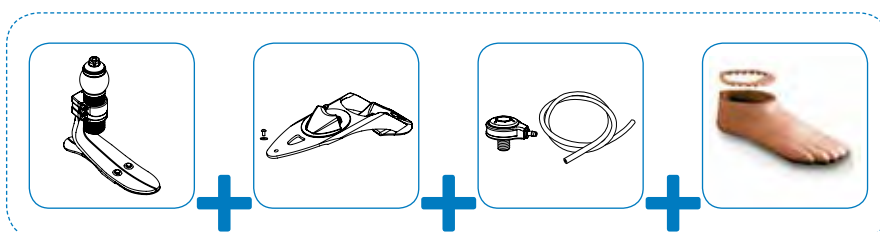

Complete foot: Unity for Re-flex Shock with male pyramid

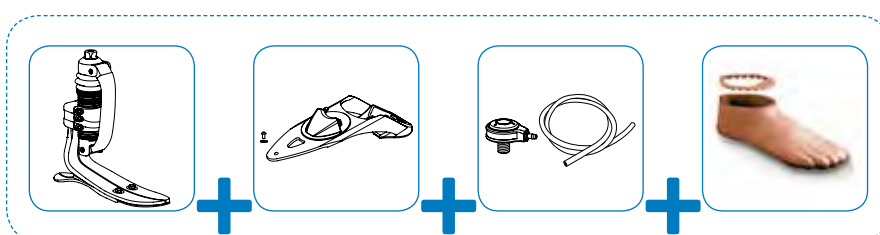

Unity® for Re-Flex™

Category Selection Guide



Weight Kg

Impact Level	44-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147	148-166
Low	1	1	2	3	4	5	6	7	8	9
Moderate	1	2	3	4	5	6	7	8	9	N/A
High	2	3	4	5	6	7	8	9	N/A	N/A

Kits

UNITY KITS

Picture	Part#	Description	Applicable Foot Kits
1	URS0024	Unity Pump Module RSX 22-24	Re-Flex Std Size 22-24
1	URS0027	Unity Pump Module RSX 25-27	Re-Flex Std Size 25-27
1	URS0030	Unity Pump Module RSX 28-30	Re-Flex Std Size 28-30
1	UTT0001	Unity TT Valve Kit	All

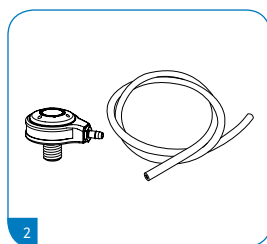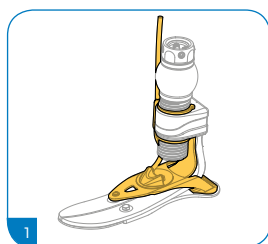


Spare Parts

PUMPS

Part#	Description
URS00011	Unity RS Spare Pump 22-24
URS00012	Unity RS Spare Pump 25-27
URS00013	Unity RS Spare Pump 28-30

Includes

- Frame
- Membrane
- Check valves
- Barb fittings
- Tube
- Tube fastener
- Heel pad

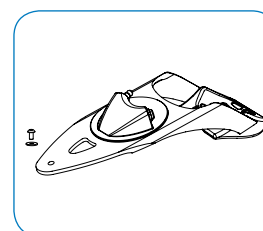

Unity® for Re-Flex®


Spare Parts

SPARE COMPONENTS

Part#	Description
URS00021	Unity RS Components 22-24
URS00022	Unity RS Components 25-30

Includes



- Membrane
- Check valves
- Barb fittings
- Tube
- Tube fastener
- Heel pad

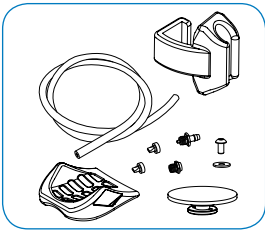

Foot Kits Available With Unity



RE-FLEX ROTATE™ WITH EVO™ STANDARD PROFILE

MALE PYRAMID ATTACHMENT

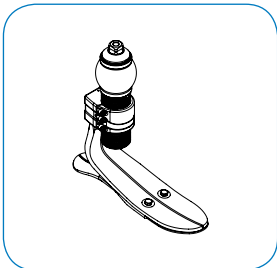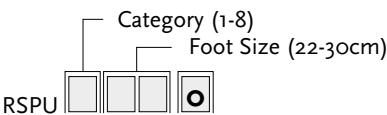



RE-FLEX ROTATE™ WITH EVO™ STANDARD PROFILE

FEMALE PYRAMID ATTACHMENT

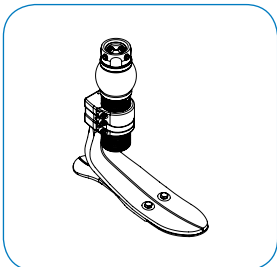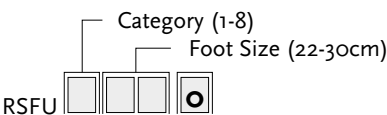

RE-FLEX SHOCK™ WITH EVO™ STANDARD PROFILE


MALE PYRAMID ATTACHMENT

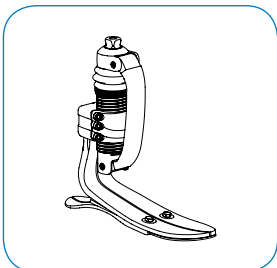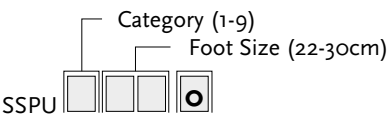

Caution: Due to restricted space with the Re-Flex Shock with EVO shock module and Unity, SSPE* kits allow for only 2 full turns of downwards adjustment from the tallest setting. Re-Flex Rotate with EVO modules can be set to any height within their designed limits.

Unity® for Re-Flex™

Cosmetic Finishing

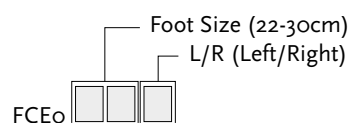

Please specify when ordering brown covers. Add "BR" to the part#.

Includes

- Foot Cover
- Flex-Foot EVO Sock
- Foot Cover Attachment

Part#	Description	Size
FCE63006	Flex-Foot EVO Sock	22-25cm
FCE63007	Flex-Foot EVO Sock	26-30cm



Unity® for LP Rotate™

Introduction



Unity for LP Rotate is suitable for people with long residual limbs which have insufficient clearance for Re-Flex Rotate or Re-Flex Shock but still prefer the benefits of a rotational shock absorption as well as the Unity sleeveless vacuum suspension system.

For the moderate to high active user, Unity for LP Rotate and Seal-In V® combined provide secure suspension, volume management, good foot dynamics and reduced rotational forces without the restricting effect of an external suspension sleeve. Unity adds minimal weight to the prosthesis and does not interfere with the foot module dynamics.

USER INFORMATION

Amputation Level	Transtibial
Impact Level	All


SYSTEM INFORMATION

Foot Sizes	22-30
Maximum Vacuum	-23inHg (600mmHg)
Normal Operating Vacuum:	-16-20 inHg
Unity System Weight	130g w/ Pump, Socket Valve and Tube

Ordering Information

LP Rotate can be ordered with Unity in the standard heel height configuration. Please select the appropriate size Unity Pump Module kit, LP Rotate shock module and foot module when placing an order and Unity will come preassembled to the LP Rotate foot. For extra convenience and system compatibility, Unity TT Valve Kit is highly recommended as a means of connecting to the socket. For liner selection and item numbers please refer to Icross Seal-In V pages.

Complete foot: Unity for LP Rotate with male/female pyramid

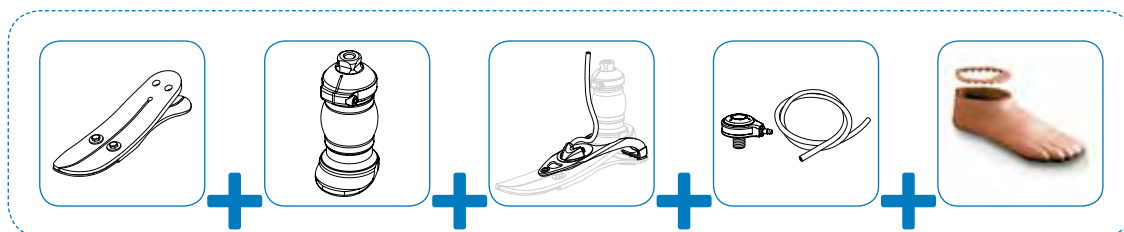

Unity® for LP Rotate™

Category Selection Guide




Weight Kg

Impact Level	45-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147
Low	1	1	2	3	4	5	6	7	8
Moderate	1	2	3	4	5	6	7	8	9
High	2	3	4	5	6	7	8	9	N/A

Kits

UNITY KITS

Picture	Part#	Description	Applicable Foot Kits
1	UVL0024	Unity Pump Module VLX 22-24	Vari-Flex with EVO Std Size 22-24
1	UVL0027	Unity Pump Module VLX 25-27	Vari-Flex with EVO Std Size 25-27
1	UVL0030	Unity Pump Module VLX 28-30	Vari-Flex with EVO Std Size 28-30
2	UTT0001	Unity TT Valve Kit	All
3	UTT0002	Unity Distal Vacuum Plate Kit	All

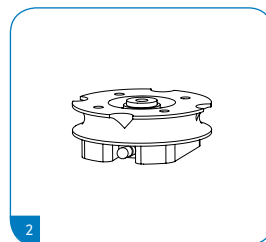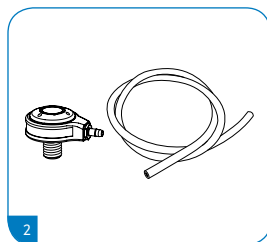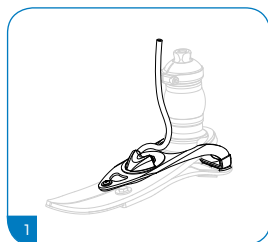

Spare Parts

SPARE COMPONENTS

Part#	Description
UVL00021	Unity VL Components 22-24
UVL00022	Unity VL Components 25-30

Includes


- Membrane
- Check valves
- Barb fittings
- Tube
- Tube fastener
- Heel pad

Unity® for LP Rotate™


Foot Kits Applicable With Unity

LP VARI-FLEX WITH EVO FOOT MODULE

PART#

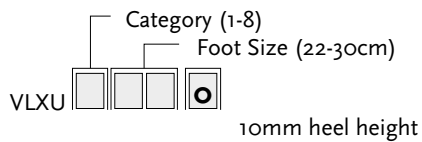

Includes:
Foot Module
Heel Module
Heel Attachment Hardware
Heel Wedges

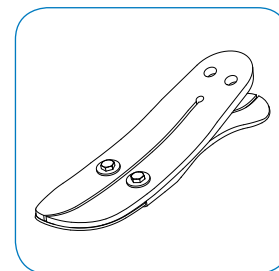


Shock Module Kits

UNITY KITS


Part#	Pyramid	Foot Size	Cat.
LRP061	Male	22-24	1-6
LRP062	Male	25-30	1-6
LRP082	Male	25-30	7-8
LRF061	Female	22-24	1-6
LRF062	Female	25-30	1-6
LRF082	Female	25-30	7-8

Includes:
• Shock module
• Male/Female adapter
• Rod Kit

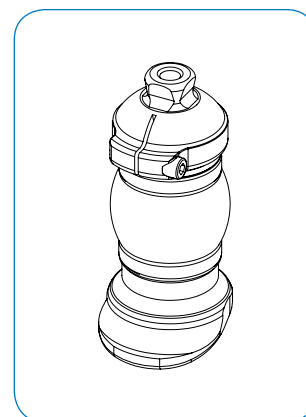

Unity® for LP Rotate™

Cosmetic Finishing

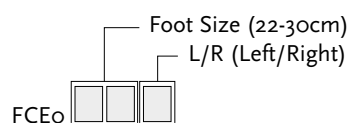


Please specify when ordering brown covers. Add "BR" to the part#.

Includes

- Foot Cover
- Flex-Foot EVO Sock
- Foot Cover Attachment

Part#	Description	Size
FCE63006	Flex-Foot EVO Sock	22-25cm
FCE63007	Flex-Foot EVO Sock	26-30cm

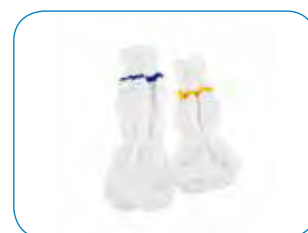

Unity® for Vari-Flex®

Introduction



Unity for Vari-Flex with EVO is an ultra lightweight foot and elevated vacuum suspension system designed to stabilise volume fluctuations and improve socket fit. For the moderate to high active user who prefers a simple, robust and lightweight prosthesis, Unity for Vari-Flex with EVO and Seal-In V® combined provide secure suspension, volume management and great foot dynamics without the restricting effect of an external suspension sleeve.

USER INFORMATION

Amputation Level	Transtibial
Impact Level	All


SYSTEM INFORMATION

Foot Sizes	22-30
Maximum Vacuum	-23inHg (600mmHg)
Normal Operating Vacuum:	-16-20 inHg
Unity System Weight	130g w/ Pump, Socket Valve and Tube

Ordering Information

Unity can be ordered with Vari-Flex with EVO. Please select the appropriate Unity Pump Module kit according to foot module size when placing an order and Unity will come preassembled on the foot. For extra convenience and system compatibility, Unity TT Valve Kit is highly recommended as a means of connecting to the socket.

For complete foot with male pyramid

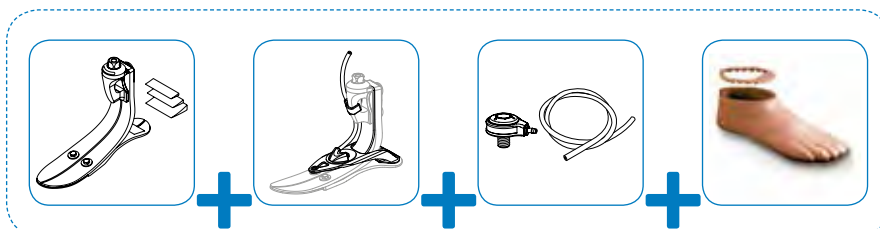

Unity® for Vari-Flex®

Category Selection Guide



Weight Kg

Impact Level	44-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147	148-166
Low	1	1	2	3	4	5	6	7	8	9
Moderate	1	2	3	4	5	6	7	8	9	Special
High	2	3	4	5	6	7	8	9	Special	Special

Kits

UNITY KITS

Picture	Part#	Description	Applicable Foot Kits
1	UVFoo24	Unity Pump Module VFX 22-24	Vari-Flex with EVO Std Size 22-24
1	UVFoo27	Unity Pump Module VFX 25-27	Vari-Flex with EVO Std Size 25-27
1	UVFoo30	Unity Pump Module VFX 28-30	Vari-Flex with EVO Std Size 28-30
1	UTT0001	Unity TT Valve Kit	All

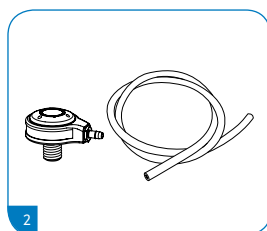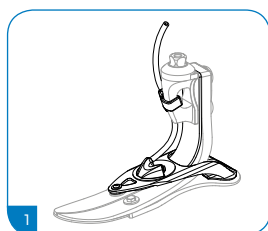


Spare Parts

PUMPS

Part#	Description
UVF00011	Unity VF Spare Pump 22-24
UVF00012	Unity VF Spare Pump 25-27
UVF00013	Unity VF Spare Pump 28-30

Includes

- Frame
- Membrane
- Check valves
- Barb fittings
- Tube
- Tube fastener
- Heel pad

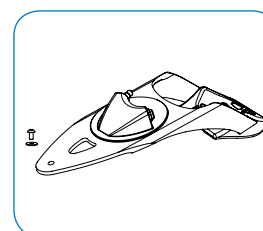

Unity® for Vari-Flex®

Spare Parts

SPARE COMPONENTS



Part#	Description
UVF00021	Unity VF Components 22-24
UVF00022	Unity VF Components 25-30

Includes

- Membrane
- Check valves
- Barb fittings
- Tube
- Tube fastener
- Heel pad



VALVE CORE

Part#	Description
UTT0005	Unity TT Valve Core

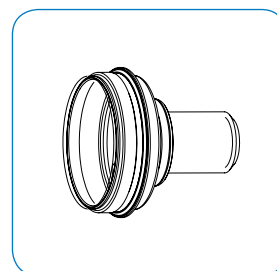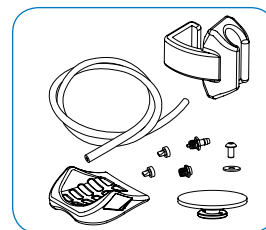


Foot Kits Applicable With Unity

VARI-FLEX WITH EVO KIT WITH MALE PYRAMID

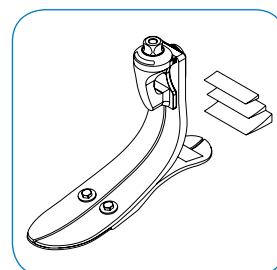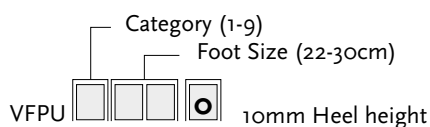

Cosmetic Finishing

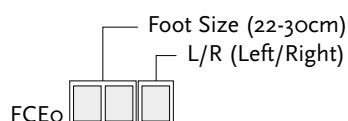

Please specify when ordering brown covers. Add "BR" to the part#.

Includes

- Foot Cover
- Flex-Foot EVO Sock
- Foot Cover Attachment

Part#	Description	Size
FCE63006	Flex-Foot EVO Sock	22-25cm
FCE63007	Flex-Foot EVO Sock	26-30cm



Unity® for LP Vari-Flex®

Introduction



Unity for LP Vari-Flex is a sleeveless vacuum suspension system designed to stabilise volume fluctuations and maintain comfortable socket fit throughout the day. For the moderate to high active user with limited build height clearance, Unity for LP Vari-Flex and Seal-In V® combined provide secure suspension, volume management and great foot dynamics without the restricting effect of an external suspension sleeve. Unity adds minimal weight to the prosthesis and does not interfere with foot function.

USER INFORMATION

Amputation Level	Transtibial
Impact Level	All


SYSTEM INFORMATION

Foot Sizes	22-30
Maximum Vacuum	-23inHg (600mmHg)
Normal Operating Vacuum:	-16-20 inHg
Unity System Weight	130g w/ Pump, Socket Valve and Tube

Ordering Information

LP Vari-Flex can be ordered with Unity and a male pyramid, standard heel height configuration. Please select the appropriate Unity Pump Module kit according to foot module size when placing an order and Unity will come preassembled on the foot. For extra convenience and system compatibility, Unity TT Valve Kit is highly recommended as a means of connecting to the socket. For liner selection charts and item numbers please refer to Iceross Seal-In V pages.

For complete foot with male pyramid

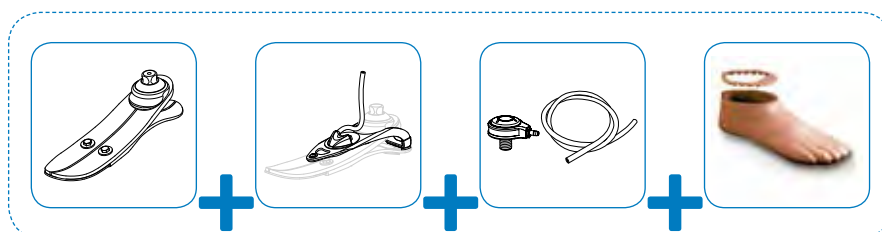

Unity® for LP Vari-Flex®

Category Selection Guide




Weight Kg

Impact Level	45-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147	148-166
Low	1	1	2	3	4	5	6	7	8	9
Moderate	1	2	3	4	5	6	7	8	9	Special
High	2	3	4	5	6	7	8	9	Special	Special

Kits

UNITY KITS

Picture	Part#	Description	Applicable Foot Kits
1	UVL0024	Unity Pump Module VLX 22-24	Vari-Flex with EVO Std Size 22-24
1	UVL0027	Unity Pump Module VLX 25-27	Vari-Flex with EVO Std Size 25-27
1	UVL0030	Unity Pump Module VLX 28-30	Vari-Flex with EVO Std Size 28-30
2	UTT0001	Unity TT Valve Kit	All
3	UTT0002	Unity Distal Vacuum Plate Kit	All

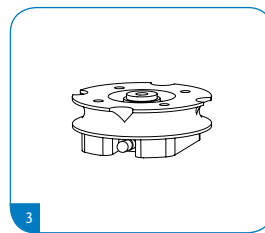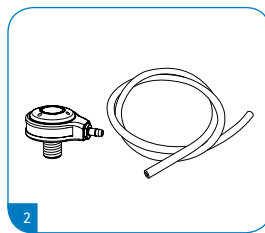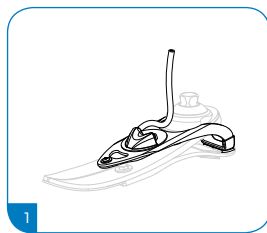

Spare Parts

SPARE COMPONENTS

Part#	Description
UVL00021	Unity VL Components 22-24
UVL00022	Unity VL Components 25-30



Includes

- Membrane
- Check valves
- Barb fittings
- Tube
- Tube fastener
- Heel pad


Unity® for LP Vari-Flex®

Foot Kits Applicable With Unity

LP VARI-FLEX WITH MALE PYRAMID

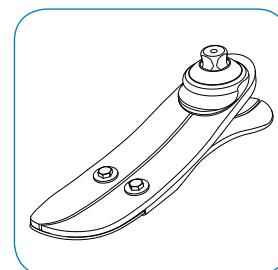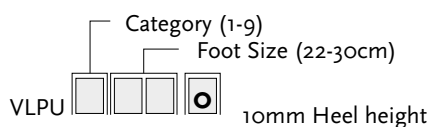

Cosmetic Finishing

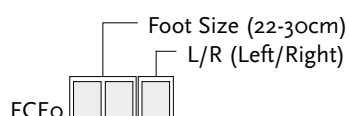

Please specify when ordering brown covers. Add "BR" to the part#.

Includes

- Foot Cover
- Flex-Foot EVO Sock
- Foot Cover Attachment



Part#	Description	Size
FCE63006	Flex-Foot EVO Sock	22-25cm
FCE63007	Flex-Foot EVO Sock	26-30cm



Sockets

Introduction

While Iceross can be used with conventional socket designs, the best outcomes are realised with a Total Surface Bearing (TSB) socket design.

The TSB design can offer several advantages over traditional socket designs including:

- Force distribution over the entire surface of the residual limb, which reduces peak pressures in the socket.
- The resulting total contact with the residual limb enhances blood circulation.
- Proprioception is improved.

Össur recommends the following casting methods for creating a transtibial TSB socket for use with Iceross:

- Hand cast or Pressure cast with plaster of Paris
- Modular Socket System direct casting using the Icecast® Anatomy Kit (for transtibial applications only)

Training courses available on request, please contact your Össur representative.

Each procedure is described in detail in the Össur Technical Manuals, which can be downloaded at www.ossur.co.uk. The success of any of the methods depends on the correct selection and application of the Iceross® liner.

Accurate, efficient socket-making

Össur works with prosthetists worldwide, developing socket-making techniques that deliver consistently good clinical outcomes. Icecast Anatomy provides additional control over soft tissue during pressure casting, reducing rectification time and ensuring repeatable, accurate results. The Modular Socket System offers a mobile socket-making system that requires only one appointment to produce the socket. Precasting adjustments for pressure points eliminate the need for rectifications.



Contents

Icecast® Anatomy	89
Icecast® Compact	90
Modular Socket System	91
Pressure Pads	93

Icecast® Anatomy

With MSS Toolkit



Icecast Anatomy represents a further development of the pressure-casting concept. Contour Cell™ shaping-chambers positioned over soft tissue areas can be inflated independently to achieve displacement of the tissue creating a more anatomically correct socket shape. The increased control over the socket shape enhances user comfort, stability and rotation control. As the residual limb is “loaded” by the pressure-casting instrument during casting, it is possible to eliminate the need for modification of the resultant socket shape.

This kit also includes all the tools required for the manufacture of MSS sockets.

Kits

Part#	Description
M-500300	Icecast MSS Tool Kit Anatomy

Includes:

- Icecast Anatomy Bladder
- Icecast Pump
- Resin Injection Tool
- Distal Attachment Tool
- Taping Ring
- Icecast Case (with attachments for trolley - sold separately)



Components

ICECAST MSS MATERIAL KIT

Part#	Description
M-500105	MSS Material Kit/5" Glass Fiber Braids
M-500107	MSS Material Kit/7" Glass Fiber Braids
M-500205	MSS Material Kit/5" Carbon Fiber Braids
M-500207	MSS Material Kit/7" Carbon Fiber Braids

Includes:

- Glass Fiber Braids (5" or 7") or Carbon Fiber Braids (5" or 7")
- Resin Cartridge
- Silicone Insulation Sheets 2x0.6m (2x2ft)
- Silicone Insulation Cap
- Molding Tube

Icecast® Compact

With MSS Toolkit



The Icecast Compact provides a simple and effective way of guaranteeing that a TSB socket fit is achieved for use with a 3mm profile Iceross silicone locking liner. The system provides a simpler and more durable socket fabrication tool than the Icecast Anatomy as it does not contain the Contour Cell technology. Socket fabrication is therefore simpler and modifications of the anatomical shape are minimal, without affecting user comfort.

This kit also includes all the tools required for the manufacture of MSS sockets.

Kits

Part#	Description
M-500301	Icecast MSS Toolkit Compact

Includes:

- Icecast Compact Bladder
- Icecast Pump
- Resin Injection Tool
- Distal Attachment Tool
- Taping Ring
- Icecast Case (with attachments for trolley)



Components

ICECAST MSS MATERIAL KIT

Part#	Description
M-500105	MSS Material Kit/5" Glass Fiber Braids
M-500107	MSS Material Kit/7" Glass Fiber Braids
M-500205	MSS Material Kit/5" Carbon Fiber Braids
M-500207	MSS Material Kit/7" Carbon Fiber Braids

Includes:

- Glass Fiber Braids (5" or 7") or Carbon Fiber Braids (5" or 7")
- Resin Cartridge
- Silicone Insulation Sheets 2x0.6m (2x2ft)
- Silicone Insulation Cap
- Molding Tube

Modular Socket System

Components

GLASS FIBRE BRAIDS WITH DISTAL ATTACHMENT

Size	1pc	10pc
5"	M-100101	M-100110
7"	M-100201	M-100210

CARBON FIBRE BRAIDS WITH DISTAL ATTACHMENT

Size	1pc	10pc
5"	M-100401	M-100410
7"	M-100501	M-100510





RESIN

Part#	Description
M-200101	Resin Cartridge (1pc)
M-200110	Resin Cartridge (10pc)
M-200150	Resin Cartridge (50pc)
M-200100	Resin Cartridge (100pc)



RELATED PRODUCTS

Picture	Part#	Description
1	L-661000	600 Modular Spacer for Icelock 600 Series*
	L-621100	Icelock Ratchet 621 Adaptation Kit
	L-631100	Icelock Lanyard 631 Adaptation Kit
	L-651100	Icelock Smooth 651 Adaptation Kit
2	A-242500	4 Hole Female Pyramid with 20/20 Offset
	L-214100	Icelock Clutch 4H 214 for Modular Socket System

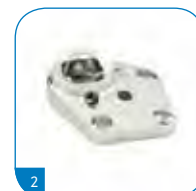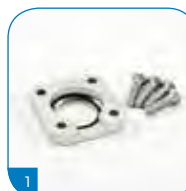


* Use with 600 series adaptation kits.

Icecast® & MSS Tooling

Components

Picture	Part#	Description
1	D-040510	Icecast Pump
2	D-040520	Icecast Anatomy Bladder
3	D-040560	Icecast Compact Bladder
4	M-400201	Resin Injection Tool
5	D-040532	Distal Attachment Tool for Icecast
6	D-040550	Distal Attachment Tool for Plaster Casting
7	X-990368	Taping Ring
8	M-500330	Trolley for Icecast Case

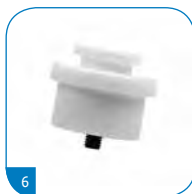

ACCESSORIES

Picture	Part#	Description
1	M-400101	Silicone Insulation Sheets (30m)
2	M-400206	Bottom Seal
3	M-400301	Silicon Insulation Sheets 2x0.6m (2x2ft)
4	M-400102	Silicone Insulation Caps (10pc)



Pressure Pads

Introduction



Pressure Pads are designed to reduce peak pressure in sockets on prominent bony areas. Made from a non-compressible silicone that optimises pressure redistribution, Pressure Pads can improve the comfort of laminated, thermoformed or directly manufactured sockets such as the Modular Socket System. The pads feature a self-adhesive backing to enable easy application.

Kits

Picture	Part#	Description
1	K-190100	Set No. 1
2	K-190101	Discs
3	K-190102	Spoon, Large
4	K-190104	Cylindrical, Large
5	K-190110	Set No. 2
6	K-190115	Oval
7	K-190150	Ergo Pads






Locks

Introduction

RATCHET PIN (Used with Iceross locking liners)

For moderate and highly active amputees with bony residual limbs who prefer the security of stepping into the prosthesis to don the socket. The ratchet pin makes an audible clicking sound when it engages, which many amputees find reassuring.

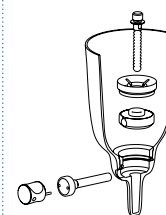


CLUTCH PIN (Used with Iceross locking liners)

For amputees of all activity levels who have fleshy residual limbs. The patient dons the prosthesis by winding into the socket rather than pushing or stepping into it.




SMOOTH PIN (Used with Iceross locking liners)

For moderate and highly active amputees with bony residual limbs who prefer the security of stepping into the prosthesis to don the socket. The lock is silent when it engages, and provides a secure suspension along the entire length of the pin.

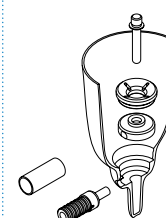


LANYARD (Used with Iceross locking liners)

For new amputees or low to moderately active amputees with fleshy residual limbs who prefer to be pulled into the socket to don the prosthesis. This suspension system allows the user to apply their prosthesis while sitting down.

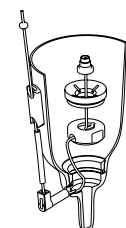

DISTAL VALVE (Designed for use with Iceross cushion liners)

For amputees of all activity levels who cannot tolerate a distal attachment or prefer suction suspension to a mechanical method. Used in conjunction with an Iceflex® suspension sleeve, this suspension system allows the user to apply their prosthesis while sitting down.

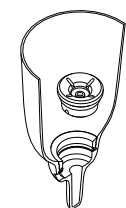

EXPULSION VALVE

For use with all Iceross Seal-In liners and for Iceross Cushion liners with a sleeve.



Contents

Icelock® 100 Series.....	97
Icelock® 200 Series.....	99
Icelock® 300 Series.....	103
Icelock® 500 Series.....	105
Icelock® 600 Series.....	106
Icelock® 700 Series.....	112


Icelock® 100 Series

Introduction



Ratchet lock with integrated pyramid.



- Interchangeable
- Easy to use - single lamination
- Can change to lanyard or valve
- Rated up to 120kg
- All impact levels

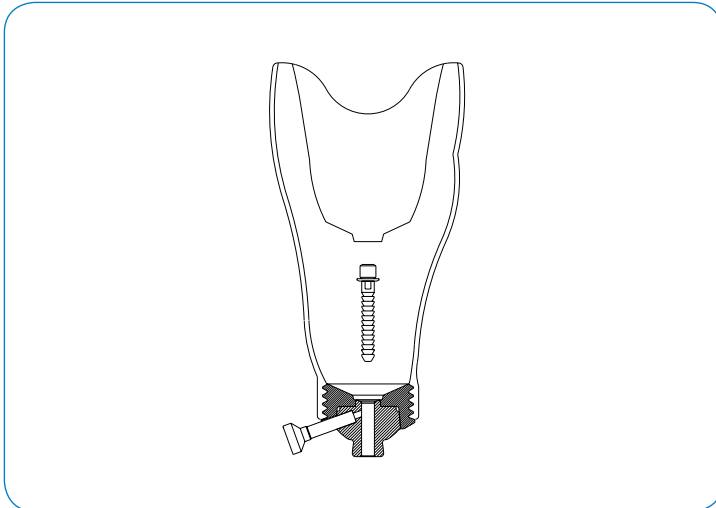

Part Number Information

ICELOCK 100 SERIES

Picture	Part#	Description	Ratings	Build Height	Weight
1	L-125000	Icelock 125 Ratchet Stainless Steel Pyramid	120kg	22mm	115g
2	L-180001	Icelock 100 Fabrication Kit without braid	N/A	22mm	99g
3	L-180002	Icelock 100 Fabrication Kit with braid	N/A	22mm	285g
4	L-192010	Lanyard Adaptation Kit			
5	L-192015	Expulsion Valve			

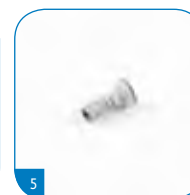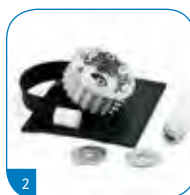

Icelock® 100 Series

Components

SPARE PARTS

Picture	Part#	Description
1	L-192000	Attachment Pin Ratchet, standard with cap, length 42mm, 13 notches
2	L-192003	Attachment Pin Ratchet, short with cap, length 22mm, 7 notches
3	L-192005	Attachment Pin Ratchet, long with cap, length 52mm, 17 notches
4	L-193144	Screw M5 x 20 for Titanium Pyramid and Stainles Steel Pyramid
5	L-192027	Release Unit Pin/Valve



FOR FABRICATION KITS

Picture	Part#	Description
1	L-192020	Puck and Ring for Puck
2	L-192031	Carbon Braid 50cm on Puck and Ring
3	L-192025	SmartCast Tool
4	L-192185	Ring for Puck



Icelock® 200 Series

Introduction




Clutch locks

- Low build height and lightweight
- Wear-resistant stainless steel pin guide
- Rated up to 100kg, 120kg and 166kg
- Easy release
- Secure clutch mechanism tested to tolerate at least a 300kg pull out load
- All impact levels

Icelock Clutch 211 (no adapter)


- Push Button can be shortened
- Pin Tunnel creates space for Attachment Pin
- Attachment Pin slides easily into Icelock Clutch 214 (aluminium 4-hole adapter)
- Push Button can be shortened
- Rotation Control on Lam/Therm Ring

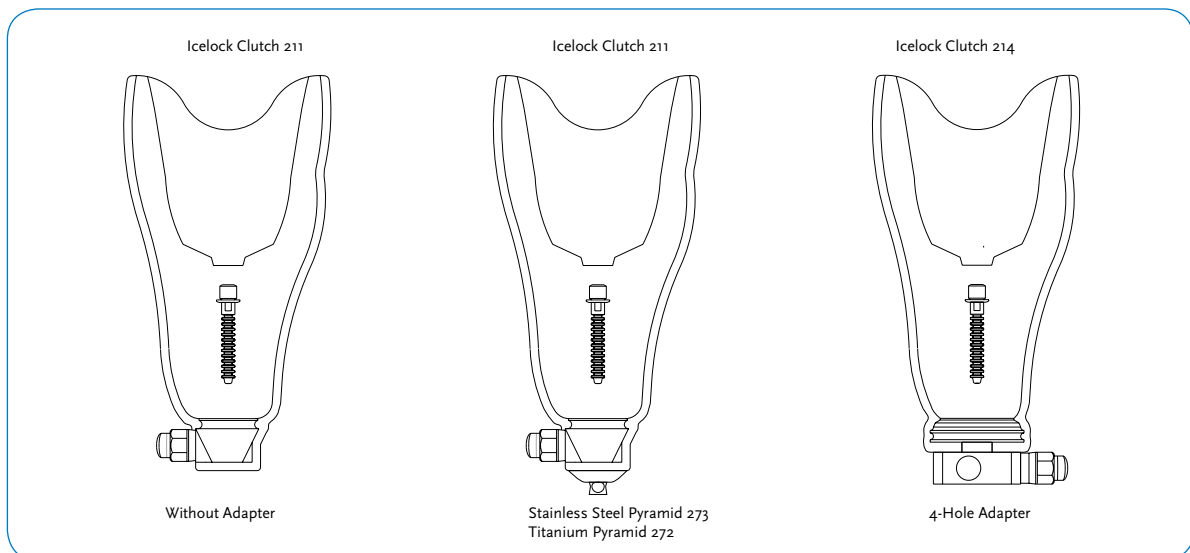

Part Number Information

ICELOCK 200 SERIES

Picture	Part#	Description	Ratings	Build Height	Weight
1	L-211000	Icelock 211 Clutch (no adapter)	N/A	21mm	78g
2	L-214000	Icelock 214 Clutch 4H	120kg	25mm	186g

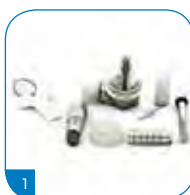

Icelock® 200 Series

Part Number Information

LAMINATION PYRAMIDS FOR ICELOCK CLUTCH 211

Picture	Part#	Description	Ratings	Build Height	Weight
1	L-273000	Icelock 273 Stainless Steel Pyramid	100kg	8mm	64g
2	L-272000	Icelock 272 Titanium Pyramid	166kg	9mm	76g




4-HOLE PYRAMIDS FOR ICELOCK CLUTCH 4H 214

Part#	Description	Weight limit
A-233100	4-Hole Male Pyramid, Aluminium	100kg
A-242100	4-Hole Female Pyramid, Aluminium	100kg
A-235100	4-Hole Male Pyramid, Titanium	166kg
A-245100	4-Hole Female Pyramid, Titanium	166kg

Components

FOR ICELOCK 211 AND 214

Picture	Part#	Description
1	L-292000	Attachment Pin Clutch, standard with cap, length 40mm, 12 notches
2	L-292003	Attachment Pin Clutch, short with cap, length 22mm, 7 notches
3	L-292005	Attachment Pin Clutch, long with cap, length 52mm, 17 notches
4	L-292020	Clutch Mechanism and T-wrench NOW WITH EASY RELEASE
5	L-293060	T-wrench
6	L-292131	Release Button

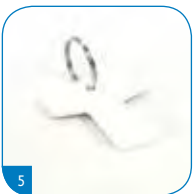


Icelock® 200 Series

Components


FOR ICELOCK 211

Picture	Part#	Description
1	L-292015	Lock Body with guide (without Dummy Screw)
2	L-292122	Dummy Screw
3	L-292126	Delrin Guide
4	L-293045	Silicone Seal
5	L-292121	Pin Tunnel with O-ring
6	L-293046	O-ring 10 x 1.5 for Pin Tunnel
7	L-693127	O-ring 8 x 2 for Pyramid Bolt
8	L-292032	Pyramid Bolt, long with O-ring

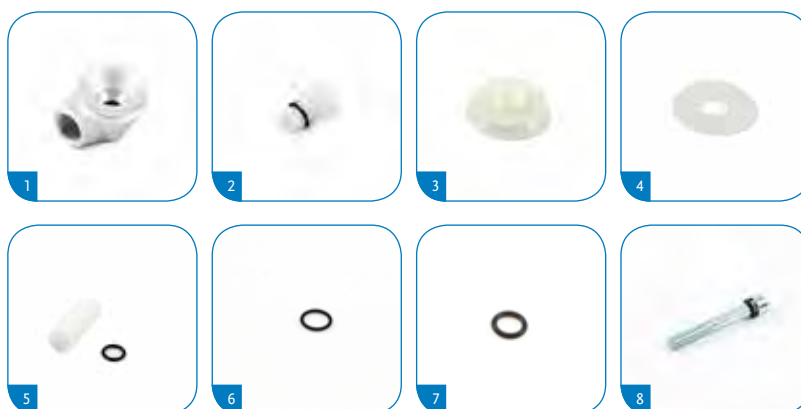

FOR PYRAMIDS 272 AND 273

Picture	Part#	Description
1	L-293046	O-ring 10 x 1.5 for Pin Tunnel
2	L-692142	Pyramid Cap
3	L-292032	Pyramid Bolt, long with O-ring (for 272 Pyramid)
4	L-692032	Pyramid Bolt with O-ring (for 273 Pyramid)
5	L-293048	Silicone Fill for 272 Pyramid
6	L-293049	Silicone Fill for 273 Pyramid
7	L-693127	O-ring 8 x 2 for Pyramid Bolts


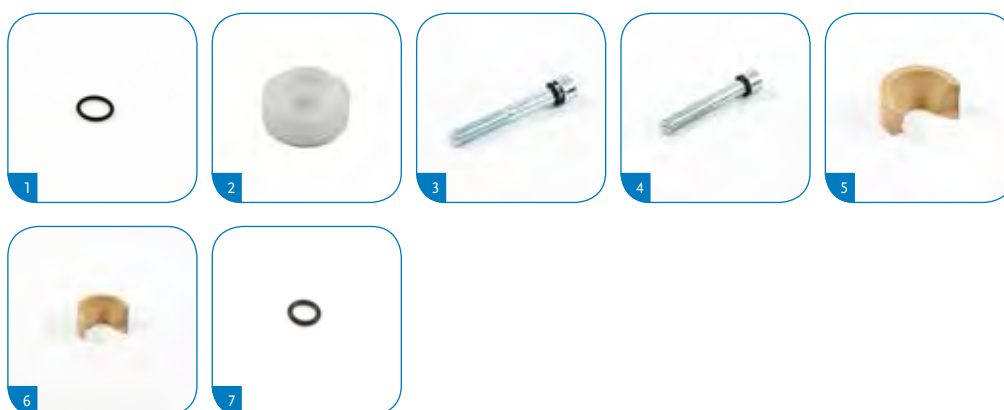


Icelock® 200 Series

Components

FOR ICELOCK 214

Picture	Part#	Description
1	L-292011	Lock Body 4-Hole with Guide (without clutch mechanism)
2	L-292040	Fabrication Kit (includes the five parts listed below)
3	L-292125	Lam/Therm Ring
4	L-292118	Delrin Tooling Body
5	L-292129	Delrin Tooling Screw
6	L-293041	Screw M6 x 20
7	L-293042	Screw M6 x 35

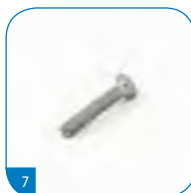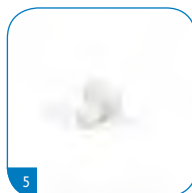

Icelock® 300 Series

Introduction



Lanyard locks without distal adapters

- Low build height
- Very lightweight
- Ideal for shower and water prostheses
- Low and moderate impact levels

Lanyard 331-L (no adapter)


- Specially designed for laminated sockets
- Silicone Dummy creates a tunnel in-between layers of laminate for the lanyard cord, for a neat cosmetic solution

Lanyard 331-T (no adapter)

- Specially designed for thermoplastic sockets

Icelock Lanyard 390 (no adapter)

- Lanyard kit with spindle locking device

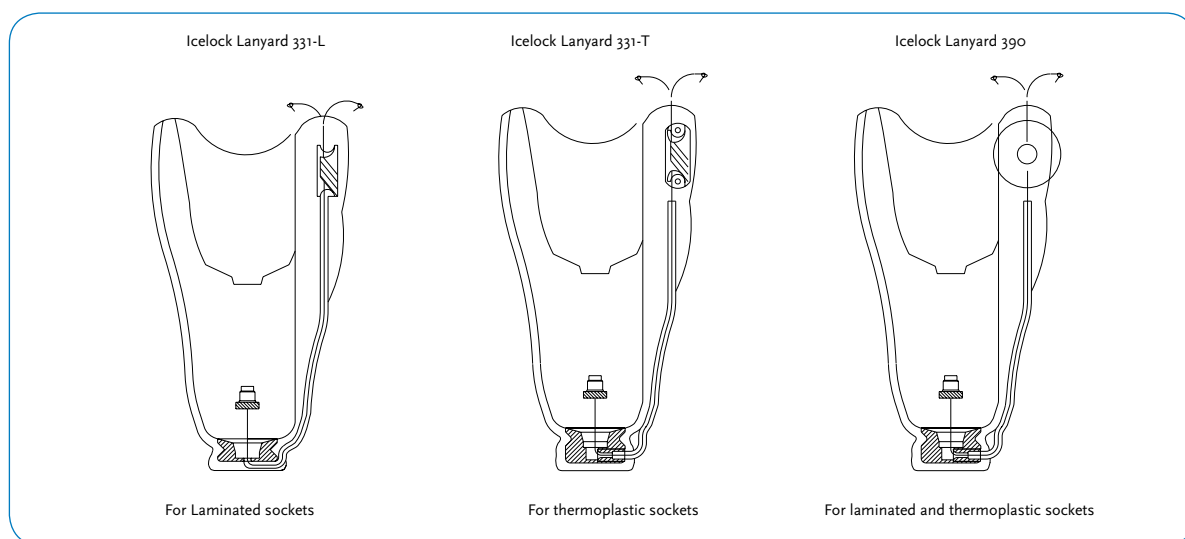

Part Number Information

ICELOCK 300 SERIES

Picture	Part#	Description	Build Height	Weight
1	L-331001	Icelock Lanyard 331-L	12mm	22g
2	L-331002	Icelock Lanyard 331-T	14mm	42g
3	K-390000	Icelock Lanyard 390	25mm	52g



Icelock® 300 Series

Components

FOR L-331001

Picture	Part#	Description
1	L-392012	Cleat with silicone dummy
2	L-392016	Pull Handle, lanyard
3	L-392018	Lanyard Cord with connecting screw
4	L-392110	Bottom Attachment, lamination




FOR L-331002

Picture	Part#	Description
1	L-392010	Cleat with cover and rivets
2	L-393012	Cord Cover 40cm, lanyard
3	L-392016	Pull Handle, lanyard
4	L-392112	Bottom Attachment, thermoplastic
5	L-392122	Delrin Guide
6	L-393018	Screw 3.5 x 30mm



FOR K-390000

Picture	Part#	Description
1	K-390001	Spindle Adapter
2	K-390002	Nut
3	K-390003	Washer
4	K-390004	Bottom Attachment
5	K-390005	Distal Connector
6	K-390007	Spindle

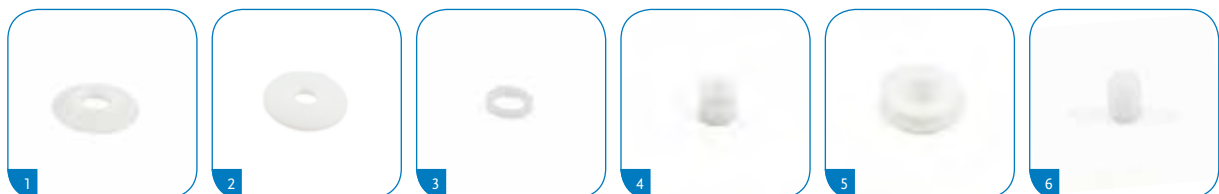

Icelock® 500 Series

Introduction



Expulsion Valve for transtibial and transfemoral applications. Recommended for use with Iceross Seal-In & Cushion liners.

Icelock Expulsion Valve 551

- Icelock 551 is an auto-expulsion valve recommended for both Transtibial and Transfemoral users.



Icelock Expulsion Valve 552

- Icelock 552 is an auto-expulsion valve recommended for Transfemoral users.
- The core of the valve can be unscrewed to facilitate rapid donning.

Part Number Information

ICELOCK 500 SERIES

Picture	Part#	Description
1	L-551002	Icelock Expulsion Valve 551
2	L-552000	Icelock Expulsion Valve 552

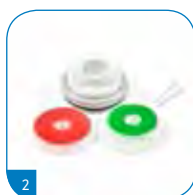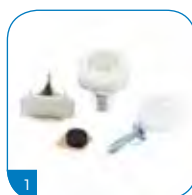

FOR ICELOCK 551

Part#	Description
L-551206	Icelock Replacement Membrane x5

FOR ICELOCK 552


Part#	Description
L-552206	Icelock Replacement Membrane x5

Icelock® 600 Series

Introduction



- Extremely low build height
- Fabrication Tools make the fabrication process very effective
- Very strong and durable
- Socket Adapters rated up to 100kg and 166kg
- Easily interchangeable
- All impact levels
- Removable from the inside of the socket

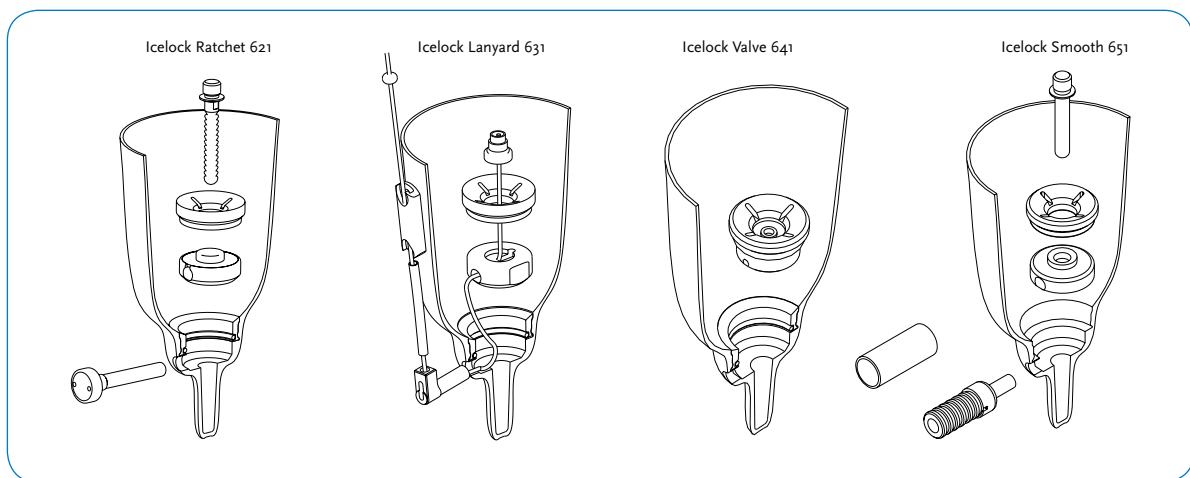


Part Number Information

ICELOCK 600 SERIES

Picture	Part#	Description	Build Height	Weight
1	L-621000	Icelock 621 Ratchet (no adapter)	15mm	59g
2	L-631000	Icelock 631 Lanyard (no adapter)	15mm	49g
3	L-641000	Icelock 641 Valve (no adapter)	15mm	40g
4	L-651000	Icelock 651 Smooth (no adapter)	15mm	68g

Use with L-680000 Fabrication Tooling Kit

Picture	Part#	Description
5	L-680000	Icelock 680 Fabrication Tools (reusable)

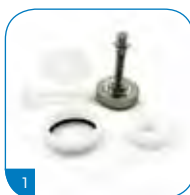

Icelock® 600 Series

Part Number Information

LAMINATION PYRAMIDS FOR ICELOCK 600


Picture	Part#	Description	Weight limit
1	L-672000	Icelock 672 Titanium Pyramid	166kg
2	L-673000	Icelock 673 Stainless Steel Pyramid	100kg
3	L-674000	Icelock 674 Aluminium 4-Hole	100kg
4	L-675000	Icelock 675 Titanium 4-Hole	166kg
5	L-692016	Release Key	



Components

FOR ICELOCK 621

Picture	Part#	Description
1	L-692022	Icelock 600 Fabrication Ring with O-ring and silicone ring
2	L-693125	O-ring 35 x 2.5mm for Fabrication Ring
3	L-692020	Ratchet Lock Body
4	L-692110	Funnel
5	L-692027	Release Unit, long
6	L-192000	Attachment Pin Ratchet, standard with cap, length 42mm, 13 notches
7	L-192003	Attachment Pin Ratchet, short with cap, length 22mm, 7 notches
8	L-192005	Attachment Pin Ratchet, long with cap, length 52mm, 17 notches
9	L-693112	Silicone Seal

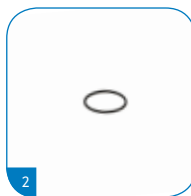

Icelock® 600 Series

Components


FOR ICELOCK 631

Picture	Part#	Description
1	L-692022	Fabrication ring with O-ring and silicone ring
2	L-693125	O-ring 35 x 2.5mm for fabrication ring
3	L-692025	Bottom Attachment, lanyard
4	L-692110	Funnel
5	L-392018	Lanyard Cord with connecting screw
6	L-393012	Cord Cover 40cm
7	L-392010	Cleat with cover and rivets
8	L-392016	Pull Handle
9	L-693112	Silicone Seal
10	L-392117	Elbow lanyard



FOR ICELOCK 641

Part#	Description
L-692022	Fabrication ring with O-ring and silicone ring
L-693125	O-ring 35 x 2.5 for fabrication ring
L-641100	Expulsion Valve/Adaptation Kit 641
L-693112	Silicone Seal

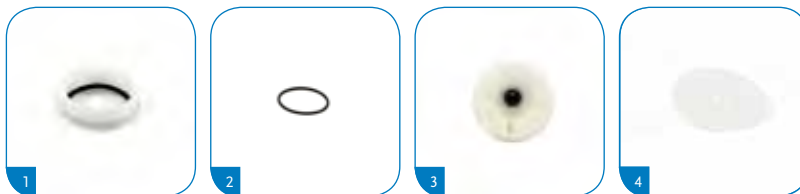

Icelock® 600 Series

Components


FOR ICELOCK 651

Picture	Part#	Description
1	L-692022	Fabrication ring with O-ring and silicone ring
2	L-692110	Funnel
3	L-693506	Smooth Pin, 43mm
4	L-693507	Smooth Pin, standard 53mm
5	L-693508	Smooth Pin, 63mm
6	L-693509	Smooth Pin, 75mm
7	L-692505	Smooth Lock
8	L-693504	Protective sleeve for push button
9	L-692501	Release Unit
10	L-693502	Push Button
11	L-693503	Dust Cap
12	L-693112	Silicone Seal



FOR ICELOCK TITANIUM PYRAMID 672

Picture	Part#	Description
1	L-692135	Titanium Pyramid
2	L-692142	Pyramid Cap
3	L-692032	Pyramid Bolt with O-ring
4	L-693127	O-ring 8 x 2 for Pyramid Bolt
5	L-693115	Silicone Ring
6	L-693112	Silicone Seal



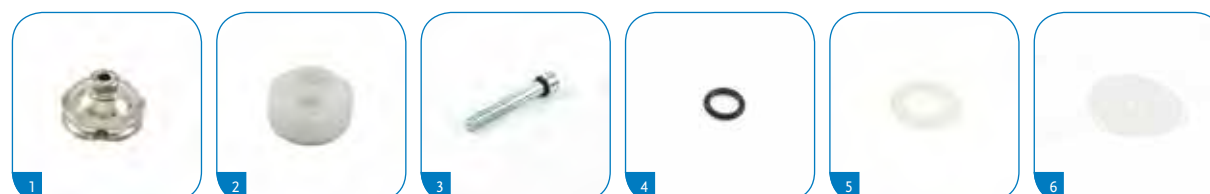


Icelock® 600 Series

Components


FOR ICELOCK STAINLESS STEEL PYRAMID 673

Picture	Part#	Description
1	L-692030	Stainless Steel Pyramid
2	L-692142	Pyramid Cap
3	L-692032	Pyramid Bolt with O-ring
4	L-693127	O-ring 8 x 2 for Pyramid Bolt
5	L-693115	Silicone Ring
6	L-693112	Silicone Seal

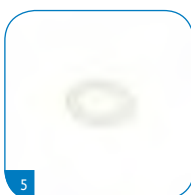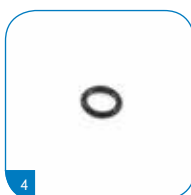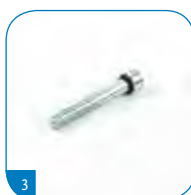



FOR ICELOCK ALUMINIUM PYRAMID 674

Picture	Part#	Description
1	L-293041	Screw M6 x 20mm
2	L-693115	Silicone Ring
3	L-693112	Silicone Seal

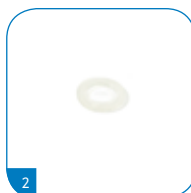

FOR ICELOCK TITANIUM PYRAMID 675

Picture	Part#	Description
1	L-693120	Screw M6 x 20mm
2	L-693115	Silicone Ring
3	L-693112	Silicone Seal



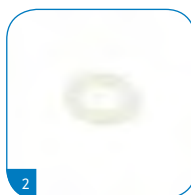
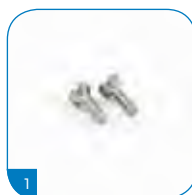


Icelock® 600 Series

Components

FOR ICELOCK FABRICATION TOOLS 680 (ALL HARDENED STAINLESS STEEL)

Picture	Part#	Description
1	L-692010	Lock Dummy with O-ring
2	L-693127	O-ring 8 x 2mm for Lock Dummy
3	L-692119	Steel Guide
4	L-692120	Release Dummy, standard
5	L-692125	Release Dummy, long
6	L-692130	Pin Dummy, standard
7	L-692133	Pin Dummy, short
8	L-693110	Pin Dummy Bolt
9	L-692015	Hex Driver with O-ring
10	L-692017	Cross Driver with O-ring
11	L-693126	O-ring 4 x 2mm for Hex and Cross Drivers

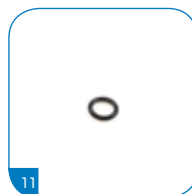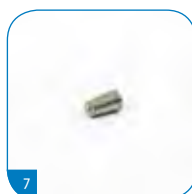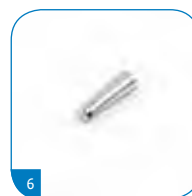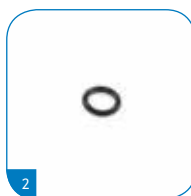

Icelock® 700 Series


Introduction



Ratchet lock designed to be used in conjunction with Iceross Upper-X® and Iceross® Junior.

- Extremely low build height
- Very strong and durable
- Removable from inside of socket
- Reusable fabrication tools

The Attachment Pin features M6 threads that are unique for Iceross Upper-X and Junior.

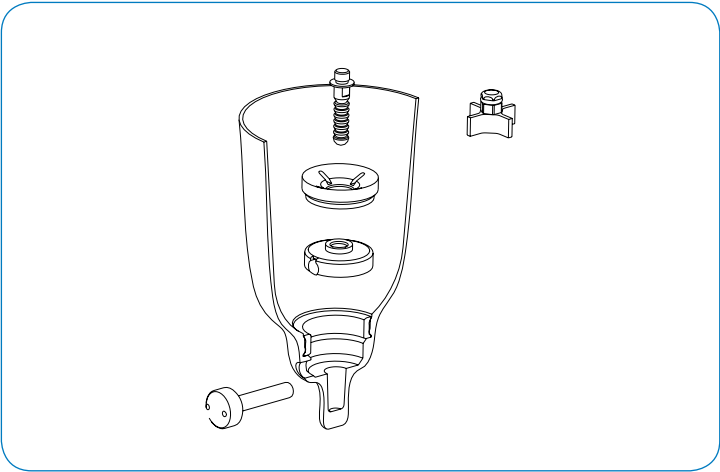

Part Number Information

ICELOCK 700 SERIES

Picture	Part#	Description	Build Height	Weight
1	L-721000	Icelock 721 Ratchet, (no adapter)	13mm	20g
2	L-780000	Icelock Fabrication Tools 780, (reusable)		




Icelock® 700 Series

Part Number Information

SPECIALLY DESIGNED PYRAMIDS FOR ICELOCK 700

Picture	Part#	Description	Weight limit
1	L-773000	Icelock 773 Stainless Pyramid	45kg
2	L-774000	Icelock 774 Aluminium 4-hole	45kg
3	L-792016	Release key to fit or remove the release unit	

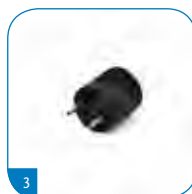

Components

FOR RATCHET LOCK 721

Picture	Part#	Description
1	L-792150	Fabrication Ring, small
2	L-792020	Ratchet Lock without Adapter, small
3	L-792110	Funnel, small
4	L-792027	Release Unit, small
5	L-792031	Pin Adapter, M10/M6
6	L-792100	Attachment Pin Ratchet, 32mm
7	L-792105	Attachment Pin Ratchet, 56mm
8	L-793112	Silicone Seal, small



Icelock® 700 Series

Components

FOR FABRICATION TOOLING (ALL HARDENED STAINLESS STEEL)

Picture	Part#	Description
1	L-792119	Steel Guide, small
2	L-792010	Lock Dummy with O-ring, small
3	L-792120	Release Dummy, small
4	L-792130	Pin Dummy, small
5	L-793110	Pin Dummy Bolt, small
6	L-692015	Hex Driver with O-ring
7	L-792017	Cross Driver with O-ring, small





Adapters

Introduction

For sockets, knees and feet, we make the strongest, lightweight adapters to help you assemble the best solution every time. Guaranteed to connect all Össur products.

Contents

Standard Adapters, Rated to 100kg	117
Standard Adapters, Rated to 166kg	120
Junior Adapters, Rated to 45kg	123
Total Shock™	124

Össur Standard Adapters

Rated to 100kg

SOCKET ADAPTERS

Picture	Part#	Description	Material	Height	Build Ht.	Weight
1	A-122100	4-Hole Socket Adapter	Aluminium	18.5mm	18.5mm	70g
2	A-134300	Male Pyramid Insert	Stainless Steel	31mm	7.5mm	94g
3	A-144300	Female Pyramid Insert	Stainless Steel	23mm	14mm	38g
4	A-114040	4-Prong Socket Adapter (Rated to 166kg)	Titanium		8mm	154g
5	A-114030	3-Prong Socket Adapter (Rated to 166kg)	Titanium		8mm	124g



4-HOLE ADAPTERS

Picture	Part#	Description	Material	Height	Build Ht.	Weight
1	A-233100	4-Hole Male Pyramid, with screws	Aluminium/Stainless Steel	23.5mm	9mm	62g
2	A-242100	4-Hole Female Pyramid, with screws	Aluminium	16mm	16mm	46g



TUBE CLAMPS

Picture	Part#	Description	Material	Height	Build Ht.	Weight
1	A-333100	Male Pyramid Tube Clamp	Aluminium/Stainless Steel	56mm	41.5mm	62g
2	A-342100	Female Pyramid Tube Clamp	Aluminium	50mm	16mm	64g
3	A-322100	4-Hole Tube Clamp, with screws	Aluminium	38mm	4mm	52g



Össur Standard Adapters

Rated to 100kg

SINGLE ADAPTERS

Picture	Part#	Description	Material	Height	Build Ht.	Weight
1	A-433110	Male Single Adapter - Short	Aluminium/Stainless Steel	50mm	35.5mm	60g
2	A-433120	Male Single Adapter - Long	Aluminium/Stainless Steel	77mm*	63mm	86g
3	A-442100	Female Single Adapter	Aluminium	100mm**	100mm	86g

*Can be cut to 70mm

**Can be cut to 60mm



DOUBLE ADAPTER


Picture	Part#	Description	Material	Height	Build Ht.	Weight
1	A-542110	Female Double Adapter	Aluminium	32mm	32mm	66g
2	A-542120	Female Double Adapter	Aluminium	45mm	45mm	76g
3	A-542130	Female Double Adapter	Aluminium	60mm	60mm	88g
4	A-542140	Female Double Adapter	Aluminium	75mm	75mm	100g
5	A-533100	Male Double Adapter	Aluminium/Stainless Steel	41mm	12mm	84g



PYLONS AND TUBES

Picture	Part#	Description	Material	Length
1	A-712010	Tube - Short	Aluminium	250mm
1	A-712020	Tube - Long	Aluminium	400mm
2	A-742010	Female Pylon - Short	Aluminium/Stainless Steel	250mm
2	A-742020	Female Pylon - Long	Aluminium/Stainless Steel	400mm
3	A-762111	Female Pylon Kit 1- Short*	Aluminium/Stainless Steel	250mm
3	A-762121	Female Pylon Kit 1- Long*	Aluminium/Stainless Steel	400mm

*Tube and female pyramid adapter

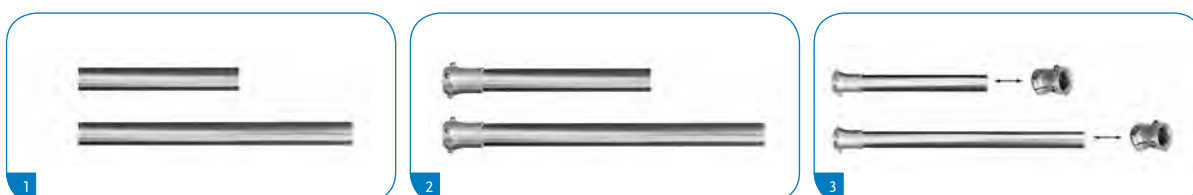

Össur Standard Adapters

Rated to 100kg

FOOT ADAPTER

Picture	Part#	Description	Material	Height	Build Ht.	Weight
1	A-833100	SACH Foot Adapter	Aluminium/Stainless Steel	25mm	9.5mm	105g



SHIFT ADAPTERS

Picture	Part#	Description	Material	Height	Build Ht.	Weight
1	A-342430	Female Pyramid Tube Clamp 14mm	Aluminium	75mm	34mm	135g
2	A-645210	Shift Adapter - Short	Titanium	60mm	60mm	208g



4-HOLE ADAPTERS

Picture	Part#	Description	Material	Height	Weight
1	A-242500	4-Hole Female Pyramid with 20/20 Offset	Aluminium	18mm	118g



Össur Standard Adapters

Rated to 166kg

SOCKET ADAPTERS

Picture	Part#	Description	Material	Height	Build Ht.	Weight
1	A-135100	Male Pyramid Socket Adapter	Titanium	23mm	8.5mm	58g
2	A-125100	4-Hole Socket Adapter	Titanium	18.5mm	18.5mm	112g
3	A-114040	4-Prong Socket Adapter	Titanium		8mm	154g
4	A-114030	3-Prong Socket Adapter	Titanium		8mm	124g
5	A-135300	Male Pyramid Insert Prong	Titanium	31mm	7.5mm	56g
6	A-145300	Female Pyramid Insert Prong	Titanium	23mm	14mm	53g
7	A-145310	Female Pyramid Insert Prong (for long residual limbs)	Titanium	9mm	9mm	33g



4-HOLE ADAPTERS, WITH SCREWS

Picture	Part#	Description	Material	Height	Build Ht.	Weight
1	A-235100	4-Hole Male Pyramid	Titanium	23mm	8.5mm	54g
2	A-235300	4-Hole Male Pyramid With Rotation	Titanium	25mm	10.5mm	72g
3	A-235400	4-Hole Male Pyramid With 8mm Offset	Titanium	23mm	8.5mm	65g
4	A-245100	4-Hole Female Pyramid	Titanium	15.5mm	15.5mm	64g
5	A-245300	4-Hole Female Pyramid With Rotation	Titanium	17.5mm	17.5mm	76g



Össur Standard Adapters

Rated to 166kg

TUBE CLAMPS

Picture	Part#	Description	Material	Height	Build Ht.	Weight
1	A-325100	4-Hole Tube Clamp, with screws	Titanium	38mm	2mm	68g
2	A-335100	Male Pyramid Tube Clamp	Titanium	56mm	5.5mm	82g
3	A-345100	Female Pyramid Tube Clamp	Titanium	50mm	14mm	78g



SINGLE ADAPTERS

Picture	Part#	Description	Material	Length	Build Ht.	Weight
1	A-435110	Male Single Adapter - Short	Titanium	50mm	35.5mm	80g
2	A-435120	Male Single Adapter - Long	Titanium	77mm	63mm	103g
3	A-445100	Female Single Adapter	Titanium	100mm	100mm	118g



DOUBLE ADAPTERS

Picture	Part#	Description	Material	Height	Build Ht.	Weight
1	A-545110	Female Double Adapter	Titanium	32mm	32mm	90g
2	A-545120	Female Double Adapter	Titanium	45mm	45mm	102g
3	A-545130	Female Double Adapter	Titanium	60mm	60mm	119g
4	A-545140	Female Double Adapter	Titanium	75mm	75mm	135g
5	A-535100	Male Double Adapter	Titanium	41mm	12mm	98g
6	A-554700	Male-Female Double Adapter	Stainless Steel	40mm	24mm	166g




Össur Standard Adapters

Rated to 166kg

PYLONS AND TUBES


Picture	Part#	Description	Material	Length
1	A-717010	Tube - Short	Carbon	250mm
1	A-717020	Tube - Long	Carbon	400mm
2	A-746010	Female Pylon - Short	Titanium/Carbon	250mm
2	A-746020	Female Pylon - Long	Titanium/Carbon	400mm

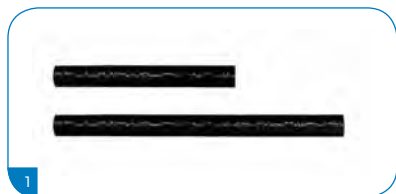

PYLONS AND TUBES

Picture	Part#	Description	Material	Length
1	A-766111	Female Pylon Kit 1 - Short*	Titanium/Carbon	250mm
1	A-766121	Female Pylon Kit 1 - Long*	Titanium/Carbon	400mm
2	A-766112	Female Pylon Kit 2 - Short**	Carbon	250mm
2	A-766122	Female Pylon Kit 2 - Long**	Carbon	400mm
3	A-718009	Tube Insert for Category 7-9	Delrin	70mm

*Kit contains Female Pylon & Female Tube Clamp. **Also contains 4-Prong Male Socket Adapter.

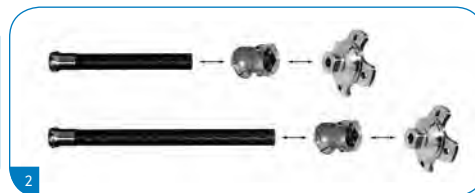

FOOT ADAPTER

Picture	Part#	Description	Material	Height	Build Ht.	Weight
1	A-835100	SACH Foot Adapter	Titanium	25mm	9.5mm	110g



Össur Junior Adapters

Rated to 45kg

SOCKET ADAPTERS

Picture	Part#	Description	Material	Build Ht.	Weight
1	AJ-114040	4-Prong Socket Adapter	Stainless Steel	9.5mm	79g
2	AJ-122100	4-Hole Socket Adapter	Aluminium	14mm	46g
3	AJ-133300	Male Pyramid Insert for Prong	Aluminium	5.5mm	36g
4	AJ-142300	Female Pyramid Insert for Prong	Stainless Steel	12mm	36g



4-HOLE ADAPTERS, WITH SCREWS

Picture	Part#	Description	Material	Build Ht.	Weight
1	AJ-233100	4-Hole Male Pyramid	Aluminium/Stainless Steel	6mm	39g
2	AJ-242100	4-Hole Female Pyramid	Aluminium	14mm	34g



TUBE CLAMPS AND FOOT ADAPTERS

Picture	Part#	Description	Material	Build Ht.	Weight
1	AJ-342100	Female Pyramid Tube Clamp	Aluminium	14mm	48g
2	1158	4-Hole Tube Clamp, with screws	Aluminium		56g
3	AJ-833100	SACH Foot Adapter	Aluminium/Stainless Steel	6mm	54g



TUBES

Picture	Part#	Description	Material	Length	Weight
1	AJ-712010	Tube Junior, 22mm	Aluminium	250mm	104g



Total Shock™

Introduction



- Vertical compression and axial rotation (10° each direction)
- Independent adjustments of compression and rotation
- Female Pyramid or tube clamp distal adapters
- Maximum patient weight 100kg



Kits

Picture	Part#	Weight limit	Height	Weight
1	4400	100kg	139mm	270g
2	4410	100kg	145mm	300g

Includes

- Shock
- Lubricant
- Rods



Components

ROD KIT

Picture	Part#	Colour	Stiffness
1	4402	Black	Normal
2	4403	Red	Firm
3	4404	Yellow	X-Firm
4	4405	Orange	XX-Firm

Includes

- Vertical Compression Rod
- Rotation Rod
- O-rings
- Retaining Ring

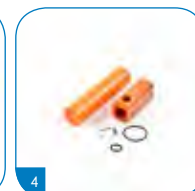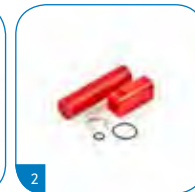

Vertical Compression Rods Part#	Rotation Rods Part#	Colour	Stiffness
4402-V	4402-R	Black	Normal
4403-V	4403-R	Red	Firm
4404-V	4404-R	Yellow	X-Firm
4405-V	4405-R	Orange	XX-Firm

Part#	Description
4400-013	Seal for Total Shock Housing
4401	Centra Lube™ Oil
4401-G	Centra Lube™ Grease*

* Use only when Total Shock is used “upside-down” (opposite to illustrations above)

All Total Shock units are lubricated with Centra Lube Oil.

Note: Never mix oil and grease.



Knees

Introduction

Total Knee®

The uniqueness of the Total Knee is in the “locking moment”, which keeps the knee stable when extended. Its polycentric design imitates natural knee motion, while the adjustable stance-flex bumper acts as a shock absorber, simulating the flexing action that occurs in the normal knee while walking or running.

Balance™ Knee

At the start of the rehabilitation process, most transfemoral amputees require maximum stability. Later, as they gain confidence, the goal is to recover a more natural gait. That's why Össur's latest mechanical knee can be easily adjusted to optimise the required balance between stability and walking dynamics for each individual.

Simple to fit and adjust, the Balance Knee is the ideal companion on the road to rehabilitation for the less active user with limited voluntary control.

Mauch® Knee

Mauch's progressively designed hydraulic systems provide amputees with a less strenuous and more controlled means of ambulation on varied terrains and in sports activities. From running, to walking step over step down stairs, Mauch offers a smooth natural gait and a high level of flexibility for more active lifestyles.

Product Overview



Contents

Mauch® Knee.....	127
Mauch® Knee Plus.....	129
Balance® Knee	131
Locking Knee.....	133
Total Knee® 1900.....	134
Total Knee® 2000.....	137
Total Knee® 2100.....	140
Total Knee® Junior	143

Mauch® Knee


Introduction



- Single axis hydraulic knee system with swing and stance control (SNS®)
- Designed for multi-speed ambulation
- Increased durability and more secure motion
- Well suited for low to high activity users
- Flexion Stop at 115°
- Ultra durable aluminium frame
- Weight limit of 136kg
- Roller bearings at knee axis and cylinder attachment points for smooth action
- Low resistance knee kits and swing only knee kits are available by special order
- Yielding stance control
- Mode selector switch allows manual locking and free swing functions
- Smooth profile assists with cosmetic finishing
- Weight of unit is 1140g



Dimensions

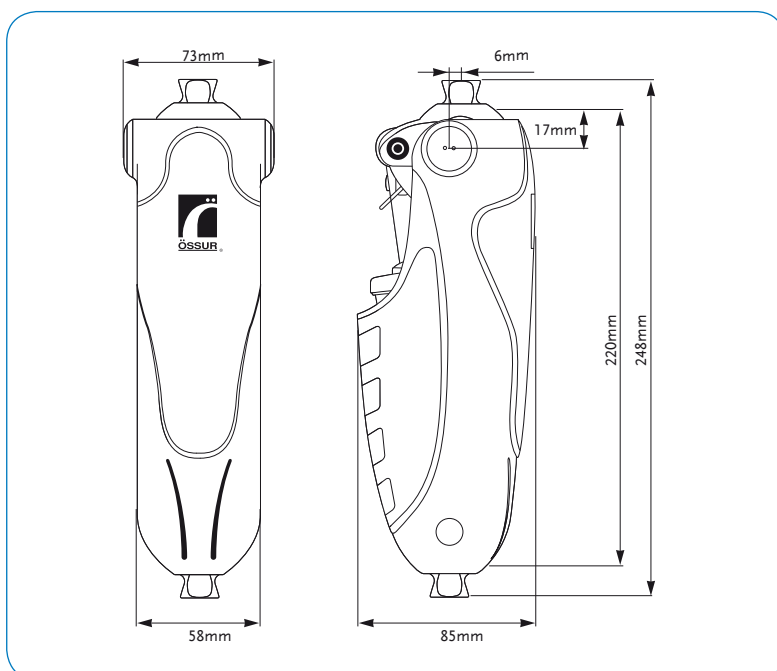

Mauch® Knee

Kits

Part#	Description
MKNo1360	Mauch Knee kit SNS
MKNo1361	Mauch Knee kit SNS - Low Resistance (Special Order)
MKNo1362	Mauch Knee kit Swing Only (Special Order)
MKNo1363	Mauch Knee kit SNS - High Stance Resistance

Components

Picture	Part#	Description
1	MKAo1361	Knee Ball, Screw
2	MKAo1362	Front Cover
3	MKAo1363	Rear Cover
4	MKAo1364	Side Plug
	MKAo1365	Soft Extension Stop
	MKAo1366	Hard Extension Stop




Mauch® Knee Plus

Introduction



- Single axis hydraulic knee system with swing and stance control (SNS®)
- Designed for multi-speed ambulation
- Yielding stance control
- Mode selector switch allows manual locking and free swing functions
- Ultra durable aluminium frame
- Increased durability and more secure motion
- Well suited for low to high activity users
- Weight limit of 166kg
- Flexion Stop at 125°
- Roller bearings at knee axis and cylinder attachment points for smooth action
- Low resistance knee kits and swing only knee kits are available by special order
- Universal 4-hole mounting patterns for versatile adapter options
- Option for through hole or threaded hole top bracket
- Weight of unit is 1210g

Clearance

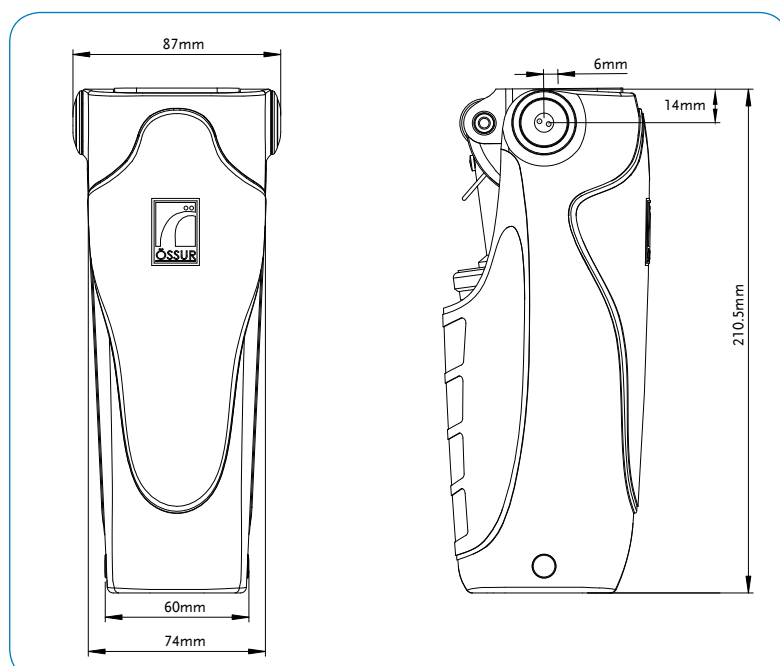

Mauch® Knee Plus

Kits

Part#	Description
MKNo1660	Mauch Knee Plus kit SNS
MKNo1661	Mauch Knee Plus kit SNS - Low Resistance (Special Order)
MKNo1662	Mauch Knee Plus kit Swing Only (Special Order)
MKNo1663	Mauch Knee Plus kit SNS - Through Hole (Special Order)
MKNo1664	Mauch Knee Plus kit SNS - Low Resistance - Through Hole (Special Order)
MKNo1665	Mauch Knee Plus kit Swing Only - Through Hole (Special Order)
MKNo1666	Mauch Knee Plus kit SNS - High Stance Resistance
MKNo1667	Mauch Knee Plus kit SNS - High Stance Resistance - Through Hole (Special Order)


Components

Picture	Part#	Description
1	MKAo1661	Knee Ball, Screw
2	MKAo1662	Front Cover
3	MKAo1663	Rear Cover
4	MKAo1664	Side Plug
	MKAo1365	Soft Extension Stop
	MKAo1366	Hard Extension Stop




Balance™ Knee

Introduction

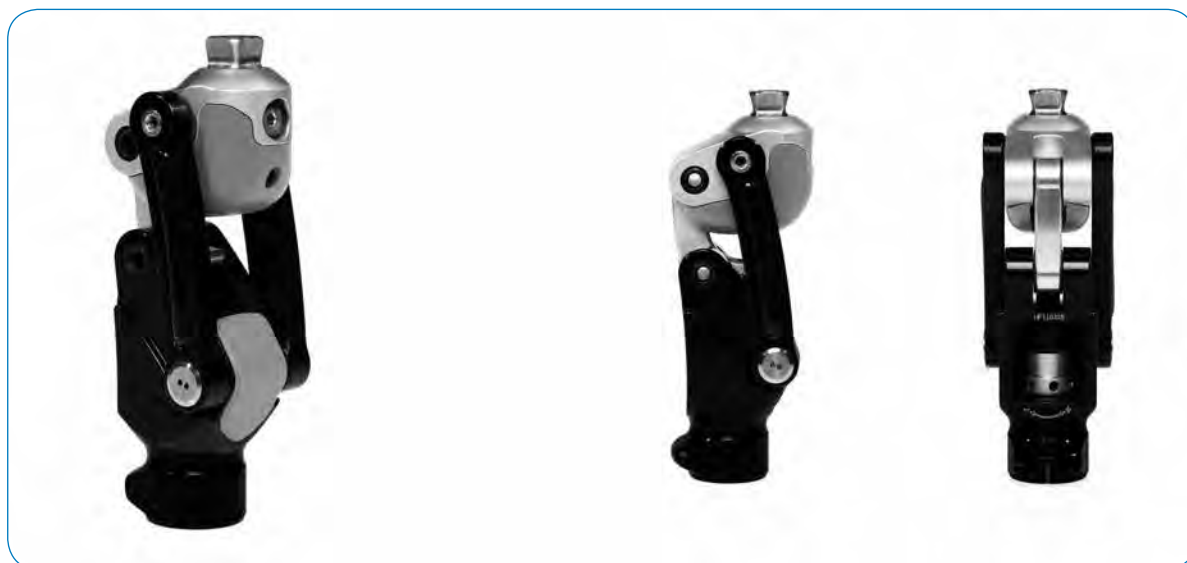

- Offers security and stability with 'four-bar' geometry.
- Designed to optimise the balance between knee stability and dynamic behaviour for each user.
- Adjustable extension assist and constant friction allowing a very easy setup and adjustment of the knee.
- All adjustment features accessible without prosthesis disassembly.
- The polycentric geometry provides mid-swing shortening.
- Designed for household and limited community ambulators, single slow speed ambulators, low impact users and low voluntary control users.
- Simple to set up, the knee's adjustment features (including the adjustable extension assist) can all be accessed externally with no disassembly of the prosthesis required.
- The design is maintenance free.
- Balance Knee is available in locking and non-locking options.
- Weight limit 125kg
- Flexion Stop at 180°
- Height of unit is 154mm
- Weight of unit is 590g

Clearance

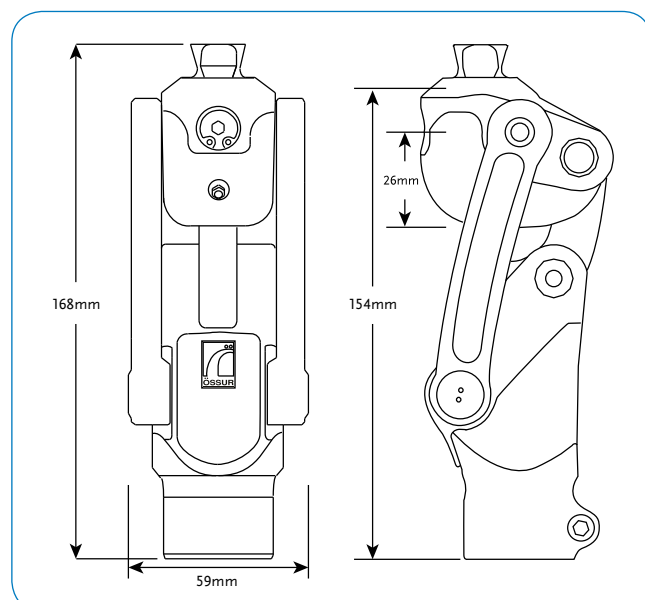

Balance™ Knee

Kits

Part#	Description
BKN12500	Balance Knee
BKN12511	Balance Knee Locking*

Includes:

· Balance Knee

· Shin Ferrule

*Lanyard Kit



Cosmetic Finishing

Part#	Description
BKN12550	Shin Ferrule for Balance Knee



Locking Knee


Introduction



Locking Knee from Össur is a simple, reliable and lightweight knee joint.

Locking Knee is intended to be used during rehabilitation and as a permanent solution for the less active amputee.

- Weight limit 125kg
- Flexion Stop at 140°
- Knee Axis height is 28mm
- Weight of unit is 248g

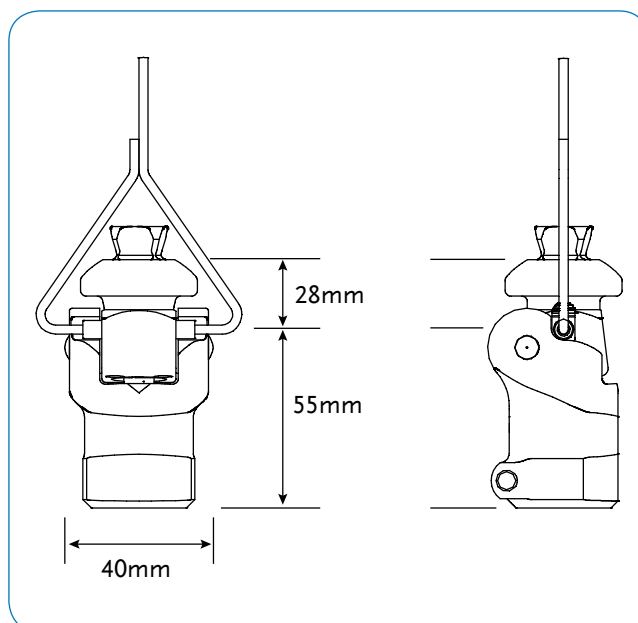

Kits

Part#	Description
LKN12500	Locking Knee


Includes:

- Lanyard Kit


Part#	Description
LKN12590	Lanyard Release Kit

Total Knee® 1900

Introduction

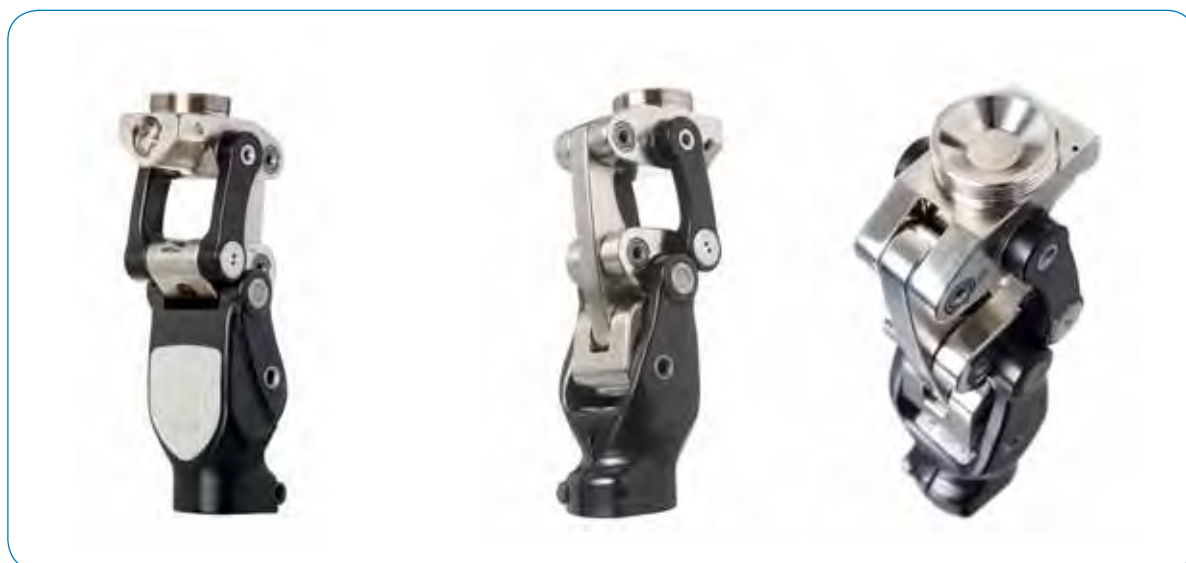

- Polycentric knee with geometric locking system
- Unique swing phase control through use of friction within an elastic polymer medium
- Used for single speed ambulation
- Activity level: Low
- Weight limit 100kg
- Adjustable stance flexion
- Adjustable extension promoter
- Weight of unit 620g
- 160° of knee flexion available

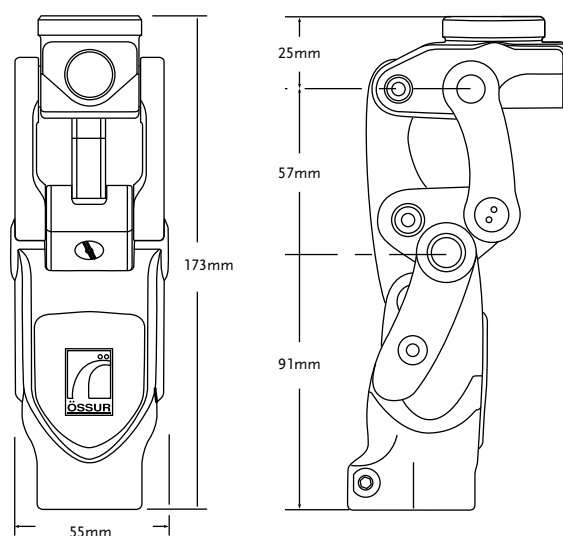

Kits

Part#	Description
1900	Total Knee 1900

Cosmetic Finishing

Picture	Part#	Description
1	2047	Shin ferrule, black
2	2052	Knee Cap, small, black
3	2053	Knee Cap, large, beige



Total Knee® 1900

Components

BUMPER SELECTION CHART

Colour	Durometer	Patient Weight
Green	Soft	Approx. 35-50kg
Black	Regular	Approx. 50-70kg
Red	Firm	Approx. 70-90kg
Yellow	Extra firm	Approx. 90-100kg

Note: The knee comes with the red bumper installed.

BUMPER KIT*

Part#	Description
2000-Kit	A complete set of bumpers and shims

Includes

- 1 Bumper of each durometer
- 1 Shim 0.38mm
- 1 Shim 0.76mm

BUMPER*

Part#	Description
2000-204G	Bumper - Green, soft durometer
2000-204B	Bumper - Black, regular durometer
2000-204R	Bumper - Red, firm durometer
2000-204Y	Bumper - Yellow, extra firm durometer

SHIMS*

Part#	Description
2000-205A	0.38mm Plastic Shim
2000-205B	0.76mm Plastic Shim

*For units with serial numbers 8000 and higher

EXTENSION PROMOTER



Part#	Description
1905	Assists in limiting heel rise and moving the knee into a fully extended position

EXTENSION ASSIST

Part#	Description
2070S	Fully assembled Extension Assist with yellow (Regular) spring inserted

Includes

- Blue Spring (Firm)
- Red Spring (Extra Firm)

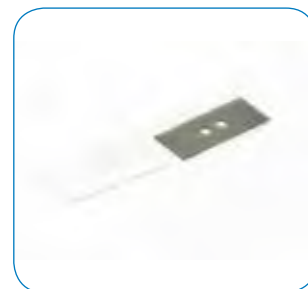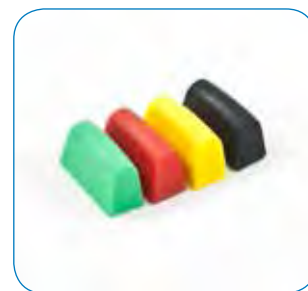

Total Knee® 1900

Components

EXTENSION PROMOTER


Part#	Description
1905	Assists in limiting heel rise and moving the knee into a fully extended position



BACK LINK (POSTERIOR LINKAGE)

Part#	Description
1900-SSB	A high load back link can be factory fitted to the Total Knee 1900.

Note: There is an additional charge for these special back links.

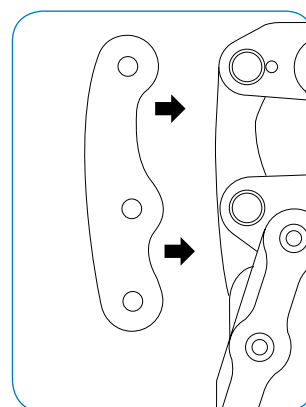

TOTAL KNEE ADAPTERS

Picture	Part#	Description	Material	Weight limit	Weight
1	A-834300	Male Pyramid Adapter	Aluminium	100kg	83g
2	A-845300	Female Pyramid Adapter	Titanium	166kg	75g
3	A-122300	4-Hole Adapter	Aluminium	125kg	65g
4	A-835300	Male Pyramid Adapter	Titanium	166kg	83g




Total Knee® 2000

Introduction



- Polycentric knee with geometric locking system
- Swing control adjustments through 3-phase and 3-valve hydraulic system
- Used for multi-speed ambulation
- Activity level: moderate
- Weight limit 100kg
- Adjustable stance flexion
- Weight of unit 690g
- 160° of knee flexion available

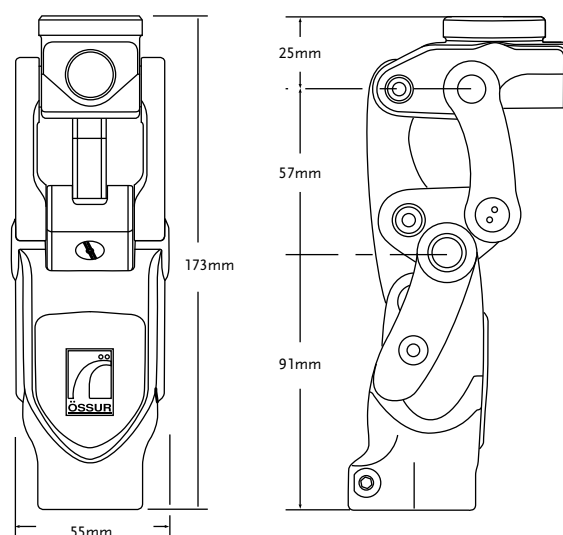

Kits

Part#	Description
2000	Total Knee 2000

Cosmetic Finishing

Picture	Part#	Description
1	2047	Shin ferrule, black
2	2052	Knee Cap, small, black
3	2053	Knee Cap, large, beige



Total Knee® 2000

Components

BUMPER SELECTION CHART

Colour	Durometer	Patient Weight kg
Green	Soft	Approx. 35-50kg
Black	Regular	Approx. 50-70kg
Red	Firm	Approx. 70-90kg
Yellow	Extra firm	Approx. 90-100kg

Note: The product comes with the red bumper installed.

BUMPER KIT*

Part#	Description
2000-Kit	A complete set of bumpers and shims

Includes

- 1 Bumper of each durometer
- 1 Shim 0.38mm
- 1 Shim 0.76mm

BUMPER*

Part#	Description
2000-204G	Bumper - Green, soft durometer
2000-204B	Bumper - Black, regular durometer
2000-204R	Bumper - Red, firm durometer
2000-204Y	Bumper - Yellow, extra firm durometer

SHIMS*

Part#	Description
2000-205A	0.38mm Plastic Shim
2000-205B	0.76mm Plastic Shim

*For units with serial numbers 14000 and higher

EXTENSION ASSIST

Part#	Description
2070S	Fully assembled Extension Assist with yellow (Regular) spring inserted

Includes

- Blue Spring (Firm)
- Red Spring (Extra Firm)



Total Knee® 2000

Components

EXTENSION PROMOTER


Part#	Description
1905	Assists in limiting heel rise and moving the knee into a fully extended position



BACK LINK (POSTERIOR LINKAGE)

Part#	Description
2000-SSB	A high load back link can be factory fitted to the Total Knee 2000.

Note: There is an additional charge for these special back links.

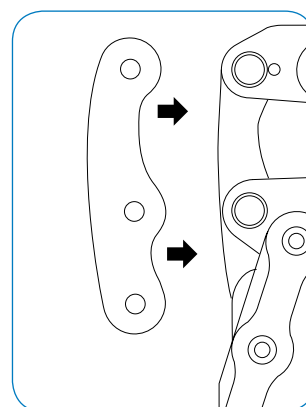

TOTAL KNEE ADAPTERS

Picture	Part#	Description	Material	Weight limit	Weight
1	A-834300	Male Pyramid Adapter	Aluminium	100kg	83g
2	A-845300	Female Pyramid Adapter	Titanium	166kg	75g
3	A-122300	4-Hole Adapter, with screws	Aluminium	125kg	65g
4	A-835300	Male Pyramid Adapter	Titanium	166kg	83g




Total Knee® 2100

Introduction



- Polycentric knee designed for multi-speed ambulation, higher impact levels and heavier loads
- Activity level: moderate to high
- 25% increased fluid capacity over 2000 model and cooling fins for greater reliability
- Geometric locking system
- 30mm pylon adaptability
- Weight limit 125kg
- Adjustable stance flexion
- Adjustable extension promoter
- Weight of unit 920g
- 160° of knee flexion available

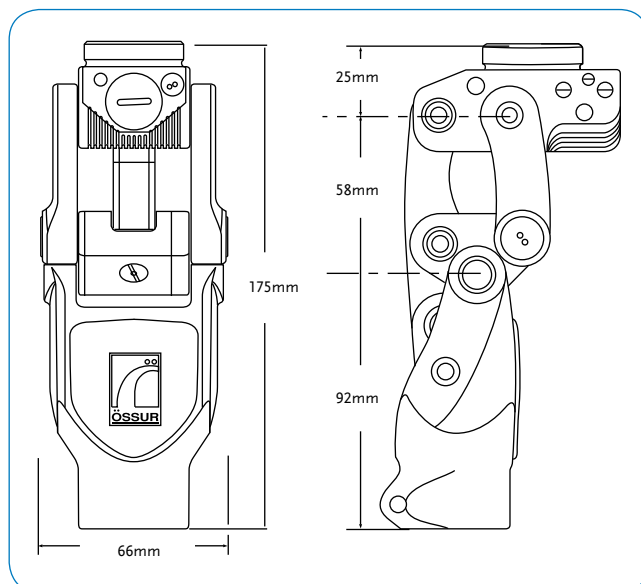


Kits

Part#	Description
2100	Total Knee 2100



Cosmetic Finishing

Picture	Part#	Description
1	2145	Shin ferrule, black
2	2150	Knee Cap, small, black
3	2151	Knee Cap, large, beige

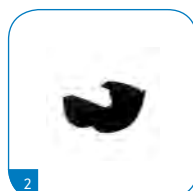

Total Knee® 2100

Components

BUMPER SELECTION CHART

Colour	Durometer	Patient Weight kg
Black	Soft	Approx. 50-70kg
Red	Regular	Approx. 70-90kg
Yellow	Firm	Approx. 90-110kg
Orange	Extra firm	Approx. 110-125kg

BUMPER KIT

Part#	Description
2100-Kit	A complete set of bumpers and shims

Includes

- 1 Bumper of each durometer
- 1 Shim 0.38mm
- 1 Shim 0.76mm

BUMPER

Part#	Description
2100-317B	Bumper - Black, soft durometer
2100-317R	Bumper - Red, regular durometer
2100-317Y	Bumper - Yellow, firm durometer
2100-317O	Bumper - Orange, extra firm durometer

SHIMS

Part#	Description
2100-312A	0.38mm Plastic Shim
2100-312B	0.76mm Plastic Shim

EXTENSION ASSIST

Part#	Description
2070S	Fully assembled Extension Assist with yellow (Regular) spring inserted

Includes

- Blue Spring (Firm)
- Red Spring (Extra Firm)

EXTENSION PROMOTER

Part#	Description
1905	Assists in limiting heel rise and moving the knee into a fully extended position



Extension promoter is found in the balance block

Total Knee® 2100

Components

GREASE GUN

Part#	Description
2105	Grease for Total Knee 2100



TOTAL KNEE ADAPTERS

Picture	Part#	Description	Material	Weight limit	Weight
1	A-834300	Male Pyramid Adapter	Aluminium	100kg	83g
2	A-845300	Female Pyramid Adapter	Titanium	166kg	75g
3	A-122300	4-Hole Adapter, with screws	Aluminium	125kg	65g
4	A-835300	Male Pyramid Adapter	Titanium	166kg	83g




Total Knee® Junior

Introduction



- Paediatric polycentric knee with geometric locking system
- Weight limit 45kg
- Adjustable stance flexion feature
- Adjustable extension promoter
- 160° of knee flexion available
- Weight of unit 380g

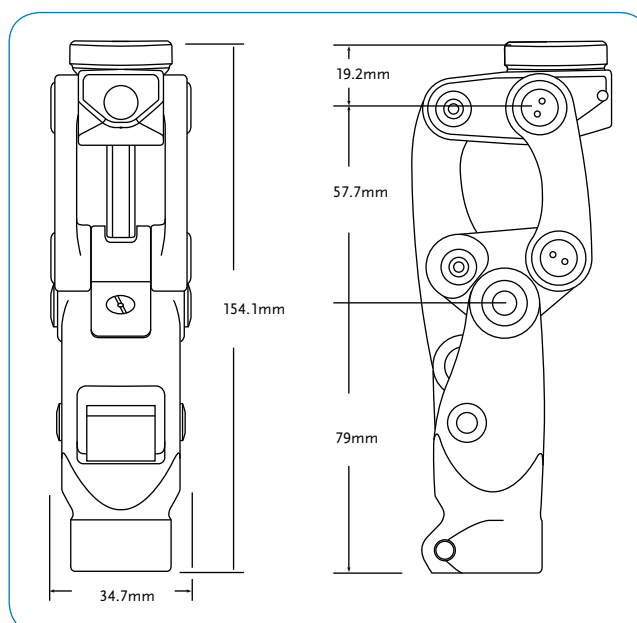

Kits

Part#	Description
1100	Total Knee Junior

Cosmetic Finishing

Picture	Part#	Description
1	1145	Shin ferrule, black
2	1150	Knee Cap, small, black
3	1151	Knee Cap, large, beige



Total Knee® Junior

Components

BUMPER SELECTION CHART

Colour	Durometer	Patient Weight kg
Green	Soft	Approx. 14-20kg
Black	Regular	Approx. 20-29kg
Red	Firm	Approx. 29-36kg
Yellow	Extra Firm	Approx. 36-45kg

BUMPER KIT

Part#	Description
1100-Kit	A complete set of bumpers and shims

Includes

- 1 Bumper of each durometer
- 1 Shim 0.38mm
- 1 Shim 0.76mm

BUMPER

Part#	Description
1100-119G	Bumper - Green, soft durometer
1100-119B	Bumper - Black, regular durometer
1100-119R	Bumper - Red, firm durometer
1100-119Y	Bumper - Yellow, extra firm durometer

SHIMS

Part#	Description
1100-121A	0.38mm Plastic Shim
1100-121B	0.76mm Plastic Shim

EXTENSION ASSIST



Part#	Description
1170	Fully assembled Extension Assist with yellow spring inserted, blue spring included

EXTENSION PROMOTER

Part#	Description
1105-A	Assists in limiting heel rise and moving the knee into a fully extended position


TOTAL KNEE JUNIOR ADAPTERS

Picture	Part#	Description	Material	Weight
1	AJ-831300	Paediatric Pyramid Adapter (Rated to 45kg)	Titanium/Aluminium	41g
2	AJ-122300	Paediatric 4-Hole Adapter, with screws	Aluminium	36g

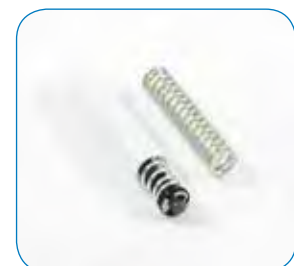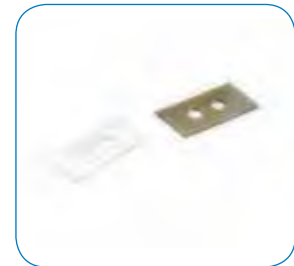

Extension promoter is found in the balance block


Cosmetic Finishing (not for Junior)

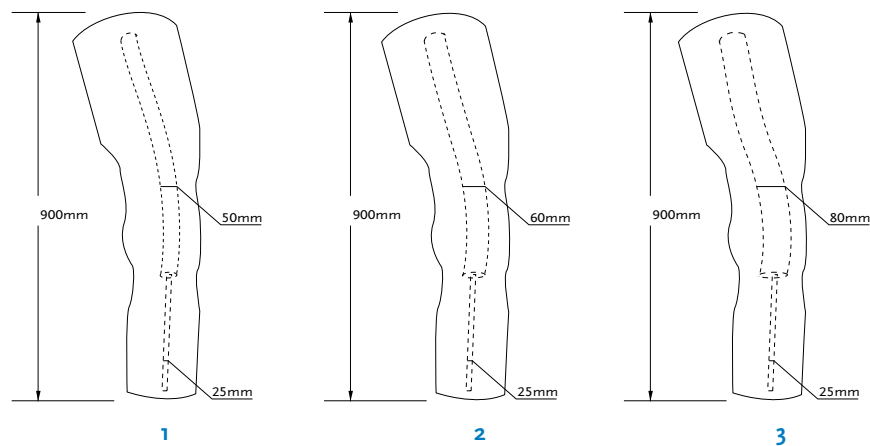

CONTINUOUS COSMETIC FOAM

Picture	Part#	Bore Size
1	C-880877	50-25
2	C-880878	60-25
3	C-880879	80-25

- Designed to fit left or right
- Open Cell Polyethylene Foam (60g)
- Fire retardant

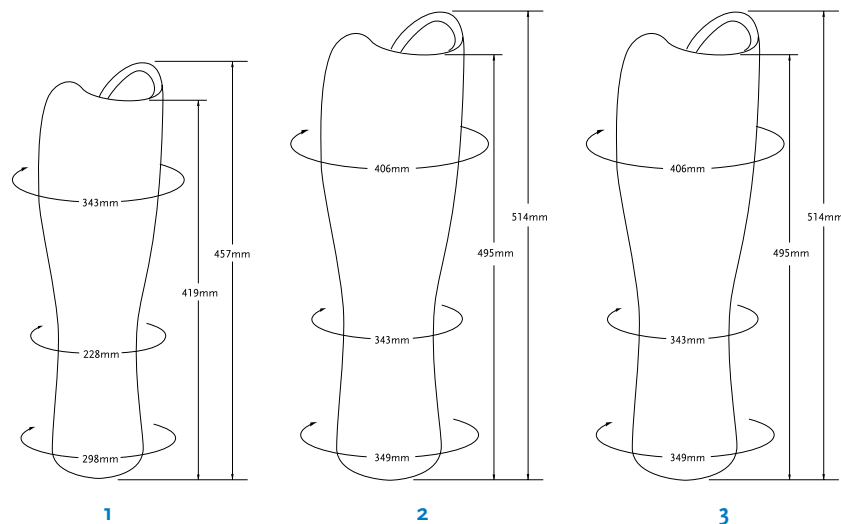

DISCONTINUOUS COSMETIC FOAM

Picture	Part#	Description
1	2075R	Cosmetic Foam Cover, Regular - for use with Total Knee 1900 and 2000
2	2075L	Cosmetic Foam Cover, Large - for use with Total Knee 1900 and 2000
3	2175L	Cosmetic Foam Cover, High activity, Large - for use with Total Knee 2000 and 2100

- Beige toned
- Designed to fit over Total Knee® Range
- Designed to fit left or right
- Internally shaped to accommodate knee cap, shin ferrule and 30mm tube
- Fire retardant

Cosmetic Finishing



SINGLE DENSITY SHIN FAIRING BLOCK

Part#	Size
C-860977	One Size

- Closed Cell Polyethylene Foam (25g)
- Fire retardant



DUAL DENSITY SHIN FAIRING BLOCK

Part#	Size
FCX0002	One Size

- Flexible section at ankle
- Fire retardant



Feet

Introduction

Flex-Foot® was first introduced in 1984 by prosthetic user, inventor and research prosthetist Van Phillips. Three critical breakthroughs made this product unique and revolutionised the everyday aspirations of amputees:

- The first, energy storage and release, is a function inherent in the patented carbon fibre design of Flex-Foot.
- The second, Active Tibial Progression transfers the stored energy from heel contact and initial loading into forward propulsion towards toe off and swing phase.
- The third, vertical shock absorption, enables a more natural gait and protects the sound limb and remaining joints of the amputated limb from excessive shock.

Today, different functions are available within the Flex-Foot range to suit individual needs. All models offer significant advantages over conventional prosthetic feet. More than 25 years of scientific material research have led to this unparalleled design, making Flex-Foot the ideal choice for durability, comfort and performance.

Product Overview



Contents

Flex-Foot® Balance	150	Elation®	207
Balance™ Foot J.....	152	Flex-Symes™	211
Flex-Foot® Assure.....	154	Flex-Foot® Junior	217
Talux® ST.....	158	Chopart™	220
Vari-Flex®	164	Flex-Run™ with Nike Sole	222
LP Vari-Flex®	170	Flex-Run™ Junior.....	225
Vari-Flex® Junior	174	Flex-Sprint™	228
Vari-Flex® XC	177	Cheetah®	230
Vari-Flex® XC Rotate™	182	Cheetah® Xtreme with Nike Spike Pad	233
Re-Flex Rotate™.....	185	Cheetah® Xtend with Nike Spike Pad	236
LP Rotate™	191	Cheetah® Xplore Junior	239
Vari-Flex® Modular	196	Cheetah® Junior.....	242
Re-Flex Shock™	202		

Product Characteristics

PROPORTIONAL RESPONSE


The layering of carbon fibre, optimised through extensive computer analysis and mechanical testing, ensures that the deflection of the forefoot from mid-stance to toe-off is proportional to the user's weight and impact level.

Benefit: Reduced fatigue by energy-efficient walking

**CARBONX® ACTIVE HEEL**


The CarbonX Active Heel absorbs the energy created during initial contact through the loading response as the amputee transfers his or her body weight onto the prosthetic foot in the early stance phase.

Benefit: Comfort and protection by active shock absorption

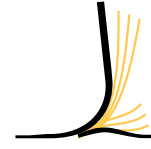**ACTIVE TIBIAL PROGRESSION**


Vertical forces generated at heel contact are stored and translated into a linear motion described as Active Tibial Progression, from the foot-flat to toe-off phase. This action reduces the need to actively push the body forward using the contralateral foot and also equalises stride length. The foot's unique design and technology produces all movement; there are no mechanical joints.

Benefit: Normal gait through natural mobility

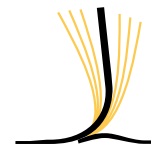**FULL LENGTH TOE LEVER**


The full length toe lever matches the length of the sound foot, giving a smoother, more natural gait. It also ensures that users spend equal time on the prosthetic and natural limb.

Benefit: Stability and security in stance phase

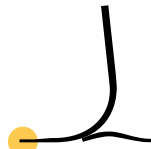**VERTICAL SHOCK PYLON**

The vertical shock pylon reduces trauma to the residual limb, joints and lower back during daily activities. It also ensures a more energy-efficient gait by minimising vertical displacement of the centre of gravity.

Benefit: Protection of joints, the spine and residual limb by reduction of the shock loading

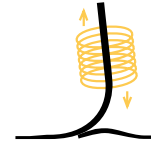**ROTATION**

Torsion/rotation is a key factor in the human gait. The trunk rotation affects every joint in the lower limbs. If joints are missing, or the rotation is limited by an interface (socket), natural gait is compromised. The full range torsion/rotation brings the prosthetic user much closer to this all-important natural gait and is also very comfortable.

Benefit: Residual limb and joint protection by natural freedom of movement

**HEEL HEIGHT ADJUSTABILITY**

Heel height can be adjusted for various types of shoes to ensure correct alignment. Users can easily change shoes to fit most occasions as the heel height can be adjusted within the given range.

Benefits: Adjustable heel height allows correct alignment of the prosthesis for a variety of footwear and therefore decreases the load on the body and the contralateral limb resulting from misalignment.

**MULTI-AXIAL GROUND COMPLIANCE**

Multi-axial compliance of the foot provides improved ground contact during ambulation on uneven surfaces.

Benefit: Improved stability in stance.

**SPLIT TOE**

A split toe feature allows the foot to adjust to the surface below, allowing amputees to walk naturally on uneven ground.

Benefit: Provides stability on uneven ground

**SANDAL TOE**

Feet designed with a sandal toe are used with an anatomically correct foot cover with a sandal toe.

Benefits: Sandal toe is cosmetically appealing, making the foot particularly suitable for sandal type shoes.

**EVO™**

Energy Vector Optimisation utilises the unique integration between the foot module and foot cover.

Benefit: Optimised energy transfer during stance phase resulting in a fluid roll-over and smoother gait for amputees.




Flex-Foot® Balance


Introduction




- Flex-Foot Balance is a comfortable multiaxial foot that combines proven durability and stability with comfort leading to optimum patient compliance for low activity amputees.
- The full length toe lever and integrated multiaxial design offer a smooth transition from heel strike to toe-off. This means added comfort and relief for those that may be struggling with dysvascular complications.

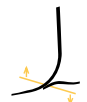USER INFORMATION		
Maximum Patient Weight	136kg	
Impact Level	Low	
FOOT INFORMATION		
Categories	0-5	
Sizes	21-28cm	
Weight of the standard foot module (Size 27, Cat 3) with male pyramid and foot cover	490g	
Features		
		
Full toe lever	Multi-axial	Sandal toe

Flex-Foot Balance with Male Pyramid

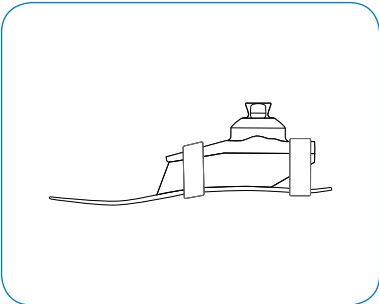

Flex-Foot Balance kit with male pyramid

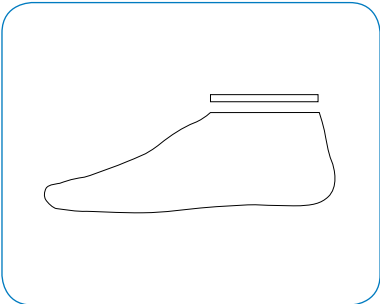


Foot cover kit

Flex-Foot® Balance

Category Selection Guide


Weight	Foot Size							
	21	22	23	24	25	26	27	28
36-44kg	0	0	0	0				
45-61kg	1	1	1	1	1	1		
62-80kg			2	2	2	2	2	2
81-97kg					3	3	3	3
98-114kg					4	4	4	4
115-136kg					5	5	5	5

Clearance

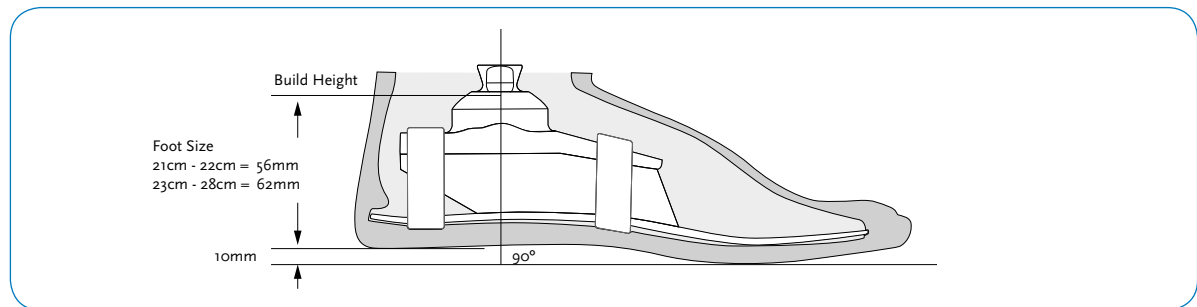

Kits

FLEX-FOOT BALANCE KIT WITH MALE PYRAMID

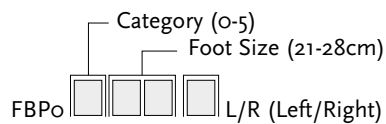

Includes

- Foot Module with Male Pyramid
- Foot cover kit

Cosmetic Finishing

See page 146 for details on cosmetic finishing options.



Please specify when ordering brown covers. Add "BR" to the part#.

Includes

- Foot Cover
- Attachment Plate
- Flex-Foot Sock

Part#	Description	Size
FCX63006	Flex-Foot Sock	21-25cm
FCX63007	Flex-Foot Sock	26-28cm
Part#	Description	
CM100068	Shoe Horn	



Balance™ Foot J

Introduction



- Balance Foot J is designed to meet the needs of low impact walkers. It is lightweight with a heel that provides stability and cushioning at heel strike. Together with the full-length keel the Balance J provides excellent stability and support.
- The foot cover that comes with Balance Foot J incorporates our successful EVO™ technology which promotes a more fluid and energetic forward progression by mimicking the natural movement of the human foot from heel-strike to toe-off.


USER INFORMATION

Maximum Patient Weight	136kg
Impact Level	Low


FOOT INFORMATION

Categories	1-5
Sizes	22-30cm
Weight of the standard foot module (Size 27, Cat 2) with male pyramid and foot cover	576g


Features

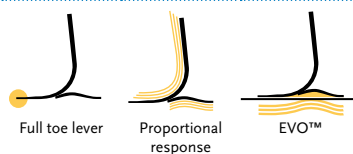

Balance Foot J with Male Pyramid

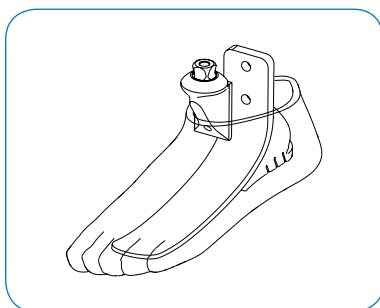

Balance Foot J with male pyramid

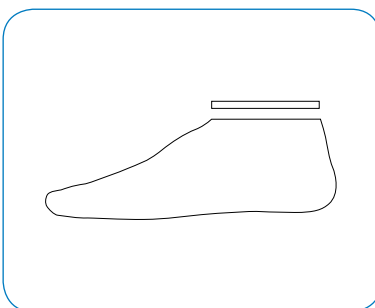


Foot cover kit

Balance™ Foot J

Category Selection Guide


Weight	Foot Size								
	22	23	24	25	26	27	28	29	30
45-61kg	1	1	1	1	1				
62- 80kg		2	2	2	2	2	2	2	2
81-97kg				3	3	3	3	3	3
98-114kg				4	4	4	4	4	4
115-136kg				5	5	5	5	5	5

Clearance

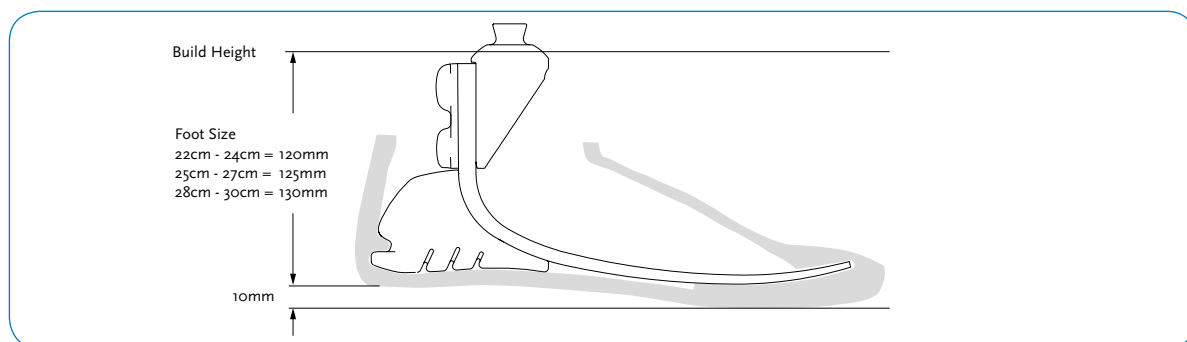

Kits

BALANCE FOOT J WITH MALE PYRAMID

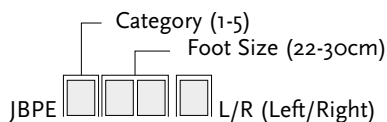

Includes

- Foot Module with Male Pyramid
- Foot cover kit*

*For Balance Foot J with brown foot cover add a "BR" to the part#

Components

MALE PYRAMID ADAPTER (METRIC) KIT

Part#	Weight	Cat.
SFX00165	104g	1-5

Includes

- Pyramid
- Hardware (M8 bolts)
- Loctite



Cosmetic Finishing

See page 146 for details on cosmetic finishing options.




Please specify when ordering brown covers. Add "BR" to the part#.

Includes

- Foot Cover
- Attachment Plate
- Flex-Foot Sock

Part#	Description	Size
FCX63006	Flex-Foot Sock	22-25cm
FCX63007	Flex-Foot Sock	26-30cm
Part#	Description	
CM100068	Shoe Horn	

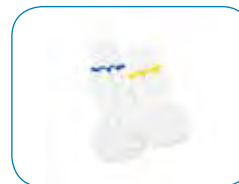

Flex-Foot® Assure™

Introduction



- Lightweight, carbon fibre foot with moderate energy storage and release
- Active heel with cushion wedge
- Provides the user with optimal roll over, safety and stability during stance phase
- Reduces impact load on the contralateral side
- Male pyramid or 30mm tube option

USER INFORMATION

Maximum Patient Weight	136kg
Impact Level	Low

FOOT INFORMATION

Categories	1-5
Sizes	19-30cm
Weight of the standard foot module (Size 27) with male pyramid	575g

Features



The EVO™ feature enables amputees to experience natural gait progression by utilising the foot plate to its fullest potential. A more natural gait reduces fatigue and puts less strain on the lower back and sound side



Full toe lever



Proportional response



Tibial progression




Active heel




EVO™

Flex-Foot Assure with Male Pyramid

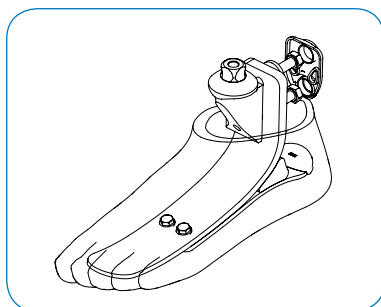


Flex-Foot Assure kit with male pyramid

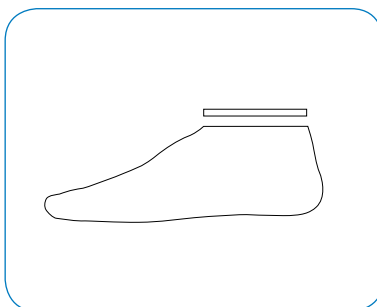


Foot cover kit


Flex-Foot Assure with Glass Fibre Tube

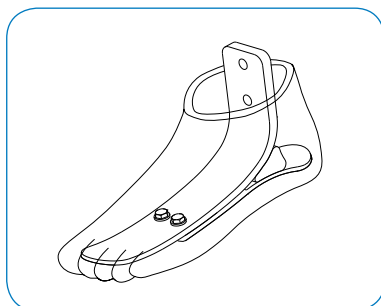

Flex-Foot Assure kit without adapter

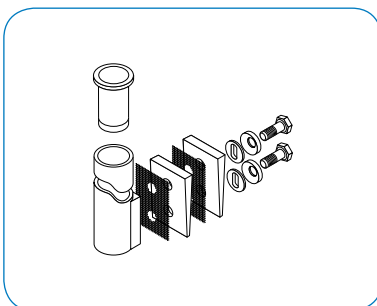

Tube kit

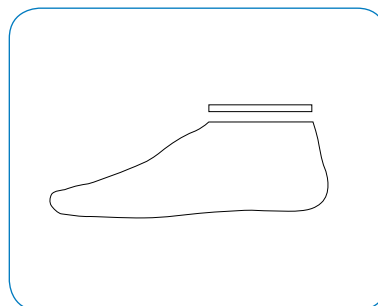


Foot cover kit

Flex-Foot® Assure

Category Selection Guide

Weight	Foot Size											
	19	20	21	22	23	24	25	26	27	28	29	30
45-61kg	1	1	1	1	1	1	1	1	1			
62-80kg	2	2	2	2	2	2	2	2	2	2		
81- 97kg					3	3	3	3	3	3	3	
98-114kg						4	4	4	4	4	4	4
115-136kg								5	5	5	5	5

Dimensions

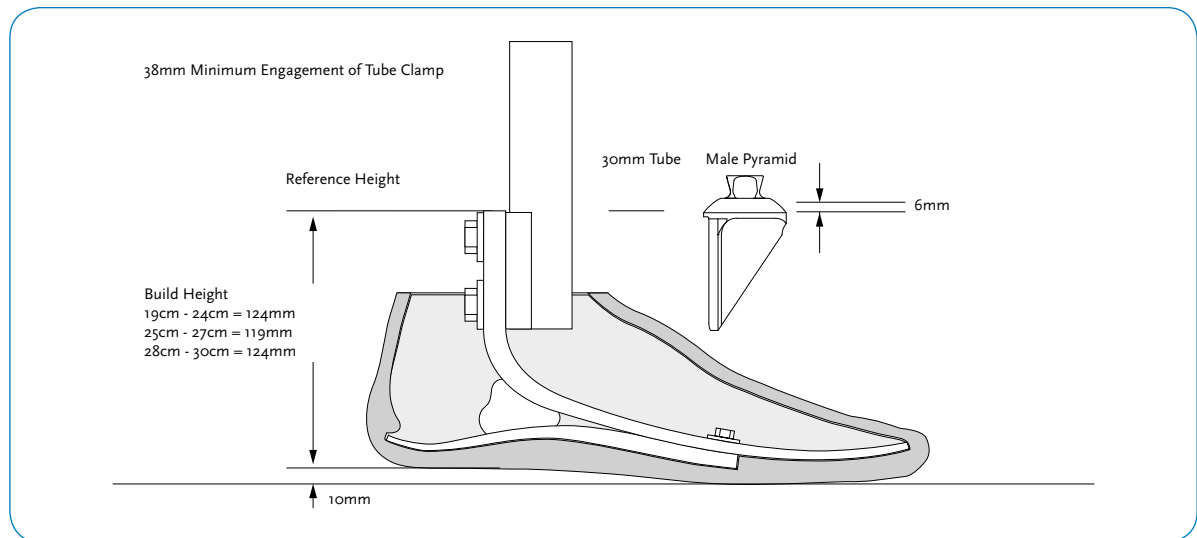

All measurements using 10mm heel height.

Kits

Part#


FLEX-FOOT ASSURE KIT WITH MALE PYRAMID

Category (1-5)
 Foot Size (19-30cm)
 FAPE  L/R (Left/Right)

Includes

- Foot Module with Male Pyramid
- Foot cover kit*



FLEX-FOOT ASSURE KIT WITHOUT ADAPTER

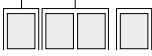Category (1-5)
 Foot Size (19-30cm)
 FAXE  L/R (Left/Right)

Includes

- Foot Module
- Foot cover kit*

*For Flex-Foot Assure with brown foot cover add a "BR" to the part#

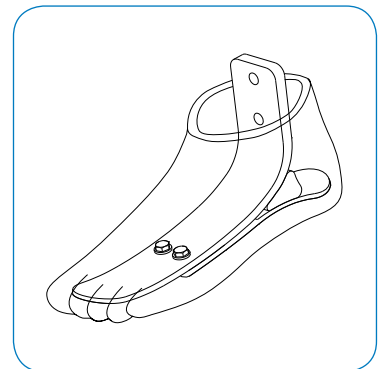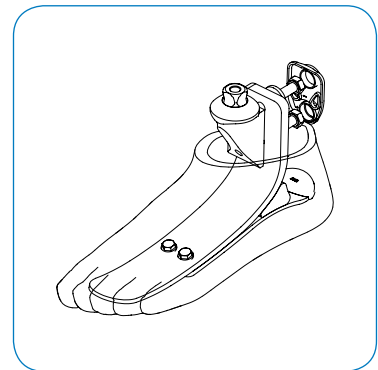

Flex-Foot® Assure

Components

MALE PYRAMID ADAPTER (METRIC) KIT

Part#	Weight	Cat.
SFX00165	104g	1-5

Includes

- Pyramid
- Loctite
- Hardware (M8 bolts)

30MM GLASS FIBRE TUBE KIT

Part#	Weight	Length	Wall Thickness	Cat.
SFX00001	292g	279mm	3mm	1-3 Red
SFX00002	333g	381mm	3mm	1-3 Red
SFX00003	331g	279mm	4mm	4-5 Green
SFX00004	401g	381mm	4mm	4-5 Green

Includes



- Tube
- Bolts (5/16" UNF)
- Spherical Washers
- Friction Pads
- Loctite 410
- Alignment Wedges

TUBE INSERT

Part#	Weight	Diameter
A-718009	18g	30mm

Part#	Description
SFX00123	30mm Alignment Fitting Tool (Metric)

Cosmetic Finishing

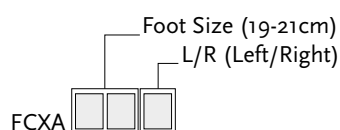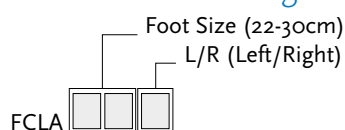

Please specify when ordering brown covers. Add "BR" to the part#.

Includes

- Foot Cover
- Flex-Foot EVO Sock
- Attachment Plate

Part#	Description	Size
FCE63006	Flex-Foot EVO Sock	22-25cm
FCE63007	Flex-Foot EVO Sock	26-30cm


Part#	Description
CM100068	Shoe Horn

See page 146 for details on cosmetic finishing options.



Flex-Foot® Assure

Diagram

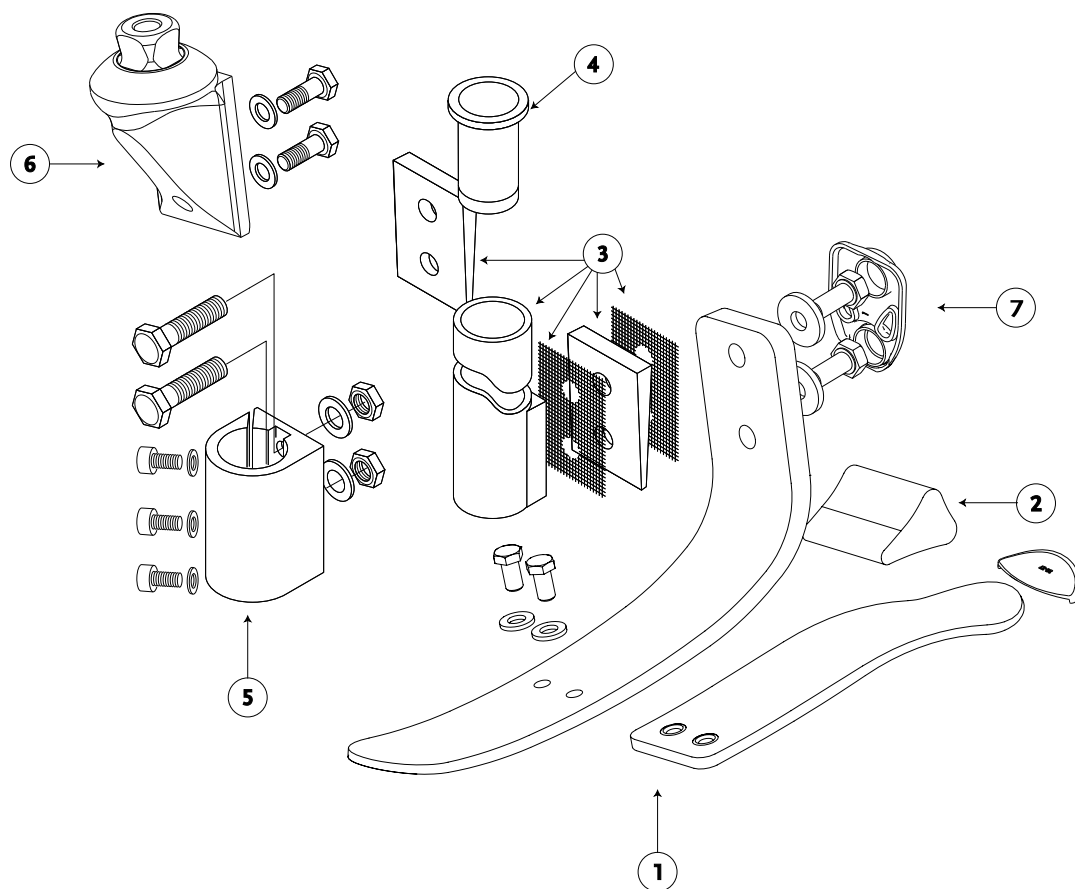

Picture	Description
1	Foot Module
2	Foam Wedge
3	30mm Fibreglass Tube Kit
4	Tube insert for 5-6
5	30mm Alignment Fitting Tool
6	Male Pyramid
7	Attachment Bolts

Talux® ST

Introduction



- Provides fluid, natural walking motion over a variety of terrains
- Dual “J” springs and carbon fibre keel ensure energy storage and release
- Tarsal Core offers multi-axial function and Achilles Strap enhances forward propulsion
- Male pyramid or 30mm tube option
- Sandal Toe


USER INFORMATION

Maximum Patient Weight	147kg
Impact Level	Low and Moderate


FOOT INFORMATION

Categories	1-8
Sizes	23-30cm
Weight of the foot (Size 27) with male pyramid and foot cover	740g


Features

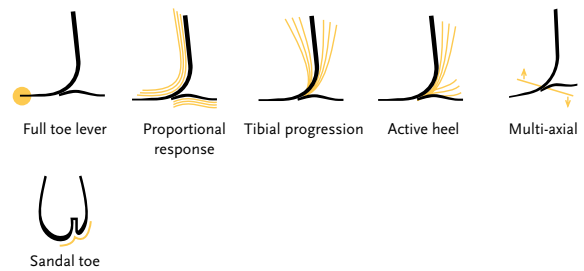

Talux with male pyramid

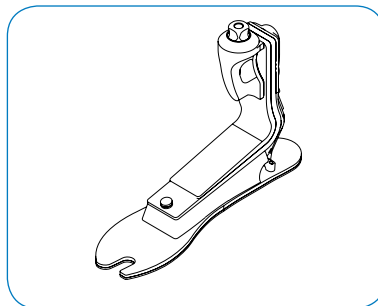


Talux kit with male pyramid

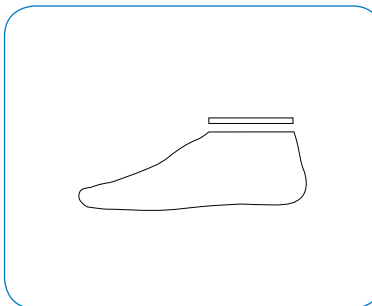


Foot cover kit


Talux with carbon fibre tube

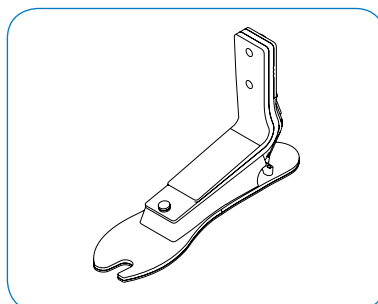

Talux kit without pyramid

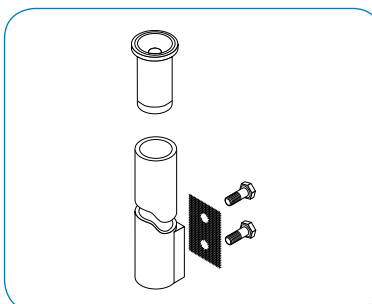

Tube kit

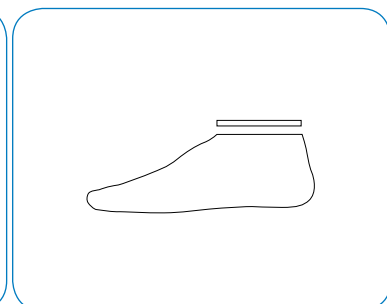

Foot cover kit

Talux® ST

Category Selection Guide

Low Impact Level

Weight


Foot Size	45-52kg	53-59kg	60-68kg	69-77kg	78-88kg	89-100kg	101-116kg	117-130kg	131-147kg
23/24	1	1	2	3	4	5			
25/26			2	3	4	5	6	7	
27				3	4	5	6	7	
28				3	4	5	6	7	8
29/30					4	5	6	7	8

Moderate Impact Level

Weight

Foot Size	45-52kg	53-59kg	60-68kg	69-77kg	78-88kg	89-100kg	101-116kg	117-130kg	131-147kg
23/24	1	2	3	4	5				
25/26		2	3	4	5	6	7		
27			3	4	5	6	7		
28			3	4	5	6	7	8	
29/30				4	5	6	7	8	

Dimensions

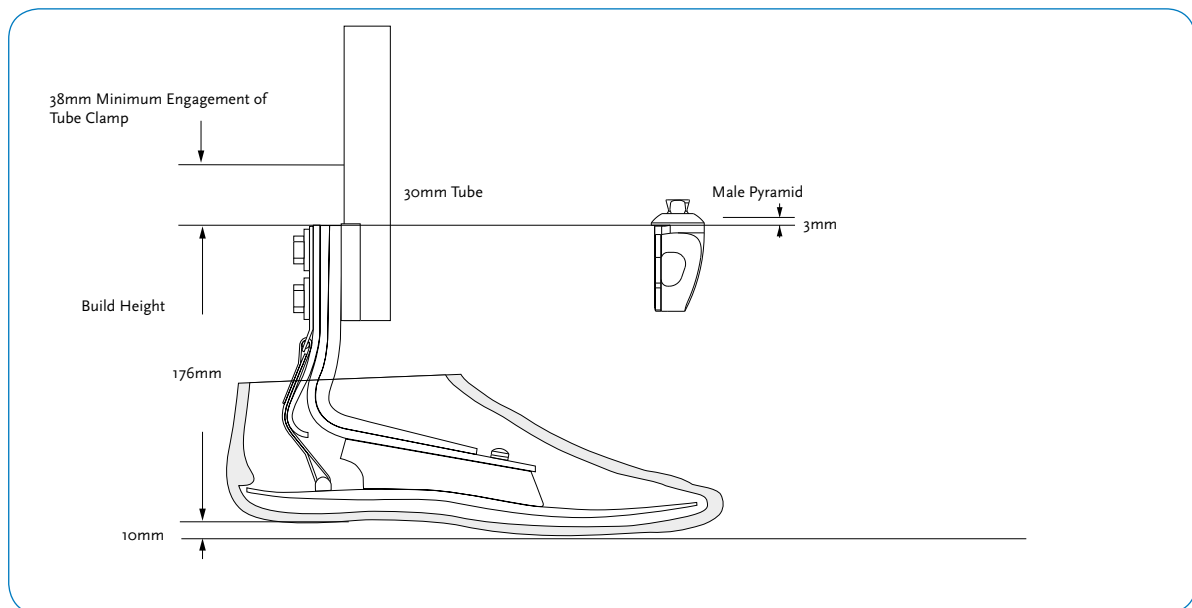

All measurements using 10mm heel height.


Kits

TALUX KIT WITH MALE PYRAMID

Category (1-8)
Foot Size (23-30cm)
TLPO  L/R (Left/Right)

Includes


- Foot Module
- Achilles Strap
- Male Pyramid
- Tarsal Core

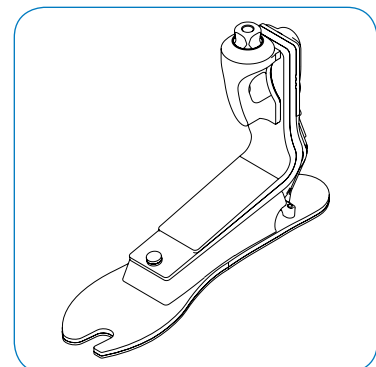

Talux® ST


Kits

TALUX KIT WITHOUT ADAPTER

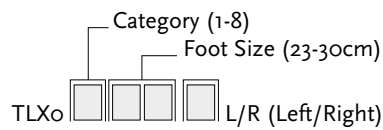

Includes

- Foot Module
- Achilles Strap
- Tarsal Core

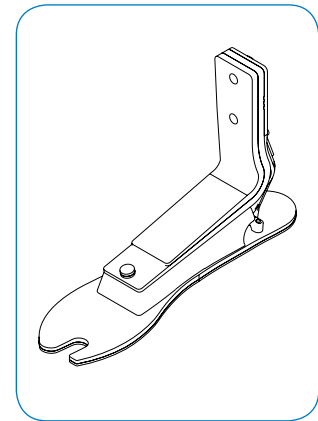

Components

BOLT-ON MALE PYRAMID ADAPTER (METRIC)

Part#	Weight	Cat.
FFC2210003	175g	1-6
FFC2210004	205g	7-8

Includes

- Pyramid
- Hardware (M8 or M10 bolts)
- Loctite



30MM CARBON FIBRE TUBE KIT

Part#	Weight	Length	Cat.
TLX01011	216g	229mm	1-6 Black
TLX01012	275g	229mm	7-8 Black
TLX01013	293g	381mm	1-6 Black
TLX01014	369g	381mm	7-8 Black

Includes

- Carbon Fibre Tube
- Attachment Hardware (5/16" UNF bolts)
- Friction Pads
- Loctite
- Tube Insert (cat. 7-8)



Talux® ST

Components

TUBE INSERT

Part#	Weight	Diameter	Cat.
A-718009	18g	30mm	7-8




ACHILLES STRAP KITS

Part#	Foot Size	Cat.
TLX01002	23-26cm	1-6
TLX01003	27-30cm	1-6
TLX01004	23-26cm	7-8
TLX01005	27-30cm	7-8

Includes

- Strap
- Hanger
- Pin
- Bolts

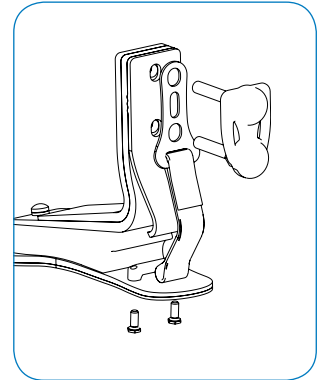


PLANTAR FLEXION CONTROL BOLTS

Part#	Foot Size	Cat.
TLX01006	23-24cm	1-6
TLX01007	25-26cm	1-6
TLX01008	27-28cm	7-8
TLX01009	29-30cm	7-8

Includes

- Bolt
- Urethane Washer
- Nut
- Loctite

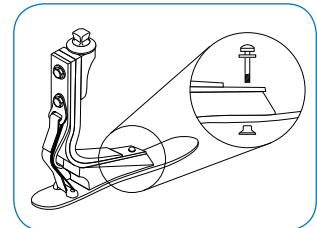


Part#	Description
CM160016	Loctite 410 black



Part#	Description
TLX01001	Heel Wedge Kit

Includes

- 3 inserts of different durometer
 - Red - Soft
 - Blue - Medium
 - White - Firm
- Adhesive


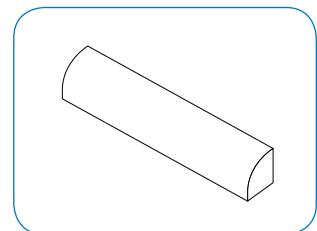


Talux® ST

Components

BOLTS

Picture	Part#	Description	Length	Size	Cat.
1	CM130219	Attachment Bolts for 30mm Tube	1-3/4"	5/16"-24	1-6
	CM130213	Attachment Bolts for 30mm Tube	1-3/4"	3/8"-24	7-9
2	CM130234	Attachment Bolts for Pyramid	30mm	M8	1-6
	CM130223	Attachment Bolts for Pyramid	35mm	M10	7-9

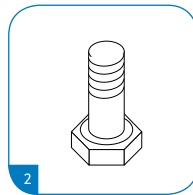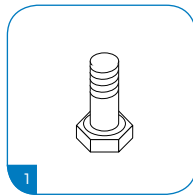

Cosmetic Finishing

See page 146 for details on cosmetic finishing options.

Foot Size (23-30cm)
FCTo  L/R (Left/Right)

Please specify when ordering brown covers. Add "BR" to the part#.

Includes

- Foot Cover
- Attachment Plate
- Flex-Foot Sock


Part#	Description	Size
FCX63006	Flex-Foot Sock	22-25cm
FCX63007	Flex-Foot Sock	26-30cm

Part#	Description
CM100068	Shoe Horn



Talux® ST

Diagram

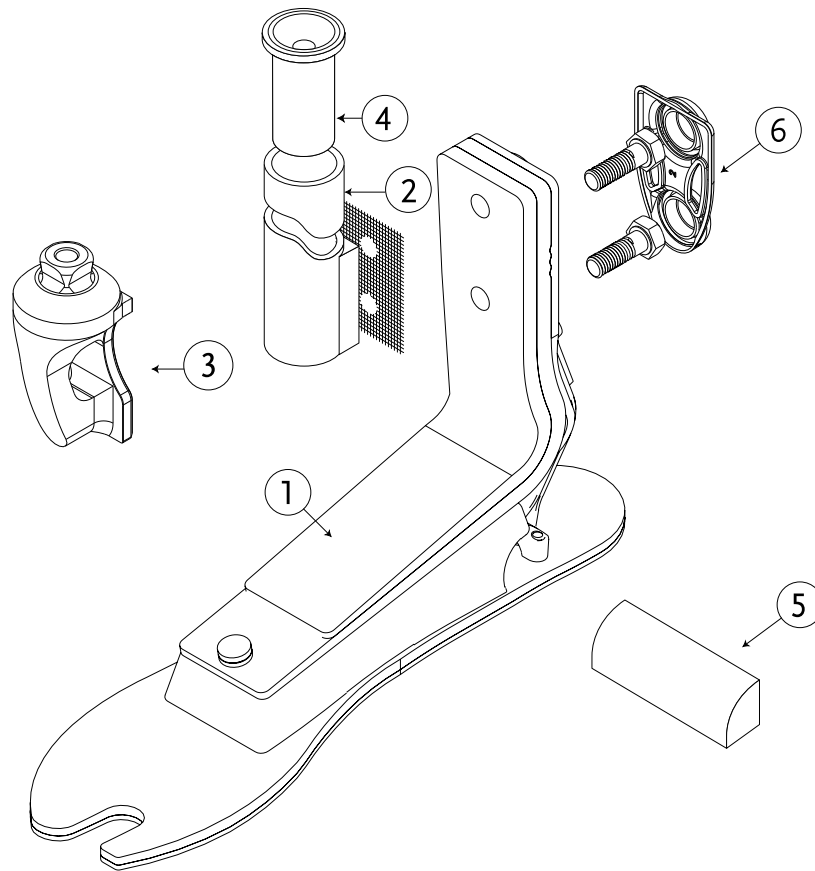

Picture	Description
1	Foot Module
2	30mm Carbon Fibre Tube Kits
3	Bolt-on Male Pyramid
4	Tube Insert
5	Talux Wedge Kit
6	Plastic cap

Vari-Flex®

Introduction



- High energy return
- Lightweight yet extremely strong for maximum performance
- Active CarbonX® heel
- Choice of male pyramid or 30mm carbon tube
- Available in low (10mm) and high (19mm) heel
- For users who walk with a varied cadence and require a reliable, lightweight foot with fast reaction

USER INFORMATION

Maximum Patient Weight	166kg
Impact Level	All

FOOT INFORMATION

Categories	1-9
Sizes	22-30cm
Weight of the foot (Size 27) with male pyramid, sock and foot cover	720g
Features	



The EVO™ feature enables amputees to experience natural gait progression by utilising the foot plate to its fullest potential. A more natural gait reduces fatigue and puts less strain on the lower back and sound side




Full toe lever



Proportional response



Tibial progression

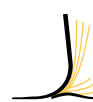

Active heel




Split toe




EVO™

Vari-Flex with male pyramid

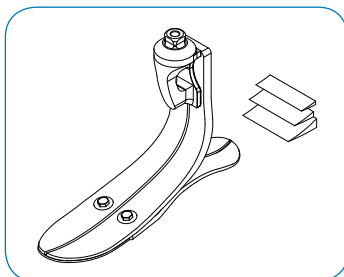


Vari-Flex with EVO kit with male pyramid

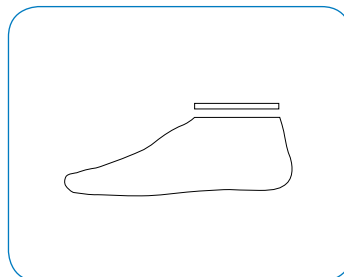


Foot cover kit


Vari-Flex with carbon fiber tube

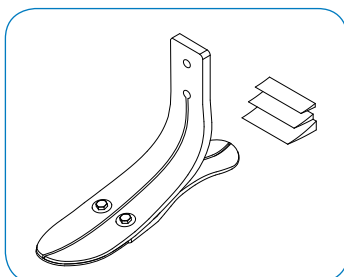

Vari-Flex with EVO kit without pyramid

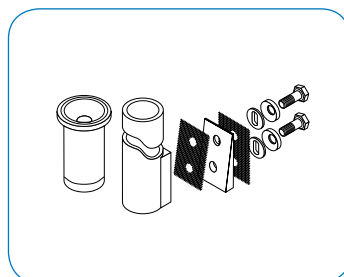

Tube kit

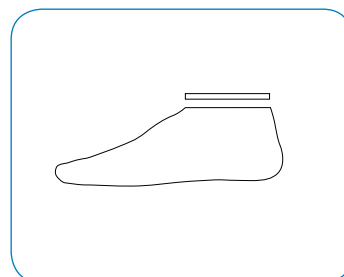


Foot cover kit

Vari-Flex®

Category Selection Guide

Impact Level	Weight Kg									
	44-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147	148-166
Low	1	1	2	3	4	5	6	7	8	9
Moderate	1	2	3	4	5	6	7	8	9	Special
High	2	3	4	5	6	7	8	9	Special	Special

Dimensions

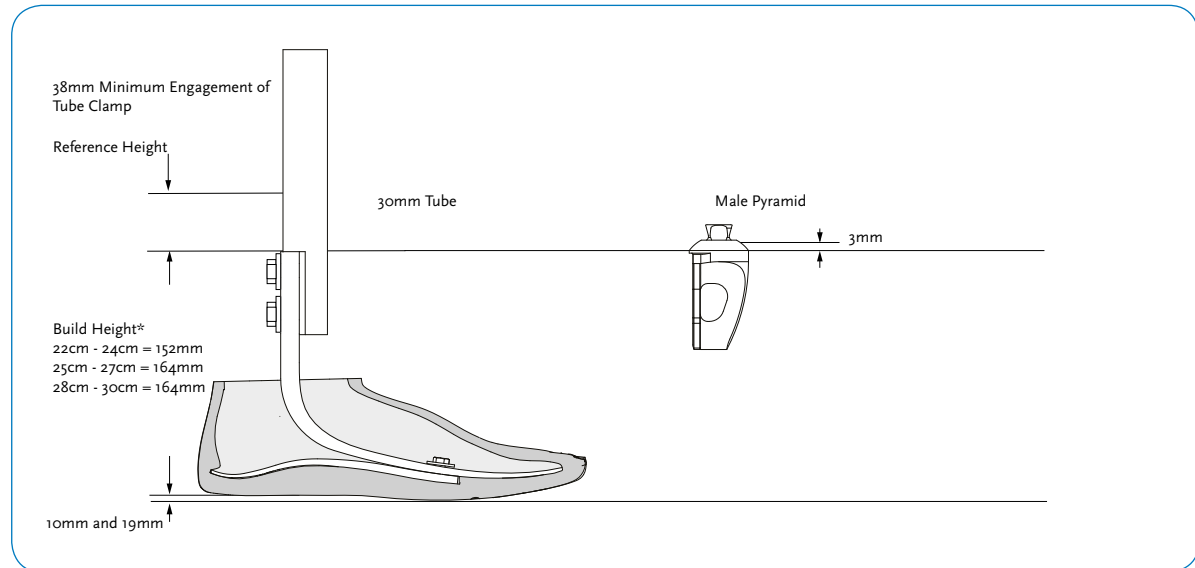

All measurements using 10mm heel height.

Kits

Part#

VARI-FLEX KIT WITH MALE PYRAMID



Category (1-9)
Foot Size (22-30cm)
VFPE [] [] [] [0] 10mm Heel height

Category (1-9)
Foot Size (22-30cm)
VFPE [] [] [] [2] 19mm Heel height

VARI-FLEX KIT WITHOUT ADAPTER

Category (1-9)
Foot Size (22-30cm)
VFXE [] [] [] [0] 10mm Heel height

Category (1-9)
Foot Size (22-30cm)
VFXE [] [] [] [2] 19mm Heel height

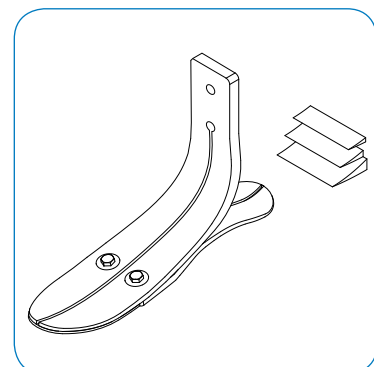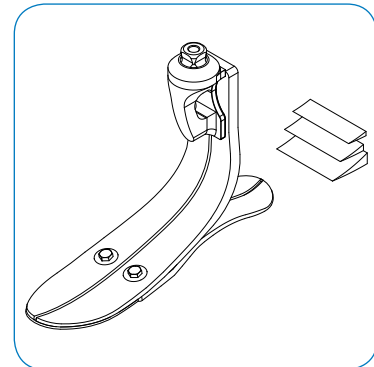

For solid toe special orders please call your Customer Care Representative on 08450 065 065

30MM CARBON FIBRE TUBE KIT

Picture	Part#	Weight	Length	Hole Spacing	Cat.
1	VFX00003	216g	229mm	38mm	1-6
1	VFX00004	275g	229mm	38mm	7-9
1	VFX00006	293g	381mm	38mm	1-6
1	VFX00007	369g	381mm	38mm	7-9

Includes

- Alignment Wedges
- Tube
- Spherical Washers
- Attachment Bolts and Hardware (5/16" UNF)
- Friction Pads
- Tube Insert (Cat. 7-9)

TUBE INSERT

Picture	Part#	Weight	Diameter	Cat.
2	A-718009	18g	30mm	7-9

ALIGNMENT WEDGES FOR TUBE KIT

Picture	Part#	Cat.
3	CM150041	1-3
3	CM150042	4-9

Includes

- 2 Wedges (2° and 4°)
- Friction Pad
- Loctite

BOLT-ON MALE PYRAMID ADAPTERS (METRIC)

Picture	Part#	Weight	Cat.
4	FFC2210001	175g	1-6
4	FFC2210002	205g	7-9

Includes

- Pyramid
- Hardware (M8 or M10 bolts)
- Plastic Bolt Cap
- Loctite

KIT TEMPORARY SLOTTED WEDGE

Picture	Part#	Cat.
5	CM100204	1-9

30MM ALIGNMENT FITTING TOOL KIT

Picture	Part#	Cat.
6	VFX00008	1-9

Includes

- Attachment Hardware
- Alignment Tool



Vari-Flex®

Components

HEEL WEDGE KITS

Part#	Cat.	Foot Size
CM150006	1-9	22-27cm
CM150009	1-9	28-30cm



Includes

- 3 Wedges
- Adhesive




BOLTS


Picture	Part#	Description	Length	Size	Cat.
1	CM150004	Heel Fastener Hardware Kit (Bolts, Washer)	18mm	1/4"	1-6
	CM150005	Heel Fastener Hardware Kit (Bolts, Washer)	18mm	5/16"	7-9
2	CM130206	Heel Fastener Bolts	18mm	1/4"	1-6
	CM130220	Heel Fastener Bolts	18mm	5/16"	7-9
3	CM130219	Attachment Bolts for 30mm Tube	1-3/4"	5/16"-24	1-6
	CM130213	Attachment Bolts for 30mm Tube	1-3/4"	3/8"-24	7-9
4	CM130234	Attachment Bolts for Pyramid	35mm	M8	1-6
	CM130223	Attachment Bolts for Pyramid	35mm	M10	7-9

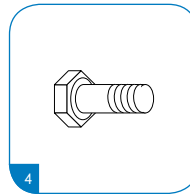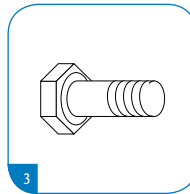

Picture	Part#	Description
1	CM160016	Loctite 410 black
2	CM140003	Heel Divider

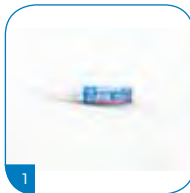

Cosmetic Finishing

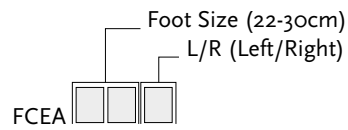

Please specify when ordering brown covers. Add "BR" to the part#.

Includes

- Foot Cover
- Flex-Foot EVO Sock
- Foot Cover Attachment

Part#	Description	Size
FCE63006	Flex-Foot EVO Sock	22-25cm
FCE63007	Flex-Foot EVO Sock	26-30cm


Part#	Description
CM100068	Shoe Horn

See page 146 for details on cosmetic finishing options.



Vari-Flex®

Diagram

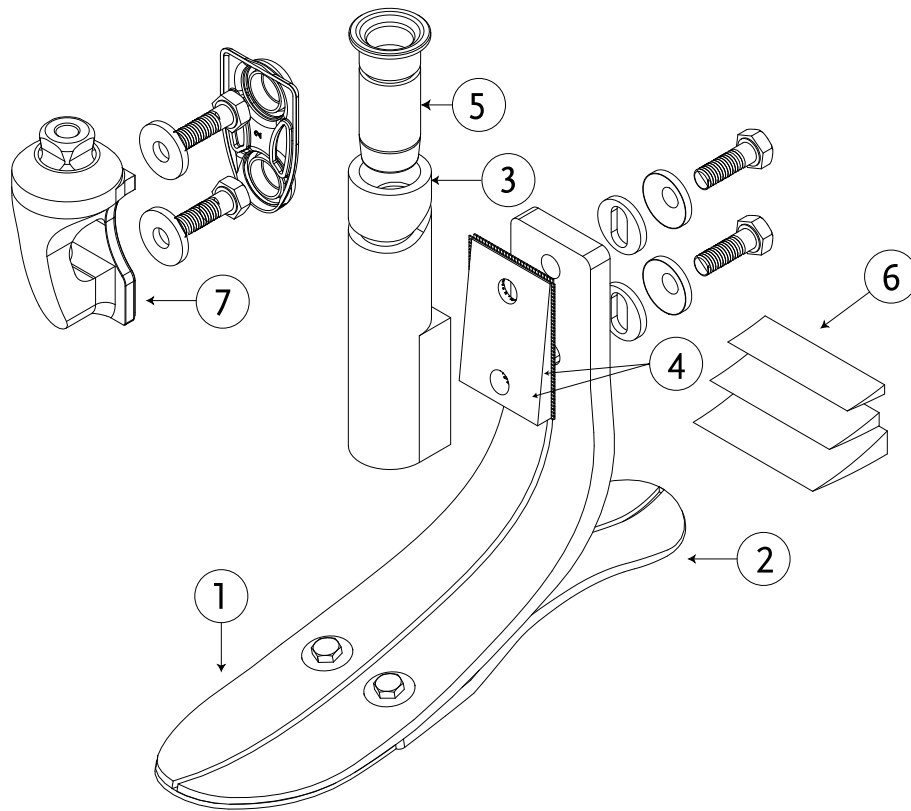

Picture	Description
1	Foot Module
2	Heel
3	30mm Carbon Fibre Tube Kit
4	Alignment Wedges for Tube Kit
5	Tube Insert
6	Heel Wedges
7	Bolt-on Male Pyramid Adapter

LP Vari-Flex®

Introduction



- A low profile option to Vari-Flex for amputees with long residual limbs
- High energy return and superior flexibility
- Active CarbonX® heel
- Lightweight yet extremely strong for maximum performance
- Easy to cosmetically finish and align
- Integrated titanium male pyramid
- Available in low (10mm) and high (19mm) heel
- For users who prefer a variable mobility and require a reliable, lightweight foot with fast reaction


USER INFORMATION

Maximum Patient Weight	166kg
Impact Level	All

FOOT INFORMATION


Categories	1-5	1-9
Sizes	22-24cm	25-30cm
Weight of the the foot (Size 27) with male pyramid, sock and foot cover	590g	

Features

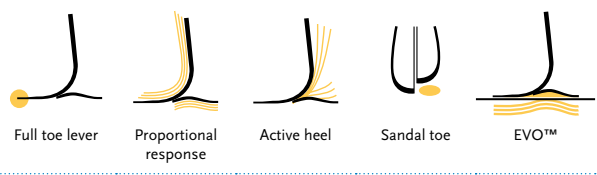


The EVO™ feature enables amputees to experience natural gait progression by utilising the foot plate to its fullest potential. A more natural gait reduces fatigue and puts less strain on the lower back and sound side

LP Vari-Flex with male pyramid

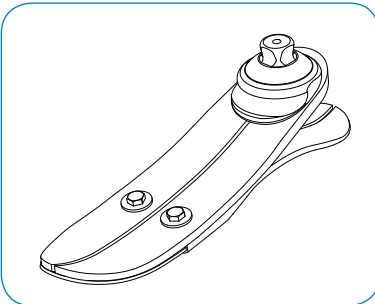


LP Vari-Flex with EVO kit with male pyramid

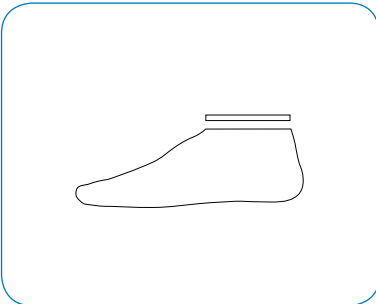


Foot cover kit

LP Vari-Flex with lamination adapter

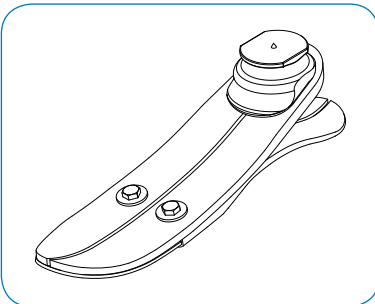

LP Vari-Flex with EVO kit with lamination adapter

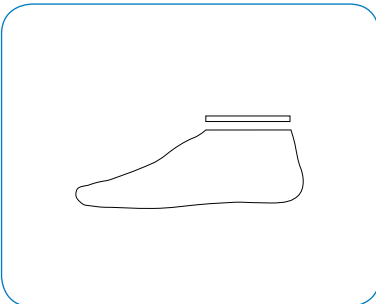


Foot cover kit

LP Vari-Flex®

Category Selection Guide

Impact Level	Weight Kg									
	44-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147	148-166
Low	1	1	2	3	4	5	6	7	8	9
Moderate	1	2	3	4	5	6	7	8	9	N/A
High	2	3	4	5	6	7	8	9	N/A	N/A

Dimensions

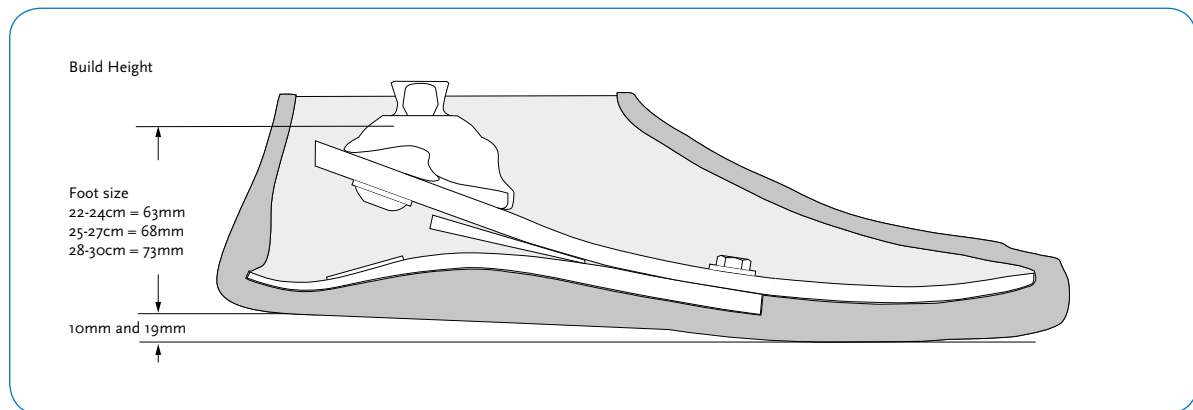

Add 20mm to the build height for LP Vari-Flex with Lamination Adapter.
All measurements using 10mm heel height.


Kits

Part#

LP VARI-FLEX KIT WITH MALE PYRAMID

Category 1-5 (Size 22-24), 1-9 (Size 25-30)
Foot Size (22-30cm)


VLPE 10mm Heel height

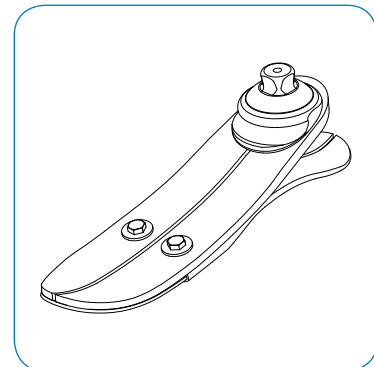

LP VARI-FLEX KIT (foot module only - no adapter included) FOR EXO-SKELETAL APPLICATIONS

Category 1-5 (Size 22-24), 1-9 (Size 25-30)
Foot Size (22-30cm)

VLXE 10mm Heel height

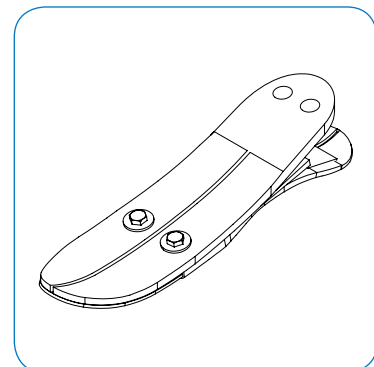

Note: Solid toe is a special order. Please call your Customer Care Representative on 08450 065 065.

LP Vari-Flex®

Components

Part#	Description	Cat.
VLX00001	Lamination Adapter Foot size 22-24	1-5
VLX00002	Lamination Adapter Foot size 25-30	1-6
VLX00003	Lamination Adapter Foot size 25-30	7-9

Includes

- Bolts

HEEL WEDGE KIT

Part#	Cat.	Foot Size
CM150011	1-9	22-30cm

Includes

- 3 Wedges
- Adhesive



HEEL FASTENER HARDWARE KIT

Part#	Cat.	Size	Diameter
CM150004	1-6	22-30cm	1/4"
CM150005	7-9	28-30cm	5/16"


Includes

- Bolts
- Washer

Picture	Part#	Description
1	CM160016	Loctite 410 black
2	CM140003	Heel Divider

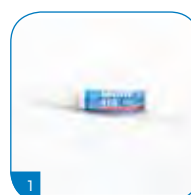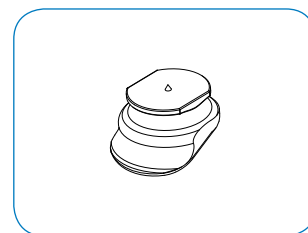

Cosmetic Finishing

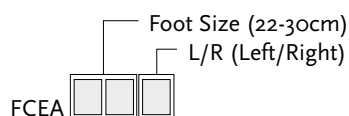

Please specify when ordering brown covers. Add "BR" to the part number.

Includes

- Foot Cover
- Flex-Foot EVO Sock
- Foot Cover Attachment

Part#	Description	Size
FCE63006	Flex-Foot EVO Sock	22-25cm
FCE63007	Flex-Foot EVO Sock	26-30cm


Part#	Description
CM100068	Shoe Horn

See page 146 for details on cosmetic finishing options.



LP Vari-Flex®

Diagram

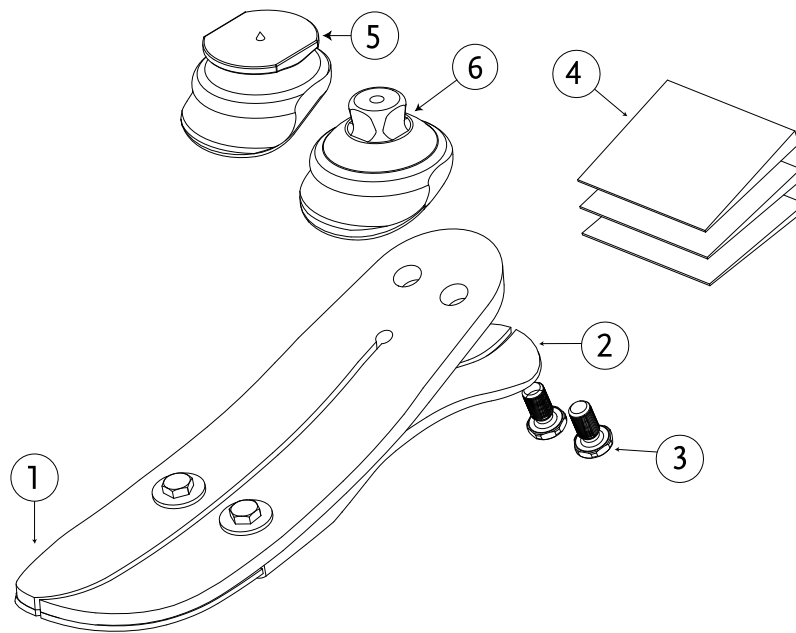

Picture	Description
1	LP Vari-Flex with EVO, Foot Module
2	Heel
3	Bolts
4	Heel Wedges
5	Lamination Adapter
6	Male Pyramid

Vari-Flex® Junior

Introduction



Vari-Flex Junior provides an exclusive combination of comfort and dynamics. It is specially made for those children who require the additional energy response.


Vari-Flex Junior is lightweight, easy to assemble and has a slender profile which makes it easy to cosmetically cover. It ensures the highest levels of user confidence and security. Promoting a natural gait with less fatigue and strain on the lower back and sound side, it is the ideal foot for any activity.

USER INFORMATION

Maximum Patient Weight	55kg
Impact Level	Low to High

FOOT INFORMATION

Categories	1-4
Sizes	16-24cm
Weight of the foot (Size 18) with male pyramid, sock and foot cover	314g
Features	

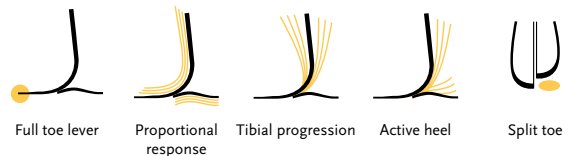

Category Selection Guide


Weight Kg

Sizes	15-24	25-34	35-45	46-55
16	1	N/A	N/A	N/A
17	1	N/A	N/A	N/A
18	1	2	N/A	N/A
19	1	2	N/A	N/A
20	N/A	2	3	N/A
21	N/A	2	3	N/A
22	N/A	2	3	4
23	N/A	2	3	4
24	N/A	N/A	N/A	4

NOTE: The above categories apply only to Vari-Flex Junior

Vari-Flex® Junior

Dimensions

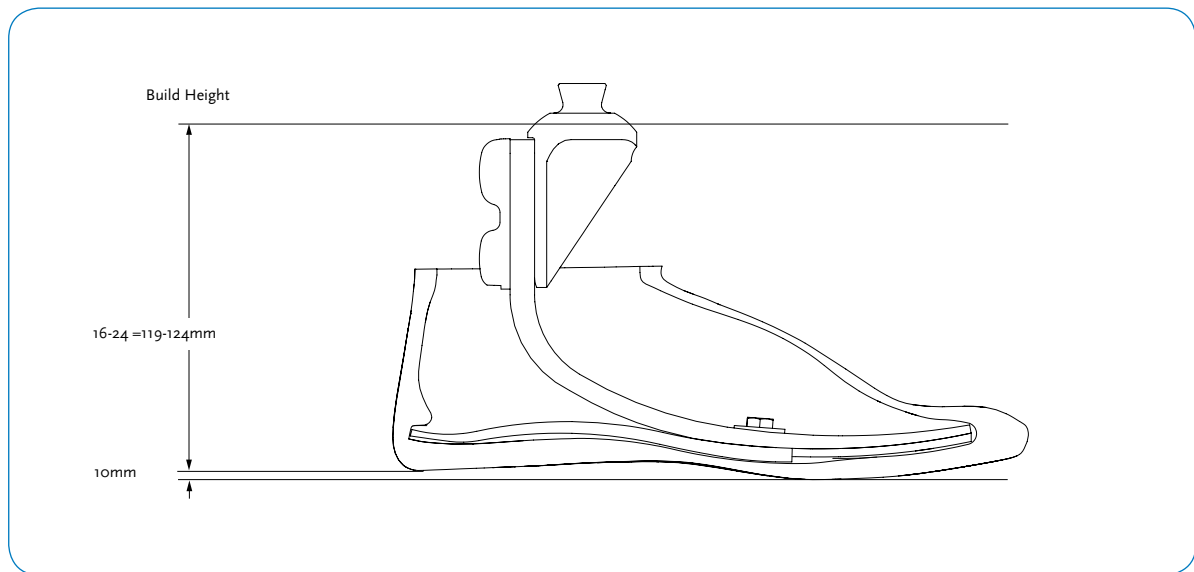


175

Kits

Part#

VARI-FLEX JUNIOR KIT WITH MALE PYRAMID

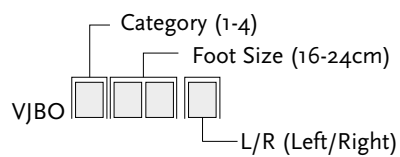

Includes

- Foot Module
- Heel Wedge
- Male Pyramid


Components

JUNIOR HEEL WEDGE KIT

Part#	Cat.	Foot Size
CM150017	All	16-24cm

Includes


- 1 Wedge
- Adhesive

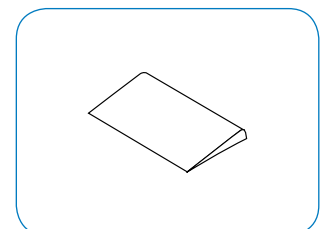


Vari-Flex® Junior

Cosmetic Finishing

SIZES 16-21

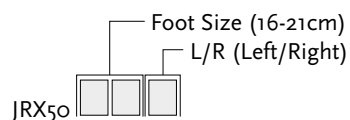

SIZES 22-24

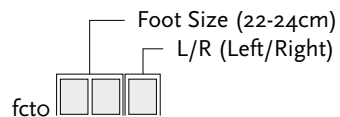

Please specify when ordering brown covers. Add "BR" to the part number.

Includes

- Foot Cover
- Flex-Foot Sock

Part#	Description	Size
FCE63006	Flex-Foot Sock	16-24cm

See page 146 for details on cosmetic finishing options.



Vari-Flex® XC

Introduction



- Vari-Flex XC is a high performance, lightweight foot
- Suitable for a wide range of activities
- C-shape dual blades
- Shock absorbing
- Lightweight bonded pyramid
- EVO™ feature

USER INFORMATION

Maximum Patient Weight	166kg
Impact Level	Low to High

FOOT INFORMATION

Categories	1-9
Sizes	22-30cm
Weight of the foot (Size 27) with male pyramid, sock and foot cover	712g

Features




Full toe lever



Proportional response




Tibial progression

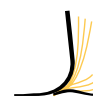

Active heel



Split toe

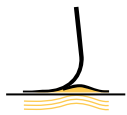

The EVO™ feature enables amputees to experience natural gait progression by utilising the foot plate to its fullest potential. A more natural gait reduces fatigue and puts less strain on the lower back and sound side



EVO™




Katy Farr

Vari-Flex® XC

Category Selection Guide

Impact Level	Weight Kg									
	45-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147	148-166
Low	1	1	2	3	4	5	6	7	8	9
Moderate	1	2	3	4	5	6	7	8	9	N/A
High	2	3	4	5	6	7	8	9	N/A	N/A

Dimensions

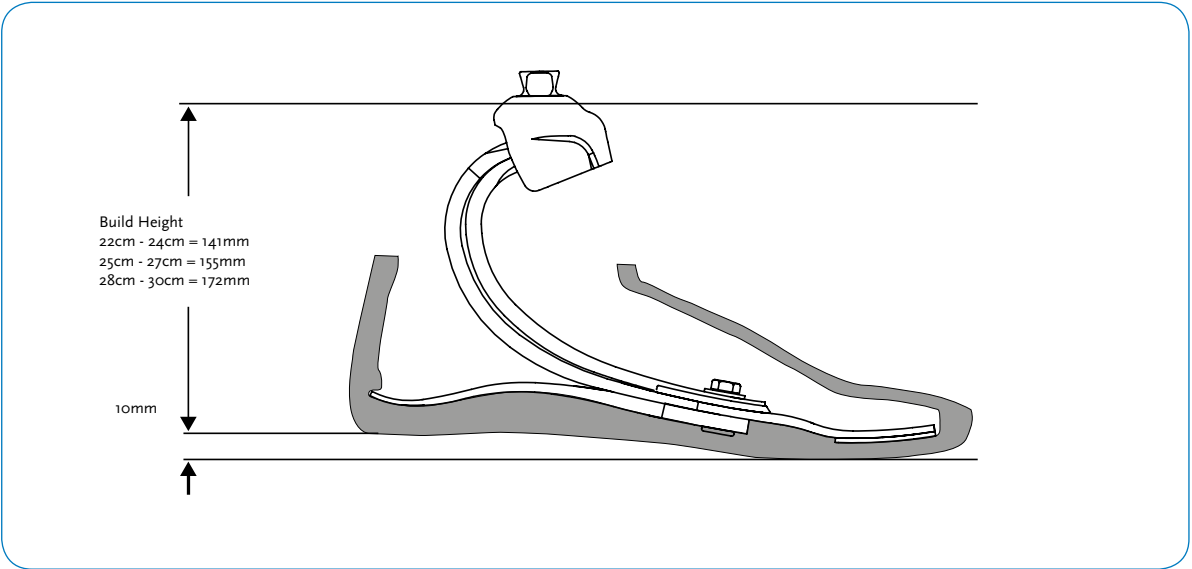

All measurements using 10mm heel height.

Kits

Part#

VARI-FLEX XC

Category (1-9)
 Foot Size (22-30cm)
 XCPE

0


 10mm Heel height

Vari-Flex® XC

Components

Part#

HEEL FOOT MODULE

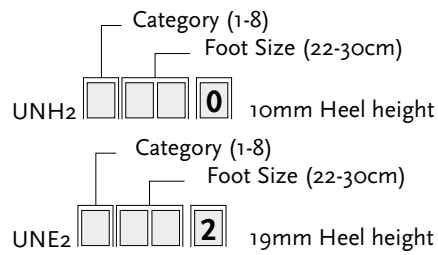

HEEL WEDGE KIT


Part#	Cat.	Foot Size
CM150011	1-9	22-30cm

Includes

- 3 Wedges
- Adhesive




Picture	Part#	Description
1	CM160016	Loctite 410 black
2	CM140003	Heel Divider

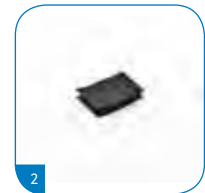

Vari-Flex® XC

Cosmetic Finishing

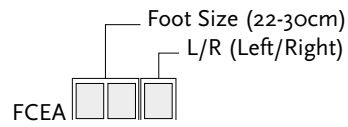

Please specify when ordering brown covers. Add "BR" to the part#.

Includes

- Foot Cover
- Flex-Foot EVO™ Sock
- Foot Cover Attachment


Part#	Description	Size
FCE63006	Flex-Foot EVO Sock	22-25cm
FCE63007	Flex-Foot EVO Sock	26-30cm



See page 146 for details on cosmetic finishing options.

Vari-Flex® XC

Diagram

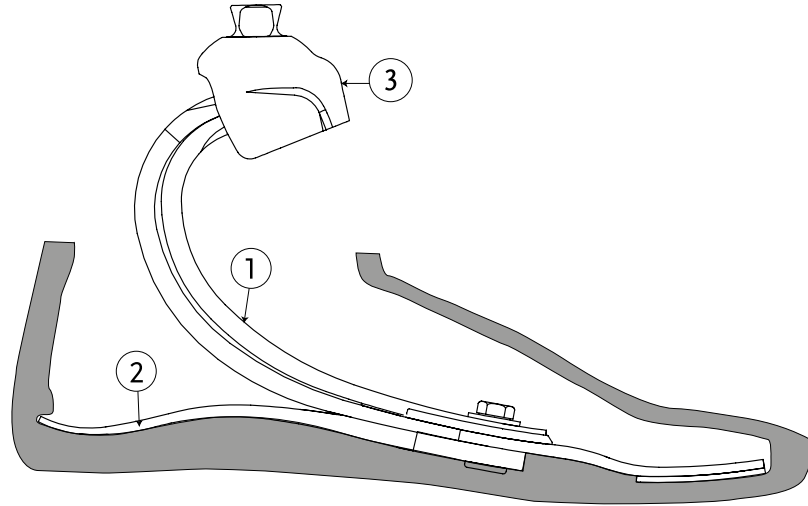

Picture	Description
1	Foot Module
2	Heel
3	Male Pyramid

Vari-Flex® XC Rotate™

Introduction



The EVO™ feature enables amputees to experience natural gait progression by utilising the foot plate to its fullest potential. A more natural gait reduces fatigue and puts less strain on the lower back and sound side

- Vari-Flex XC is a high performance, lightweight foot
- Suitable for a wide range of activities
- C-shape dual blades
- Shock absorbing
- Lightweight bonded pyramid
- EVO™ feature


USER INFORMATION

Maximum Patient Weight	147kg
Impact Level	Low to High

FOOT INFORMATION

Categories	1-8
Sizes	22-30cm
Weight of the foot (Size 27) with male pyramid, sock and foot cover	972g

Features

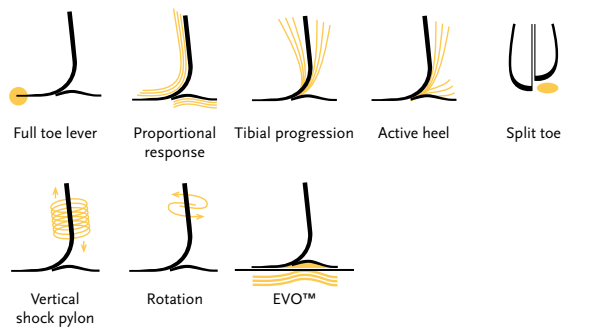


Vari-Flex® XC Rotate™

Category Selection Guide

Weight Kg

Impact Level	45-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147
Low	1	1	2	3	4	5	6	7	8
Moderate	1	2	3	4	5	6	7	8	N/A
High	2	3	4	5	6	7	8	N/A	N/A


Dimensions

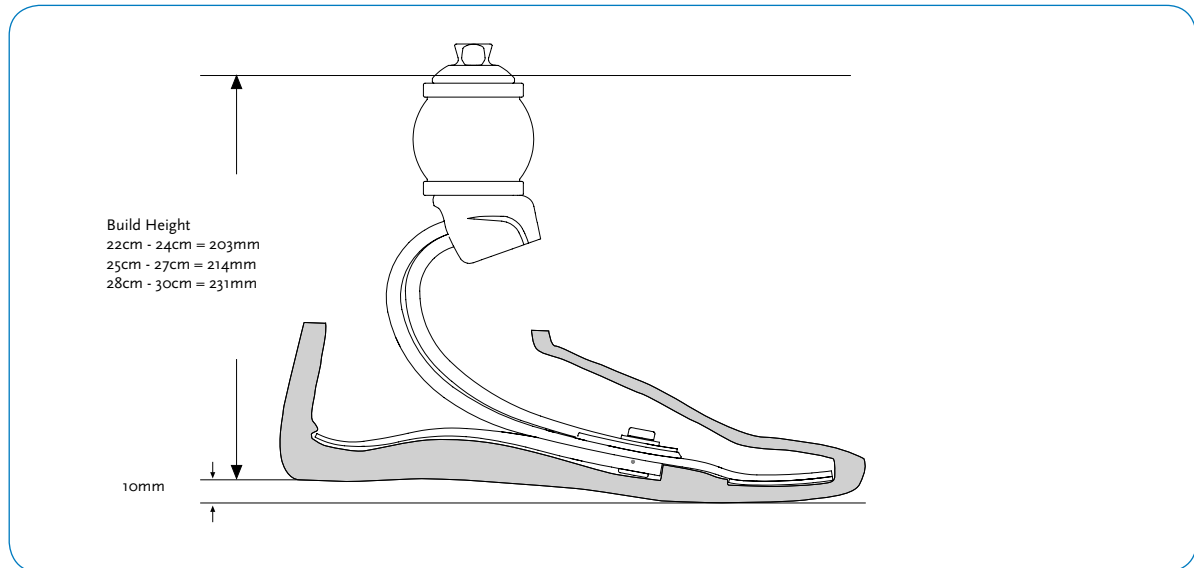

All measurements using 10mm heel height.

Kits

VARI-FLEX XC ROTATE

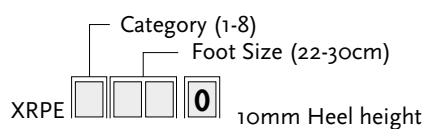

Includes:

- Foot module
- Heel wedges

Components

HEEL WEDGE KIT

Part#	Cat.	Foot Size
CM150011	1-8	22-30cm

Includes


- 3 Wedges
- Adhesive

Picture	Part#	Description
1	CM160016	Loctite 410 black
2	CM140003	Heel Divider



Vari-Flex® XC Rotate™

Cosmetic Finishing

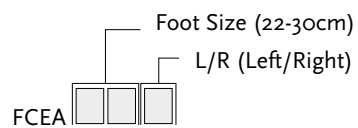

Please specify when ordering brown covers. Add "BR" to the part#.

Includes

- Foot Cover
- Flex-Foot EVO™ Sock
- Foot Cover Attachment

Part#	Description	Size
FCE63006	Flex-Foot EVO Sock	22-25cm
FCE63007	Flex-Foot EVO Sock	26-30cm



See page 146 for details on cosmetic finishing options.

Re-Flex Rotate™

Introduction



- High performance foot for active users who value excellent comfort without compromising on energy return
- Sleek profile that is easy to cosmetically finish
- Active CarbonX® heel
- Male or Female Pyramid
- Torsion Cell designed to reduce both rotary forces on the knee joint and shear forces on the residual limb
- An integrated Torsion Cell allows rotational yield designed to protect the residual limb and knee joint. Shock absorption by titanium coil spring reduces impact on the body while enhancing both control and comfort during walking and other activities



USER INFORMATION

Maximum Patient Weight	147kg
Impact Level	Low to High

FOOT INFORMATION


Categories	1-8
Sizes	22-30cm
Weight of the short assembly (Size 27) with pyramid adapter and foot cover	1136g

Features

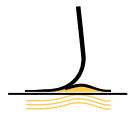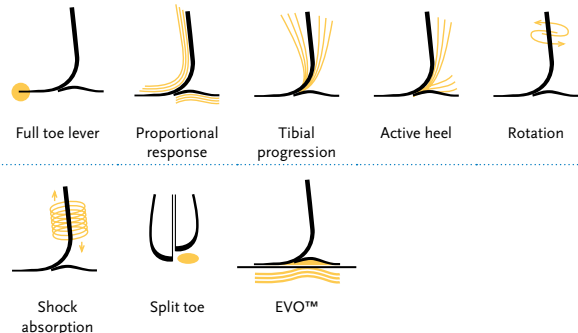


The EVO™ feature enables amputees to experience natural gait progression by utilising the foot plate to its fullest potential. A more natural gait reduces fatigue and puts less strain on the lower back and sound side

Re-Flex Rotate

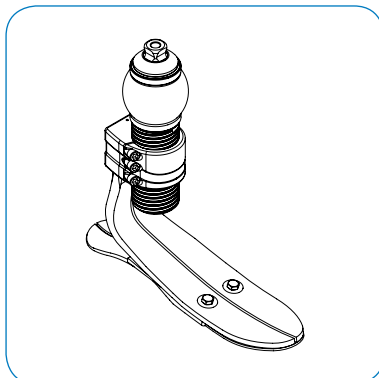

Foot module (short or tall profile)

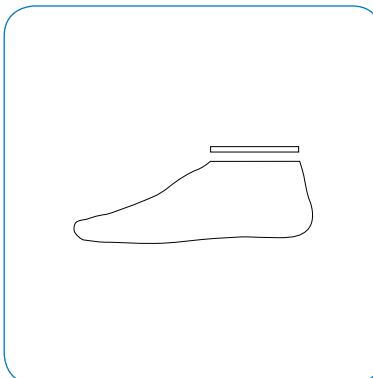



Foot cover kit

Re-Flex Rotate™

Category Selection Guide

Impact Level	Weight Kg								
	44-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147
Low	1	1	2	3	4	5	6	7	8
Moderate	1	2	3	4	5	6	7	8	N/A
High	2	3	4	5	6	7	8	N/A	N/A


Dimensions

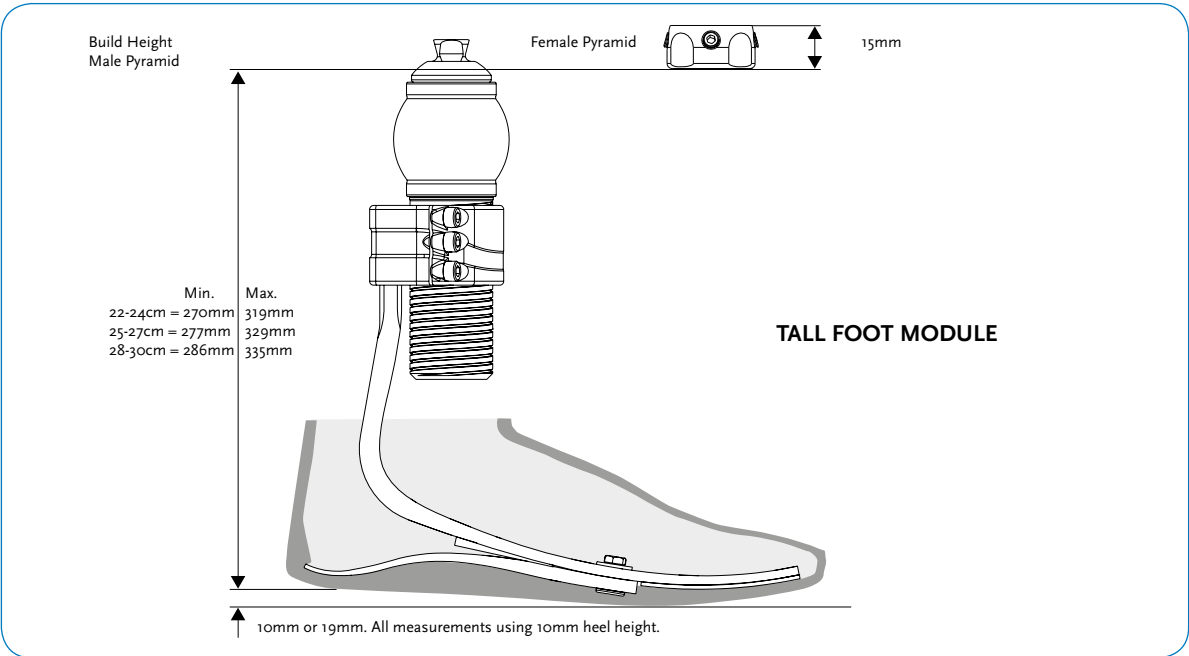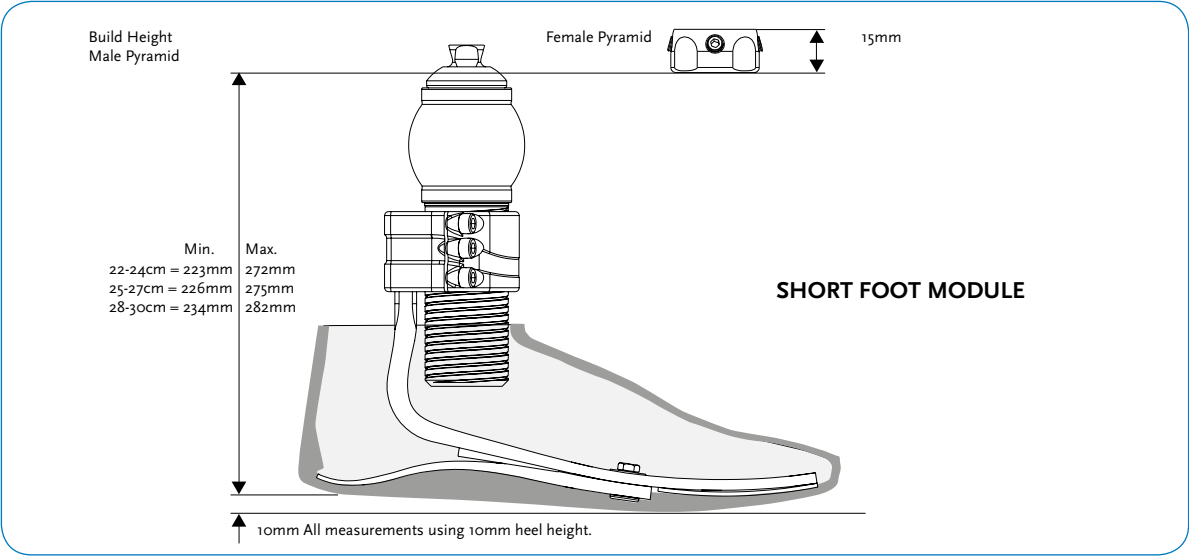


Re-Flex Rotate™

Kits

SHORT PROFILE WITH SHOCK AND MALE PYRAMID


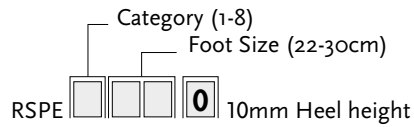


SHORT PROFILE WITH SHOCK AND FEMALE PYRAMID


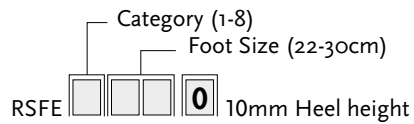


Note: 19mm heel height for short profile, foot size 25-30, available upon special order.

TALL PROFILE WITH SHOCK AND MALE PYRAMID

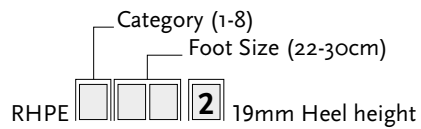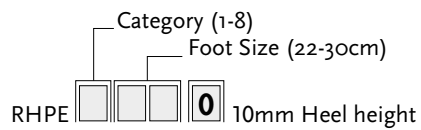





TALL PROFILE WITH SHOCK AND FEMALE PYRAMID

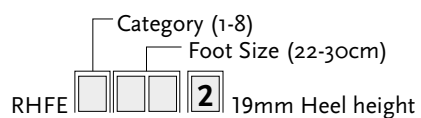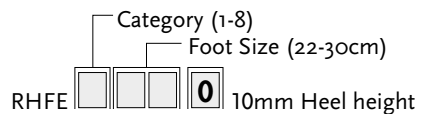

Includes

- Foot Kit
- Shock Kit
- Female/Male Attachment
- Heel Wedges


Note: Solid toe is a special order. Please call your Customer Care Representative on 08450 065 065

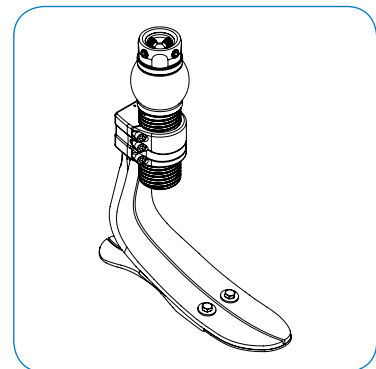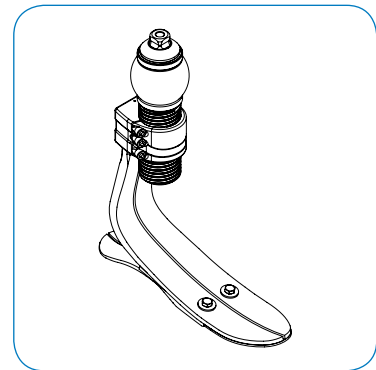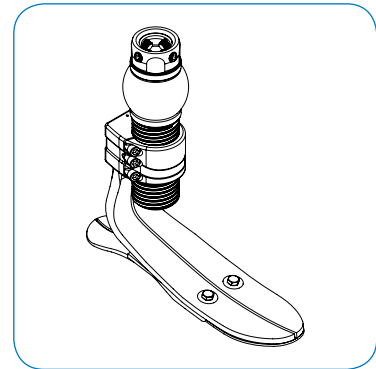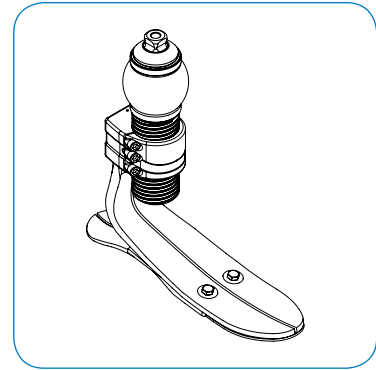

Re-Flex Rotate™


Kits

SHORT PROFILE WITHOUT SHOCK

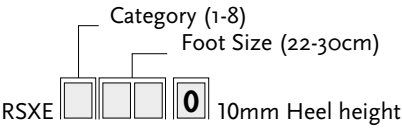

Note: 19mm heel height for standard profile, foot size 25-30, available upon special order.

TALL PROFILE WITHOUT SHOCK

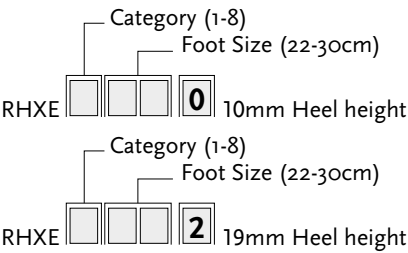


Includes

- Foot Module
- Heel Wedges
- Loctite

Note: Solid toe is a special order. Please call your Customer Care Representative on 08450 065 065

SHOCK KITS

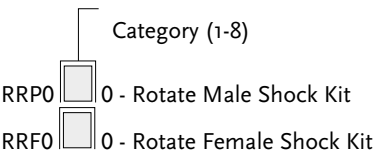

Includes

- Male/Female Pyramid Shock Module
- Loctite (Female)

Components

HEEL WEDGE KITS





Part#	Size	Cat.
CM150011	22-30cm	1-8

Includes

- 3 Wedges
- Adhesive


Tools

Part#	Description
FFC7980001	Torque Wrench 3Nm Preset

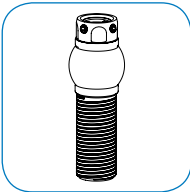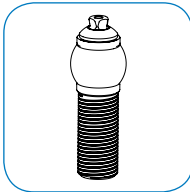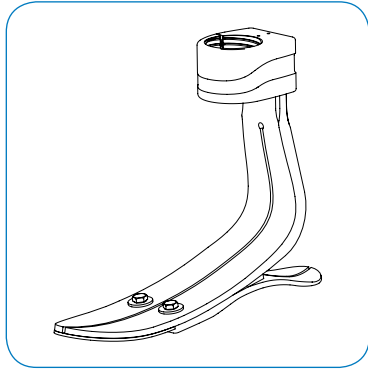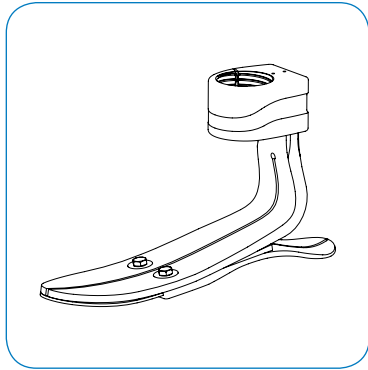

Re-Flex Rotate™

Cosmetic Finishing

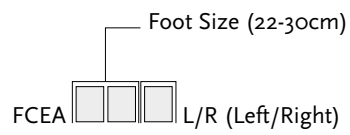

Please specify when ordering brown covers. Add "BR" to the part#.

Includes

- Foot Cover
- Attachment Plate
- Flex-Foot EVO Sock

Part#	Description	Size
FCE63006	Flex-Foot EVO Sock	22-25cm
FCE63007	Flex-Foot EVO Sock	26-30cm


Part#	Description
CM100068	Shoe Horn

See page 146 for details on cosmetic finishing options.



Re-Flex Rotate™

Dimensions

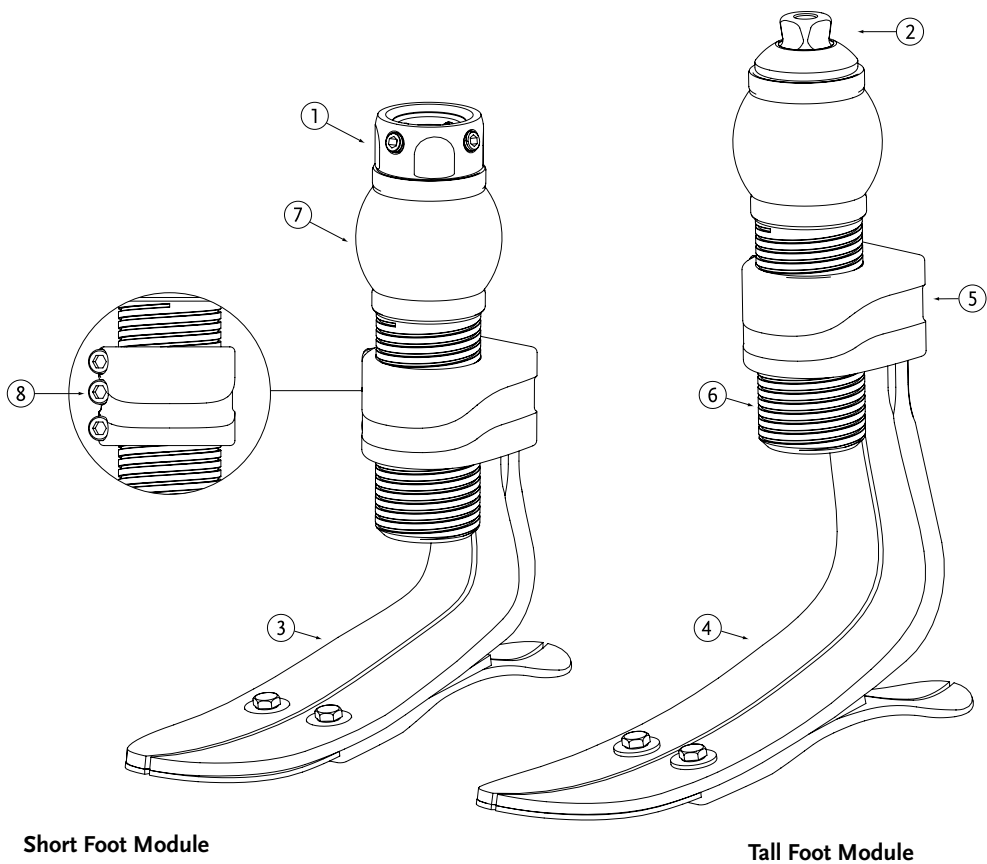

Picture	Description
1	Female Pyramid
2	Male Pyramid
3	Short Foot Kit
4	Tall Foot Kit
5	Fixed Attachment
6	Threaded Outer Tube
7	Polyurethane Torsion Cell
8	Attachment Clamp Screws

LP Rotate™

Introduction



- The LP Rotate has been specially designed for people with long residual limbs, who have insufficient clearance for the Re-Flex Rotate.
- This compact energy-storing foot offers rotational shock absorption designed to reduce shear forces on the residual limb. The Rotate shock module is combined with the LP Vari-Flex foot module. Össur's EVO™ feature* means the full potential of the foot plate is realised with every step, helping to enhance natural gait progression.


USER INFORMATION

Maximum Patient Weight	147kg
Impact Level	Low to High

FOOT INFORMATION

Categories	1-8
Sizes	22-30cm
Weight of the foot (Size 27) with male pyramid	870g

Features

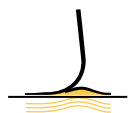


The EVO™ feature enables amputees to experience natural gait progression by utilising the foot plate to its fullest potential. A more natural gait reduces fatigue and puts less strain on the lower back and sound side



Full toe lever



Proportional response

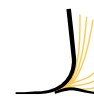

Active heel



Split toe



Rotation




Shock absorption




EVO™


LP Rotate

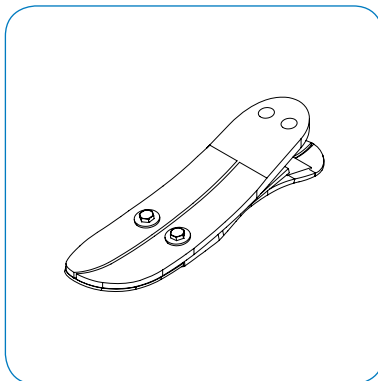

Foot module

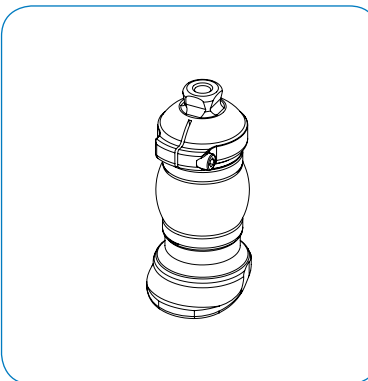

Shock module

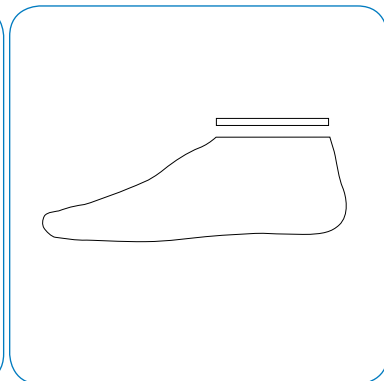


Foot cover kit

LP Rotate™

Category Selection Guide

Impact Level	Weight Kg								
	44-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147
Low	1	1	2	3	4	5	6	7	8
Moderate	1	2	3	4	5	6	7	8	N/A
High	2	3	4	5	6	7	8	N/A	N/A

Dimensions

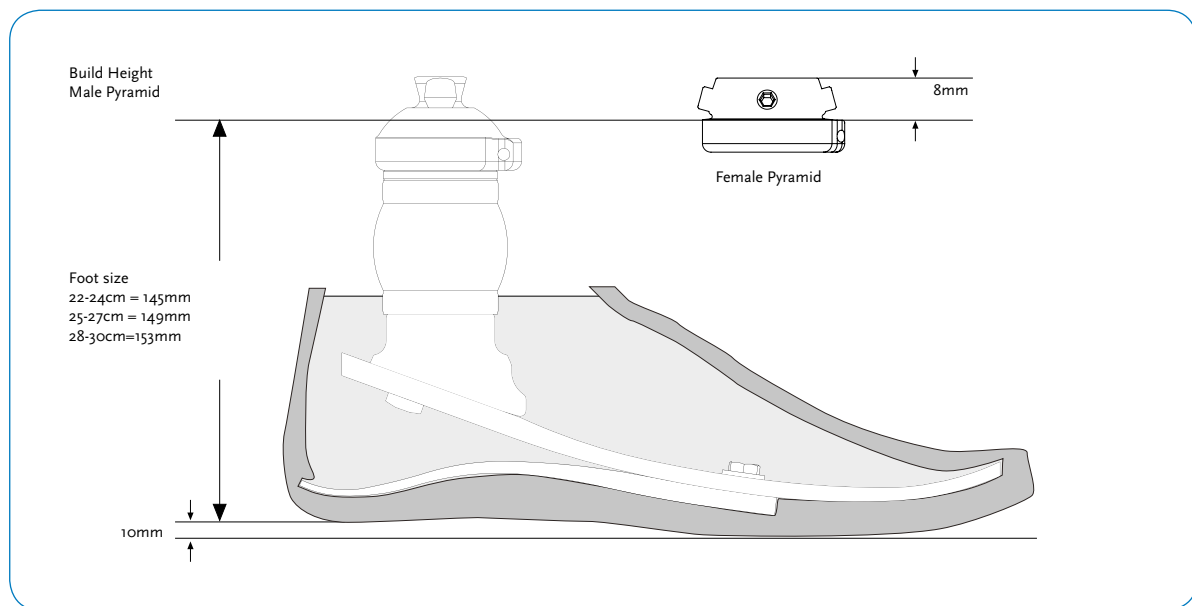

All measurements using 10mm heel height.

Kits


Part#

FOOT MODULE

Category (1-8)
Foot Size (22-30cm)
VLXE 0 10mm Heel height

Includes

- Foot Module
- Heel Module
- Heel Attachment Hardware
- Heel Wedges

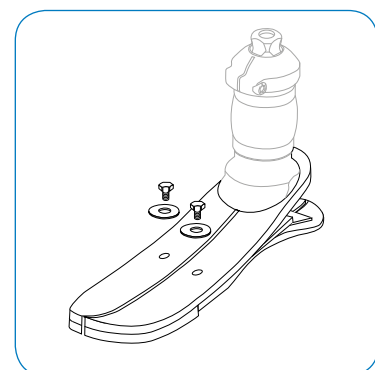

Note: Solid toe is a special order. Please call your Customer Care Representative on 08450 065 065.


LP Rotate™

Shock Modules and Rod Kits

SHOCK MODULE




Part#	Pyramid	Foot Size	Cat.
LRP061	Male	22-24	1-6
LRP062	Male	25-30	1-6
LRP082	Male	25-30	7-8
LRF061	Female	22-24	1-6
LRF062	Female	25-30	1-6
LRF082	Female	25-30	7-8

Includes: Shock module, Male/Female adapter, Rod Kit

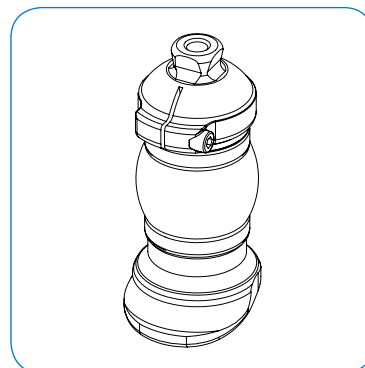

ROD KIT

Part#	Cat.	Colour
LRX50005	1-3	Yellow
	4-6	Blue
	7-8	Red

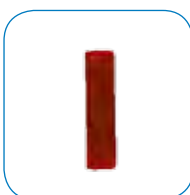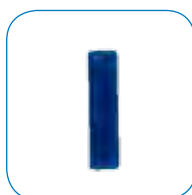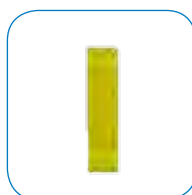

* NB 3 rods in one pack

Components

PYRAMID ADAPTERS

Picture	Part#	Description	Material	Weight limit	Weight
1	A-835300	Male Pyramid	Titanium	147kg	83g
	A-900005	Pinch Bolt	Titanium	N/A	1g
2	A-845300	Female Pyramid	Titanium	166kg	75g



LP Rotate™

Components

HEEL WEDGE KIT

Part#	Foot Size	Cat.
CM150011	22-30cm	1-8


Includes

- 3 Wedges
- Adhesive




BOLTS

Picture	Part#	Description
1	CTL30012	Adjustment Plug
2	CM140003	Heel Divider
3	CM160016	Loctite 410 black

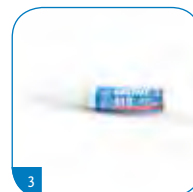

Cosmetic Finishing

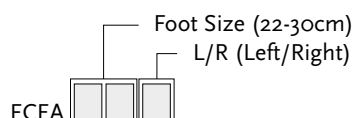

Please specify when ordering brown covers. Add "BR" to the part#.

Includes

- Foot Cover
- Attachment Plate
- Flex-Foot EVO Sock

Part#	Description	Size
FCE63006	Flex-Foot EVO Sock	22-25cm
FCE63007	Flex-Foot EVO Sock	26-30cm


Part#	Description
CM100068	Shoe Horn

See page 146 for details on cosmetic finishing options.



LP Rotate™

Diagram

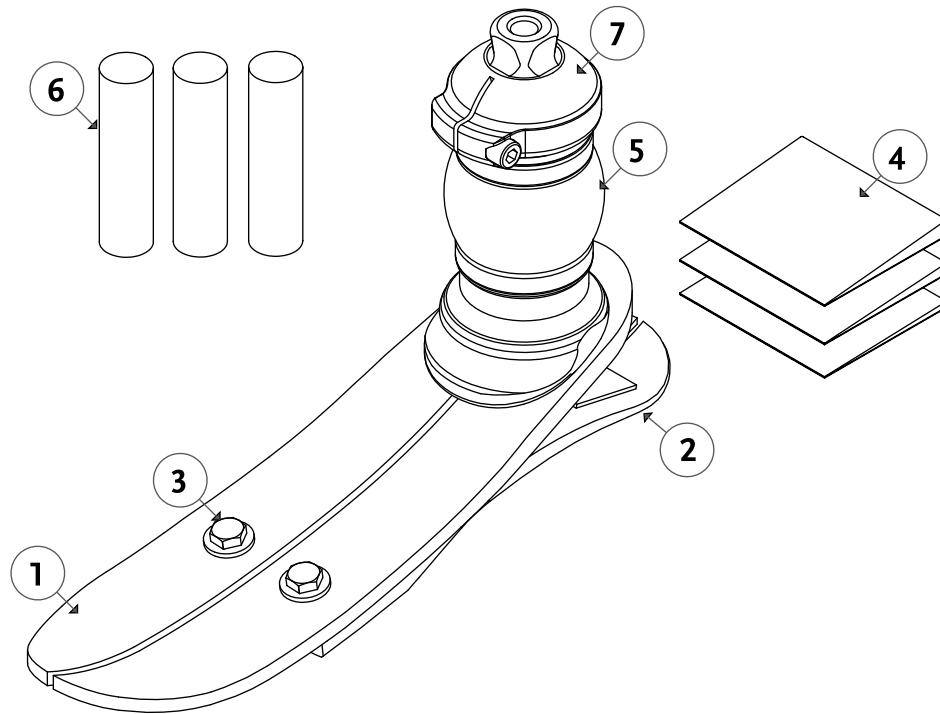

Picture	Description
1	Foot Module
2	Heel
3	Bolts
4	Heel Wedges
5	Shock Module
6	Rod Kit
7	Proximal Adapter

Vari-Flex® Modular™

Introduction



- High energy storage and release
- Lightweight yet extremely strong for maximum performance
- Available in low and high heel options
- Active CarbonX® heel

USER INFORMATION

Maximum Patient Weight	166kg +
Impact Level	All

FOOT INFORMATION

Categories	1-9 +
Sizes	22-30cm
Weight of the foot (Size 27)	403g

Features



Full toe lever Proportional response Active tibial progression Active heel Split toe




EVO™




The EVO™ feature enables amputees to experience natural gait progression by utilising the foot plate to its fullest potential. A more natural gait reduces fatigue and puts less strain on the lower back and sound side


Vari-Flex Modular

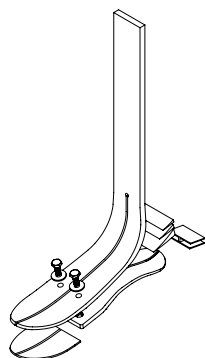

Modular III kit

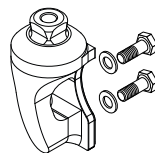

Bolt-on male pyramid

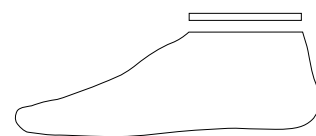


Foot cover kit

Vari-Flex® Modular™

Category Selection Guide

Impact Level	Weight Kg									
	44-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147	148-166
Low	1	1	2	3	4	5	6	7	8	9
Moderate	1	2	3	4	5	6	7	8	9	Special order
High	2	3	4	5	6	7	8	9	Special order	Special order

Dimensions

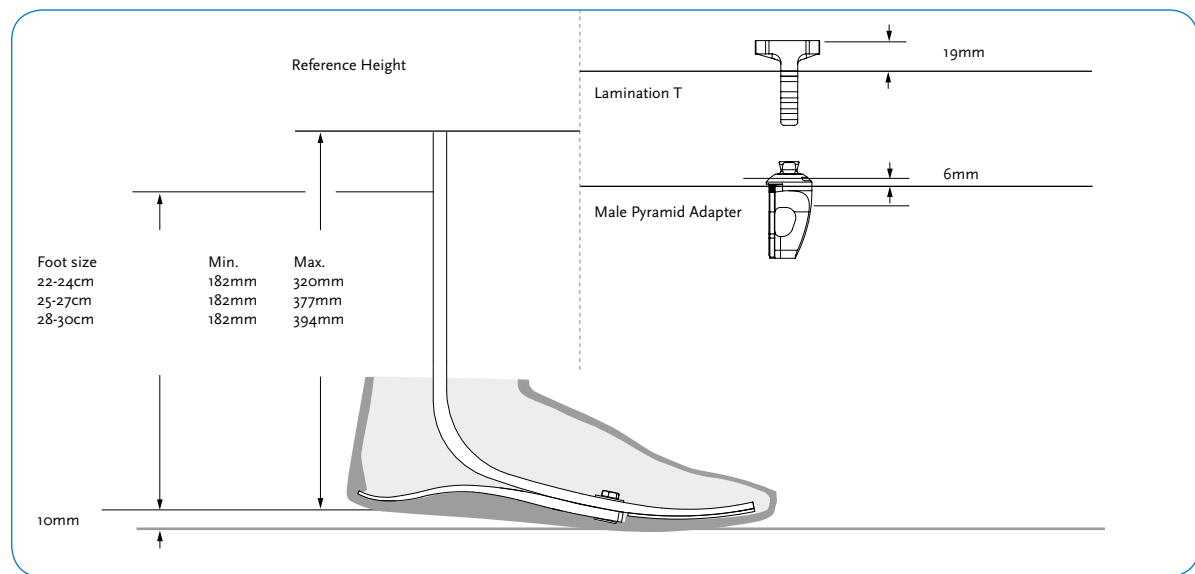

Kits


Part#

FOOT MODULE

Category (1-9)
Foot Size (22-30cm)
MDXE 10mm Heel height

Includes

- Foot Module
- Heel Module
- Heel Attachment Hardware
- Heel Wedges

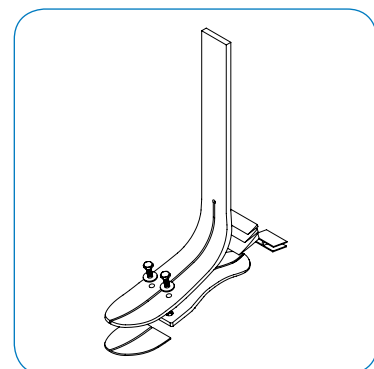

Note: Solid toe is a special order. Please call your Customer Care Representative on 08450 065 065.

Vari-Flex® Modular™

Standard Adapters

BOLT-ON MALE PYRAMID ADAPTER (METRIC)

Part#	Weight	Cat.
FFC2210001	175g	1-6
FFC2210002	205g	7-9

Includes

- Pyramid (bolt centres 38mm)
- Hardware (M8 or M10 bolts)
- Loctite



ALUMINIUM T KITS

Part#	Weight	Cat.
MDX00020	360g	1-6
MDX00021	360g	7-9

Includes

- Aluminium T (bolt centres 25.4mm)
- Alignment Wedges
- Attachment Bolts
- M6 Bolts, Washers and Nuts
- Spherical and Round Washers
- Abrasive Pads
- Hex Nuts



ALUMINIUM T CONNECTORS


Part#	Weight	Cat.
CM140326	360g	1-6
CM140327	360g	7-9



Components

Part#

HEEL FOOT MODULE

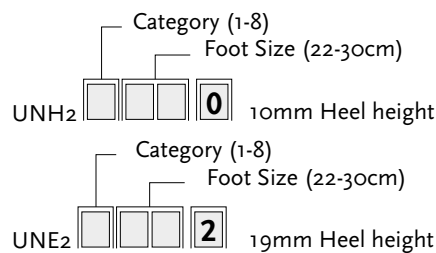

Part#	Description
MDX50007	Direct Lamination Kit Cat 10+

Includes

- Fibreglass Cloth
- Potting Epoxy





Vari-Flex® Modular™

Standard Adapters



BOLTS

Picture	Part#	Description	Length	Size	Cat.
1	CM130223	Attachment Bolt For Pyramid	35mm	M10	7-9
2	CM130237	Male Pyramid Connector Bolts	30mm	M8	1-6

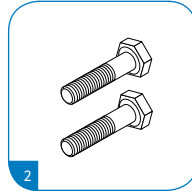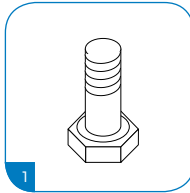

WEDGES

Picture	Part#	Description	Cat.
1	CM100203	Kit Temporary Slotted Wedge, Cat 1-9	1-9
2	CM150051	Alignment Wedges/Spacer Block Kit	1-6
2	CM150052	Alignment Wedges/Spacer Block Kit	7-9

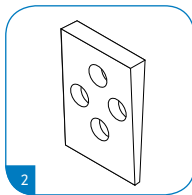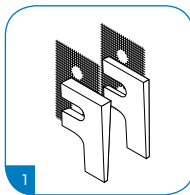

HEEL WEDGE KITS


Part#	Foot Size	Cat.
CM150006	22-27cm	1-9
CM150009	28-30cm	1-9

Includes

- 3 Wedges
- Adhesive




Picture	Part#	Description
1	CM160016	Loctite 410 black
2	CM140003	Heel Divider
3	CM160009	Kit 3oz Potting Epoxy

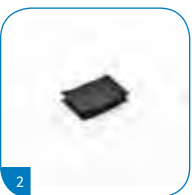

Vari-Flex® Modular™

Cosmetic Finishing

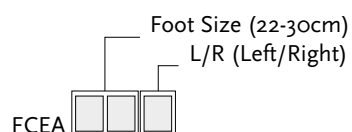

Please specify when ordering brown covers. Add "BR" to the part#.

Includes


- Foot Cover
- Attachment Plate
- Flex-Foot Sock

Part#	Description	Size
FCX63006	Flex-Foot Sock	22-25cm
FCX63007	Flex-Foot Sock	26-30cm

Part#	Description
CM100068	Shoe Horn

See page 146 for details on cosmetic finishing options.



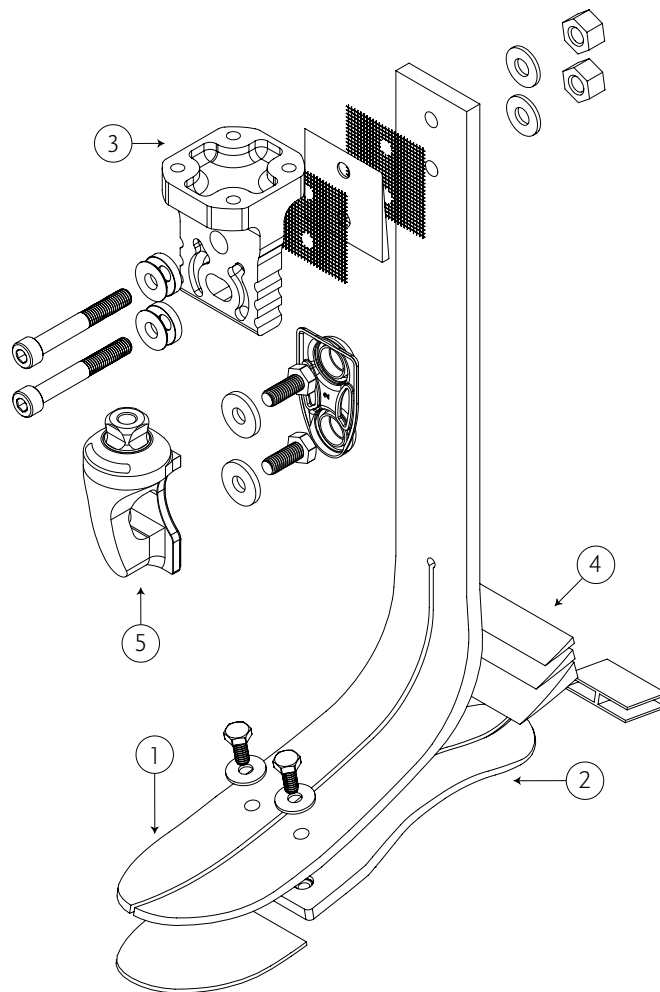

Picture	Description
1	Foot Module Kit
2	Low/High Heel
3	Aluminium T Connector
4	Wedge Kit
5	Bolt-on Male Pyramid


Re-Flex Shock™

Introduction




The EVO™ feature enables amputees to experience natural gait progression by utilising the foot plate to its fullest potential. A more natural gait reduces fatigue and puts less strain on the lower back and sound side


- Re-Flex Shock offers the best vertical shock absorption of all the Össur mechanical feet. With dynamic energy return and all around responsiveness the foot is a great choice for active users.
- The composite spring at the front is dedicated to reduce shear forces on the residual limb by offering optimal shock absorption that reduces impact subjected to the body. Enhancing both control and comfort during walking and other activities.
- Incorporation of the EVO™ feature enables users to experience natural gait progression by utilising the foot plate to its fullest potential. A more natural gait reduces fatigue and puts less strain on the lower back and sound side.

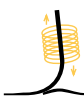USER INFORMATION				
Maximum Patient Weight	166kg			
Impact Level	Low to High			
FOOT INFORMATION				
Categories	1-9			
Sizes	22-30cm			
Weight of the standard foot module (Size 27) with pyramid adapter and foot cover	1048g			
Features				
				
Full toe lever	Proportional response	Active tibial progression	Active heel	Vertical shock pylon
				
Split toe	EVO™			

Re-Flex Shock

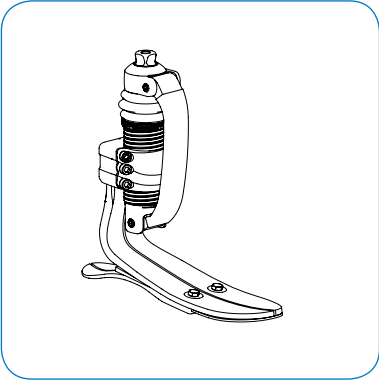

Foot module (standard or tall profile)

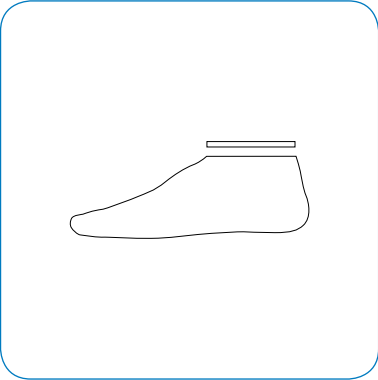



and Foot cover kit

Re-Flex Shock™

Category Selection Guide

Impact Level	Weight (kg)									
	45-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147	148-166
Low	1	1	2	3	4	5	6	7	8	9
Moderate	1	2	3	4	5	6	7	8	N/A	N/A
High	2	3	4	5	6	7	8	N/A	N/A	N/A

Clearance

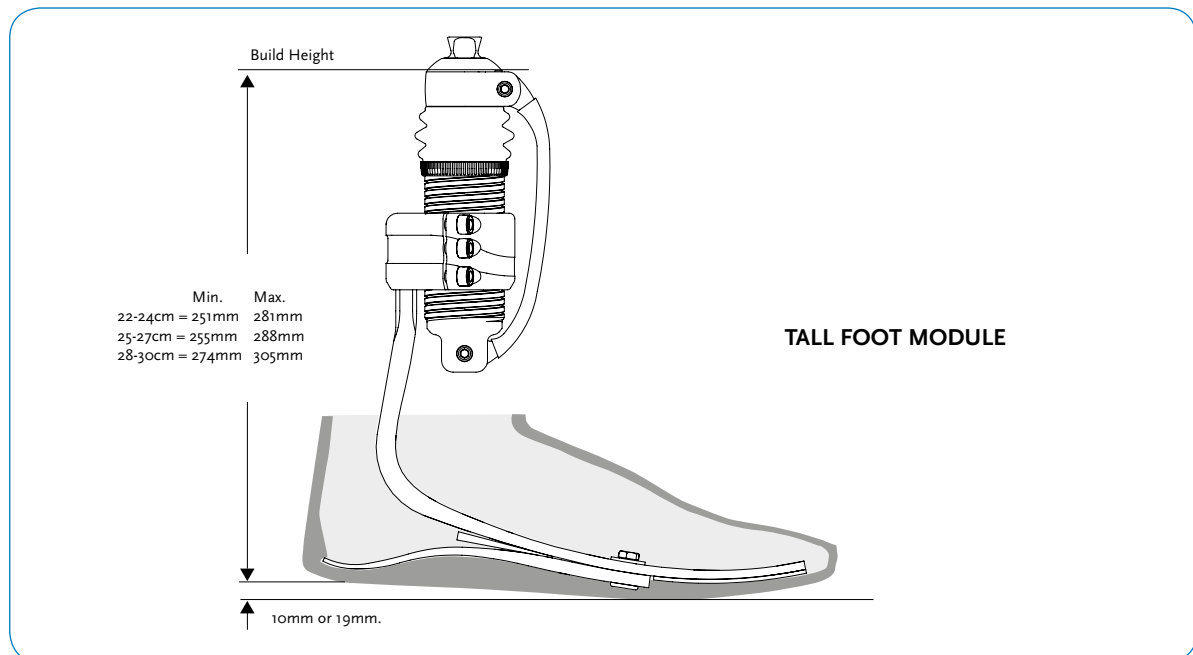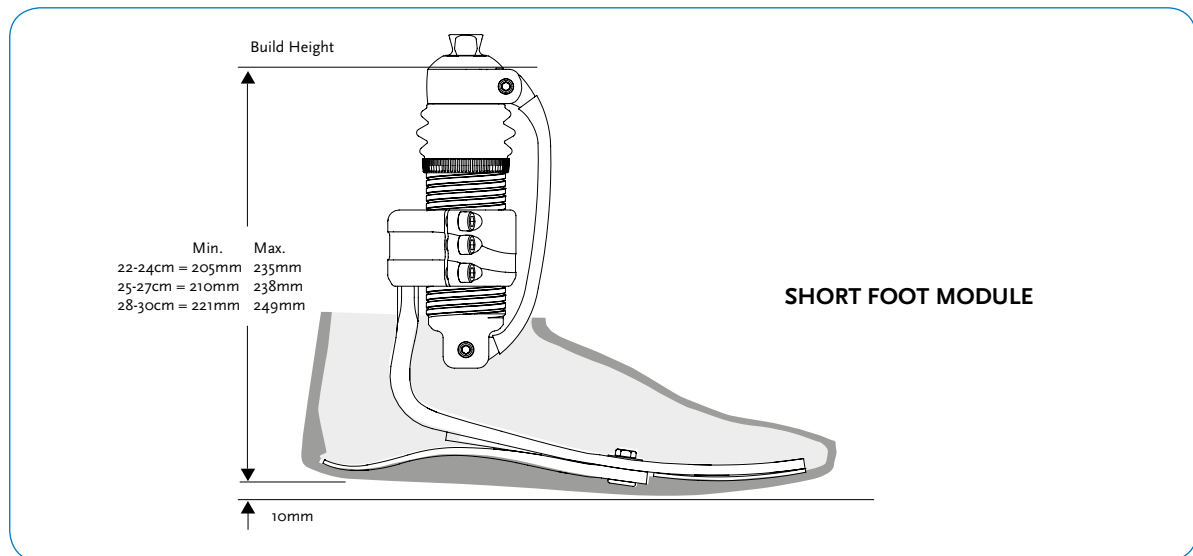


All measurements using 10mm heel height.


Re-Flex Shock™

Kits


SHORT PROFILE WITH SHOCK AND MALE PYRAMID

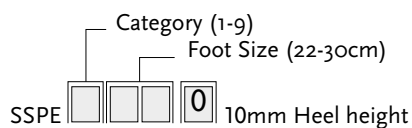

Note: 19mm heel height for standard profile, foot size 25-30, available upon special order

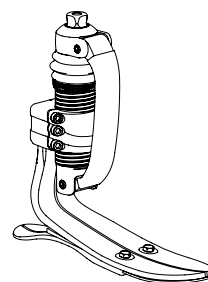


TALL PROFILE WITH SHOCK AND MALE PYRAMID

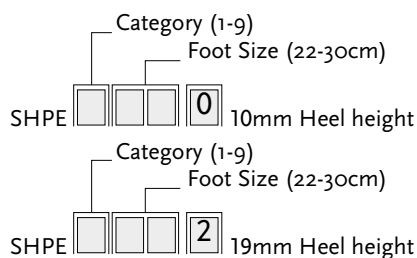


Includes:

- Foot Kit
- Shock Kit
- Heel Wedges
- Torque Wrench


Note: Solid toe is a special order

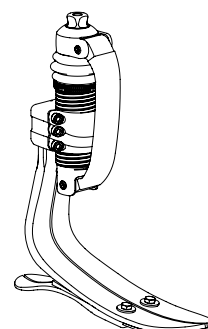


SHORT PROFILE WITHOUT SHOCK

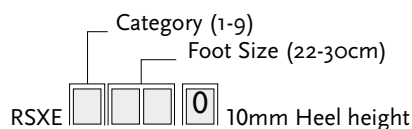

Note: 19mm heel height for standard profile, foot size 25-30, available upon special order.

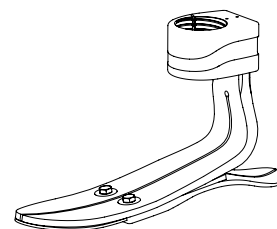


TALL PROFILE WITHOUT SHOCK

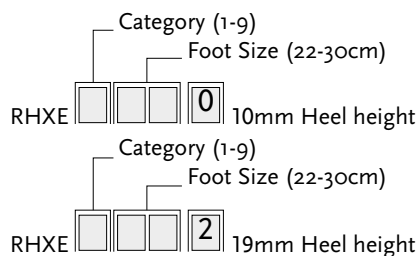

Includes

- Foot Kit
- Torque Wrench


Note: Solid toe is a special order. Please call your Customer Care Representative on 08450 065 065.

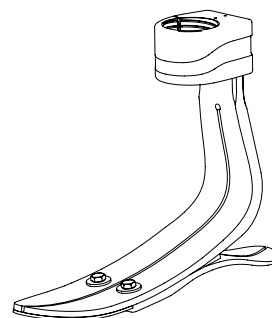

Re-Flex Shock™


Kits

SHOCK KITS

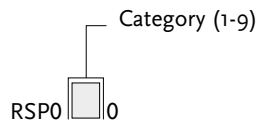

Includes:

- Male pyramid shock module, front spring
- Loctite

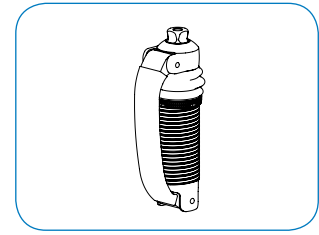

Components

HEEL WEDGE KITS

Part#	Size	Cat.
CM150011	22-30cm	1-9

Includes

- 3 Wedges
- Adhesive



Part#	Description
CM160025	Loctite 243
RSS00001	Pin Kit
RSS40011	Bellow
RSS40012	Collar

Tools

Part#	Description
FFC7980001	Torque Wrench 3Nm Preset



Cosmetic Finishing



Please specify when ordering brown covers. Add "BR" to the part#.

Includes

- Foot Cover
- Attachment Plate
- Flex-Foot EVO Sock

Part#	Description	Size
FCE63006	Flex-Foot EVO Sock	22-25cm
FCE63007	Flex-Foot EVO Sock	26-30cm


Part#	Description
CM100068	Shoe Horn

See page 146 for details on cosmetic finishing options.



Re-Flex Shock™

Diagram

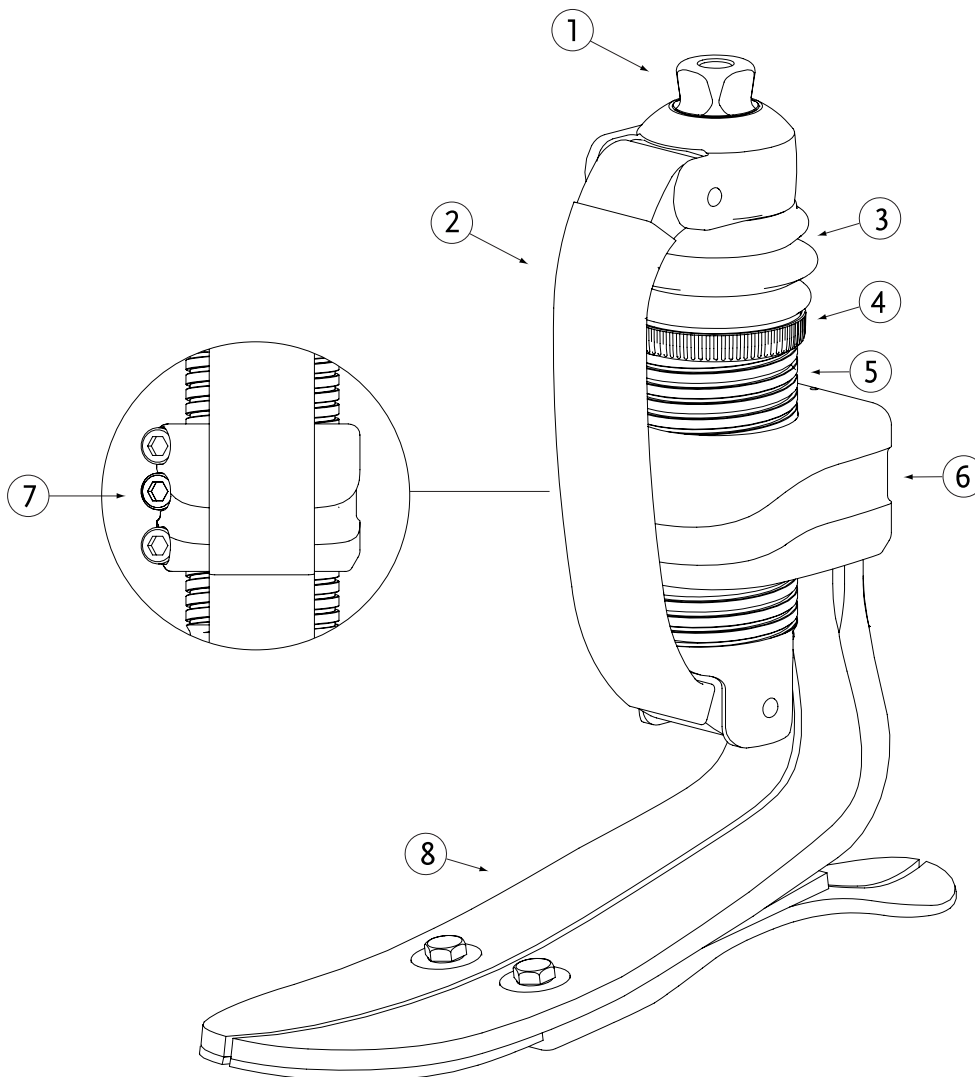

Picture	Description
1	Male Pyramid
2	Front Spring
3	Bellow
4	Collar
5	Threaded Outer Tube
6	Fixed Attachment
7	Attachment Clamp Screws
8	Foot Module

Elation®

Introduction



Elation offers precision heel height adjustment to match the foot to the chosen shoe - over a range from 0 - 5cm. Adjustment is simple and secure.

- Heel height adjustable foot: 0 - 5cm
- Energy storing with progressive stiffening
- Cosmetically appealing cover with sandal toe
- Easy to use

USER INFORMATION


Maximum Patient Weight	100kg
Impact Level	Low and Moderate

FOOT INFORMATION


Sizes	22-28cm
-------	---------

Weight of the foot (Size 25) with male pyramid	625g
---	------

Features

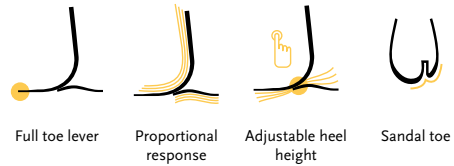


Elation with male pyramid

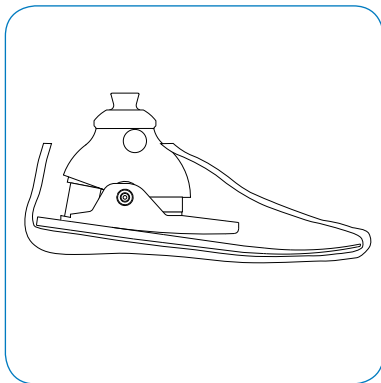

Elation kit with male pyramid

Dimensions

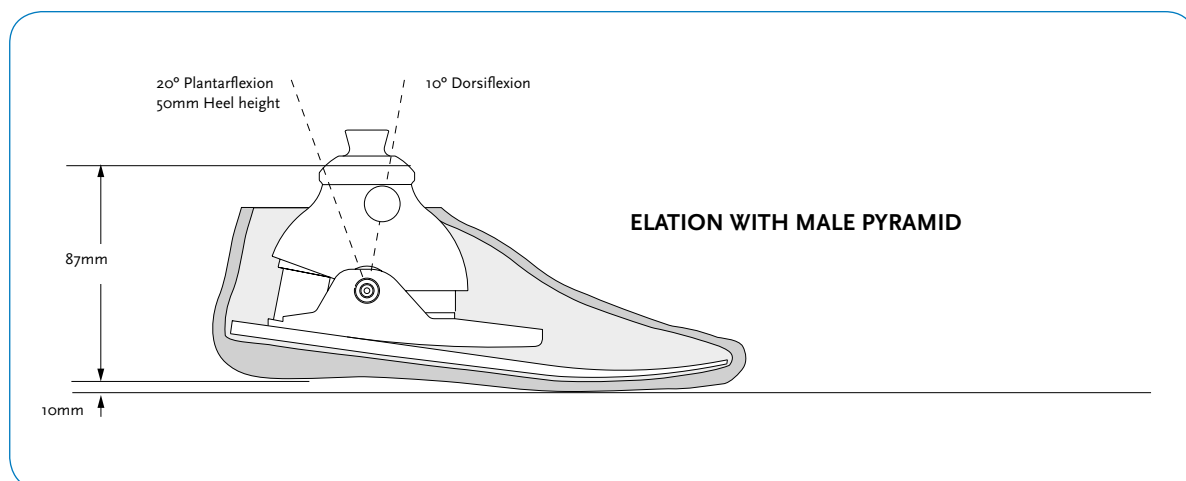

All measurements using 10mm heel height.

Kits

Part#

ELATION KIT WITH MALE PYRAMID

Foot Size (22-30cm)
ELPoo  L/R (Left/Right)

Includes

- Foot/Ankle module fully assembled with front and rear bumpers (black rear bumper) with bonded foot cover
- 3 extra rear bumpers (green, red and orange) included

*For Elation with brown foot cover add "BR" to the part#.

Elation®

Components

MALE PYRAMID ANKLE KIT

Left Foot	Right Foot	Foot Size
ELP00001	ELP00002	22-23cm
ELP00003	ELP00004	24-25cm
ELP00005	ELP00006	26-28cm



Includes

- Ankle for Elation foot module with front and rear bumpers (black rear bumper)
- 3 extra bumpers (green, red and orange)
- Attachment Hardware

FOOT MODULE



Includes

- Foot Module integrated into foot cover
- *For Elation with brown foot cover add "BR" to the part#.

BUMPERS

Part#	Description
ELP0000F	Front Bumper (one double durometer bumper)
ELP0000K	Rear Bumper Kit (Green, Black, Red and Orange bumpers, one of each colour)
ELP0000G	Rear Bumper Green, Soft
ELP0000B	Rear Bumper Black, Medium
ELP0000R	Rear Bumper Red, Firm
ELP0000A	Rear Bumper Orange, Extra Firm



HARDWARE KITS

Part#	Foot Size
ELX50001	22-25cm
ELX50002	26-28cm


Includes

- 2 x Bolts
- 2 x Washer
- Loctite



Elation®

Diagram

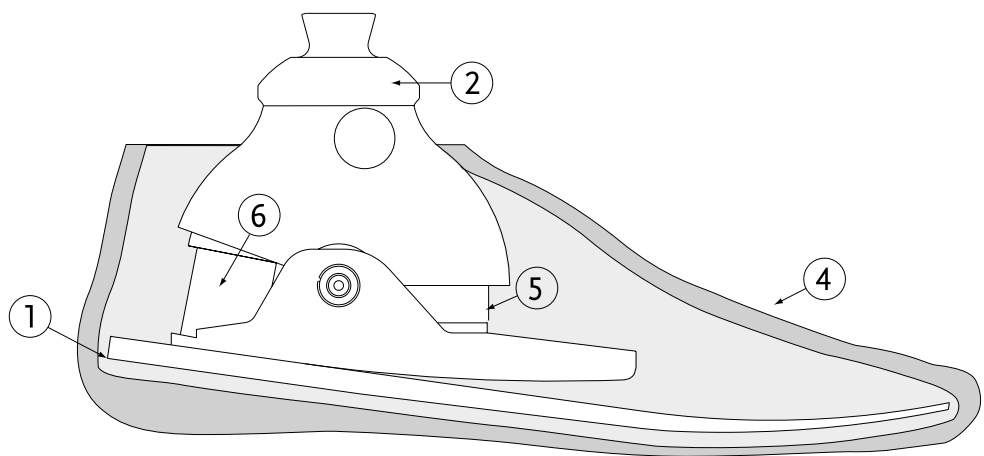

Picture	Description
1	Foot module
2	Elation with pyramid
3	Elation with tube clamp
4	Foot cosmesis
5	Front bumper
6	Rear bumper

Flex-Symes™

Introduction



- Specifically designed for Symes users
- CarbonX® Active Heel
- Lightweight and very low build height
- Circular alignment wedges allow 6° angulation


USER INFORMATION

Maximum Patient Weight	166kg
Impact Level	All


FOOT INFORMATION

Categories	1-6	1-9
Sizes	22-24cm	25-30cm
Weight of the the foot (Size 25) with socket adapter	350g	


Features

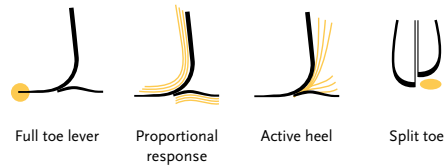

Flex-Symes with Socket Adapter

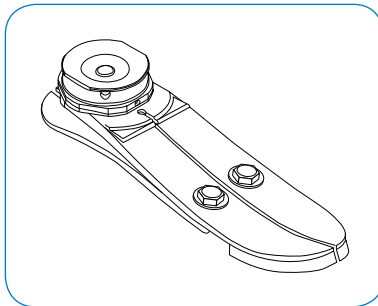


Flex-Symes kit with socket adapter

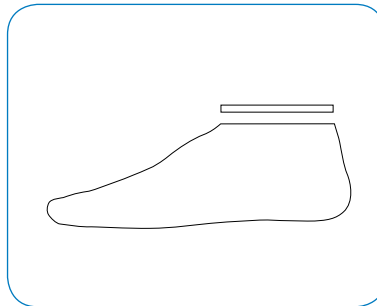


Foot cover kit

Flex-Symes without adapter

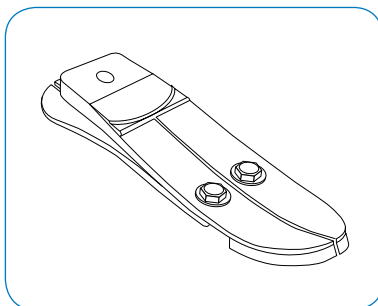


Flex-Symes kit without adapter

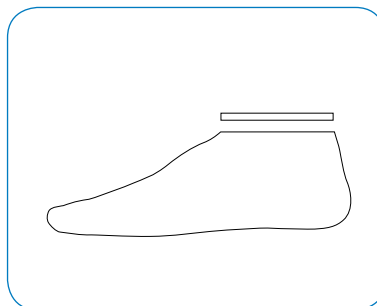

Foot cover kit

Impact Level	Weight Kg									
	44-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147	148-166
Low	1	1	2	3	4	5	6	7	8	9
Moderate	1	2	3	4	5	6	7	8	9	Special order
High	2	3	4	5	6	7	8	9	Special order	Special order

Dimensions

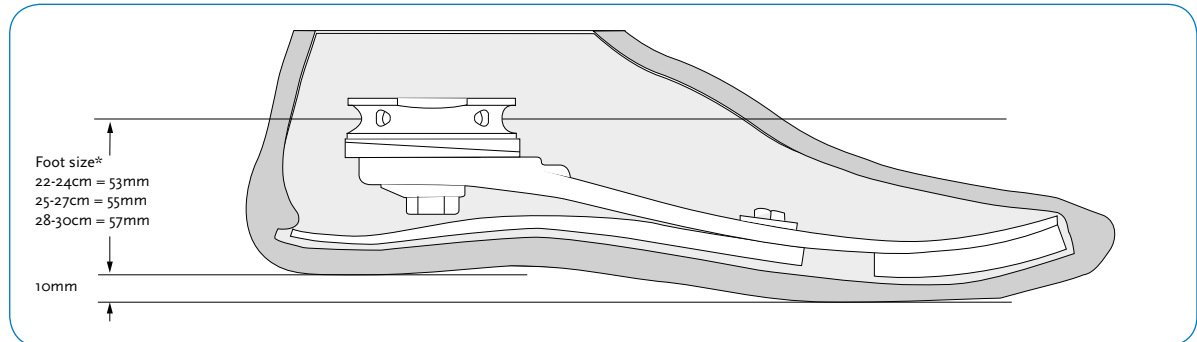

All measurements using 10mm heel height.

Kits


Part#

FLEX-SYMES KIT WITH SOCKET ADAPTER

Category 1-5 (Size 22-24), 1-9 (Size 25-30)
Foot Size (22-30cm)
SLPo 10mm Heel height

Includes

- Carbon foot module with heel
- 2 Circular alignment wedges
- 2 Alignment pins
- 1 Anti rotation bolt
- Washer
- Loctite 410
- Potting Epoxy
- Socket adapter

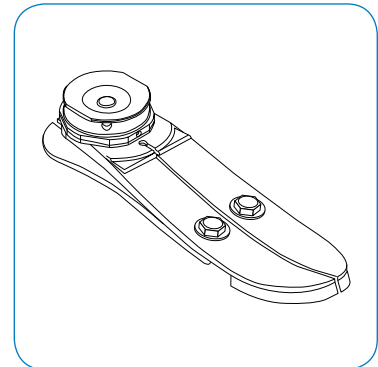


FLEX-SYMES KIT WITHOUT ADAPTER

Category 1-5 (Size 22-24), 1-9 (Size 25-30)
Foot Size (22-30cm)
SLXo 10mm Heel height

Includes

- Carbon foot module with heel
- 2 Alignment pins
- 1 Anti rotation bolt
- Washer
- Loctite 410
- Potting Epoxy

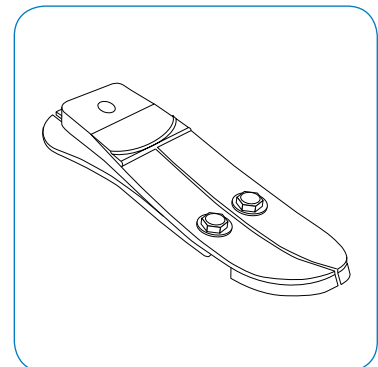

Note: Solid toe is a special order. Please call your Customer Care Representative on 08450 065 065.

Flex-Symes™


Components

ALIGNMENT DISC KIT

Part#	Foot Size	Cat.
SLP30001	22-24cm	1-6
SLP30002	25-27cm	1-6
SLP30003	28-30cm	1-6
SLP30004	25-27cm	7-9
SLP30005	28-30cm	7-9

Includes

- 2 Circular alignment wedges
- 2 Alignment pins
- 1 Anti rotation bolt
- Washer
- Loctite 410

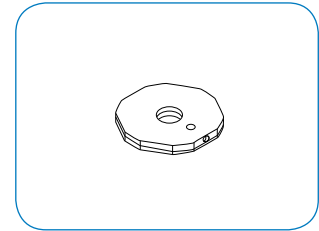

LAMINATION KIT

Part#	Foot Size	Cat.
SLP30011	22-24cm	1-6
SLP30012	25-27cm	1-6
SLP30013	28-30cm	1-6
SLP30014	25-27cm	7-9
SLP30015	28-30cm	7-9

Includes

- Socket adapter with insert
- Potting Epoxy



SOCKET ADAPTER WITH INSERT

Part#	Cat.	Foot Size
SLX50001	(1-6)	22-24cm
SLX50002	(1-6)	25-27cm
SLX50003	(1-6)	28-30cm
SLX50004	(7-9)	25-27cm
SLX50005	(7-9)	28-30cm



Flex-Symes™

Components

HEEL FASTENER HARDWARE KIT

Part#	Length	Size	Cat.
CM150004	18mm	1/4"	1-6
CM150005	18mm	5/16"	7-9

Includes

- Bolt
- Washer



CENTRAL BOLTS

Part#	Foot Size	Bolt Size	Cat.
SLX30021	22-27cm	M12	1-6
SLX30022	25-27cm	M14	7-9
SLX30022	28-30cm	M14	1-6
SLX30023	28-30cm	M16	7-9



ANTI ROTATION BOLT KIT

Part#	Foot Size	Bolt Size
SLX30005	22-27cm	M6
SLX30007	22-30cm	M8

Includes

- Bolt
- Washer



HEEL WEDGE KIT


Part#	Foot Size	Cat.
CM150011	22-30cm	1-9

Includes

- 3 Wedges
- Adhesive




Picture	Part#	Description
1	SLX30006	Alignment Pin
2	CM140003	Heel Divider
3	CM160009	Potting Epoxy

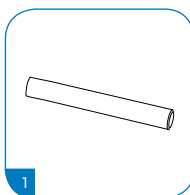

Flex-Symes™

Cosmetic Finishing

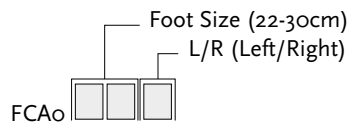

Please specify when ordering brown covers. Add "BR" to the part#.

Includes


- Foot Cover
- Attachment Plate
- Flex-Foot Sock

Part#	Description	Size
FCX63006	Flex-Foot Sock	22-25cm
FCX63007	Flex-Foot Sock	26-30cm

Part#	Description
CM100068	Shoe Horn

See page 146 for details on cosmetic finishing options.



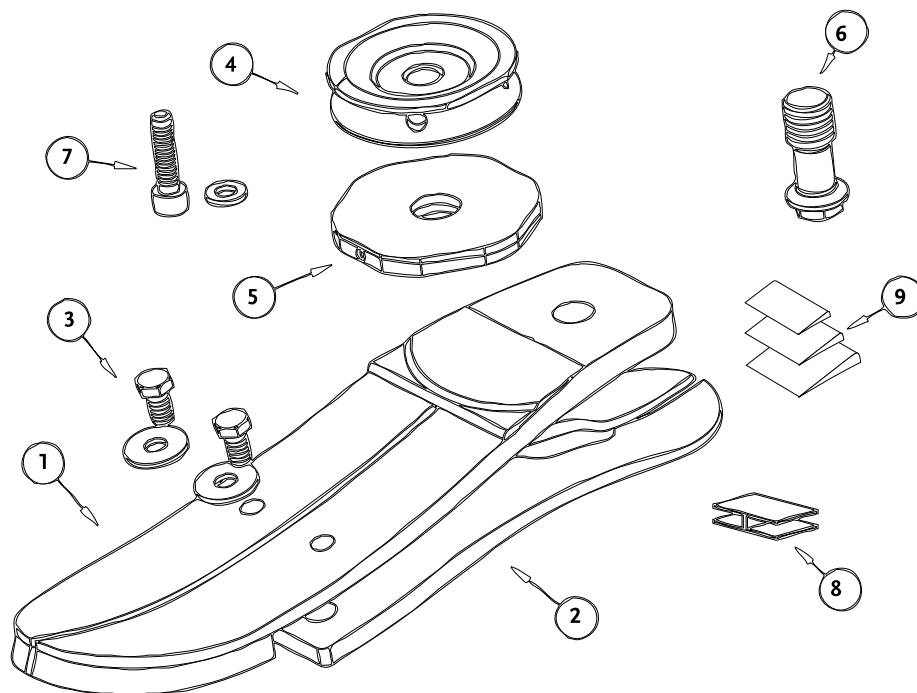

Picture	Description
1	Foot module
2	Heel
3	Heel Fastener Hardware Kit
4	Socket adapter
5	Alignment disk kit
6	Central bolt
7	Anti rotation bolt kit
8	Heel divider
9	Heel wedge kit

Flex-Foot® Junior

Introduction



Flex-Foot Junior is a low profile and energy storing foot module with integrated male pyramid. The shape of the male pyramid contributes to progressive stiffening of the foot and the layering of the carbon fibre is specifically designed to offer the support and flexibility needed for the active lifestyle of children.

- Foot module deflection proportional to childrens gait patterns
- Full length keel
- Sandal toe
- 10mm heel height

USER INFORMATION

Maximum Patient Weight	45kg
------------------------	------


FOOT INFORMATION

Categories	1-3
------------	-----


Sizes	16-21cm
-------	---------

Weight of the the foot (Size 21) with male pyramid	180g
---	------


Features

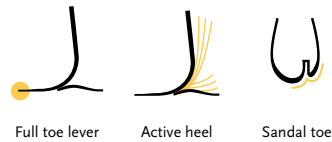

Flex-Foot Junior

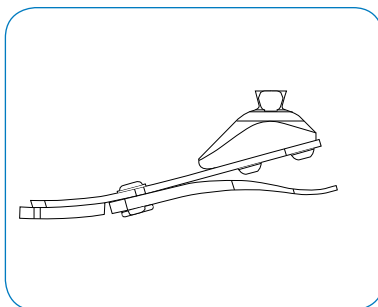

Flex-Foot Junior kit

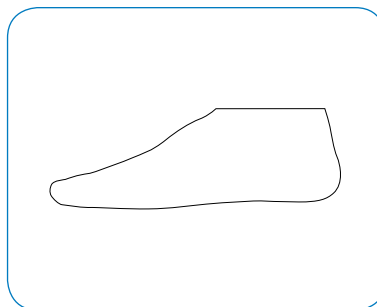

Foot cover kit


Flex-Foot® Junior

Category Selection Guide

Weight Kg	Foot Size					
	16	17	18	19	20	21
15-24kg	1	1	1	1		
25-34kg			2	2	2	2
35-45kg					3	3

The foot is recommended for daily activities involving rigorous and repetitive actions that are typical for a child and at moderate to high levels of impact on the prosthesis. For example: Jumping, climbing and running.

Dimensions

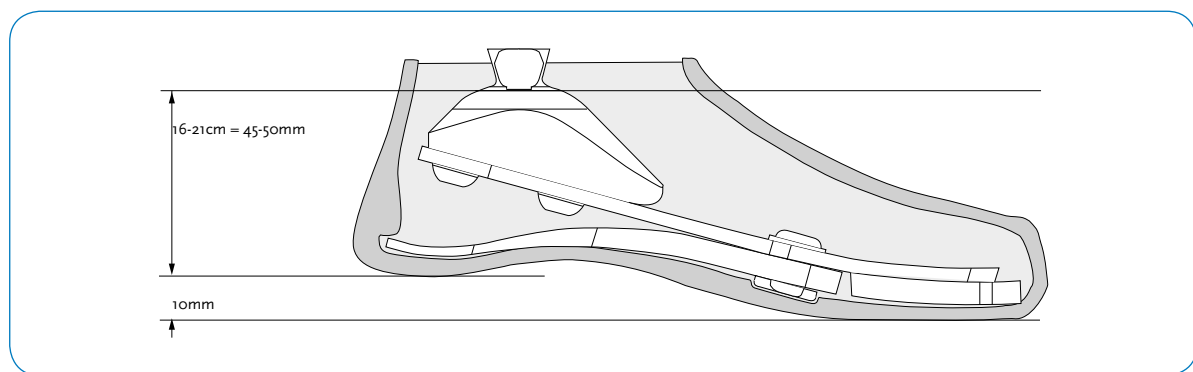


All measurements using 10mm heel height.

Kits

Part#

FLEX-FOOT JUNIOR KIT WITH MALE PYRAMID

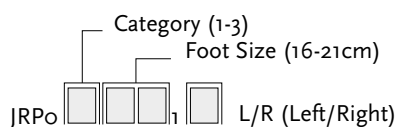

Includes

- Toe Cushion
- Foot Module
- Heel Module
- Heel Fasteners
- Heel Wedges
- Male Pyramid

Components

JUNIOR HEEL WEDGE KIT

Part#	Foot Size	Cat.
CM150017	16-21cm	1-3


Includes

- 3 Wedges
- Adhesive



Flex-Foot® Junior

Cosmetic Finishing

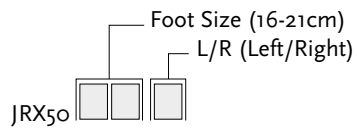

Please specify when ordering brown covers. Add "BR" to the part#.
Brown Junior foot covers are special orders.

Includes

- Foot Cover
- Flex-Foot Sock

Part#	Description	Size
FCX63008	Flex-Foot Sock	16-21cm



Chopart™

Introduction



- 100% Carbon Fibre, biaxial weave keel
- Fits into standard foot covers
- Suitable for Chopart, Lisfranc or metatarsal amputations

USER INFORMATION

Maximum Patient Weight	147kg
Impact Level	All

FOOT INFORMATION

Categories	1-8
Sizes	22-30cm

Category Selection Guide

Impact Level	Weight Kg								
	44-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147
Low	1	1	2	3	4	5	6	7	8
Moderate	1	2	3	4	5	6	7	8	Special order
High	2	3	4	5	6	7	8	Special order	Special order

Chopart foot plates are a standard order in the following sizes and categories:

Sizes 22-24 in categories 1-5


Sizes 25-26 in categories 2-6

Sizes 27-28 in categories 3-7

Sizes 29-30 in categories 4-8

For information on special or custom orders, please contact Customer Care on 08450 065 065.


Dimensions

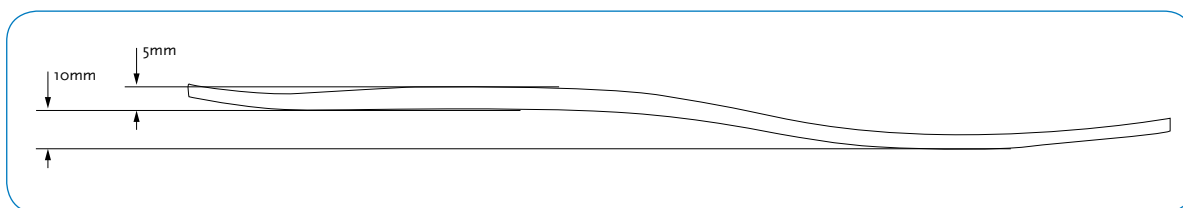

Kits

Part#

CHOPART KIT

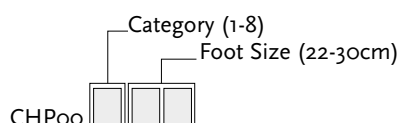


Includes

- Graphite Foot Plate
- Cloth Kit
- Potting Epoxy



Chopart™

Cosmetic Finishing

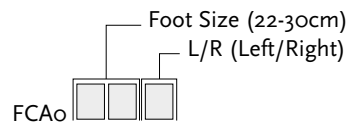

Please specify when ordering brown covers. Add "BR" to the part#.

Includes

- Foot Cover
- Flex-Foot Sock

Part#	Description	Size
FCX63006	Flex-Foot Sock	22-25cm
FCX63007	Flex-Foot Sock	26-30cm

Part#	Description
CM100068	Shoe Horn

See page 146 for details on cosmetic finishing options.



Flex-Run™ with Nike Sole

Introduction



- The Flex-Run with Nike Sole offers unmatched performance and efficiency for amputees and prosthetists. Össur and Nike have partnered to design the ideal running foot from the ground up, leveraging Össur's unmatched expertise in running prosthetics with Nike's world renowned traction and sole technologies.
- The improved Flex-Run features a longer toe lever and efficient energy return. The Nike Sole offers great traction to help improve performance for the runner while improving efficiency and results for the prosthetist.
- Flex-Run with Nike Sole is appropriate for Transfemoral and Transtibial users participating in high impact activities such as recreational jogging, trail running, distance running and triathlons.
- Custom foot module for long distance running and recreational sports
- Lightweight, vertical compliance, efficient energy return
- Bolt-on male or female pyramid


USER INFORMATION

Maximum Patient Weight	130kg
Impact Level	High to Extreme


FOOT INFORMATION

Categories	1-9
Weight of the the foot (Size 21) with sole and male pyramid	694g


Features

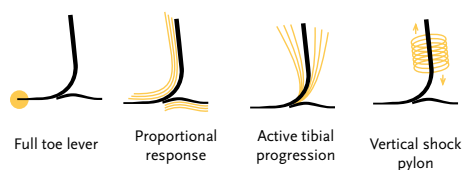

Flex-Run

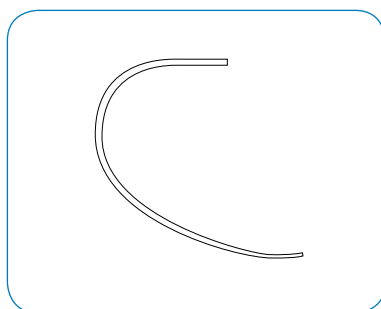

Flex-Run kit

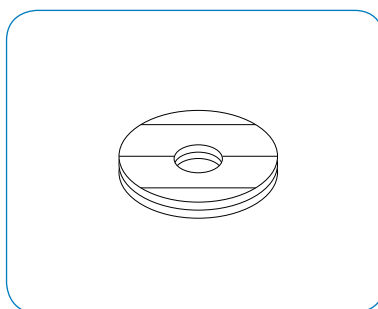

Attachment plate


Flex-Run™ with Nike Sole

Category Selection Guide

Weight Kg

	37-44	45-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130
Category	1	2	3	4	5	6	7	8	9

Dimensions

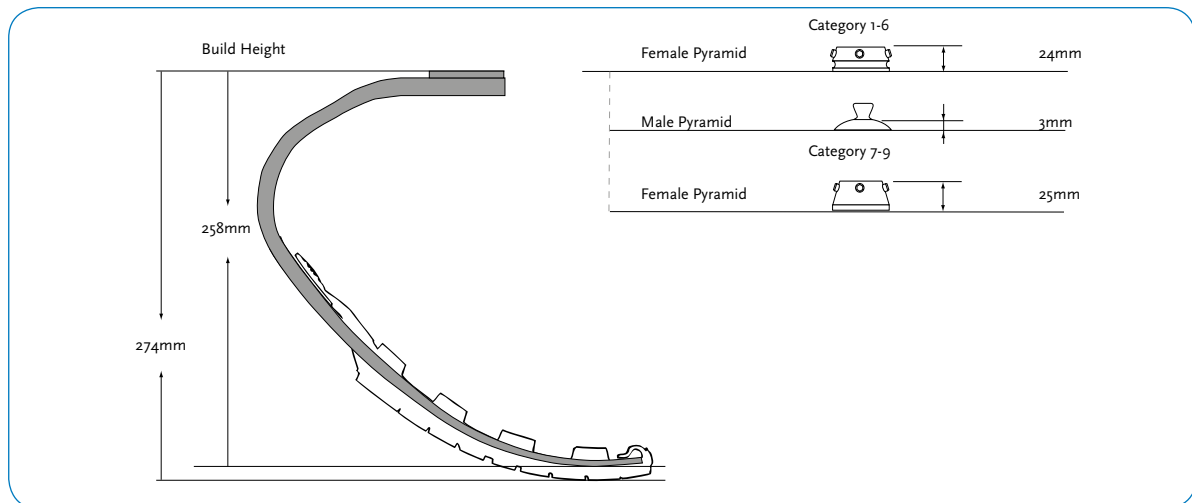

Kits

Part#

FLEX-RUN KIT

Category (1-9)
 FSX008
 HI / LO (High/Low)

Flex-Run is a special order and is excluded from the trial period. Trial feet are available on request, please call your Customer Care Representative on 08450 065 065.

Includes:



- Foot Module
- Attachment Plate (can be supplied bonded on the foot upon request)
- Nike Sole
- Bag for Nike Sole

NIKE SOLE KIT

Part#	Description	Cat.
CM150163	Nike Sole Kit	1-9

Includes:

- Nike Sole
- Medallion
- Loctite
- Bag for Nike Sole

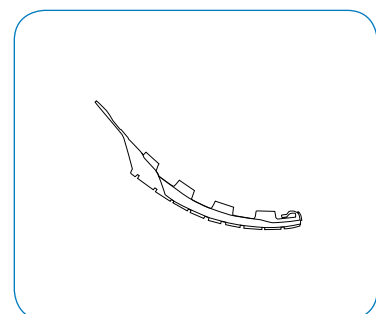
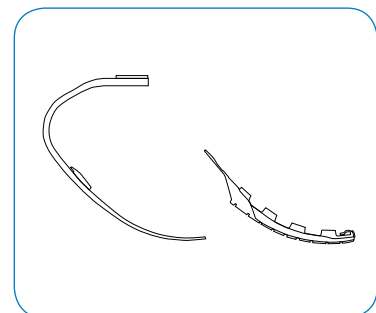


Flex-Run™ with Nike Sole

Kits


FOOT ADAPTER SETS

Picture	Part#	Description	Weight	Bolt Size	Cat.
1	FWX00030	Male Pyramid Foot Adapter Set	75g	M10 x 30mm	1-6
2	CM150161	Female Pyramid Foot Adapter Set	100g	M10 x 35mm	1-6
	CM150162	Female Pyramid Foot Adapter Set	186g	M12 x 50mm	7-9

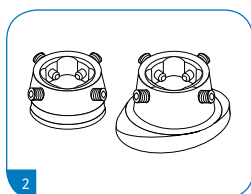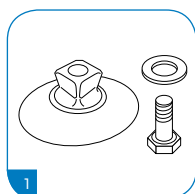


REPLACEMENT BOLT

Part#	Size	Cat.
CM130221	M10	1-6
CM130225	M12	7-9

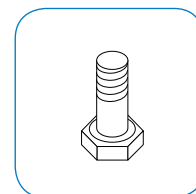

Components

Part#	Description	Cat.
FWX40001	Attachment Plate	2-6
FWX40002	Attachment Plate	7-9

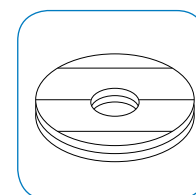

Flex-Run™ Junior

Introduction



- Lower build height than standard Flex-Run
- Custom foot module for long distance running and recreational sports
- Lightweight, vertical compliance, efficient energy return
- Bolt-on male or female pyramid

USER INFORMATION

Maximum Patient Weight	84kg
Impact Level	High to Extreme

FOOT INFORMATION

Categories	1-6
Weight of the the foot (Size 21) with male pyramid	336g

Features




Proportional
response




Active tibial
progression

Flex-Run Junior

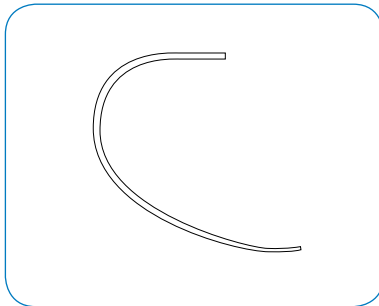

Flex-Run kit

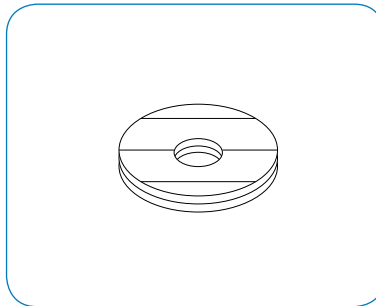


Attachment plate

Flex-Run™ Junior

Category Selection Guide

	Weight Kg					
	15-24	25-34	35-45	46-56	57-67	68-84
Category	1	2	3	4	5	6

Dimensions

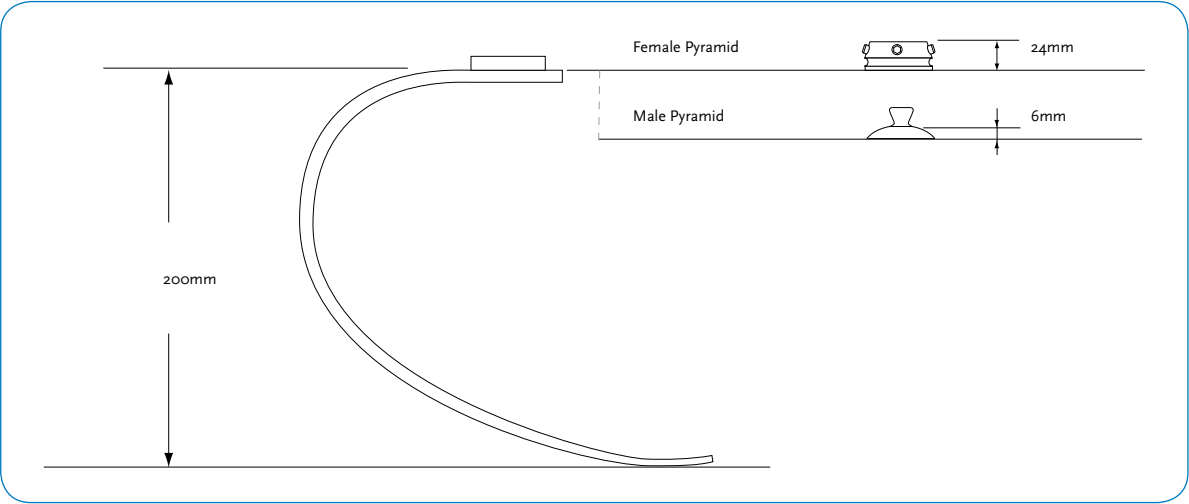


Kits

Part#

FLEX-RUN KIT


FSX007  Category (1-6)

Flex-Run is a special order and is excluded from the trial period.

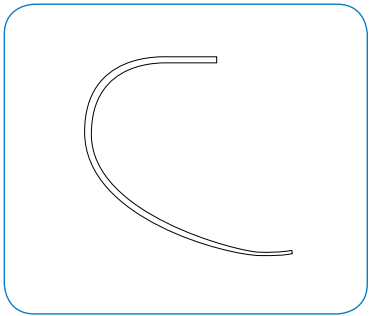



Components

Part#	Description	Cat.
FWX40001	Attachment Plate	1-6

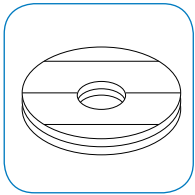

FOOT ADAPTER SETS

Picture	Part#	Description	Weight	Bolt Size	Cat.
1	FWX00030	Male Pyramid Foot Adapter Set	75g	M10 x 30mm	1-6
2	CM150161	Female Pyramid Foot Adapter Set	100g	M10 x 35mm	1-6

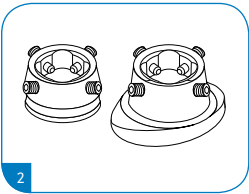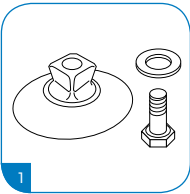


Flex-Run™ Junior

Components

REPLACEMENT BOLT

Part#	Size	Cat.
CM130221	M10	1-6

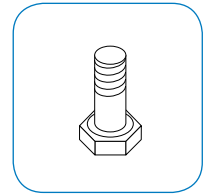

Flex-Sprint™

Introduction



- Custom foot module for track sprinting
- Ideal for transfemoral amputees
- One adult foot size
- L-bracket connector with 4-hole pattern for easy knee connection

USER INFORMATION

Maximum Patient Weight	147kg
Impact Level	High


FOOT INFORMATION

Categories	1-9
------------	-----

Category Selection Guide

Weight Kg									
	44-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147
Category	1	2	3	4	5	6	7	8	9

Dimensions

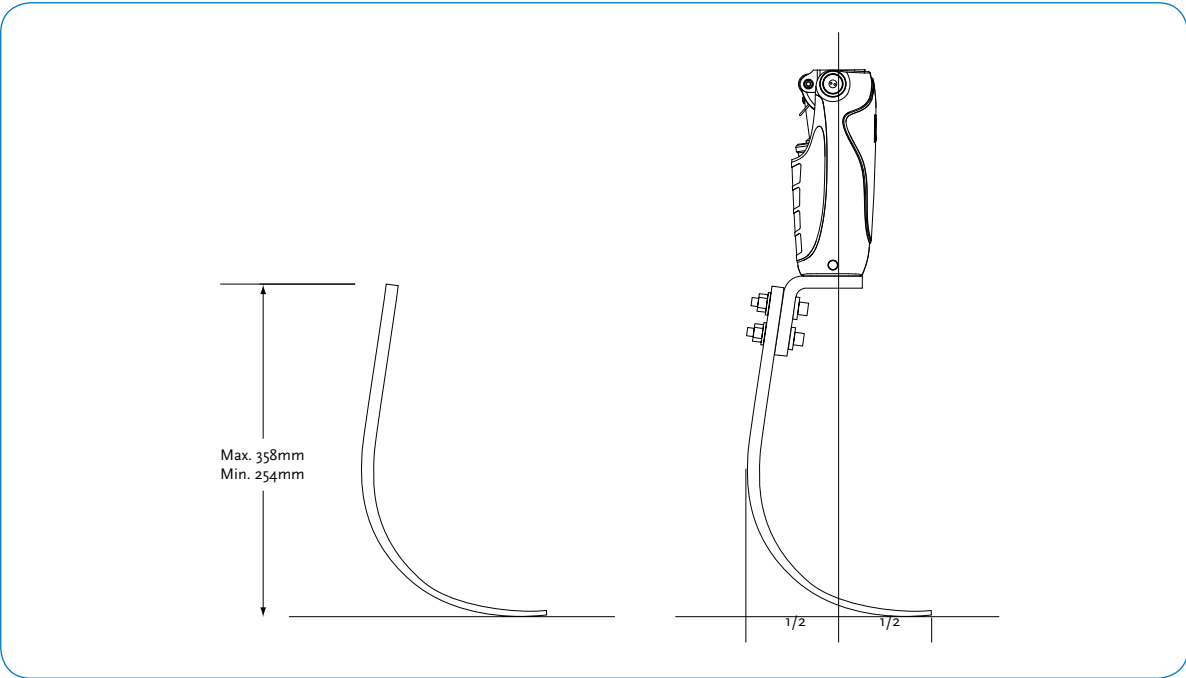

Flex-Sprint™


Kits

Part#

FLEX-SPRINT KIT

FSX002  Category (1-9)

Flex-Sprint is a special order and is excluded from the 60-day trial period.

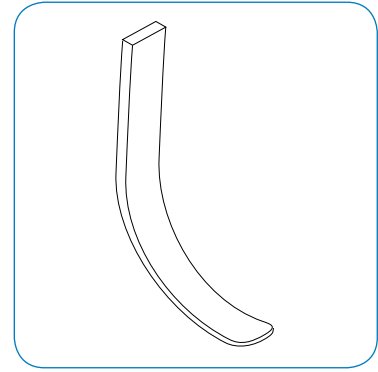



Components

L-BRACKET WITHOUT INSERTS

Picture	Part#	Description	Cat.
1	FSX50005	5° L-Bracket without Inserts	1-6
2	FSX50006	20° L-Bracket without Inserts	1-6

Includes

- L-Bracket without inserts
- Pyramid Attachment Hardware
- Alignment Wedges
- Alignment Hardware

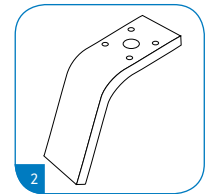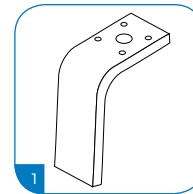


HARDWARE KIT FOR T CONNECTORS

Part#	Weight	Size	Cat.
CM100054	108g	5/16" x 2-1/2"	1-6
CM100055	144g	3/8" x 2-1/2"	7-9

Includes

- Bolts
- Nuts
- Spherical Washer
- Flat Washer

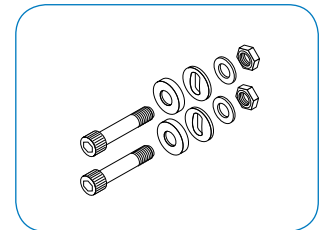

Cheetah®

Introduction




- Custom foot module for track sprinting
- Ideal for transtibial amputees
- Attaches to posterior of the socket
- Lamination connector for clearance less than 250mm
- Pylon connector for clearance greater than 250mm
- One adult foot size

USER INFORMATION	
Maximum Patient Weight	147kg
Impact Level	High
FOOT INFORMATION	
Categories	2-9

Category Selection Guide

Weight Kg								
	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147
Category	2	3	4	5	6	7	8	9

Dimensions

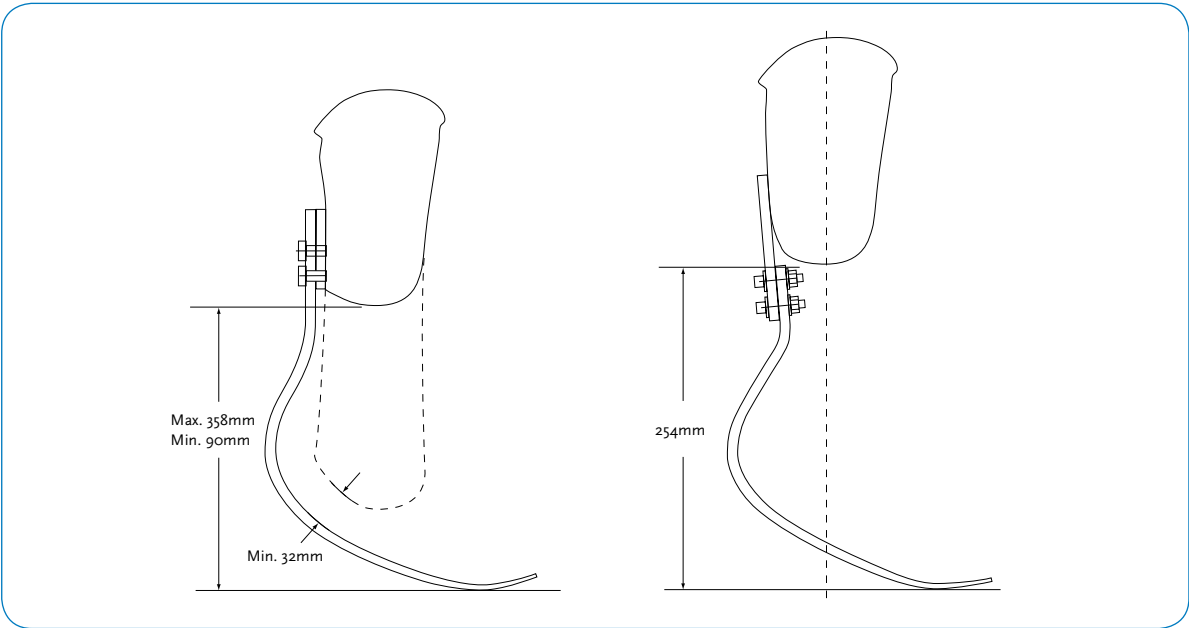

Cheetah®


Kits

Part#

CHEETAH KIT

FSX004  Category (2-9)

Cheetah is a special order and is excluded from the 60-day trial period.

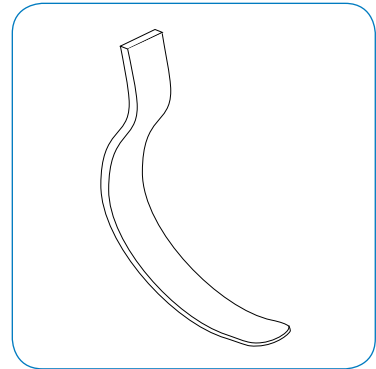

Components

CONNECTOR KIT

Picture	Part#	Description	Cat.
1	FSX50001	Lamination Connector Kit	2-6
1	FSX50002	Lamination Connector Kit	7-9
2	FSX50003	Pylon Connector Kit	2-6
2	FSX50004	Pylon Connector Kit	7-9





Includes

- Connector
- Alignment Wedges
- Attachment Hardware (4)
- Lamination Kit

ATTACHMENT BOLT

Picture	Part#	Description	Size	Cat.
1	CM130218	Attachment Bolt for Lamination Connector	5/16"-24 x 1-1/2"	2-6
1	CM130213	Attachment Bolt for Lamination Connector	3/8"-24 x 1-3/4"	7-9
2	CM130305	Attachment Bolt for Pylon Connector	5/16"-24 x 1-1/2"	2-6
2	CM130304	Attachment Bolt for Pylon Connector	3/8"-24 x 1-3/4"	7-9

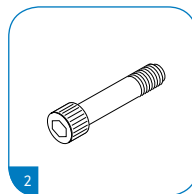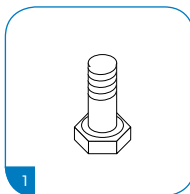


HARDWARE KIT FOR T CONNECTORS

Part#	Weight	Size	Cat.
CM100054	108g	5/16" x 2-1/2"	2-6
CM100055	144g	3/8" x 2-1/2"	7-9

Includes

- Bolts
- Nuts
- Spherical Washer
- Flat Washer

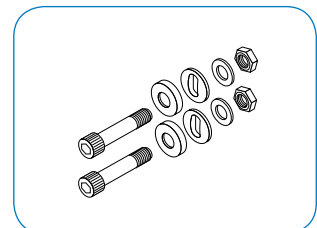

ALIGNMENT WEDGES

Part#	Description	Cat.
CM140502	Alignment Wedges 2°	2-9
CM140505	Alignment Wedges 4°	2-9
CM140508	Alignment Wedges 6°	2-9
CM140511	Alignment Wedges 8°	2-9



Cheetah Xtreme® with Nike Spike Pad

Introduction



The blade of champions – inspired by the world's fastest. Cheetah Xtreme is the optimal short-distance sprinting foot, a holder of multiple world records.

Designed specifically for the fast, short-distance sprints, the carbon blade features a more extreme curve than the other Cheetahs. This dynamic shape allows the foot to flex more offering a powerful energy kick. A longer, flatter toe enhances push off, while the plantar-flexed pylon supports better forward progression.

The Nike Spike Pad has been designed exclusively for the Cheetah. It provides specialized traction that has been optimized for the blade.

A new alignment adapter for the Cheetah is now available for quicker set up and easier alignment of the blade.

USER INFORMATION

Maximum Patient Weight	147kg
Impact Level	Extreme

FOOT INFORMATION

Categories	1-9
Weight of the the foot	918g

Features




Proportional response



Active tibial progression

For complete Cheetah Xtreme kit

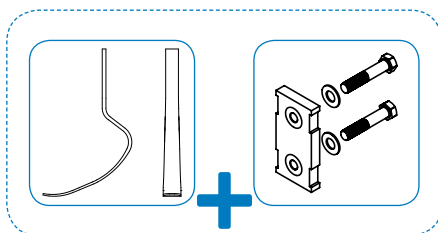


Cheetah Xtreme® with Nike Spike Pad

Category Selection Guide

	Weight Kg								
	42-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147
Category	1	2	3	4	5	6	7	8	9

Dimensions

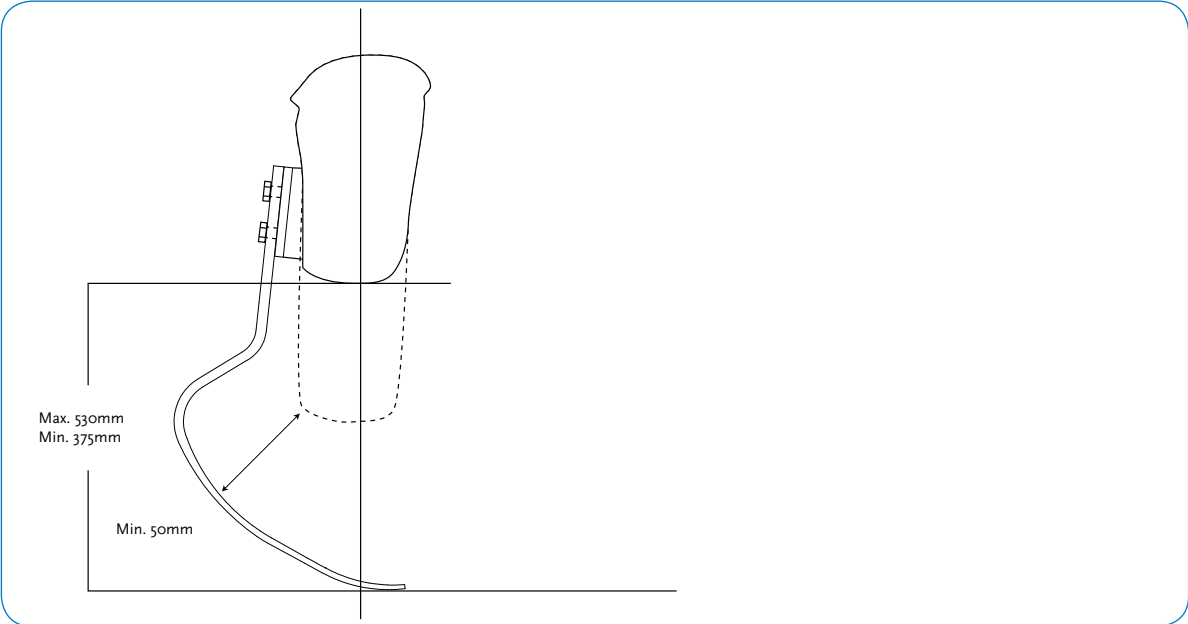

Kits

Part#


CHEETAH XTREME KIT

FSX004  Category (1-9)

Includes

- Cheetah Xtreme Blade
- Nike Spike Pad Kit

Cheetah Xtreme is a special order and is excluded from the 60-day trial period.

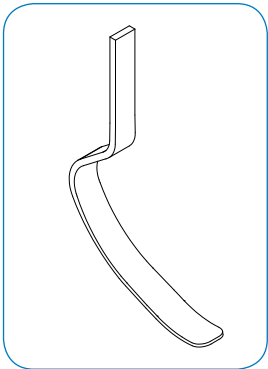

NIKE SPIKE PAD KIT

Part#	Description	Cat.
FSX50008	Nike Spike Pad	1-9

Includes

- Nike Spike Pad
- 7mm Pyramid Spikes
- Wrench



Note: Nike Spike Pad has same spikes as Nike track and field shoes, additional spikes can be purchased at Nike.

Cheetah Xtreme® with Nike Spike Pad


Components

LAMINATION CONNECTOR KIT

Part#	Description	Cat.
FSX50009	Lamination Connector Kit	1-9

Includes

- Connector
- Alignment Wedges (2°, 4°, 6°, 8°)
- Attachment Hardware
- Lamination Kit

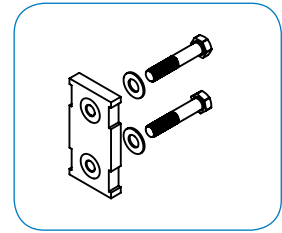


HARDWARE KIT FOR T CONNECTORS

Part#	Weight	Size	Cat.
CM100055	144g	x 2-1/2"	1-9


Includes

- Bolts
- Nuts
- Spherical Washer
- Flat Washer

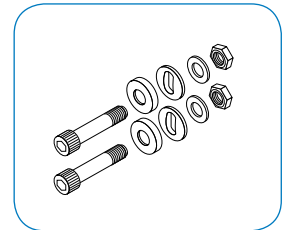


ALIGNMENT ADAPTER KIT

Part#	Description	Cat.
FSX50007	Alignment Adapter	All

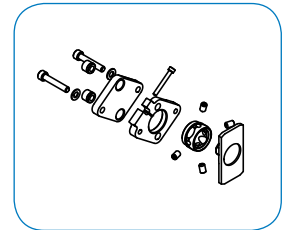


ATTACHMENT BOLT

Part#	Description	Size	Cat.
CM130211	Attachment Bolt for Lamination Connector	3/8"-24 x 1 1/2"	1-9
CM130213	Attachment Bolt for Lamination Connector	3/8"-24 x 1 3/4"	1-9
CM130214	Attachment Bolt for Lamination Connector	3/8"-24 x 1 3/8"	1-9

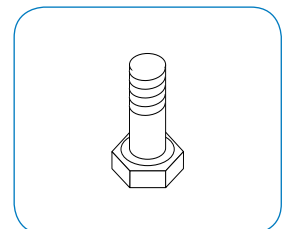

ALIGNMENT WEDGES

Part#	Description	Cat.
CM140515	Alignment Wedges 2°	All
CM140516	Alignment Wedges 4°	All
CM140517	Alignment Wedges 6°	All
CM140518	Alignment Wedges 8°	All

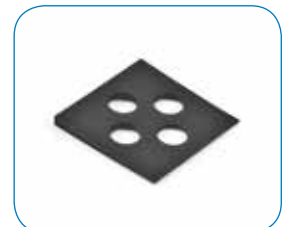

Cheetah Xtend® with Nike Spike Pad

Introduction



Kick into high gear with a responsive and comfortable ride from start to finish.

The Cheetah Xtend is the perfect companion for longer sprints and short-distance running (400-5000m). This highly effective carbon fiber blade has a continuous curve toward the toe offering a smooth rollover. The longer, flatter toe enhances push off, while the plantar flexed pylon supports better forward progression.

The Nike Spike Pad has been designed exclusively for the Cheetah. It provides specialized traction that has been optimized for the blade.

A new alignment adapter for the Cheetah is now available for quicker set up and easier alignment off the blade.



USER INFORMATION

Maximum Patient Weight	147kg
Impact Level	Extreme

FOOT INFORMATION

Categories	1-9
Weight of the the foot	750g

Features




Proportional response



Active tibial progression

For complete Cheetah Xtreme kit

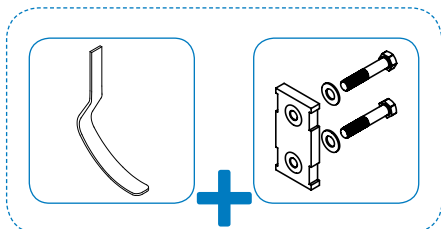


Cheetah Xtend® with Nike Spike Pad

Category Selection Guide

Weight Kg

	42-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147
Category	1	2	3	4	5	6	7	8	9

Dimensions

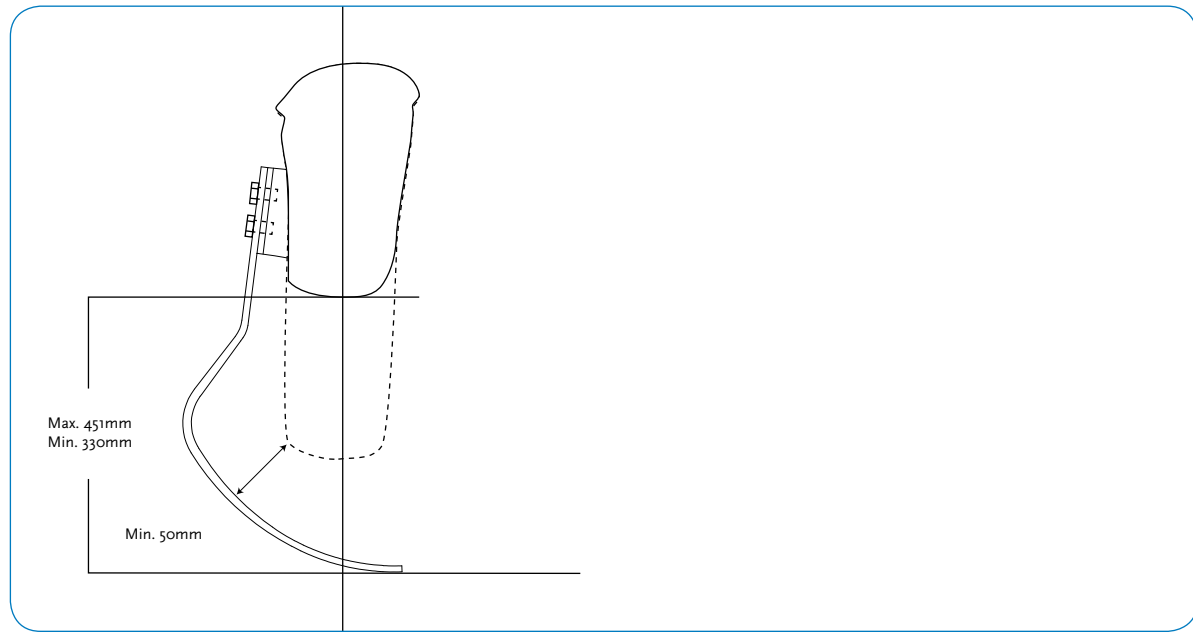

Kits

Part#

CHEETAH XTEND KIT

FSX006  Category (1-9)

Includes

- Cheetah Xtend Blade
- Nike Spike Pad Kit

Cheetah Xtend is a special order and is excluded from the 60-day trial period.


NIKE SPIKE PAD KIT

Part#	Description	Cat.
FSX50008	Nike Spike Pad	1-9

Includes

- Nike Spike Pad
- 7mm Pyramid Spikes
- Wrench

Note: Nike Spike Pad has same spikes as Nike track and field shoes, additional spikes can be purchased at Nike.

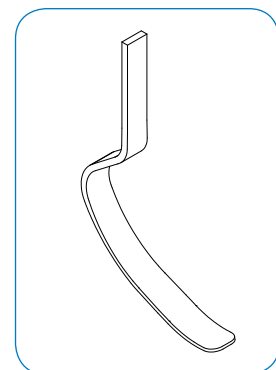

Cheetah Xtend® with Nike Spike Pad


Components

LAMINATION CONNECTOR KIT

Part#	Description	Cat.
FSX50009	Lamination Connector Kit	1-9

Includes

- Connector
- Alignment Wedges (2°, 4°, 6°, 8°)
- Attachment Hardware
- Lamination Kit

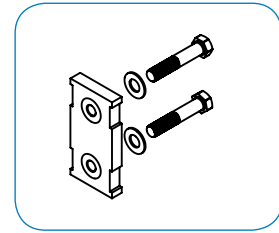


HARDWARE KIT FOR T CONNECTORS

Part#	Weight	Size	Cat.
CM100055	144g	x 2-1/2"	1-9


Includes

- Bolts
- Nuts
- Spherical Washer
- Flat Washer

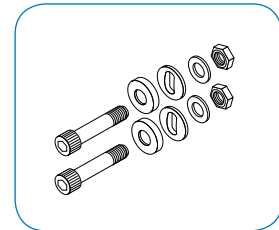


ALIGNMENT ADAPTER KIT

Part#	Description	Cat.
FSX50007	Alignment Adapter	All

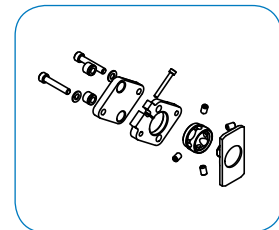


ATTACHMENT BOLT

Part#	Description	Size	Cat.
CM130211	Attachment Bolt for Lamination Connector	3/8"-24 x 1 1/2"	1-9
CM130213	Attachment Bolt for Lamination Connector	3/8"-24 x 1 3/4"	1-9
CM130214	Attachment Bolt for Lamination Connector	3/8"-24 x 1 3/8"	1-9

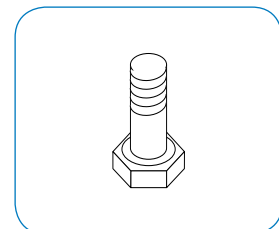

ALIGNMENT WEDGES

Part#	Description	Cat.
CM140515	Alignment Wedges 2°	All
CM140516	Alignment Wedges 4°	All
CM140517	Alignment Wedges 6°	All
CM140518	Alignment Wedges 8°	All

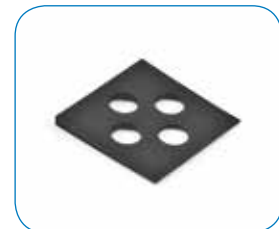

Cheetah® Xplore Junior

Introduction



A custom-built, hybrid carbon fiber foot designed to function as both an everyday and sports foot. For kids that require the flexibility of an everyday walking foot and a running foot in a single prosthesis, Cheetah Xplore Junior is the answer. Cheetah Xplore couples the proven performance of our Cheetah sprinting foot with a heel for everyday walking/standing.


USER INFORMATION

Maximum Patient Weight	55Kg
Impact Level	Low to High

FOOT INFORMATION

Categories	1-4
Sizes	One Size
Weight of the the foot	918g

Features

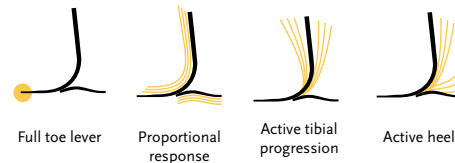


Category Selection Guide

Weight Kg

	15-24	25-34	35-45	46-55
Category	1	2	3	4

Cheetah® Xplore Junior

Dimensions

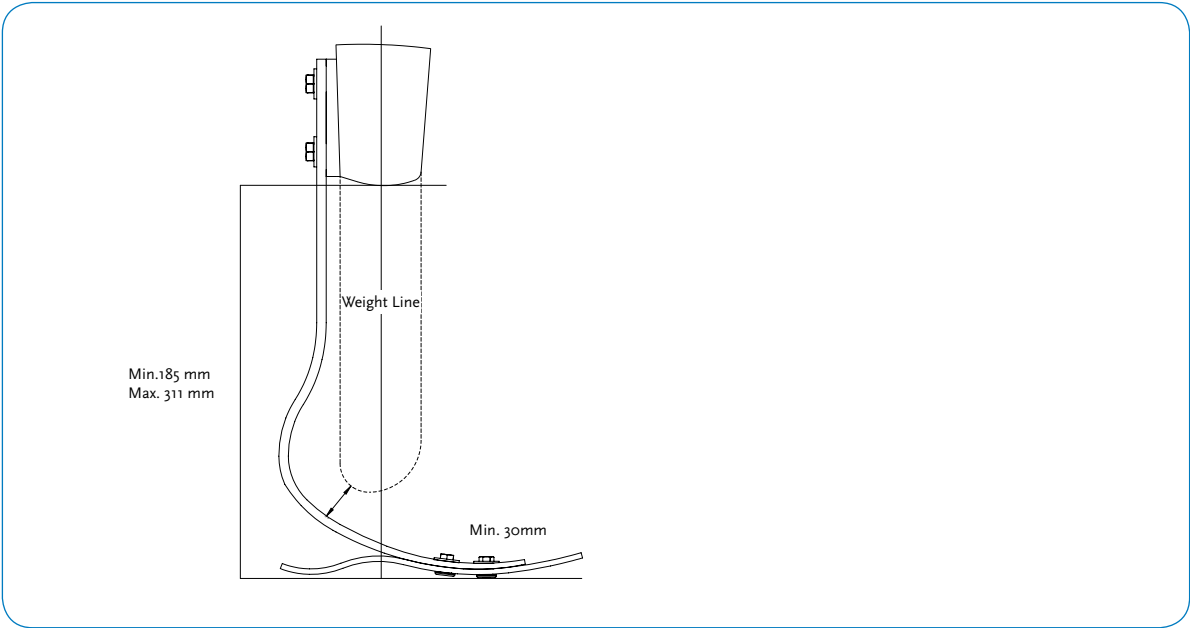


Kits

Part#

CHEETAH XPLORE JUNIOR KIT


CHXo  24 Category (1-4)

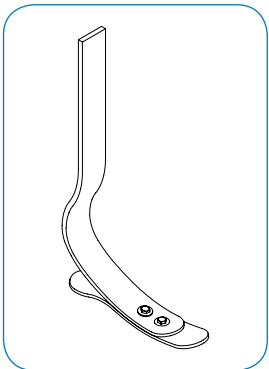

Components

LAMINATION CONNECTOR KIT


Part#	Description
CM140129	Lamination Connector

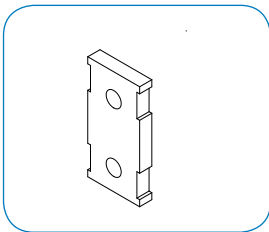


Cheetah® Xplore Junior

Cosmetic Finishing

SIZES 19-21

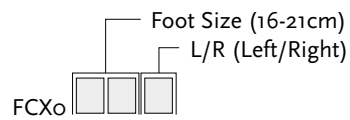

SIZES 22-24

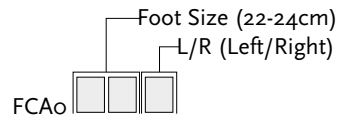

Please specify when ordering brown covers. Add "BR" to the part number.

Includes

- Foot Cover
- Flex-Foot Sock

Part#	Description	Size
FCE63008	Flex-Foot Sock	16-24cm

See page 146 for details on cosmetic finishing options.



Cheetah® Junior

Introduction



Custom-built, high-performance carbon fiber foot designed primarily for sporting activities. This is the optimal run and play foot for both transtibial and transfemoral users. It attaches posterior to the socket, making it agile, strong and a proven performer for those children who want to practice sports.

USER INFORMATION

Maximum Patient Weight	55Kg
Impact Level	High to Extreme

FOOT INFORMATION

Categories	1-4
Sizes	One Size
Weight of the the foot	232g

Features



Proportional response



Active tibial progression


Category Selection Guide

Weight Kg

	15-24	25-34	35-45	46-55
Category	1	2	3	4

Cheetah® Junior

Dimensions

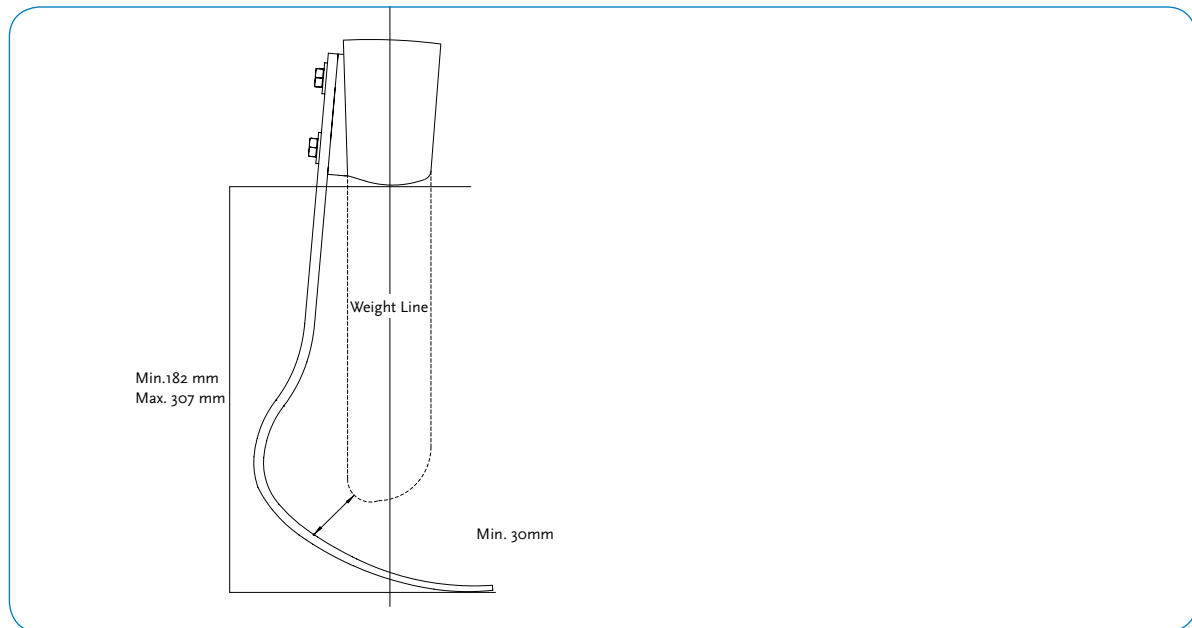


Kits

Part#

CHEETAH JUNIOR KIT


CJXo  Category (1-4)

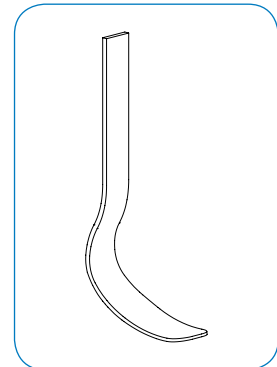

Components

LAMINATION CONNECTOR KIT

Part#	Description
CM140129	Lamination Connector

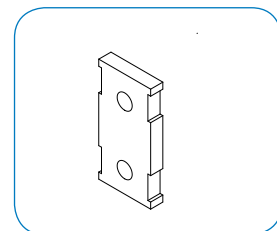



Rehabilitation

Introduction

Good prosthetic outcomes depend on integrated rehabilitation that starts before the first prosthesis is used. Rehabilitation tools like the Össur Rigid Dressing (ÖRD) and Femurett ensure that when prosthetic rehabilitation starts, the user is best prepared in terms of their residual limb and general health to ensure they get the best from their prosthesis.

The ÖRD immobilises, protects and controls post-operative oedema in the transtibial residual limb. Used as part of a controlled rehabilitation programme including compression treatment, the transition to prosthesis use can be seamless and effective.

The Post-Op liner can be used in conjunction with the ÖRD to promote early mobilisation of the patient and a speedy progression into a prosthesis. Using compression therapy, the Post-Op liner moulds the residual limb into the best possible shape for prosthetic socket fitting.

The Femurett is an early walking device for transfemoral amputees, used to assist therapists to assess and prepare potential prosthesis users for prosthetic rehabilitation. The adjustable pylons and sockets ensure that the therapy department use the Femurett for the majority of their clients.

Contents

Össur Rigid Dressing	247
Iceross Post-Op	248
Femurett	249

Össur Rigid Dressing (ÖRD)

Introduction

The Össur Rigid Dressing (ÖRD) immobilises, protects and controls post-operative oedema in the transtibial residual limb immediately following transtibial amputation surgery. Immobilisation and oedema control are achieved by providing a firm dressing formed to match the volume and shape of the residual limb. The immobilisation will assist in prevention of contractures by restricting knee flexion. Effectively constraining and matching the volume of the residual limb will prevent an increase in oedema following surgery. Containing the residual limb provides protection from discomfort, caused by bed trauma or possible falls during early rehabilitation.

The ÖRD is made out of soft and pliable material and thus conforms to the shape of the residual limb. The dressing is made rigid by creating a vacuum in the internal chambers of the product.

The ÖRD allows inspection of the residual limb during early rehabilitation.

Easy to adjust for volume matching to the residual limb during post-operative use. After 5 - 10 days we recommend compression treatment using the Icross® post-operative liner.

The lightweight ÖRD facilitates early mobilisation.



Kits

Part#	Description
K-239100	ÖRD Small
K-239101	ÖRD Large

Size Selection

Two measurements are required:

- Circumference 4cm from distal end of residual limb
- Circumference 50cm from distal end of residual limb

Size	Circumference at 4cm	Circumference at 50cm
Small	Up to 42cm	Up to 58cm
Large	42cm - 65cm	58cm - 70cm

Note: All circumference measurements are taken over dressings and padded material.

Components

Part#	Description
K-239105	Vacuum Pump for Rigid Dressing




Iceross® Post-Op Liner

Introduction

Iceross Post-Op liner is used for compression therapy. Best results are obtained with an Iceross Post-Op liner if a rigid dressing (ÖRD) has been applied to the residual limb immediately following amputation surgery and used for 5–7 days.

Early treatment with Iceross Post-Op liner has been shown to facilitate early mobilisation and speed up the maturation process of the residual limb. This makes it possible to provide the amputee with a prosthesis at an earlier stage leading to less dependence on medical staff.

Compression therapy with the Iceross Post-Op liner gives the residual limb an optimal shape for prosthetic socket fitting, thereby reducing the risk of prosthetic complications.

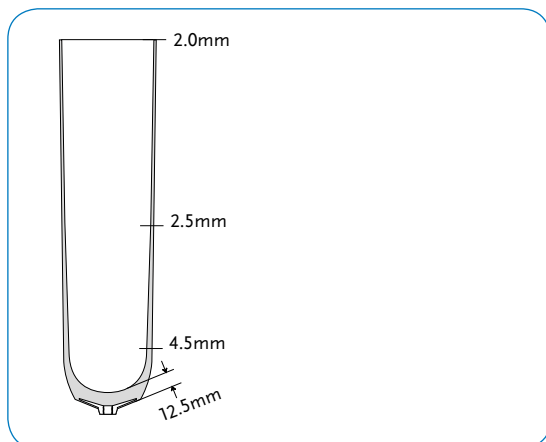

The main objectives of Iceross Post-Op are to:

- Diminish oedema
- Shape the residual limb
- Contribute to pain relief
- Accelerate the rehabilitation process

USER INFORMATION

Amputation Level	Transtibial
User	New amputee

LINER INFORMATION

Size	18, 20, 22, 23.5, 25, 26.5, 28, 30, 32, 34, 36, 38, 40, 42, 45
------	--

For sizing and selection advice please refer to pages 24 - 28.

Profile	2mm
---------	-----

Matrix	10cm
--------	------

Suspension



Locking

Features



Matrix

Ordering Information

ICEROSS POST-OP

Part#	Profile	Suspension method	Matrix
I-PO20XX	2mm	Locking	10cm

Femurett

Introduction



Femurett is a training and assessment device. Rehabilitation can proceed more smoothly when patient assessment and training is started and carried out using this prosthesis.

With simple adjustments, Femurett can be set to meet the exact dimensional requirements of the patient.

- Markings on the adjustable section simplify accurate setting
- Adjustable to varus and valgus alignment

Femurett, complete kit for transfemoral amputees consists of:

- Single articulated knee joint with extension spring
- Foot section (left and right with shin)
- 6 sockets, one of each size, left and right
- Weight limit 105kg

Components

Part#	Description
6674007	Femurett Locking Knee
6674021	Femurett Shin (pair)

SOCKET

Size	Left	Right	Thigh circumference (cm)
Small	6674306	6674313	38 - 46
Medium	6674320	6674337	44 - 54
Large	6674344	6674351	52 - 60

Includes shoulder strap, frame and connecting tube



General Information

Torque Values For Össur Products

	Torque	Loctite	Tool Required
ICEROSS®			
Iceross Ratchet/Smooth/Clutch Pin	4Nm	410/411	7mm Spanner 
M10-M6 Adapter for use with Standard Liners	6Nm	425/412	10mm Wrench 
ICELOCK®			
100 Series			
Mount Screws, pyramid	10Nm	425	4mm Hex 
Mount Screws, 4-hole	8Nm	242	4mm Hex 
200 Series			
Mount Screws (214 only)	12Nm	410/411	4mm Hex 
Clutch Release Mechanism	15Nm	None	17mm Socket 
600 Series			
Funnel	8Nm	None	Tool Provided in Kit
700 Series			
Funnel	8Nm	None	Tool Provided in Kit
ADAPTERS			
Angular set screws (pyramids, al & ti)	15Nm	242	4mm Hex 
Mount screws (4-hole, al & ti)	10Nm	242	4mm Hex 
Pinch bolts (tube clamps, al & ti)	10Nm x 2	242	4mm Hex 
TOTAL KNEE®			
TK 2100			
Pinch Bolt - 30mm Tube Adapter	15Nm	No	5mm 
TK 2000			
Pinch Bolt - 30mm Tube Clamp	15Nm	No	5mm 
TK 1900			
Pinch Bolt - 30mm Tube Clamp	15Nm	No	5mm 
Polymer Friction Set Screw	2Nm	No	2.5mm 
TK 1100			
Pinch Bolt - 22mm Tube Clamp	8Nm	No	4mm 
BALANCE KNEE®			
Balance Knee Locking			
Pinch Bolt - 30mm Tube Adapter	15Nm	No	5mm 
MAUCH®			
Mauch Knee			
Knee Ball Attachment Screw	1Nm	No	2mm 
Mauch Knee Plus			
Knee Ball Attachment Screw	1Nm	No	2mm 
4-Hole Distal Attachment Bolts (SHCS)	13Nm	243	5mm 
4-Hole Distal Attachment Bolts (FHCS)	13Nm	243	4mm 
Mauch (Common Components)			
Set Screw (upper pivot pin to piston rod)	3Nm	243	2.5mm 
Distal Cylinder Studs	10Nm	243	6mm 
Mauch XG			
Knee Cap Attachment Screw	1Nm	No	2mm 
Flexion Stop Attachment Screws	3Nm	243	2.5mm 
Mauch Gaitmaster Systems			
Anterior Knee Ball Attachment Screw			
(Inner Split Shaft Attachment Screw)	0.7Nm	243	2.5mm 
Medial/Lateral Knee Ball Attachment Screws (Split Shaft)	N/A	243	2.5mm 
Medial/Lateral Knee Ball Attachment Screws (Solid Shaft)	N/A	243	1.5mm 
Distal Knee Ball Attachment Screws	N/A	243	2mm 
Knee Bracket Set Screws	1Nm	243	2mm 
Flexion Stop Attachment Screws	1Nm	243	2mm 
4-Hole Distal Attachment Bolts (SHCS)	13Nm	243	5mm 
4-Hole Distal Attachment Bolts (FHCS)	13Nm	243	4mm 
Distal Bearing Assembly Attachment Screws	3Nm	243	2.5mm 
Solid Shaft End Screw	0.7Nm	243	19mm 

General Information

Torque Values For Össur Products

		Torque	Loctite	Tool Required	
TOTAL SHOCKS™					
Total Shock 4400					
Pinch bolts		8Nm	No	5mm	●
Total Shock 4410					
Set screws		11Nm	410	6mm	●
FLEX-FOOT®					
Flex-Foot Assure					
Fibreglass tube kit attachment bolts		22Nm	410	5/16"	○
Male Pyramid adapter		22Nm	410	13mm	○
Vari-Flex / Talux					
Tube kit attachment bolts	Cat. 1-6	22Nm	410	1/2"	○
	Cat. 7-9	34Nm	410	9/16"	○
LP Vari-Flex					
Lamination Adapter Bolts	sizes 22-24	34Nm	410	5mm	●
	sizes 25-30	34Nm	410	6mm	●
Alignment fitting tool	Cat. 1-6	11Nm	410	3/8"	○
	Cat. 7-9	22Nm	410	1/2"	○
Re-Flex Rotate					
Clamp		3Nm	243	5mm Hex	●
Re-Flex Shock					
Clamp		3Nm	243	5mm Hex	●
Ceterus					
Fabric Clamp		9Nm	410	3/16"	●
Re-Flex VSP					
Fabric Clamp		5.5Nm	410	3/16"	●
Modular III / Cheetah / Flex-Sprint					
T-Connector Bolts	Cat. 1-6	33Nm	410	3/16"	●
	Cat. 7-9	45Nm	410	5/16"	●
Elation					
Attachment bolt	sizes 22-25	15Nm	242	5mm Hex	●
	sizes 26-28	21Nm	242	6mm Hex	●
Tube clamp bolt		11Nm	N/A	5mm Hex	●
Cheetah					
Lamination Connectors	Cat. 4-6	22Nm	410	1/2"	○
	Cat. 7-9	34Nm	410	9/16"	○
Flex-Run™ with Nike Sole					
Centre Bolt	Cat. 1-6	41Nm	410	17mm	○
	Cat. 7-9	43Nm	410	19mm	○
Carbon-X Active Heels					
all products	Cat. 1-6	14Nm	410	3/8"	○
	Cat. 7-9	24Nm	410	7/16"	○
Bolt-on Male Pyramid					
Modular III, Vari-Flex, Talux	Cat. 1-6	22Nm	410	13mm	○
Modular III, Vari-Flex, Talux	Cat. 7-9	34Nm	410	17mm	○
Junior Adapter					
Alignment set screws		4Nm	243	3mm Hex	●
4-Hole screws		6Nm	243	3mm Hex	●
Pinch bolt		6Nm	243	3mm Hex	●

General Information

Terms Of Sale And Policies

HOW TO ORDER

Össur part numbers and item descriptions must be included on all orders.

PRICES AND PAYMENT TERMS

All prices are subject to change without prior notification. All prices exclude V.A.T.

RETURNS, ADJUSTMENTS AND CREDITS

Please inspect all orders immediately upon receipt. In the event of items being ordered or shipped in error, please contact our Customer Care Department for a Return Authorisation number (RA#). Össur must be contacted within five business days of the receipt of goods. Prior approval and an RA# must be obtained for credit. The RA# must be printed on the outside of the package and a copy of the original invoice or packing slip must be included. All returns are subject to inspection.

TRIAL PERIOD

Össur provides a trial period with certain products. If the amputee is dissatisfied with any such Össur product for any reason, the prosthetist can return the product to Össur for a full refund. If a product needs to be returned within the Trial Period, the following must occur:

Prior approval and a Return Authorisation number (RA#) must be obtained from the Össur Customer Care Department for goods to be returned for credit. No items will be accepted for credit after the Trial Period has ended. (Trial Period begins the date the order is invoiced.)

- Any items that are received by Össur after the Trial Period has ended will be sent back to the customer at their cost and no credit will be issued.
- RA#s are valid for thirty days, starting the day the request was made. After thirty days, the RA# is invalid.
- Any items that are returned in addition to those listed and/or any items that are returned which differ from those listed on the RA# will be treated as an unauthorised return.
- The customer is responsible for the returned package prior to receipt at Össur.
- Returns must be accompanied by a copy of the original invoice or packing slip.

For Trial Periods, please refer to Warranty section in the Product Catalogue, page 208.

General Information

Warranty

All Össur products are guaranteed to be free of defects in materials and workmanship upon purchase. Any product found to be defective under normal and correct use will be replaced at no charge with the same item. If an item needs to be replaced under warranty, the customer should order this in the normal way. The replacement item will ship as soon as possible and the full retail cost will be invoiced. Upon receipt of the warranty item, a credit will be given under the terms that apply.

KNEES		LIMITED WARRANTY
RHEO® KNEE3		Variable Warranty Available
PROPRIO FOOT™		24 Months
POWER KNEE™		36 Months
SYMBIONIC® LEG		Variable Warranty Available

LINERS AND SLEEVES

Iceross® liners and Iceflex® sleeves can be returned for full credit within 30 days only if the sealed bag has not been opened. Any items that have been removed from this bag will not be accepted for credit and will be returned to the customer at the customer's expense.

LINERS		LIMITED WARRANTY	LINERS		LIMITED WARRANTY
Iceross Sport		9 Months	Iceross Seal-In X5 (TF)		6 Months
Iceross Synergy		9 Months	Iceross Seal-In Transfemoral		6 Months
Iceross Synergy Wave		9 Months	Iceross Transfemoral		6 Months
Iceross Dermo		6 Months	Iceross Upper-X		6 Months
Iceross Dermo Wave		6 Months	Iceross Junior		6 Months
Iceross Dermo Seal-In		6 Months	Distal Cup		12 Months
Iceross Comfort		6 Months	Iceross Sleeve		3 Months
Iceross Comfort Wave		6 Months	Iceflex Endurance		3 Months
Iceross Original		6 Months	Iceflex Balance		1 Month
Iceross Seal-In X5 (TT)		6 Months			
Iceross Seal-In X5 Wave (TT)		6 Months			
Iceross Dermo Cushion		6 Months			
Iceross Dermo Cushion Wave		6 Months			
Iceross Activa		9 Months			
Iceross Synergy Cushion		9 Months			
Iceross Synergy Cushion Wave		9 Months			
Iceross Comfort Cushion		6 Months			

General Information

Warranty

SOCKET SYSTEMS	LIMITED WARRANTY	MAUCH KNEE	LIMITED WARRANTY
Icecast Anatomy	12 Months	Mauch Knee	24 Months
		Mauch Knee Plus	24 Months
ICELOCK	LIMITED WARRANTY	TOTAL KNEE	LIMITED WARRANTY
Icelock 100 Series	12 Months	Total Knee 1900	24 Months
Icelock 200 Series	12 Months	Total Knee 2100	24 Months
Icelock 300 Series	12 Months	Total Knee 2000	24 Months
Icelock 500 Series	6 Months	Total Knee Junior	24 Months
Icelock 600 Series	12 Months		
Icelock 700 Series	12 Months		
STANDARD ADAPTERS	LIMITED WARRANTY	BALANCE KNEE	LIMITED WARRANTY
Össur Standard Adapters	24 Months	Balance Knee Locking	24 Months
Össur Junior Adapters	24 Months	Locking Knee	24 Months
Fixed Offset Adapters	12 Months		
Total Shock 4400 and 4410	12 Months		

FEET

Össur provides each prosthetic foot with a 60-Day Trial Period. If the amputee is dissatisfied with the prosthetic foot for any reason, the prosthetist can return it to Össur for a full refund within 60 days from the invoice date. For a detailed description of how to return products see Terms of Sale and Policies.

FEET	LIMITED WARRANTY	FEET	LIMITED WARRANTY
Flex-Foot Balance	24 Months	Elation	12 Months
Balance Foot J	36 Months	Flex-Symes	36 Months
Flex-Foot Assure	36 Months	Flex-Foot Junior	18 Months
Talux ST	36 Months	Chopart	24 Months
Vari-Flex	36 Months	Flex-Run with Nike Sole	12 Months
LP Vari-Flex	36 Months	Flex-Run Junior	12 Months
Vari-Flex Junior	36 Months	Flex-Sprint*	12 Months
Vari-Flex XC	36 Months	Cheetah*	12 Months
Vari-Flex XC Rotate	36 Months	Cheetah Xtreme with Nike Spike Pad*	12 Months
Re-Flex Rotate	36 Months	Cheetah Xtend with Nike Spike Pad*	12 Months
LP Rotate	36 Months	Cheetah Xplore Junior	12 Months
Vari-Flex Modular	36 Months	Cheetah Junior	12 Months
Re-Flex Shock	36 Months	Foot Cover	2 Months

* Excluded from the trial period guarantee, for information on trialling these products please contact your Sales Representative

General Information

Order Form - Custom Made Iceross

Date _____

Customer Data

Customer Name _____

Customer Number _____ Address _____

Contact Person _____ Postal code _____ City _____

User Data

Name _____

Circumference 4cm from distal end of residual limb (cm) _____

Product Information

Type of Liner/Silicone	Part#	Size	Profile	Umbrella	Matrix	Cover
☐ Seal-In X5 (TT) Special	I-290400	_____	☐ 3mm	☐ None	☐ Standard	☐ Standard ☐ Synergy ☐ Dermo ☐ Original ☐ Sport ☐ Comfort
☐ Dermo (Wave) Special	I-290400	_____	☐ 3mm	☐ Standard ☐ None	☐ None ☐ 10cm ☐ Other _____	☐ Standard ☐ Synergy ☐ Sport ☐ Comfort
☐ Dermo Seal-In Special	I-290450	_____	☐ 3mm	☐ None	☐ None ☐ 10cm ☐ Other _____	☐ Standard ☐ Synergy ☐ Sport ☐ Comfort
☐ Synergy (Wave) Special	I-290750	_____	☐ 3mm	☐ Standard ☐ None	☐ None ☐ 10cm ☐ Other _____	☐ Standard ☐ Dermo ☐ Sport ☐ Comfort
☐ Sport Special	I-290700	_____	☐ 3mm	☐ Standard ☐ None	☐ None ☐ 10cm ☐ Other _____	☐ Standard ☐ Synergy ☐ Dermo ☐ Comfort
☐ Original Special	I-290250	_____	☐ 2mm	☐ Standard ☐ None	☐ None ☐ 10cm ☐ Other _____	☐ Standard ☐ Synergy ☐ Dermo ☐ Comfort
☐ Comfort Special	I-290600	_____	☐ 3mm ☐ 6mm	☐ Standard ☐ None	☐ None ☐ 10cm ☐ Other _____	☐ Standard ☐ Synergy ☐ Sport ☐ Dermo
☐ Transfemoral Special	I-290800	_____	☐ 2mm	☐ Standard ☐ None	☐ None ☐ 10cm ☐ Other _____	☐ Standard ☐ Original ☐ Comfort
☐ Transfemoral Seal-In Special	I-290850	_____	☐ 2mm	☐ Standard ☐ None	☐ None ☐ 10cm ☐ Other _____	☐ Standard ☐ Original ☐ Comfort ☐ None
☐ Seal-In X5 (TF) Special	I-290850	_____	☐ 2mm	☐ None	☐ Standard	☐ Standard ☐ Original ☐ Comfort

Time of Delivery

Please contact your Customer Care representative for a confirmed time of delivery for your order.

General Information

Returns / Complaints Form

Please copy and complete this form and include it with any return. Please also apply to Customer Care for a Returns Authorisation number (RA#), as we will not be able to process credit claims otherwise. Please also put RA# on the outside of the box, since it will help us to process your return faster

Send returns to:

Össur UK Ltd.
Unit No 1,
S:Park,
Hamilton Road,
Stockport
SK1 2AE

Thank you for your co-operation,
Össur UK

Date

RA number

Customer data

Customer

Customer code Street

Technician/Prosthetist Name

Postal code City

Product code

Product code Description

Reason for reclaim

☐ on trial ☐ ordered by mistake ☐ sent by mistake ☐ defect or failure during guarantee period

Patient information

This information is needed for all credit claims for defective goods.

Without it we will not be able to process any guarantee claim.

Name or patient code

Weight..... Age Activity level

Description of defect or failure.....

.....

.....

Signature

